


**Gluco/Rx**

*Nexus*

TD-4277

**BLOOD GLUCOSE  
MONITORING SYSTEM**

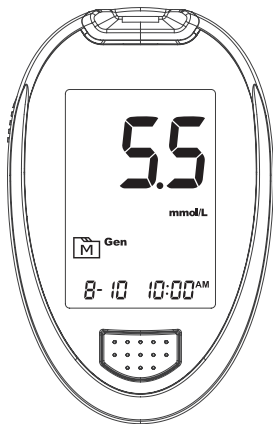

**OWNER'S  
MANUAL**



## Dear GlucoRx Nexus System Owner:

Thank you for purchasing the **GlucoRx Nexus** Blood Glucose Monitoring System. This manual provides important information to help you to use the system properly. Before using this product, please read the following contents thoroughly and carefully.

Regular monitoring of your blood glucose levels can help you and your doctor gain better control of your diabetes. Due to its compact size and easy operation, you can use the **GlucoRx Nexus** Blood Glucose Monitoring System to easily monitor your blood glucose levels by yourself anywhere, any time.

If you have other questions regarding this product, please contact your place of purchase or call the local customer service.

# **IMPORTANT SAFETY PRECAUTIONS**

## **READ BEFORE USE**

1. Use the device only for the intended use described in this manual.
2. Do not use accessories which are not supplied by the manufacturer.
3. Do NOT use the device if it is not working properly or if it is damaged.
4. Do NOT under any circumstances use on newborns or infants.
5. This device does NOT serve as a cure for any symptoms or diseases. The data measured is for reference only. Always consult your doctor to have the results interpreted.
6. Before using this device to test blood glucose, read all instructions thoroughly and practice the test. Carry out all quality control checks as directed.
7. Keep the device and testing supplies away from young children. Small items such as the battery cover, batteries, test strips, lancets, and vial caps are choking hazards.
8. Use of this instrument in a dry environment, especially if synthetic materials are present (synthetic clothing, carpets etc.) may give damaging static discharges that may cause erroneous results.
9. Do not use this instrument in close proximity to sources of strong electromagnetic radiation as these may interfere with its accurate operation.
10. Proper maintenance and periodically control solution test are essential to the longevity of your device. If you are concerned about your accuracy of measurement, please contact the local customer service or place of purchase for help.

**KEEP THESE INSTRUCTIONS IN A SAFE PLACE**

# TABLE OF CONTENTS

<b>IMPORTANT SAFETY PRECAUTIONS</b>	<b>4</b>
<b>BEFORE YOU BEGIN</b>	<b>7</b>
Important Information	7
About Alternative Site Testion (AST)	8
Intended Use	9
Test Principle	9
Contents of System	10
Meter Overview	11
Display Screen	12
Test Strip	13
<b>SETTING THE METER</b>	<b>14</b>
<b>THE FOUR MEASURING MODES</b>	<b>17</b>
<b>BEFORE TESTING</b>	<b>18</b>
Control Solution Testing	18
Performing a Control Solution Test	18
<b>TESTING WITH BLOOD SAMPLE</b>	<b>21</b>
Preparing the Lancing Device for Blood Testing	21
Preparing the Puncture Site	21
Performing a Blood Glucose Test	22
<b>METER MEMORY</b>	<b>25</b>
Reviewing Test Results	25
Reviewing Blood Glucose Day Average Results	26
<b>DOWNLOADING RESULTS ONTO A COMPUTER</b>	<b>27</b>

<b>MAINTENANCE</b>	<b>28</b>
Battery	28
Replacing the Battery	28
Caring for Your Meter	30
Caring for Your Test Strips	31
Caring for Your Carry Case	32
Important Control Solution Information	33
<b>SYSTEM TROUBLESHOOTING</b>	<b>34</b>
Result Readings	34
Error Message	35
Troubleshooting	36
<b>DETAILED INFORMATION</b>	<b>37</b>
<b>SYMBOL INFORMATION</b>	<b>37</b>
<b>SPECIFICATIONS</b>	<b>38</b>

# BEFORE YOU BEGIN

## Important Information

- ▶ Severe dehydration and excessive water loss may cause readings which are lower than actual values. If you believe you are suffering from severe dehydration, consult a healthcare professional immediately.
- ▶ If your blood glucose results are lower or higher than usual, and you do not have any symptoms of illness, first repeat the test. If you have symptoms or continue to get results which are higher or lower than usual, follow the treatment advice of your healthcare professional.
- ▶ Use only capillary whole blood sample to test your blood glucose. Using other substances will give incorrect results.
- ▶ If you are experiencing symptoms that are inconsistent with your blood glucose test results and you have followed all instructions described in this owner's manual, contact your healthcare professional.
- ▶ We do not recommend using this product on severely hypotensive individuals or patients in shock. Please consult your healthcare professional before use.
- ▶ The measurement unit used for indicating the concentration of blood or plasma glucose can either have a weight dimension (mg/dL) or a molarity (mmol/L). The approximate calculation rule for conversion of mg/dL in mmol/L is:

mg/dL	Divided by 18	= mmol/L
mmol/L	Times 18	= mg/dL

For example:

1)  $120 \text{ mg/dL} \div 18 = 6.6 \text{ mmol/L}$

2)  $7.2 \text{ mmol/L} \times 18 = 129 \text{ mg/dL}$  approximately.

## About Alternative Site Testing (AST)

**Important: There are limitations for doing AST. Please consult your healthcare professional before you do AST.**

### What is AST?


Alternative site testing (AST) means that people use parts of the body other than fingertips to check their blood glucose levels. This system allows you to test on the palm, the forearm, and the upper arm with the equivalent results to fingertip testing.

### What's the advantage?

Fingertips feel pain more readily because they are full of nerve endings (receptors). At other body sites, since nerve endings are not so condensed, you will not feel as much pain as at the fingertip.

### When to use AST?

Food, medication, illness, stress and exercise can affect blood glucose levels. Capillary blood at fingertip reflects these changes faster than capillary blood at other sites. Therefore when testing blood glucose during or immediately after meal, physical exercise, or any other events, **take blood sample from your finger only.**

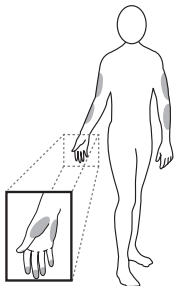

We strongly recommend you do AST **ONLY** in the following intervals:

- ▶ In a pre-meal or fasting state (more than 2 hours since the last meal).
- ▶ Two hours or more after taking insulin.
- ▶ Two hours or more after exercise.

Do **NOT** use AST if:

- ▶ You think your blood glucose is low.
- ▶ You are unaware of hypoglycaemia.
- ▶ Your AST results do not match the way you feel.
- ▶ You are testing for hyperglycaemia.
- ▶ Your routine glucose results are often fluctuating.



## Intended Use

This system is intended for use outside the body (*in vitro* diagnostic use) by people with diabetes at home and by health care professionals in clinical settings as an aid to monitor the effectiveness of diabetes control. It is intended to be used for the quantitative measurement of glucose (sugar) in fresh whole blood samples (from the finger, palm, forearm and upper arm).

It should not be used for the diagnosis of diabetes, or testing on newborns.

Professionals may test with capillary and venous whole blood; home use is limited to capillary whole blood testing.


## Test Principle

Your system measures the amount of sugar (glucose) in whole blood. The glucose testing is based on the measurement of electrical current generated by the reaction of glucose with the reagent of the strip. The meter measures the current, calculates the blood glucose level, and displays the result. The strength of the current produced by the reaction depends on the amount of glucose in the blood sample

## Contents of System

Your new GlucoRx Nexus system kit includes:

- |                          |                                   |
|--------------------------|-----------------------------------|
| 1 Meter                  | 7 Lancing Device with AST cap     |
| 2 Owner's Manual         | 8 2 x 1.5V AAA Alkaline Batteries |
| 3 Warranty Card          | 9 Control Solution                |
| 4 Quick Start User Guide | 10 Test Strips                    |
| 5 Daily Log Book         | 11 Lancets                        |
| 6 Protective Wallet      |                                   |

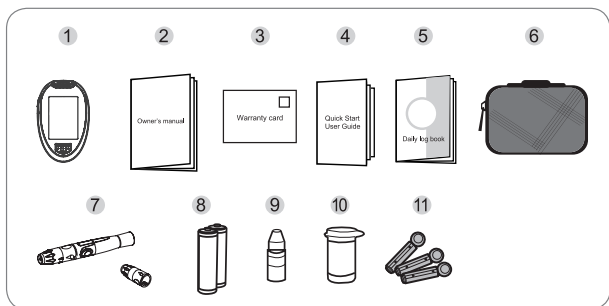


Test strips, control solutions, or sterile lancets may not be included in the kit (please check the contents on your product box). They can be purchased separately.

### NOTE

If any items are missing from your kit or opened prior to use, please contact your local customer services or place of purchase for assistance.

## Meter Overview

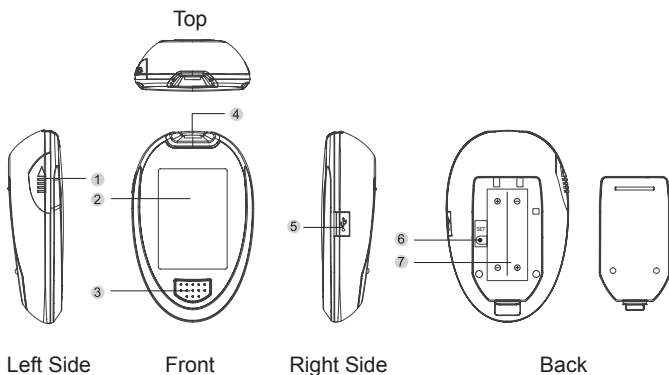

### 1 TEST STRIP EJECTOR

Eject the used strip by pushing this button up.

### 2 DISPLAY SCREEN

### 3 M BUTTON

Enter the meter memory and silence a reminder alarm.

### 4 TEST STRIP SLOT WITH STRIP INDICATION LIGHT

Insert test strip here to turn the meter on for testing.

### 5 DATA PORT


Download test results with the cable connection.

### 6 SET BUTTON

Enter and confirm the meter settings.

### 7 BATTERY COMPARTMENT

## Display Screen

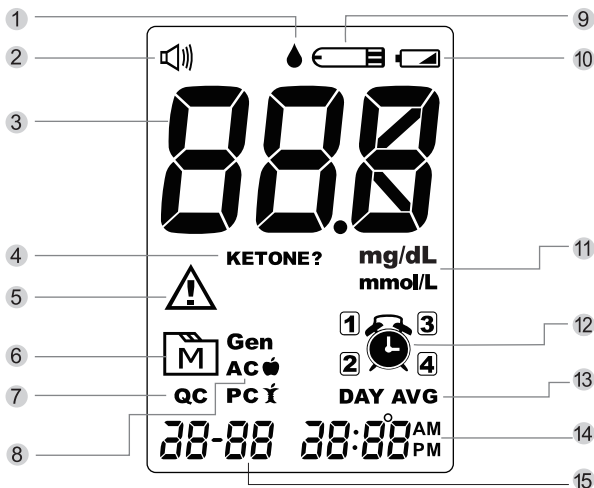

1 Blood Drop Symbol

2 Buzzer

3 Test Result

4 Ketone Warning

5 Error Message

6 Memory Mode Symbol

7 Control Solution Mode

8 Measuring Mode

9 Test Strip Symbol

10 Low Battery Symbol

11 Measurement Unit


12 Reminder Alarm

13 Day Average

14 Time


15 Date

## Test Strip

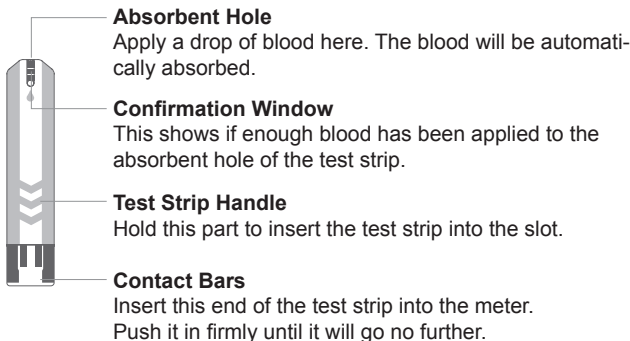


### Attention!

Test results might be wrong if the contact bar is not fully inserted into the test slot.

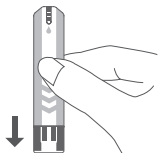

The front side of the test strip should face up when inserting the test strip.

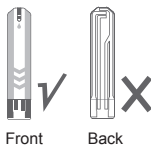

### NOTE

The GlucoRx Nexus meter should only be used with GlucoRx Nexus Test Strips. Using other test strips with this meter can produce inaccurate results.

# SETTING THE METER

Before using your meter for the first time or if you change the meter battery, you should check and update these settings. Make sure you complete the steps below and have your desired settings saved.

## To Enter the Setting Mode

Start with the meter off (no test strip inserted). Press SET to turn on the meter.



### 1. Setting the date

With the year flashing, press and release M button to choose the correct year. Press SET.



With the month flashing, press and release M button to choose the correct month. Press SET.




With the date flashing, press and release M button to choose the correct date. Press SET.



## 2. Setting the time format

Press M button to select the desired time format --- 12h or 24h.  
Press SET.

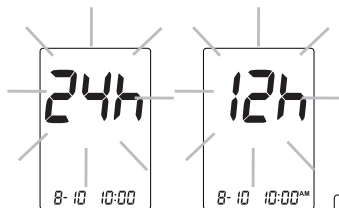

## 3. Setting the time

With the hour flashing, press and release M button to choose the correct hour. Press SET.




With the minute flashing, press and release M button to choose the correct minute. Press SET.




## 4. Setting the Buzzer

With the buzzer displays, press M to switch between On and Off. Press SET.

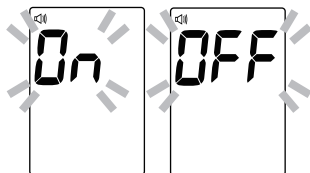

## 5. Setting the reminder alarm

You may set up any or all of the reminder alarms


(1-4). The meter displays “On” or “OFF” and  , press M button to turn on or turn off to set the first reminder alarm.

Press M button to select “On”, then press SET to set the hour. When the hour is flashing, press M button to add an hour. Press SET to confirm and go to minutes, press M button to add one minute. Hold M button longer to add faster. Press SET to confirm and go to the next alarm setting.



If you do not want to set an alarm, press SET to skip this step.

If you want to turn off an alarm, find the alarm number by pressing SET in the setting mode, press M button to change from “ON” to “OFF”.

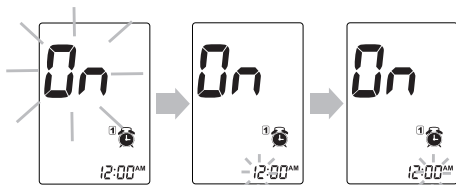

At the time of your alarm, the meter will beep and automatically turn on. You can press M button to silence the alarm and insert a test strip to begin testing. If you do not press M button, the meter will beep for 2 minutes then switch off. If you do not want to test at this time, press M button to switch off the meter.

**Congratulations! You have completed all settings!**

### NOTE

- These parameters can **ONLY be changed** in the setting mode.
- If the meter is idle for 3 minutes during the setting mode, it will switch off automatically.




# THE FOUR MEASURING MODES

The meter provides you with four modes for measuring General, AC, PC and QC.


Modes	Use when
General (displays as "Gen")	any time of day, regardless of time since last meal
AC	no food intake for at least 8 hours
PC	2 hours after a meal
QC	testing with the control solution

You can switch between each mode by:

1. Start with the meter switched off. Insert a test strip to turn on the meter, the screen will display flashing blood drop and "Gen".

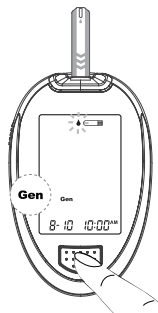

2. Press M button to switch between General, AC, PC and QC.

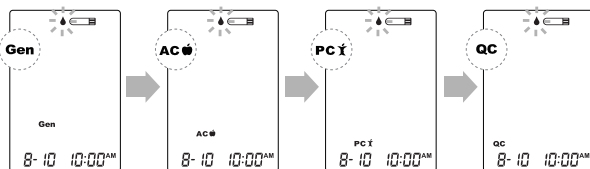

# BEFORE TESTING

## Control Solution Testing

TaiDoc Control Solution contains a known amount of glucose that reacts with test strips and is used to ensure your meter and test strips are working together correctly.


### Do a control solution test when:

- you first receive the meter,
- at least once a week to routinely check the meter and test strips,
- you begin using a new vial of test strips,
- you suspect the meter or test strips are not working properly,
- your blood glucose test results are not consistent with how you feel, or if you think the results are not accurate,
- practicing the testing process, or
- you have dropped or think you may have damaged the meter.

## Performing a Control Solution Test


### 1. Insert the Test Strip to Turn on the Meter

Insert the test strip into the meter. Wait for the meter to display the test strip and blood drop symbols.

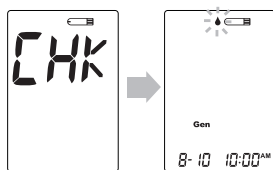

## 2. Press M button to mark this test as a control solution test.

With “QC” displayed, the meter will store your test result in memory under “QC”. If you press M button again, the “QC” will disappear and this test is no longer a control solution test.

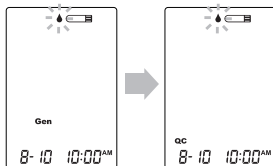


### WARNING

When doing the control solution test, you have to mark it so that the test result will NOT mix with the blood glucose TEST RESULTS stored in the memory. Failure to do so will mix up the blood glucose test results with the control solution test results in memory.

## 3. Apply Control Solution.

Shake the control solution vial thoroughly before use. Squeeze out a drop and wipe it off, then squeeze another drop and place it on the tip of the vial cap.

Hold the meter to move the absorbent hole of test strip to touch the drop. Once the confirmation window fills completely, the meter will begin counting down.

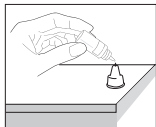


To avoid contaminating the control solution, do not directly apply control solution onto a strip.



#### 4. Read and Compare the Result

After counting down to 0, the test result of control solution will appear on the display. Compare this result with the range printed on the test strip vial and it should fall within this range. If not, please read instructions again and repeat the control solution test.

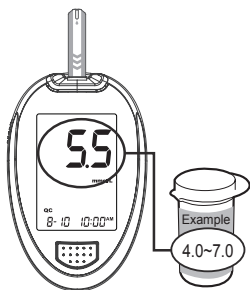

#### NOTE

Example only! Please read exact range from your test strip vial

#### Out-of-range results

If you continue to have test results fall outside the range printed on the test strip vial, the meter and strips may not be working properly. Do NOT test your blood. Contact your local customer service or place of purchase for help.

#### NOTE

- The control solution range printed on the test strip vial is for control solution use only. It is not a recommended range for your blood glucose level.
- See the **Maintenance** section for important information about your control solutions.

# TESTING WITH BLOOD SAMPLE

## NOTE

To reduce the chance of infection:

- Never share a lancet of the lancing device. GlucoRx Lancing device is for **self-use only**.
- Always use a new, sterile lancet. Lancets are for single use only.
- Avoid getting hand lotion, oils, dirt, or debris in or on the lancets and the lancing device.

## Preparing the Lancing Device for Blood Testing

Please follow the instructions in the lancing device insert for collecting a blood sample.

## Preparing the Puncture Site.

Stimulating blood perfusion by rubbing the puncture site before blood extraction has a significant influence on the glucose value obtained.

Blood from a site that has not been rubbed exhibits a measurably different glucose concentration than blood from the finger. When the puncture site was rubbed prior to blood extraction, the difference was significantly reduced.

**Please follow the suggestions below before obtaining a drop of blood:**

- \* Select the puncture site either at fingertips or another body parts (please see section “Alternative Site Testing” (AST) on how to select the appropriate sites).
- \* Hand wash the puncture site with warm water and dry thoroughly before testing.
- \* Rub the puncture site for about 20 seconds before penetration.
- \* Use a clear cap (included in the kit) while setting up the lancing device.


## ■ Fingertip testing

Press the lancing device's tip firmly against the side of your fingertip. Press the release button to prick your finger, then a click indicates that the puncture is complete.



## ■ Blood from sites other than the fingertip

Replace the lancing device cap with the clear cap for alternative site testing. Pull the cocking control back until it clicks. When lancing the forearm, upper arm and hand, avoid lancing the areas with obvious veins for excess bleeding.

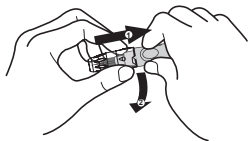




### NOTE

- Choose a different spot each time you test. Repeated punctures at the same spot may cause soreness and calluses.
- Please consult your health care professional before you begin AST.
- It is recommended that you discard the first drop of blood as it might contain tissue fluid, which may affect the test result.

## Performing a Blood Glucose Test

### 1. Insert the Test Strip to Turn on the Meter

Wait for the meter to display the test strip  and blood drop  symbols with Gen preset.

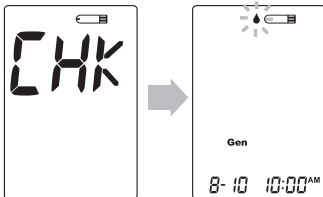

## 2. Select the appropriate measuring mode by pressing M button.

For selecting the measuring mode, please refer to the “FOUR MEASURING MODES” section.

## 3. Obtain a Blood Sample

Use the pre-set lancing device to puncture the desired site. After

penetration, discard the first drop of blood with a clean cotton swab. Gently squeeze the punctured area to obtain another drop of blood. Be careful NOT to smear the blood sample.




The volume of blood sample must be at least 0.5 microliter ( $\mu\text{L}$ ) of volume. ( . actual size).



## 4. Apply the Sample

Hold the blood drop to touch the absorbent hole of the test strip. Blood will be drawn in and after the confirmation window is completely filled, the meter begins counting down.

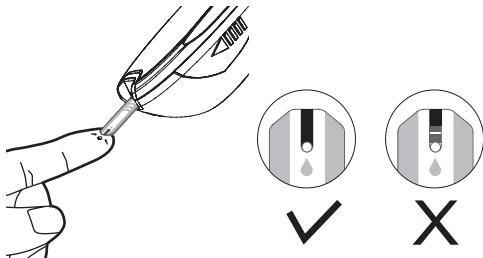

## NOTE

- Do not press the site against the test strip or try to smear the blood.
- If you do not apply a blood sample to the test strip within 3 minutes, the meter will automatically turn off. You must remove and reinsert the test strip to start a new test.
- The confirmation window should be filled with blood before the meter begins to count down. **NEVER** try to add more blood to the test strip after the drop of blood has moved away. **Discard the used test strip and retest with a new one.**
- If you have trouble filling the confirmation window, please contact your healthcare professional or the local customer service for assistance.

## 5. Read the Result


The result of your blood glucose test will appear after the meter counts to 0. This reading will automatically be stored in the memory.



## 6. Eject the Used Test Strip and Remove the Lancet

To eject the test strip, point the strip at a sharp disposal container. The meter will turn itself off automatically after the test strip is ejected.

**Always follow the instructions in the lancing device insert when removing the lancet.**

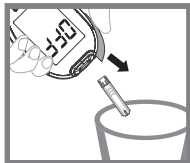

## WARNING

The used lancet and test strip may be biohazardous. Please discard them carefully according to your local regulations.




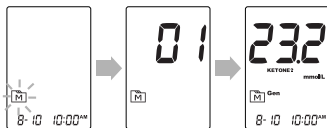
# METER MEMORY

The meter stores the 1000 most recent blood glucose test results along with respective dates and times in its memory. To enter the meter memory, **start with the meter off.**


## Reviewing Test Results

### 1. Press and release M button.

 will appear on the display. Press M button again, and the first reading you see is the last blood glucose result along with date, time and the measuring mode.




**2. Press M button** to recall the test results stored in the meter with each press.

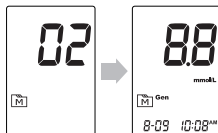


### 3. Exit the meter memory


After the last test result, press M button again and the meter will turn off.

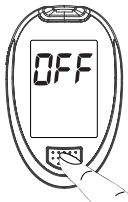

## Reviewing Blood Glucose Day Average Results

### 1. Press and release M button

When  appears on the display, keep pressing M button for 3 seconds until the flashing “DAY AVG” appears. Release M button and then your 7-day average result measured in general mode will appear on the display.

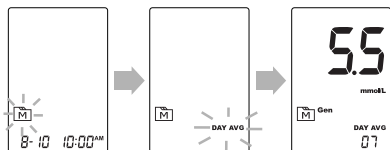

**2. Press M button to review** 14-, 21-, 28-, 60- and 90- day average results stored in each measuring mode in the order of Gen, AC, then PC.

### 3. Exit the meter memory

Keep pressing M button and the meter will turn off after displaying the last test result.



#### NOTE

- Any time you wish to exit the memory, keep pressing M button for 5 seconds or leave it without any action for 3 minutes. The meter will switch off automatically.
- Control solution results are **NOT** included in the day average.
- If using the meter for the first time, “---” displays when you recall the test results or review the average result. It indicates that there is no test result in the memory.

# DOWNLOADING RESULTS ONTO A COMPUTER


You can use the meter with an USB cable and the Health Care Software System to view your test results on your personal computer. To learn more about the Health Care Software System or to obtain an USB cable separately, please contact your local customer services or the place of purchase for assistance.

## 1. Obtaining the required cable and installing the software

For downloading Health Care Software System, please visit the GlucoRx website: <http://www.glucorx.co.uk>.

## 2. Connecting to a personal computer

Connect the cable to an USB port of your computer. With the meter turned off, connect the other end of the USB cable to the meter data port. "USB" will appear on the meter display, indicating that the meter is in communication mode.

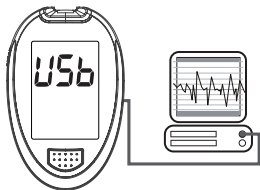

## 3. Data transmission

Follow the instructions provided in the software to transmit data. Results with date and time will be transmitted. Remove the cable and the meter will automatically turn off.

### WARNING

While the meter is connecting to the PC, it is unable to perform a blood glucose test.


# MAINTENANCE


## Battery


Your meter comes with two 1.5V AAA size alkaline batteries.

### ► Low Battery Signal

The meter will display one of the below messages to alert you when the meter power is getting low.

1. The  symbol appears along with display messages: The meter is functional and the result remains accurate, but it is time to change the batteries.

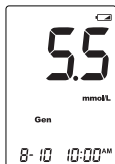

2. The  symbol appears with E-b, error message and low:


The power is not enough to do a test. Please change the batteries immediately.



## Replacing the Battery

To replace the battery, make sure that the meter is turned off.

1. Press the edge of the battery cover and lift it up to remove.
2. Remove the old batteries and replace with two 1.5V AAA alkaline batteries.
3. Close the battery cover. If the batteries are inserted correctly, you will hear a “beep” afterwards.

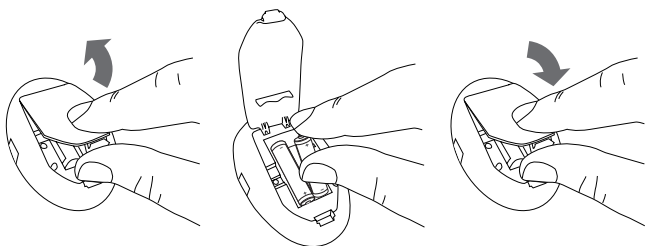

#### NOTE

- Replacing the batteries does not affect the test results stored in memory.
- As with all small batteries, these batteries should be kept away from small children. If swallowed, promptly seek medical assistance.
- Batteries might leak chemicals if unused for a long time. Remove the batteries if you are not going to use the device for an extended period (i.e., 3 months or more).
- Properly dispose of the batteries according to your local environmental regulations.

## Caring for Your Meter

To avoid the meter and test strips attracting dirt, dust or other contaminants, please wash and dry your hands thoroughly before use.

### ► Cleaning

1. To clean the meter exterior, wipe it with one disinfecting wipe (Example: Micro-Kill Plus™) to clean exposed surfaces thoroughly and remove any visible dirt, blood or any other body fluid with the wipe. Then dry the device with a soft, dry and clean cloth. Do NOT rinse the meter with water.
2. Do NOT use organic solvents to clean the meter.

### ► Meter Storage

- Storage condition: -20°C~60°C (-4°F~140°F) , below 95% relative humidity.
- Always store or transport the meter in its original storage case.
- Avoid dropping and heavy impact.
- Avoid direct sunlight and high humidity.

### ► Meter Disposal

- The used meter should be treated as contaminated that may carry a risk of infection during measurement. The batteries in this used meter should be removed and the meter should be disposed in accordance with local regulations.
- The meter falls outside the scope of the European Directive 2002/96/EC-Directive on waste electrical and electronic equipment (WEEE).



## Caring for Your Test Strips

- Storage conditions: 2°C~32°C (35.6°F~89.6°F), below 85% relative humidity. Do not freeze.
- Store your test strips in their original vial only. Do not transfer to another container.
- Store test strip packages in a cool and dry place. Keep away from direct sunlight and heat.
- After removing a test strip from the vial, immediately close the vial cap tightly.
- Touch the test strip with clean and dry hands.
- Use each test strip immediately after removing it from the vial.
- Write the opening date on the vial label when you first open it. Discard remaining test strips after 6 months.
- Do not use test strips beyond the expiry date. This may cause inaccurate results.
- Do not bend, cut, or alter a test strip in any way.
- Keep the strip vial away from children since the cap and the test strip may be a choking hazard. If swallowed, promptly see a doctor for help.

For further information, please refer to the test strip package insert.

## Caring for Your Carry Case

1. Remove protective film attached to the sticky sheet before using for first time.

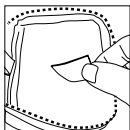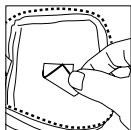


2. If the stickiness of the sticky sheet eventually worsens when too dirty, please remove the sticky sheet lightly.




### NOTE

To prevent splitting, do not forcefully pull the sticky sheet out.

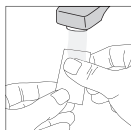

3. Rinse the sticky sheet with tap water, gently rubbing the surface with your fingers to remove dirt.

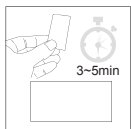


4. After rinsing with water, allow the sticky sheet to air dry for 3 to 5 minutes.



### NOTE

Never use any object to wipe off the sticky sheet.



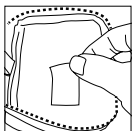

5. After drying completely, the sticky sheet can be restored. Place it back in the carry case for use.



#### NOTE

The sticky sheet will become less sticky the more you touch it, through transfer of your skin oils.

## Important Control Solution Information

- Use only **TaiDoc** control solutions with your meter.
- Do not use the control solution beyond the expiry date or 3 months after first opening. Write the opening date on the control solution vial and discard the remaining solution after 3 months.
- It is recommended that the control solution test be done at room temperature (20°C-25°C / 68°F-77°F). Make sure your control solution, meter, and test strips are at this specified temperature range before testing.
- Shake the vial before use, discard the first drop of control solution, and wipe off the dispenser tip to ensure a pure sample and an accurate result.
- Store the control solution tightly closed at temperatures between 2°C and 30°C (35.6°F and 86°F). Do NOT freeze.







# SYSTEM TROUBLESHOOTING

If you follow the recommended action but the problem persists, or error messages other than the ones below appear, please call your local customer service. Do not attempt to repair by yourself and never try to disassemble the meter under any circumstances.

## Result Readings

Appears	When Glucose
<b>Lo</b>	<1.1mmol/L (20mg/dL)
<b>KETONE?</b>	$\geq$ 15mmol/L (270mmol/L)
<b>H ,</b>	>33.3mmol/L (600mg/dL)
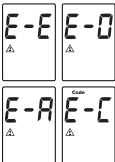
## Error Messages

MESSAGE	WHAT IT MEANS	WHAT TO DO
	Appears when the batteries cannot provide enough power for a test.	Replace the batteries immediately.
	Appears when a used test strip is inserted.	Repeat with a new test strip.
	Problem in operation.	Repeat the test with a new test strip. If the meter still does not work, please contact customer services for assistance.
	Appears when test strip is removed while counting down, or insufficient blood volume.	Review the instructions and repeat test with a new test strip. If the problem persists, please contact your local customer services for help.
	Appears when ambient temperature is below system operation range.	System operation range is 10°C to 40°C (50°F to 104°F). Repeat the test after the meter and test strip are in the above temperature range.
	Appears when ambient temperature is above system operation range.	

## Troubleshooting

1. If the meter does not display a message after inserting a test strip:

POSSIBLE CAUSE	WHAT TO DO
Batteries exhausted.	Replace the batteries.
Test strip inserted upside down or incompletely.	Insert the test strip with contact bars end first and facing up.
Defective meter or test strips.	Please contact customer services.

2. If the test does not start after applying the sample:

POSSIBLE CAUSE	WHAT TO DO
Insufficient blood sample.	Repeat the test using a new test strip with larger volume of blood sample.
Defective test strip.	Repeat the test with a new test strip.
Sample applied after automatic switch-off (3 minutes after last user action).	Repeat the test with a new test strip. Apply sample only when flashing “●” appears on the display.
Defective meter.	Please contact customer services.

3. If the control solution testing result is out of range.

POSSIBLE CAUSE	WHAT TO DO
Error in performing the test.	Read instructions thoroughly and repeat the test again.
Control solution vial was poorly shaken.	Shake the control solution vigorously and repeat the test again.
Expired or contaminated control solution.	Check the expiry date of the control solution.
Control solution that is too warm or too cold.	Control solution, meter, and test strips should be at room temperature (20°C-25°C / 68°F-77°F) before testing.
Defective test strip.	Repeat the test with a new test strip.
Meter malfunction.	Please contact customer services.

## DETAILED INFORMATION

The meter provides you with plasma equivalent results.








Time of day	Normal blood glucose range for people <b>with</b> diabetes (mmol/L)
Fasting and before meal	4 to 7 mmol/L (72 to 126 mg/dL)
2 hours after meals	< 9 mmol/L (162 mg/dL)

Source: Diabetes UK. Blood glucose target. Balance: No.234, 2010 April. P.69

**Please work with your doctor to determine a target range that works best for you.**

\*1: American Diabetes Association position statement on the Diabetes Control and Complications Trial (1993).

## SYMBOL INFORMATION

Symbol	Referent
	<i>In vitro</i> diagnostic medical device
	Consult instructions for use
	Temperature limitation
	Use by
	Batch code
	Manufacturer
<b>SN</b>	Serial number
	Authorised representative in the European Community

# SPECIFICATIONS

Model No.: TD-4277

Dimension & Weight: 96 (L) x 61 (W) x 26 (H) mm, 67.2 g

Power source: Two 1.5V AAA alkaline batteries

Display: LCD

Memory: 1000 measurement results with respective date and time

External output: USB Cable

Auto electrode inserting detection

Auto sample loading detection

Auto reaction time count-down

Auto shutdown after 3 minutes of idleness

Temperature Warning

Operating condition: 10°C to 40°C (50°F to 104°F), below 85% R.H.  
(no condensing)

Meter Storage/Transportation condition: -20°C to 60°C (-4°F to 140°F), 95%  
R.H.

Strip Storage/Transportation Conditions: 2°C to 32°C (35.6°F to 89.6°F), below  
85% R.H.

Measurement units: fixed mmol/L or mg/dL

Measurement range: 1.1 to 33.3 mmol/L (20 to 600 mg/dL)

Expected Service Life: 5 years

This device has been tested to meet the electrical and safety requirements of:  
IEC/EN 61010-1, IEC/EN 61010-2-101, EN 61326-1, IEC/EN 61326-2-6





TaiDoc Technology Corporation  
B1-7F., No. 127, Wugong 2nd Rd., Wugu Dist.,  
New Taipei City 248, Taiwan (R.O.C.)  
[www.taidoc.com](http://www.taidoc.com)



MedNet GmbH  
Borkstraße 10, 48163 Münster, Germany

Distributed by  
GlucorX Ltd.  
4 Wintonlea, Monument Way West  
Woking, Surrey, GU21 5EN, UK  
Website: [www.glucorx.co.uk](http://www.glucorx.co.uk)  
Call centre number: 01483 755133



0123



For self-testing.