
# Setting the region number

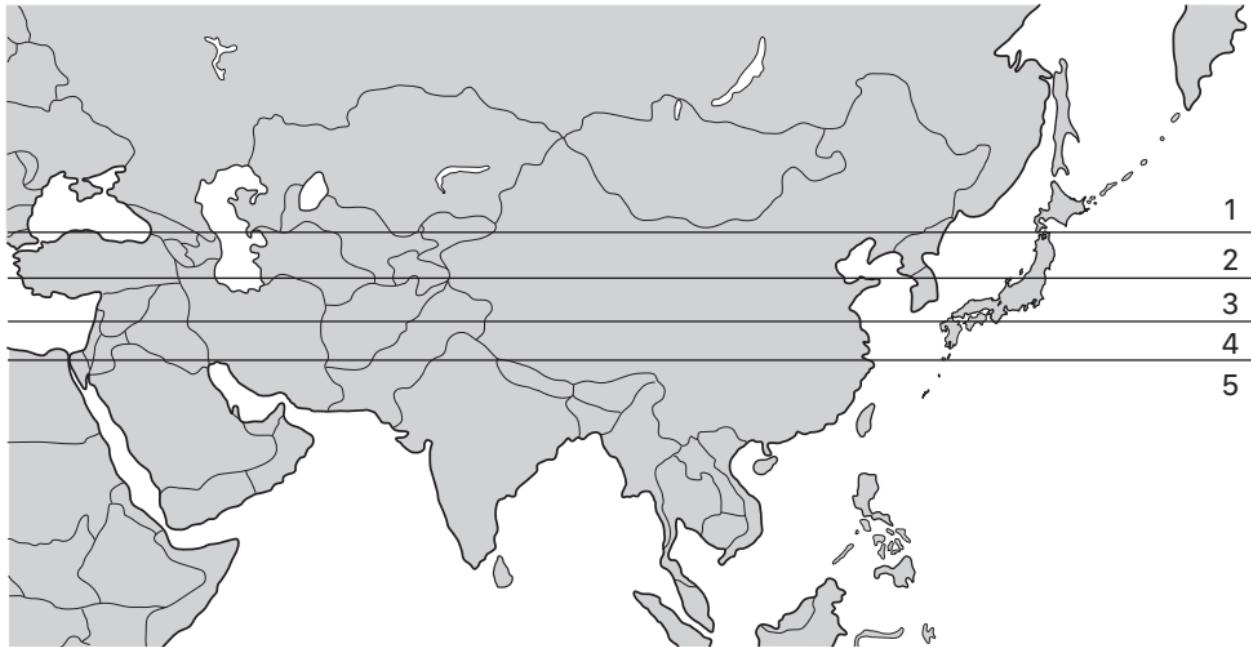
Your geographical location, longitude and latitude, affects your gravity, so you need to select the appropriate region number to get an accurate measurement.

Because this product is designed for use in Japan, specific regions are not available for every country. This can affect your measurements by up to 100g either way.

Please continue to follow the Instruction Manuals(English, Simp-Chinese, and Trad-Chinese) next page. >>>

# E C L E A R

## Body composition analyzer Instruction Manual (Quick guide) HCS-WFS01 Series

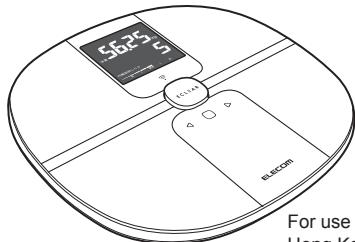For use in Singapore and Hong Kong only.

If you made your purchase in another country, please refer to manual from the relevant country.  
<https://www.elecom.co.jp/global/download-list/>

## Precautions

### DANGER


Never use this product with the following medical electronic devices.

- (1) Pacemakers and other implantable medical electronic devices
  - (2) Heart-lung machines and other artificial life-support medical electronic devices
  - (3) Electrocardiographs and other wearable medical electronic devices
- Since this unit will allow a weak electric current to flow through the body, the above devices may malfunction and cause serious damage or accidents.


## How to use \* The region must be configured during the first start-up and battery replacement.

### ■ Registering a user

- ① Make sure that the power is off, and then press the (on/setting) button. A user-number selection screen is displayed.

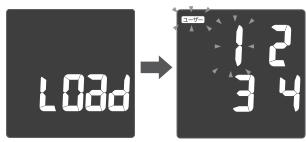- ② Select a user number (the flashing number), and then press the (on/setting) button. The user number is entered.

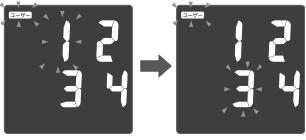- ③ Display the age, and then press the (on/setting) button. The age is entered.

- ④ Select M or F, and then press the (on/setting) button. The gender is entered.

- ⑤ Display the height of the registrant, and then press the (on/setting) button. The height is entered, and the screen changes to the setting confirmation screen.

- Height can be configured in units of 0.5 cm.  
⑥ Using the data that was entered, "gender" flashes, "Age" and "Height" are displayed in order, and then the weight management screen is displayed after.

\* If no operations are performed at the setting screens for 30 seconds, the power switches off automatically.  
\* If the power is switched off during configuration, perform the procedure again from the beginning.

### ■ Checking measurement results

- ① Visceral fat level    ② Difference with previous

- ③ BMI

- ④ Body fat %

- ⑤ Skeletal muscle %

- ⑥ Bone quantity

- ⑦ Basal metabolism

\* Measurement ranges have been established for each measurement item. If the measurement range was exceeded, L or H is displayed.

5

## Initialize Wi-Fi settings

- ① Connect your smart phone to the Wi-Fi environment.  
Please use a router that supports the 802.11b/g/n standard. 802.11a dedicated routers cannot be used.

- ② Access the "ECLEAR APP" download page and install "ECLEAR APP".

- ③ Activate the installed "ECLEAR APP" and select this product (HCS-WFS series) from "Device settings".

- ④ According to the application screen, enter the password of the Wi-Fi environment.


- ⑤ Make sure that the power is not turned on and press the (on/setting) button to turn on the body composition indicator.  
⑥ On the user selection screen, press and hold the (off/communication) button until the wireless communication display on the display section of the body composition indicator flashes.  
⑦ Tap "ECLEAR APP's" "Next".

- ⑧ When the "Pi" sound goes off, "yes" is displayed on the screen.

Wi-Fi connection of this product is completed.

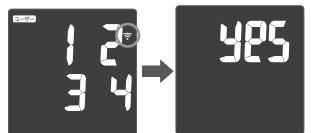
When the connection is completed, the wireless communication display on the display section of the body composition indicator changes from blinking to lit.

7

### ■ Measure after specifying the user number

- ① Press the (on/setting) button. A user-number selection screen is displayed.

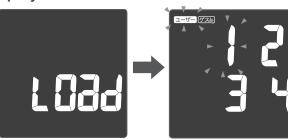- ② Select a user number (the flashing number), and then press the (on/setting) button. The selected user is entered, and then "0.00kg" is displayed.


- ③ Carefully step on this machine. Measurement is started. Stand still, and remain upright until the measurement completes.

### ■ Measure using Guest mode

- ① Press the (on/setting) button. A user-number selection screen is displayed.

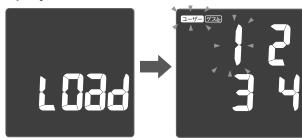- ② Press or button to select (Guest), then press (on/setting) button. The analyzer enters Guest mode.

- ③ Display the age of the guest, and then press the (on/setting) button. The age is entered.

- ④ Select M or F, and then press the (on/setting) button. The gender is entered.

- ⑤ Display the height of the registrant, and then press the (on/setting) button. The height is entered, and the screen changes to the setting confirmation screen.

Height can be configured in units of 0.5 cm.

- ⑥ Using the data that was entered, "gender" flashes, "Age" and "Height" are displayed in order, and then the weight management screen is displayed after. "0.00kg" is displayed.

- ⑦ Carefully step on the analyzer to start measurement.

6

## Reading the standards

### ■ Understanding the standards of the visceral fat level

Level	Evaluation	Understanding the evaluation
1~10	Normal	Your level of stored visceral fat is not high. Continue making a habit of balanced meals and moderate exercise.
11~14	High ① - Slightly high -	Your level of stored visceral fat is slightly high. Review your current lifestyle habits, and make a habit of balanced meals and moderate exercise.
15~18	High ② - High -	Your level of stored visceral fat is high. Firstly, be careful not to eat and drink too much. Try to eat meals that leave you 70 to 80% full. It is best to attempt exercise after your body has started to feel lighter.
19~30	High ③ - Very high -	Your level of stored visceral fat is very high. Review your lifestyle habits. In particular, don't eat and drink too much. Switch to meals that leave you 60 to 70% full. Consult with a doctor if you have any concerns.

## Troubleshooting

Problem	Check points
The display section's  lights up	<input type="checkbox"/> Switch or exchange the batteries.
(error) is displayed	
(error) is displayed	<input type="checkbox"/> Did you step on the analyzer during automatic correction ()? <input type="checkbox"/> Did you move during measurement? How was your standing position or posture? <input type="checkbox"/> Are the bottoms of your feet dirty? <input type="checkbox"/> Did you try to measure with socks or other footwear still on?
(error) is displayed during measurement	
(error) is displayed	<input type="checkbox"/> Are the router standards compatible with 802.11b/g/n? <input type="checkbox"/> Once, turn off the power and reset the Wi-Fi settings.

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Complies with  
IMDA standards  
DA106058.

Operating Frequencies: 2412-2472MHz for 11b/g/n(HT20)  
Maximum field intensity/RF Output power: 5.5mW

\* ELECOM CO., LTD. owns the copyright of this manual.  
\* Copying/reprinting any part, or the whole of this manual without written permission from ELECOM is prohibited.  
\* The specification or appearance of this product might be changed without prior notice to improve the product quality.  
\* To export the part of this product applicable to strategic goods or services, export or service trading permission according to foreign currency laws is necessary.  
\* ELECOM is not responsible for any damages to, or failure of other devices caused by using this product.

A customer who purchases outside Japan should contact the local retailer in the country of purchase for enquiries. In "ELECOM CO., LTD. (Japan)", no customer support is available for enquiries about purchases or usage in/from any countries other than Japan. Also, no language other than Japanese is available. Replacements will be made under stipulation of the Elecom warranty, but are not available from outside of Japan.

Clear Body Composition Analyzer  
Manual, Safety precautions, and Compliance with regulations  
1st edition, November 28, 2019  
**ELECOM**  
©2019 ELECOM Co., Ltd. All Rights Reserved.

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ELECOM Sales Hong Kong Ltd.
ELECOM Singapore Pte. Ltd.