



HOFFRICHTER GmbH
Mettenheimer Straße 12/14
19061 Schwerin
Germany
Telephone: +49 385 39925-0
Fax: +49 385 39925-25
E-mail: info@hoffrichter.de
www.hoffrichter.de

LAVI
Ventilator

HOFFRICHTER
SERVONA RESPIRATORY CARE



LAVI user's manual for physicians and medical professionals

User's manual

for physicians and medical professionals

Applies to device types 9LV401 and 9LV402 from device software 1.400





User's manual

for physicians and medical professionals

©2019 HOFFRICHTER GmbH
All rights reserved.

The content published in this user's manual is the sole property of HOFFRICHTER GmbH. All forms of reproduction, editing, distribution and use of any kind, even in part, require the prior written consent of HOFFRICHTER GmbH.

HOFFRICHTER GmbH reserves the right to amend or replace this user's manual without prior notice.

Please ensure that you are always working with the most current version of this user's manual. Should you have any questions, please contact the ventilation device provider, or check our information at www.hoffrichter.de

The following documents are available in addition to this user's manual:

- LAVI user's manual for patients
- LAVI quick brief instruktion
- Service manual
- Hygiene concept
- Maintenance schedule

Licensing information

The device software is in part based on freeware. You can save and read the list of software used as well as the corresponding licensing conditions by copying the data to an SD card.

Every HOFFRICHTER GmbH device is supplied with a serial number for traceability purposes.

Please enter your device’s serial number here. You will find the serial number on the rating plate on the bottom of the device.

SN

Please always provide the serial number in case of queries and complaints.

CE₀₁₂₃

CE mark and number from the notified body. The medical device complies with the applicable regulations of EU 93/42/EEC for medical products.

Table of contents

Chapter 1	11
Introduction	11
Scope of delivery	12
Symbols	14
Intended purpose	18
User qualification	20
Chapter 2	21
Safety Information	21
General safety instructions	22
Electrical safety	24
Installation requirements and transport	25
Instructions before commissioning	26
Using oxygen	26
Integration into IT networks	27
Safety-related test	27
Chapter 3	29
Description of Device	29
Front side of device	30
Left device side	31
Rear of device	32
Top of device	34
Bottom of device	36
Chapter 4	41
Commissioning	41
General information	42
Setting up the device	42
Power supply	43
System setup, non-invasive ventilation	49
System setup, invasive ventilation	53
Calibrating the tube circuit	56
Using oxygen	62
Using the functional bag	68
Switching on the device	69

Switching off the device	69
Start ventilation.....	69
Stop ventilation	69
Chapter 5	71
Ventilation Modes	71
PCV mode	72
APCV mode.....	74
PSV mode	76
PSV-S mode	78
P-SIMV mode	80
CPAP mode	82
Chapter 6	83
Operating the Device	83
Menu structure	85
Menu lock.....	86
User profiles	87
Home screen	88
Monitoring screen	90
Parameter screen.....	96
Event log screen	109
System screen	111
Statistics screen	115
Service screen	118
Chapter 7	119
Alarms and Messages.....	119
General information	120
Audible alarm output (audio alarms)	122
Visible alarm output	123
Alarm log.....	125
Forwarding alarms.....	125
Alarm overview	126
Messages	131

Chapter 8	133
Cleaning and Disinfection	133
Important information	134
Overview	134
Cleaning the device	136
Cleaning the mask	136
Cleaning the tube circuit	137
Cleaning the exhalation valve	137
Cleaning/exchanging the air filter	138
Using the device for more than one patient	140
 Chapter 9	 143
Routine Tests and Maintenance Work	143
Overview	144
Safety-related test (SRT)	145
Battery maintenance	145
 Chapter 10	 147
Annex	147
Notes for home therapy	148
Device technical specifications	148
Data management	150
Error messages	151
Technical data	153
Measured values	156
Setting ranges and control precision	157
Factory settings	159
Compliance with standards	162
Replacement parts and accessories	163
Internal battery	168
Manufacturer's declaration on electromagnetic compatibility	174
Disposal	178
Disclaimer	179

List of figures

Figure 1:	Rating plate.....	15
Figure 2:	Front side of device	30
Figure 3:	Left device side.....	31
Figure 4:	Rear of device	32
Figure 5:	Top of device.....	34
Figure 6:	Bottom of device.....	36
Figure 7:	Setting up the device	42
Figure 8:	Assembling the power supply unit and power supply unit holder	43
Figure 9:	Mains connection via power supply unit.....	44
Figure 10:	Start screen	45
Figure 11:	System setup, non-invasive ventilation with vented mask.....	49
Figure 12:	System setup, non-invasive ventilation with non-vented mask.....	49
Figure 13:	Connecting AquaTREND uni.....	51
Figure 14:	Disconnecting AquaTREND uni.....	52
Figure 15:	System setup, invasive ventilation with external humidifier.....	54
Figure 16:	System setup, invasive ventilation with HME filter.....	55
Figure 17:	Connecting the bacterial filter	58
Figure 18:	Connecting the alarm box/nurse call	59
Figure 19:	Inserting the SD card	60
Figure 20:	Removing SD card	61
Figure 21:	Connecting the oxygen source.....	62
Figure 22:	Connecting the FiO ₂ sensor.....	64
Figure 23:	Functional bag.....	68
Figure 24:	PCV mode diagram.....	73
Figure 25:	APCV mode diagram	75
Figure 26:	PSV mode diagram	77
Figure 27:	PSV-S mode diagram	79
Figure 28:	P-SIMV mode diagram	81
Figure 29:	CPAP mode diagram	82
Figure 30:	Home screen	88
Figure 31:	Measurement screen, factory setting.....	90
Figure 32:	Graph screen.....	92
Figure 33:	Parameter screen.....	96
Figure 34:	Event log screen (alarms)	109
Figure 35:	Event log screen (events).....	110
Figure 36:	System screen.....	111
Figure 37:	Statistics screen	115
Figure 38:	Alarm displays on the toolbar	122
Figure 39:	Alarm output on the toolbar	123
Figure 40:	Alarm output in the textbox.....	124

Figure 41: Alarm box 125

Figure 42: Messages on the toolbar 131

Figure 43: Filter cassette structure 138

Figure 44: Block diagram for the device..... 149

Figure 45: Pneumatic block diagram..... 150

Figure 46: Error list 152

Figure 47: Lithium-ion battery RRC2054 169

Figure 48: Exchanging the battery 172

Chapter 1

Introduction

Scope of delivery

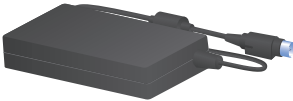
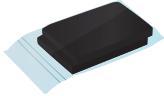
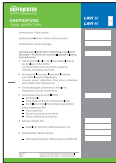
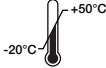
Figure	Name
	LAVI ventilator
	Switched-mode power supply FSP090-RACM1
	Power supply unit holder with integrated strain relief
	Mains cable
	Internal battery RRC2054 (Only LAVI 40)
	Leakage tube circuit (L = 180 cm, Ø 22 mm)
	Straight O ₂ connection adapter
	Carrying case
	Spare filter cassette, complete (open) with filters

Figure	Name
	Spare coarse filter, 1 pack (2 units)
	Spare fine filter, 1 pack (5 units)
	User's manual for patients
	Brief Instruction
	Final inspection report

Symbols

Symbols on the packaging

Symbol	Meaning
	European Article Number
	Article number
	Serial number
	CE mark and number from the notified body. The medical device complies with the applicable regulations of EU 93/42/EEC for medical products.
	Transport and store package with arrows pointing up at all times.
	Fragile contents
	Protect from moisture!
	Humidity range during storage and transport
	Air pressure range during storage and transport
	Temperature range during storage and transport
 CAUTION LITHIUM ION RECHARGEABLE BATTERIES INSIDE (No lithium metal) Do not damage or puncture this package. If package is damaged, batteries must be quarantined, inspected and repaired. For additional emergency information call +49 365 38925-0 VORSICHT ENTHÄLT WIEDERLADEBARE LITHIUM-IONEN-BATTERIEN (Kein Lithium-Metall) Vermeidung beschädigen des Verpackungsmaterials. Die Batterien müssen bei Beschädigung isoliert, untersucht und repariert werden. Für weitere Informationen für den Notfall: +49 365 38925-0  CAUTION! Lithium-Ion Batteries DO NOT LOAD OR TRANSPORT PACKAGE IF DAMAGED For more information, call +49 365 38925-0	CAUTION! Device contains lithium-ion batteries

Symbols on the rating plate

The rating plate is on the bottom of the device.

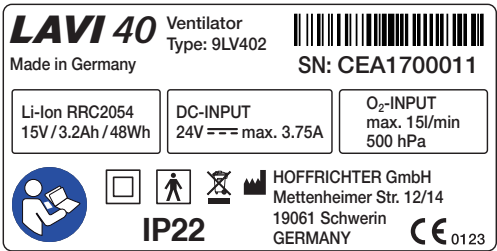


Figure 1: Rating plate

Symbol	Meaning
	Follow the instructions in the user’s manual.
	Protection class II (protective insulation)
	BF application part
IP22	Protection against: • solid foreign objects with diameters from 12.5 mm • access to hazardous parts with a finger • falling/dripping water, as long as the housing is tilted up to 15°
SN	Serial number
CE ₀₁₂₃	CE mark and number from the notified body. The medical device complies with the applicable regulations of EU 93/42/EEC for medical products.
	Manufacturer
	Do not dispose of the device with the household waste. Please contact the relevant customer services department to find out how to properly dispose of the device.

Symbols on the device

Symbol	Meaning
Housing	
	Follow the instructions in the user's manual.
Connections	
	FiO ₂ sensor connection
	DC connection
	COM interface
	Connection of remote alarm/nurse call (ESD-sensitive component – do not touch!)
	USB interface
	Oxygen connection
Control	
	Release button for the integrable humidifier
	On/off button

Symbols used in this user's manual

Important information is denoted by symbols in this user's manual. Please ensure that you follow these instructions in order to prevent accidents, personal injury and material damage.

In addition, the local accident prevention regulations and general safety regulations in force in the area of use must be observed.

DANGER

This symbol denotes hazardous situations that lead to serious injuries or death.

WARNING

This symbol denotes hazardous situations that may lead to serious injuries or death.

CAUTION

This symbol denotes hazardous situations that may lead to moderate injuries.

ATTENTION

This symbol denotes situations that may lead to material damage or damage to the device.

Please note:

Notes provide tips and information for the efficient, correct use of the device.

Intended purpose

The LAVI ventilator device may only be used for non-life-supporting respiration. It serves as intermittent respiration support as well as to provide respiration to patients with insufficient spontaneous breathing ability. The device is suitable for adults and children from a tidal volume of 100 ml and higher and is designed for home care and/or use in professional health care facilities. It must not be used for intensive ventilation.

LAVI is not intended for use in vehicles, planes or helicopters.

Description of function

The blower sucks in the ambient air via a filter and transports it to the patient at a set pressure via a leakage tube circuit (with passive exhalation valve). Ventilation can be invasive (e.g. using a tracheostoma) or non-invasive (using a mask).

Ventilation is carried out according to ventilation parameters set using the control elements. Ventilation can be monitored based on measurements and graphs on the display.

In case the custom and permanent alarm parameters are violated, visual and acoustic alarms are emitted. If the prescribed ventilation is no longer possible due to a technical error, an acoustic alarm is emitted for at least 2 minutes.

If LAVI is operated with an internal battery, ventilation can continue on without disruption in the event of a power failure.

LAVI can be connected to a low-pressure source of oxygen for ventilation with increased oxygen concentration. When oxygen is introduced but ventilation is not started, the oxygen supply is interrupted by a safety valve. Remaining oxygen can be released from the device via the oxygen outlet.

It is also possible to combine LAVI with the integrable humidifier "AquaTREND uni" or an external humidifier.

Therapy and statistical data, including alarms and events, can be copied onto an SD card and analysed using the PC software "easySET".

For service purposes, LAVI features one USB and one COM interface.

Indication

LAVI can be used for the following indications:

- Obstructive ventilation disorders (e.g. COPD)
- Restrictive ventilation disorders (e.g. scoliosis, thorax deformities)
- Neurological, muscular and neuromuscular disorders (e.g. paralysis of the diaphragm)
- Central breathing regulation disorders

Regardless of the indications named here, use of the device always depends on the doctor's individual diagnosis.

Contraindications



WARNING

Risk of injury due to contraindications!

Ventilation may be contraindicated for certain pre-existing conditions.

The following conditions may be a contraindication for non-invasive ventilation:

- Severe cardiac arrhythmia
- Severe hypotension
- Severe epistaxis
- Pneumothorax or pneumomediastinum
- Pneumoencephalus
- Cranial trauma
- Status after cranial or brain surgery
- Acute inflammation of the paranasal sinuses, middle ear infection or a perforated ear drum
- Aspiration hazard

In individual cases, the attending physician must decide on the therapy.

Side effects

The following undesired side effects may occur in connection with artificial respiration:

Invasive ventilation:

- Complications due to tube / tracheal cannula
- Gastric inflation

Mask ventilation:

- Pressure points and skin defects in the face
- Eye irritation due to leaks
- Gastric inflation
- Aspiration
- Sinusitis
- Nose bleeds

General complications of mechanical ventilation:

- Pulmonary barotrauma / volutrauma caused by ventilation
- Ventilator-associated pneumonia
- Effects on the cardio-vascular system

User qualification

The respiration and alarm parameters may only be set by trained specialist personnel under the supervision of a doctor. Care staff and patients have been instructed on how to use and handle the device. These people must be familiarised with operation of the device and must have read all of and understood this user's manual before commissioning the device. In addition, the operator must inform the users of what accessories are compatible with the device.

Maintenance and repairs may only be performed by trained and authorised service companies.

Chapter 2

Safety Information

Please note: *Heed all important information in this user's manual. Otherwise, there is a risk of accidents, injury and material damage.*

General safety instructions



WARNING

Risk of infection due to germs!

Hygienically preparing and cleaning the device must be performed according to this user's manual and the applicable regulations of the hospital or nursing home.



CAUTION

Risk of injury due to incorrect device settings!

- Only qualified, trained, specialist medical staff under the supervision of a physician may make adjustments to the ventilator.
- Please be sure to check the ventilation and alarm parameter settings after all servicing work.

Risk of injury due to incorrect accessories!

Manufacturer tested and approved accessories are recommended for the device. If other accessories are used, this may lead to insufficient ventilation or the use of hazardous materials may lead to further, secondary complications.

- Please read first this user's manual carefully and in its entirety before first using the ventilator.
- Keep the instructions in close proximity to the device for immediate reference if necessary.
- In case of problems, unexpected events or unusual device behaviour (e.g. during commissioning, use or maintenance), inform the device operator immediately and document the incidents. You will generally find the operator contact data on the device as well as in the medical device book.
- The device must only be used by persons who have fully read and understood this user's manual before commissioning and have familiarised themselves with the device. Disregarding these instructions can lead to life-threatening situations for the patient.
- In cases of emergency, an alternative ventilation option, such as a second ventilator or an emergency ventilation bag, must be available at all times and for use by the attending person.
- The device must only be used on the responsibility and prescription of the physician.
- The device must only be used on patients whose clinical record indicate it.
- LAVI is not intended for use in vehicles, planes or helicopters.
- Please take the care to ensure that the patient remains connected to the tubing circuit during ventilation.
- The device must not be used with flammable anaesthetics or ambient air that contains explosive gases. This may cause fires or explosions.
- Before being used again on another patient, all parts that come into contact with respiratory gas must be treated hygienically.

- Equipment that is not part of the ventilation system must not be connected.
- In order to ensure patient safety, the device must be operated in such a way that all adjustable alarms are activated and adjusted to the patient.
- Alarms must not be ignored. They indicate conditions that require an immediate action.
- Safety-related testing and maintenance are required every two years for the ventilator.
- In case of excessive agitation on the part of the patient, there is a risk of hyper-ventilation in all ventilation modes with inspiration triggering.
- The device must not be steam-sterilised in an autoclave.
- Filters and other parts that are connected to the device must be regularly replaced. Please dispose of the used parts according to the regulations for used medical material and/or the local environmental protection rules.
- Please ensure that the total resistance of the ventilation system does not exceed 6 hPa with a flow of 60 l/min for adults and 30 l/min for children.
- Any modification to the device poses a threat to its reliability and is thus not permitted.
- Masks may only be used on the prescription of a physician and after training by qualified medical staff.
- Only use the mask after instruction by qualified medical staff. Clarify in particular the intake of medicines and possible contraindications and side effects associated with the use of the prescribed mask.
- Please note the operating, transport and storage conditions.
- Temperatures lower than + 5°C and higher than + 40°C can impair the function of the device.
- During operation, the power supply unit can reach a surface temperature of up to 57°C. For this reason, do not touch the power supply unit for more than 1 minute to prevent burns to the skin.
- Keep small parts of the respiratory therapy system out of the reach of children and animals.



WARNING

Risk of injury due to electric shock!

- Do not try to open the device. Maintenance and repairs may only be performed by personnel authorised by HOFFRICHTER GmbH.
- Do not touch live parts of the mains cable or power supply unit if defective.
 - ⇒ Replace a defective mains cable/power supply unit.
- The device must not be used in wet rooms, as humidity penetrating the device presents a risk of electric shock.

Risk of injury due to interrupted operation!

- The device must always be located at least 30 cm away from other devices or equipment such as defibrillators, diathermy units, mobile phones, microwaves, remote controlled toys etc. Electromagnetic fields that exceed 10 V/m may adversely affect the operation of the ventilator.
- During certain examinations or treatments, mutual interference between the ventilator and other medical devices may occur. Please observe the information regarding electromagnetic compatibility and monitor the devices with regard to error-free and proper operation.
- The use of accessories or power supplies we have not approved for the ventilator can increase the emission of electromagnetic radiation, reduce interference immunity or lead to an increased patient leakage current.

- Only the supplied power supply unit may be used for operating the ventilator.
- Only LAVI and a humidifier may be connected to a power strip. Additional power strips or extension cables must not be plugged in to this power strip.
- Portable power strips the LAVI/an external humidifier are connected to must not be placed on the floor.
- Do not exceed the permitted maximum load of the power strip. Refer to the user's manual for the maximum power consumption of the LAVI/humidifier.
- Respiratory therapy may be contraindicated for certain pre-existing conditions.
- The contacts for connecting the remote alarm/nurse call and the RS232 interface must not be touched at the same time as the patient to prevent current from being discharged via the patient.
- Only accessories approved by HOFFRICHTER that are not connected to the power supply may be connected to the remote alarm/nurse call and the RS232 interface.

- To disconnect the device from the mains supply, unplug it.
- Before cleaning the device, the plug must be disconnected from the electrical outlet.
- Do not reach for the device under any circumstances if it falls into water.

Installation requirements and transport



CAUTION

Risk of injury due to the device falling down!

- ⇒ For operation, the device must be placed on a safe and level base.

Risk of injury due to unclean or insufficient air supply!

- ⇒ Please ensure the device is operated in an area where there is sufficient and clean ambient air.

Risk of injury due to overheated ventilation air!

The device must not be operated under any environmental conditions other than those stipulated. Excessive ambient temperatures can result in an increased ventilation air temperature.

- ⇒ Please note the operating, transport and storage conditions.
-

- The air inlet at the rear of the device, as well as all ventilation slots, must not be blocked.
- The display and the info LEDs must not be covered and must be visible to the user at all times.
- No objects must be placed on the device.
- The system must never be stored or transported at ambient temperatures under - 20°C and over + 50°C (with battery).
- The device must not be exposed to direct sunlight.
- Due to possible electromagnetic interference the ventilation device must not be placed directly next to other devices in which the electromagnetic radiation is not CE compliant and/or the limits values are exceeded (see page 175). If this is unavoidable, then ventilator operation must be monitored for trouble-free and correct operation.
- Do not place the device near water containers (baths).

Instructions before commissioning

- A malfunctioning device can endanger the patient or operator. Should the appliance not start properly, or if the device's automatic self-tests should fail, you must stop operating the device. In such cases, the service provider must be informed.
- Position the device so the mains plug is easily accessible and can be unplugged quickly in the event of a potential hazard.
- Do not use the device if the housing or the cable of the device or the power supply are damaged.

Using oxygen



CAUTION

Risk of injury due to increased oxygen supply!

The oxygen supplied must not exceed a pressure of 500 hPa and a flow of 15 l/min. The oxygen must be dosed using an external flow meter.

- Observe the user's manual of the manufacturer or distributor from whom you obtain your oxygen.
- If the patient is supplied with oxygen via the device, the FiO_2 should be measured.
- On the LAVI, FiO_2 measurement is possible using the optional FiO_2 sensor. We recommend using this particular sensor exclusively.
- To prevent incorrect calibration of the FiO_2 sensor, make sure fresh air is supplied when operating the device.
- The FiO_2 sensor contains a caustic liquid. Avoid skin or eye contact if there is a sensor leak! Replace the sensor.
- Keep the FiO_2 sensor out of the reach of children and animals.
- When supplying oxygen via the O_2 connection on the device, ensure that only dry gas is used. Moisture may lead to device defects. If necessary, a humidifier can be connected between the air outlet of the device and the patient.
- The connection between the O_2 connection and the external O_2 source must be absolutely airtight. Otherwise, leakage losses may occur during ventilation.
- The oxygen supply should be stopped before ventilation is interrupted. We also recommend allowing the device to run for several respiratory cycles without oxygen supply after stopping ventilation.
- In the event of an oxygen leak, the oxygen supply should be closed off immediately. The room must be ventilated immediately. Any sparks, fire or potential fire sources in the vicinity of the device must be avoided.

- Oxygen supports combustion. Therefore, observe the fire protection regulations applicable for using oxygen. Please ensure that the oxygen fittings, as well as all ports and surfaces near the oxygen lines are free of grease. Do not smoke and do not handle naked flames. When using oxygen, an increased oxygen concentration in the ambient air can occur.

Integration into IT networks

- Integrating the device into an IT network including other devices can pose unknown risks. The user is therefore responsible for assessing and evaluating the risks as well as for risk minimisation.
- Changes to the IT network can result in new, unknown risks. This includes changes to the network configuration, integrating and removing elements as well as updating or upgrading devices on the IT network.

Safety-related test

- In order to ensure the operating safety of the device, a safety-related test or maintenance must be carried out at the prescribed intervals.

Chapter 3

Description of Device

Front side of device

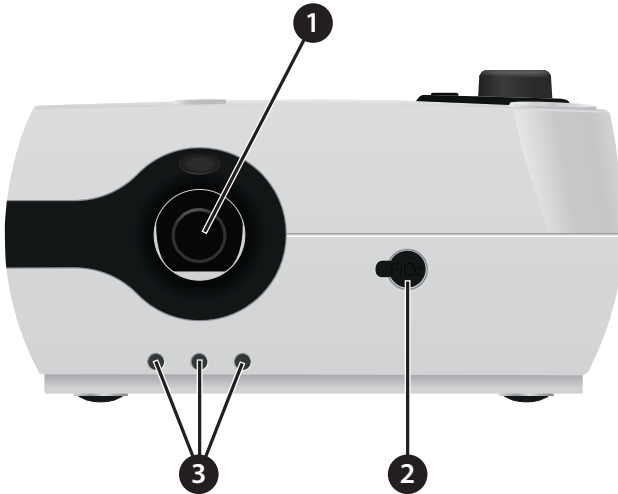


Figure 2: Front side of device

- 1 Leakage tube system/plug-in humidifier connection
The leakage tube circuit/humidifier AquaTREND uni is connected here.
- 2 Connection FiO₂ sensor cable
Connect the FiO₂ sensor cable here for measuring the oxygen concentration.
See also page 64.
- 3 Contact sockets for the plug-in humidifier AquaTREND uni

Left device side

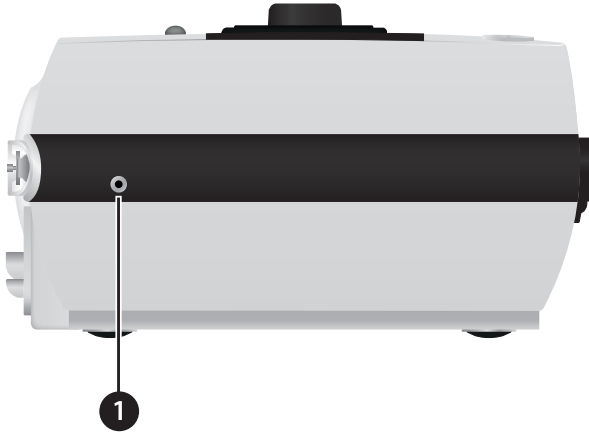


Figure 3: Left device side

- 1 Oxygen output**
This is the exit for excess oxygen from the oxygen valve of the unit when ventilation has been turned off.

Rear of device

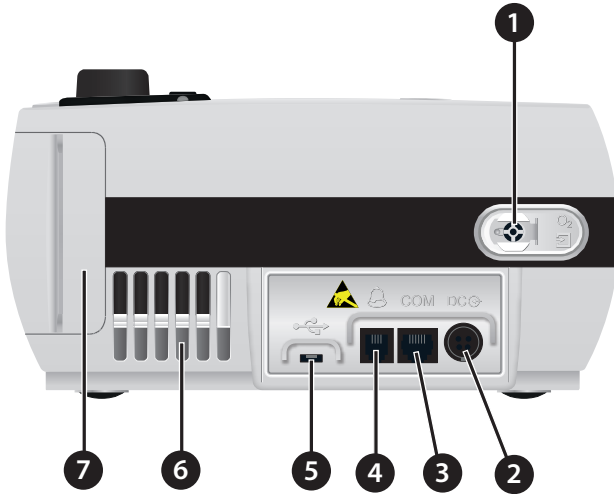


Figure 4: Rear of device

- 1 Oxygen connection
During oxygen input the oxygen source is connected here. Use the supplied oxygen connection adaptor for this purpose. See also page 62.
- 2 DC connection
The DC plug is connected here. See also page 43.
- 3 RS232 port (service interface)
- 4 Remote alarm/nurse call connection
An alarm box (optional accessory) or a nurse call system may be connected here. See also page 59.
Please note:
 - Only electrically isolated nurse call systems may be connected.
 - In rare cases, touching the contacts can cause false alarms due to electrostatic charge. Do not touch the contacts on the interface!

- 5 **Micro-USB port (PC connection/service interface)**
Here, you can connect a PC using a Micro-USB 2.0 cable. In order to be able to communicate with the device, the PC software "EASYset" must be installed on the PC.
Please note: Only devices that meet standards IEC 60601-1 and IEC 60950-1 may be connected.
- 6 **Ventilation slits**
The ambient air is sucked in here.
- 7 **Filter cassette**
The filter cassette contains the two air filters (coarse and fine filter). For information on how to replace and clean the filter, refer to page 138.

Top of device



Figure 5: Top of device

- 1 Display
- 2 Info LED
- 3 The info LED lights up/flashes in the event of an error and provides information on the alarm priority. It also indicates the operating status of the device.

Colour	State	Priority/status
Red	Flashing	HIGH
Yellow	Flashing	MEDIUM
Turquoise	Lit	LOW
Green	Lit	Ventilation is running
White	Flashing	Device is booting
	Lit	Device for operation

4 **Softkeys**
The softkeys can be assigned various functions. These functions are then shown on the display. See also page 37.

5 **Multifunction knob MFK**
The multi-function knob is used for menu navigation, parameter configuration and other operating actions.

6 **On/off button**

Function	Action
Start ventilation	Press briefly
Stop ventilation	Press briefly and confirm with MFK
Switching on the device	Press for > 4 s
Switching off the device	Press for > 4 s

7 **Release button for integrable humidifier AquaTREND uni.** See also page 52.

Bottom of device



Figure 6: Bottom of device

- 1 SD card slot
An SD card here can be inserted here.
- 2 Rating plate
- 3 Battery compartment cover
The exchangeable internal battery is located under the cover. See also page 47 and page 168.

Explanation and functions of the softkeys

Symbol	Meaning															
	Alarm key															
	<table><tr><th>Function</th><th>Condition</th><th>Action</th></tr><tr><td>Confirmation of all current alarms</td><td>Active alarms</td><td>Press briefly</td></tr><tr><td>Confirm no longer active alarms</td><td>Saved alarms</td><td>Press briefly</td></tr><tr><td>Mute the audible alarm for 2 min (audio alarm pause)</td><td></td><td>Press briefly</td></tr><tr><td>Cancel the audible alarm suppression</td><td>Audio alarm paused</td><td>Press briefly</td></tr></table>	Function	Condition	Action	Confirmation of all current alarms	Active alarms	Press briefly	Confirm no longer active alarms	Saved alarms	Press briefly	Mute the audible alarm for 2 min (audio alarm pause)		Press briefly	Cancel the audible alarm suppression	Audio alarm paused	Press briefly
	Function	Condition	Action													
	Confirmation of all current alarms	Active alarms	Press briefly													
	Confirm no longer active alarms	Saved alarms	Press briefly													
Mute the audible alarm for 2 min (audio alarm pause)		Press briefly														
Cancel the audible alarm suppression	Audio alarm paused	Press briefly														
	Escape key															
	<table><tr><th>Function</th><th>Action</th></tr><tr><td>Exit current screen</td><td>Press briefly</td></tr><tr><td>Exit selected parameter</td><td>Press briefly</td></tr><tr><td>Cancel action</td><td>Press briefly</td></tr></table>	Function	Action	Exit current screen	Press briefly	Exit selected parameter	Press briefly	Cancel action	Press briefly							
	Function	Action														
	Exit current screen	Press briefly														
Exit selected parameter	Press briefly															
Cancel action	Press briefly															
	Heater key If an integrable humidifier is connected, the heating key is available:															
	• On the home screen and • In the first level when the home screen icon is pressed on the home screen.															
	<table><tr><th>Function</th><th>Action</th></tr><tr><td>Switch on/off heating for the integrable humidifier</td><td>Press briefly</td></tr></table>	Function	Action	Switch on/off heating for the integrable humidifier	Press briefly											
Function	Action															
Switch on/off heating for the integrable humidifier	Press briefly															
	Home key															
	<table><tr><th>Function</th><th>Action</th></tr><tr><td>Return to the home screen</td><td>Press briefly</td></tr></table>	Function	Action	Return to the home screen	Press briefly											
Function	Action															
Return to the home screen	Press briefly															
	Error key															
	<table><tr><th>Function</th><th>Action</th></tr><tr><td>Error list display</td><td>Press briefly</td></tr></table>	Function	Action	Error list display	Press briefly											
Function	Action															
Error list display	Press briefly															

Explanation of the toolbar icons

Symbol	Meaning
	Clinic mode active
	Home mode active
	Alarm active  Red: high priority  Yellow: middle priority  Turquoise: low priority
	Alarm inactive  Grey: inactive, regardless of the continuation of the alarm condition
	Audio alarm paused The audible alarm has been paused for 2 min. The audible alarm of a new alarm event will also be paused for 2 min. The audible alarm may be deactivated by pressing the alarm key before an alarm event occurs. Pressing the key again reactivates the audible alarm in case an alarm event has occurred.
	Counter for paused audio alarm Indicates how much longer the audible alarm will be paused.
	Error detected On the home screen, press the error key to display the error message(s). You can find a list of all potential errors on page 151.
	PC is connected via the Micro-USB port, connection to “easySET” available
	Humidifier  Heating on  Heating in standby  Heating off  Heating deactivated due to battery operation or a fault in mains operation
	FiO ₂ sensor connected

Symbol	Meaning
	FiO ₂ sensor connected, but not calibrated
	SD card is inserted into the device
	Menu lock active The home and escape keys are deactivated and the MFK functions are restricted. It is not possible to access the menu.
	Internal battery charging
	Power level of the internal battery
	green: about 60% ... 100%
	yellow: about 20% ... 60%
	red: about 0% ... 20%
	Please note: <i>If no percentage is shown, the battery requires maintenance (see page 145).</i>
	Replace internal battery - Defective battery or - Battery capacity too low or - Battery incompatible
	Mains operation active

Explanation of symbols on the pressure bar

Symbol	Meaning
	Setting of trigger lock "on"
	Trigger lock triggered
	Spontaneous respiration detected The device has detected spontaneous breathing by the patient. This triggered the inspiration trigger. The symbol will remain visible during inspiration and will shut off with the beginning of the expiration.
	Back-up frequency is active The device is operating in PSV mode. The patient is not breathing spontaneously and is ventilated at the set frequency.

Chapter 4

Commissioning

General information

- Before commissioning the device, read the safety information from page 21.
- Before commissioning the ventilation system (ventilator, tube, humidifier etc.), check all connections for leaks, as well as the stability of the connected accessories.
- Never operate the device without the air filter.
- Only use original HOFFRICHTER filters.
- If the device was previously in an environment where the air temperature was not the same as in the new operating location, allow approximately 1 hour until the temperatures have evened out before commissioning.

Setting up the device

Place the device on a flat and stable surface. Make sure that the device is placed securely and that the air inlet at the rear of the device is not blocked. Make sure that the user is able to see the display and the info LED during ventilation. The device is designed for operation within arm's reach.

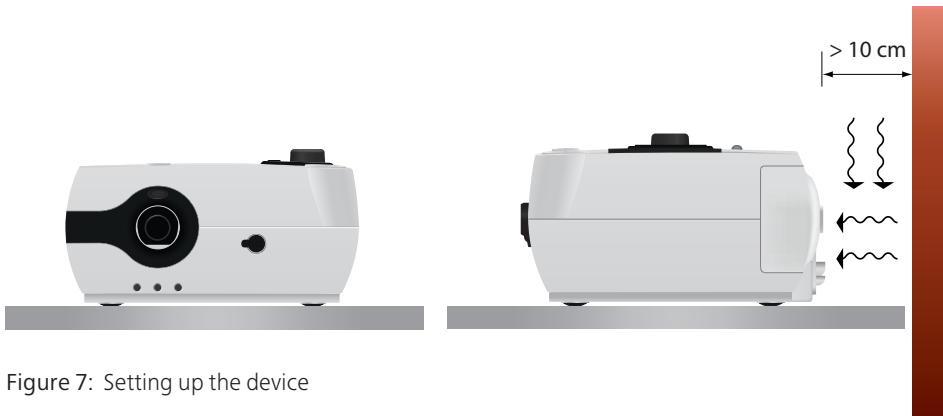


Figure 7: Setting up the device

Power supply

The ventilator can be supplied by two different power sources:


- Mains connection via power supply unit
- Internal battery

The ventilator automatically detects which power sources are available. If the device is connected to a power supply unit, the power supply unit is used as the primary source and the internal battery as the secondary source.

Mains operation

Mains operation means that the device is supplied with power by the power supply unit. The device can remain connected to the mains continuously without posing a risk. If inserted in the device, the internal battery is also charged.


Connect the device to the mains supply as follows:

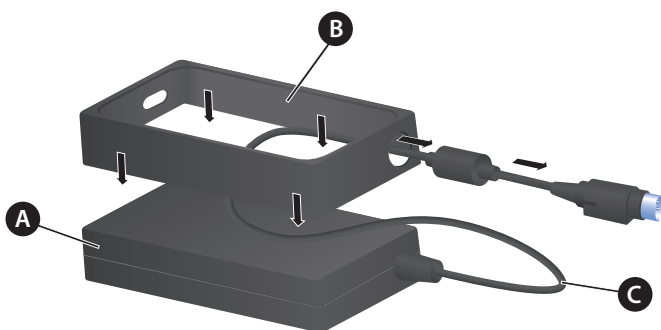

A Power supply unit B Power supply unit holder with integrated strain relief
C Power supply unit cable

Figure 8: Assembling the power supply unit and power supply unit holder

1. Pull the power supply unit cable through the round opening on the power supply unit holder.
 2. Press the power supply unit into the power supply unit holder until you feel it click into place.
-

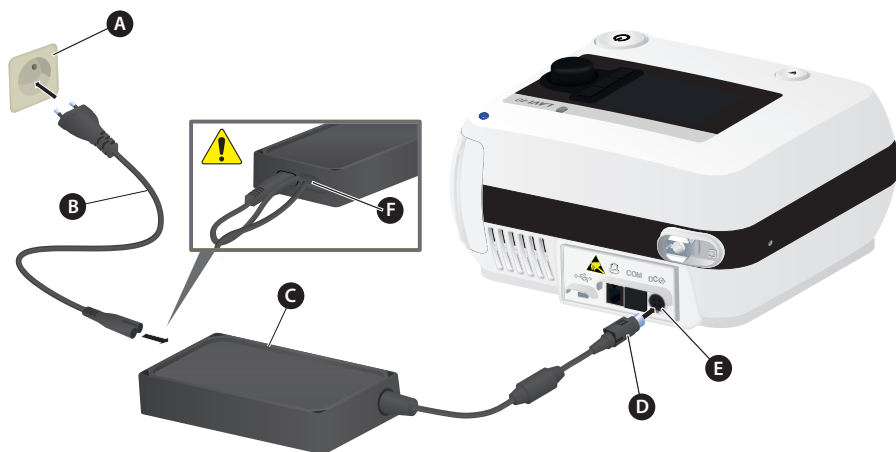


A Socket B Power cable C Power supply unit D DC plug with snap lock
E DC connector socket F Strain relief

Figure 9: Mains connection via power supply unit

3. Connect the mains cable to the power supply and hook the mains cable behind the strain relief to secure it against accidental unplugging.
4. Insert the DC plug into the DC connector socket.
Please note: The DC plug is a POWER-DIN plug with snap lock. Do not pull on the cable to disconnect the cable from the device—instead, grasp the plug and disconnect it by pulling it in a straight line away from the device.
5. Insert the mains cable plug into the power socket (100 - 240 V, 50/60 Hz).
6. Press the on/off button for more than 4 s to switch on the LAVI.
7. The device boots while performing a self-test:
 - Testing the primary and secondary alarm sounds: Both alarm sounds give a short beep one-by-one.
 - Verifying alarm LED: the alarm LED lights up white.
 - Checking of other hardware components

The booting progress is displayed by a progress bar (A). The status (B) is also displayed:

- Booting... → Device is booting
- ⚠ Error message → Error detected, press MFK to continue. All error messages are listed on page 151. The device cannot be started if certain errors are in effect. In this case, contact your service provider.

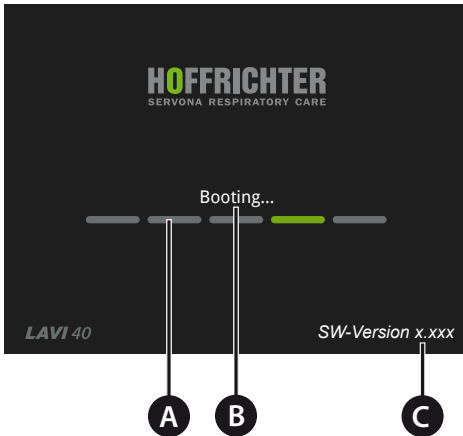

A Progress bar B Status C Software version

Figure 10: Start screen

If no errors were detected during the self-test or the errors have been confirmed, the display will switch to the home screen.

It takes up to 1 minute to boot the device, after which the device will be ready for operation.




Standby operation

If device mains operation is switched on and respiration is switched off, the device can be placed into standby operation. This means that the device automatically moves into a sleep mode after the last operation. A standby screen will be displayed. If the device has an internal battery, this will be charged and the power level will be shown.

The time at which the device will switch to standby operation can be modified using the “Standby” parameter on the system screen. The default factory setting is 5 min. If you do not want the device to switch to standby operation, set the parameter to “off”.

To exit standby operation, press a key or the MFK.

In the event of an alarm, standby stops and the home screen is displayed.

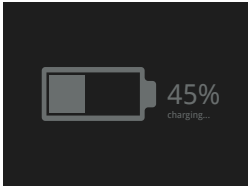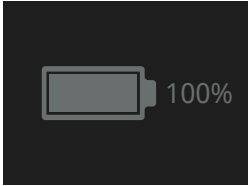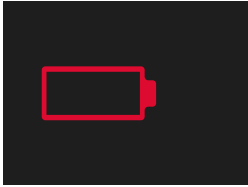Standby screen	Meaning
	No battery in the device
	Battery charging
	Battery fully charged
	Replace internal battery • Defective battery or • Battery capacity too low or • Battery incompatible

Battery operation with internal battery

Battery operation means that the device is supplied with power by the internal battery. If the device is supplied with power by the internal battery, the device automatically switches off after two minutes if ventilation is not active.

Please note:

- *The internal battery is solely intended to provide temporary power in the event of outages in the mains supply and when changing the power source. It must not be used as the primary power source for ventilation.*
- *In battery operation, keep track of the battery power level and recharge it as needed.*
- *In order to ensure the full function of the battery, the battery must be maintained in accordance with the "Battery maintenance" on page 145.*
- *In battery operation, the device cannot be operated with a humidifier.*

Device operating times with a new, fully charged battery:

Battery power level	Time	Alarm
100 – 0%	approx. 200 min	-
> 20 ¹ – 0%	approx. 35 min	Low Internal Battery
approx. 5 – 0%	approx. 5 min	Internal Battery Empty

Measuring conditions:

Output volume V_{del} = 800 ml, ventilation frequency f = 20 rpm

I:E ratio = 1:2, resistance R = 5 hPa (l/s)⁻¹ ± 10%, compliance C = 50 ml (hPa)⁻¹ ± 5%

The internal battery enables operation of at least 1 hour at maximum power consumption.

 CAUTION

Ventilation failure!

If the alarm "Internal battery empty" occurs, the device will switch off after approx. 5 minutes.

⇒ Connect the device to an alternate power supply right away.

1 The display "Internal battery low" is shown at 20% of the battery's rated value.

The device will automatically switch off the ventilation 5 minutes after the “Internal battery empty” alarm occurs. A notice window will appear “Ventilation currently not possible. Internal battery empty”. The device will switch off once another minute has passed. This prevents the battery from being fully drained and ensures the device shuts down properly.

Please note: This 5-minute period cannot be maintained if the battery was deeply discharged and was not charged up to 10%. In such a case, ventilation will end immediately.

Charge battery

To charge the battery, operate the device using the mains power. Recharging a fully discharged battery takes approximately 2.5 hours in mains operation. The device is fully functional during recharging.

Please note: The battery begins charging when the power level falls below 95%.

Power failure

Please note: During a power failure, the battery capacity display must be monitored and an alternative power source kept ready.

If the power supply is interrupted by a power failure and the device is on, the device is supplied with power via the internal battery.

Power failure and thus the switch to the internal battery is indicated by an alarm sound, as well as by the message “Power Failure”.

When the power supply returns, the device is supplied with power from the mains supply and the internal battery is charged.

Operation without internal battery

Please note: If the device is operated without a battery, it shuts off immediately in case of an interruption to the power supply. Connect the device to an alternate power supply right away.

If the power supply is interrupted during ongoing ventilation, an acoustic signal is output for approx. 2 minutes. The acoustic signal can be deactivated by pressing the on/off button.

System setup, non-invasive ventilation



WARNING

Risk of injury due to tube circuit and cable!

If laid incorrectly, the tube circuit or other cables (e.g. pulse oximetry) could cause patient strangulation.

⇒ Tubes and cables must always be positioned so that they cannot wrap around the neck or limbs of the patient.

Risk of suffocation due to closed exhalation valve!

If the opening of the exhalation valve is closed, the patient could suffocate.

⇒ Make sure that the opening of the exhalation valve remains open so that expired air can escape.

During non-invasive ventilation, the user must always be able to exhale (expire). This can either be accomplished using a vented mask with integrated exhalation valve or a non-vented mask with separate exhalation valve.

The exhalation valve is optionally available as an accessory.

Setup with vented mask

Connect the components as follows:



A Vented mask B Leakage tube

Figure 11: System setup, non-invasive ventilation with vented mask

Setup with non-vented mask

Connect the components as follows:



A Non-vented mask B Passive exhalation valve C Leakage tube

Figure 12: System setup, non-invasive ventilation with non-vented mask



WARNING

Lung damage due to insufficient humidifier output!

The AquaTREND uni humidifier must not be used for invasive ventilation, as the humidifier does not have a temperature or humidity control.

⇒ Use an external humidifier according to DIN EN ISO 8185 with a humidity output of > 33 mg/l.

Risk of injury due to overheated ventilation air!

If the humidifier does not have enough water, overheated ventilation air can dry out the airways.

⇒ Always top up the humidifier to the maximum fill level and make sure there is always enough water.

Risk of infection due to germs!

Condensation of water in the tube circuit can result in germ formation.


⇒ Increase the room temperature, reduce the heat output of the humidifier or use heated tubes.

Please note:

- *Before using AquaTREND uni, make sure to read the safety and cleaning instructions in the user's manual.*
- *In battery operation (see page 47) the humidifier heating is deactivated.*
- *The humidification output of > 10 mg/l according to DIN EN ISO 8185:2009-07 is possible for a therapy pressure of up to 30 hPa.*

During non-invasive ventilation, LAVI can be used to humidify the breathing air with the integrated humidifier AquaTREND uni. The humidifier features an integrated heater.


Connecting AquaTREND uni to the device:

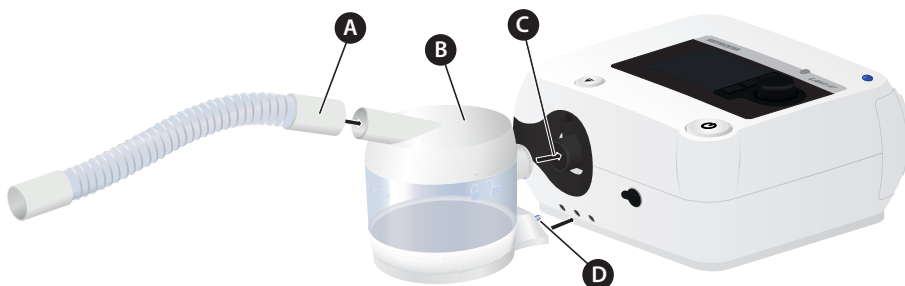

A Leakage tube B Humidifier C Air outlet D Contact pins

Figure 13: Connecting AquaTREND uni

1. Connect the humidifier to LAVI. You will feel it click into place. Make sure that you position the contact pin and the air outlet on the humidifier on the device without twisting or canting.
2. Connect the leakage tube to the tube connection on the humidifier.
3. Change the heating level on the system screen as needed. You can choose a value from 1 to 5, with 1 being the lowest heating level and 5 the maximum heating level.

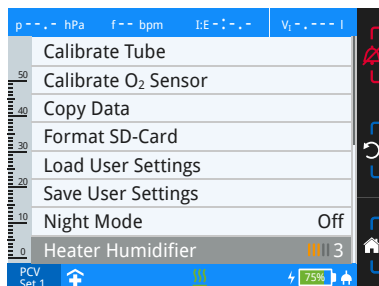

Switching on the heating

When the humidifier is in standby operation, the heating is automatically switched on when ventilation begins.

To preheat the water, you can switch on the heating before beginning with ventilation. Press the heating key.


Please note: If the humidifier heating is switched on while therapy is not in progress, the heating will be turned off after an hour for safety reasons.

Switching off the heating

The heating is automatically switched off when ventilation ends.

To switch off the heating manually, press the heating key on the home screen or in the first level of the other screens.

Disconnecting AquaTREND uni from the device:

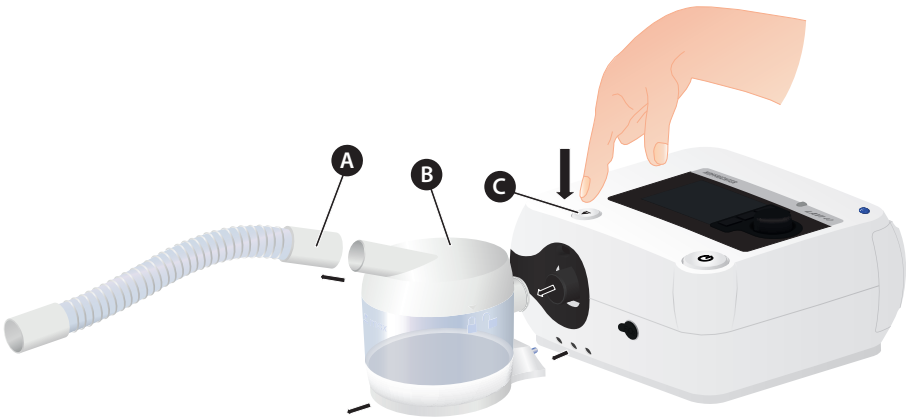

A Leakage tube B Humidifier C Release button

Figure 14: Disconnecting AquaTREND uni

1. Make sure that the ventilator has been switched off.
 2. Disconnect the leakage tube from the humidifier.
 3. Press the release button while removing the humidifier from the device.
-



WARNING

Risk of injury due to tube circuit and cable!

If laid incorrectly, the tube circuit or other cables (e.g. pulse oximetry) could cause patient strangulation.

⇒ Tubes and cables must always be positioned so that they cannot wrap around the neck or limbs of the patient.

Risk of suffocation due to missing exhalation valve!

If no option for exhaling through a separate valve is available in the invasive ventilation system, the patient could suffocate.

⇒ Use a separate exhalation valve.

Risk of suffocation due to closed exhalation valve!

If the opening of the exhalation valve is closed, the patient could suffocate.

⇒ Make sure that the opening of the exhalation valve remains open so that expired air can escape.

Risk of injury due to cold or dry breathing air!

Cold or dry breathing air presents an increased risk of collapsed lung and thus impaired gas exchange. This can also dry out the mucous membranes, resulting in an increased risk of infection. For this reason, the breathing air should be heated and humidified.


⇒ Use an external humidifier or HME filter.

Setup with external humidifier

Please note:

- Use an approved humidifier according to DIN EN ISO 8185 with a humidity output of $> 33 \text{ mg/l}$.
- Follow the manufacturer user's manual.
- The humidifier should be positioned below the patient and the device, so that no water can accumulate in the patient's lungs or in the ventilator. If water accumulates in the tube circuit, we recommend using water traps.

Connect the components as follows:

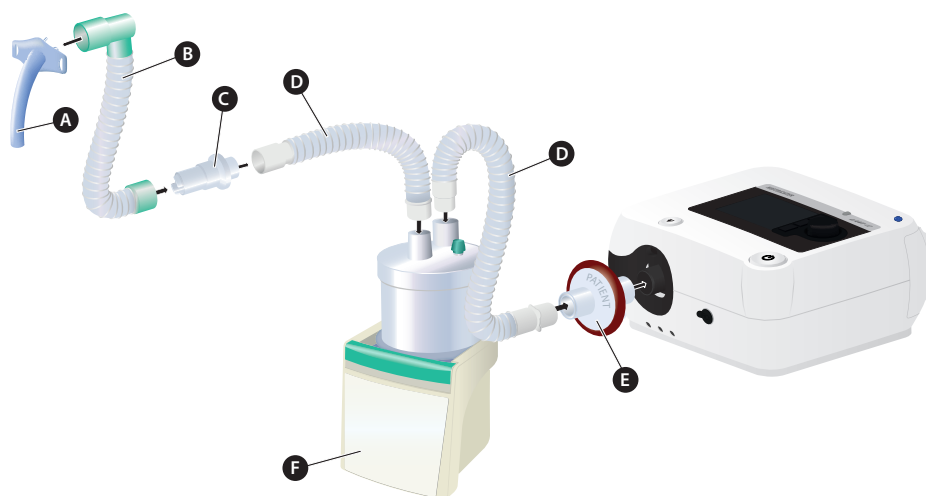

A Tracheal cannula B Catheter mount tube C Passive exhalation valve D Leakage tube
E Bacterial filter F Humidifier

Figure 15: System setup, invasive ventilation with external humidifier


1. Connect all components as per Figure 15.
 2. Calibrate the connected tube circuit (see page 56).
-

Setup with HME filter

If no humidifier is used, we recommend using an HME filter to maintain the moisture of the respiratory gas. A “combined filter” consisting of an HME filter and a bacterial filter is well-suited (e.g. Medisize Hygrovent HMEF).

Please note: When using a HME filter, read the user’s manual from the respective manufacturer. In particular, follow all information related to replacement intervals.

Connect the components as follows:

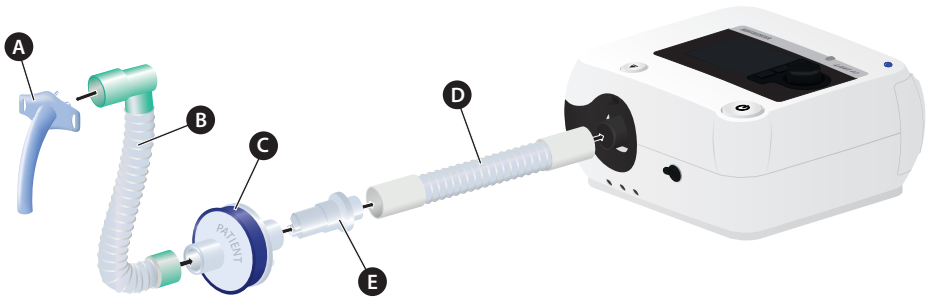

A Tracheal cannula B Catheter mount tube C HME filter D Leakage tube
E Passive exhalation valve

Figure 16: System setup, invasive ventilation with HME filter

1. Connect all components as per Figure 16.
 2. Calibrate the connected tube circuit (see page 56).
-

Calibrating the tube circuit

A tube calibration function is available in the device to enable use of different tube circuits and accessories. During tube calibration, the resistance of the circuit upstream of the air outlet is determined, which forms the basis for correct pressure measurement. Perform tube calibration when changes have been made to the system upstream of the air outlet. These may include connecting and disconnecting of the following components, for example:

- Bacterial filter, humidifier, tube circuit, FiO₂ sensor etc.

**CAUTION**


Risk of injury due to incorrect pressure measurement!
If the device is operating with incorrect calibration data, the pressure measurement may be distorted.
⇒ Before commissioning and changing the tube circuit, perform tube calibration.

To calibrate the tube circuit:

Please note: *During tube calibration, the patient tube circuit must be open and ventilation switched off. A mask can be connected.*

1. On the home screen, navigate to “System” and press the MFK.
2. Navigate to “Calibrate Tube” and press the MFK.

During calibration, the system shows “Run...”.

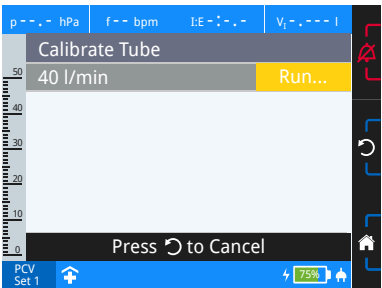


3. If the calibration was successful, “Ok” will appear after a few seconds.

To finish, press the MFK. The calibration data are then applied. Until then, you can cancel calibration at any time by pressing the Escape key.



If the calibration was not successful, "Error" will appear.

In the event of an error, check the entire system. Resistance in the overall system may be too high. You may, for example, have to exchange the bacterial filter or use another humidifier. Then rerun the calibration.

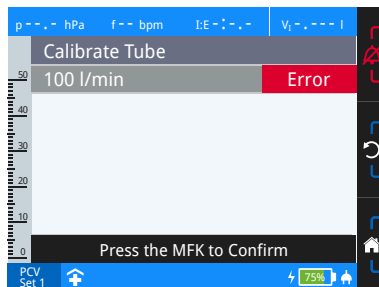

If you do not perform a tube calibration, the calibration data from the last tube calibration are used. When commissioning the device for the first time, the default calibration data¹ stored in the device are used.

1 Operation without humidifier:

Leakage tube circuit: Ø 22 mm; 1.80 m

Height: ca. 44 m NHN, air pressure: approx. 1008 hPa

Connected ventilation components and accessories: no

Operation with humidifier:

Leakage tube circuit: Ø 22 mm; 1.80 m

Height: ca. 44 m NHN, air pressure: approx. 1008 hPa

Connected ventilation components and accessories: Humidifier AquaTREND uni

Connecting additional accessories


Connecting the bacterial filter

If the device is intended for use by more than one patient (e.g., in operation in clinics), a suitable bacterial filter (e.g., Medisize Barr-vent S) should be used continuously to protect the device and the patient from contamination by human pathogens.

Please note:

- *Please be sure to replace the bacterial filter daily and follow the manufacturer's user's manual.*
- *If the optionally available humidifier AquaTREND uni is used for ventilation, do not use a bacterial filter.*

Connect the components as follows:

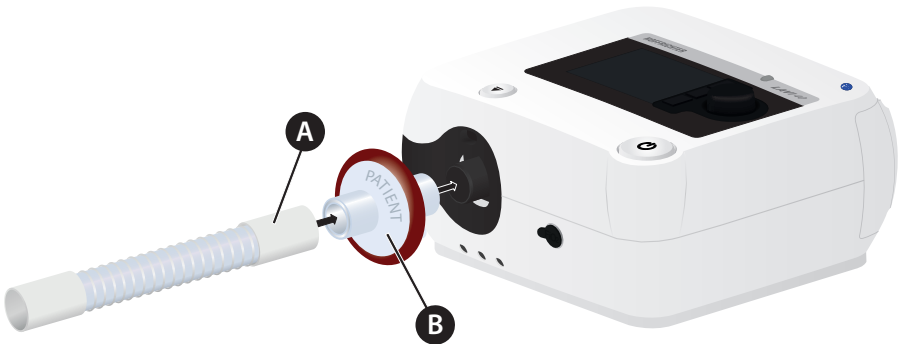


A Leakage tube B Bacteria filter

Figure 17: Connecting the bacterial filter

1. Connect all components as per Figure 17.
 2. Calibrate the connected tube circuit (see page 56).
-

Connecting the alarm box or nurse call

Connect the components as follows:

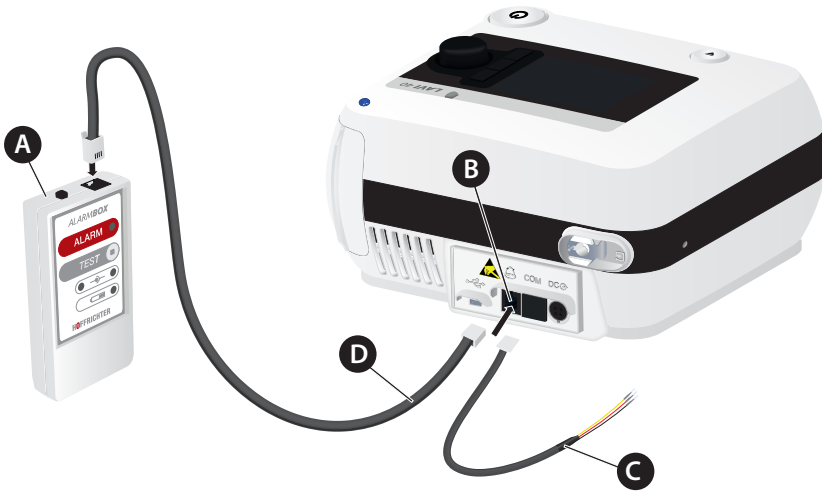

A Alarm box B Remote alarm/nurse call connection C Nurse call cable
D Alarm box cable

Figure 18: Connecting the alarm box/nurse call

Connect the alarm box to the device as shown in Figure 18.

Alarm boxes are available as an accessory (see page 167).

An on-site nurse call can also be connected to the remote alarm/nurse call connection point as well. You will need a cable with a RJ10 plug to do this. The cables are available as an accessory (see page 167).

Additional information on alarm boxes and forwarding alarms is available on page 125.

Using an SD card

To be able to pass on therapy data to a doctor, you can copy the data to an SD card. You can use SD and SDHC cards up to 32 GB. You can format the SD card on LAVI and delete any existing data. This method is described under “Copy data”.

Please note: *We recommend using SD cards with a minimum of 8 GB memory.*

Copy the data as follows:

Please note: *We only recommended copying data when therapy is not in progress. If therapy is active, the data cannot be copied in full.*

1. Insert the SD card into the SD card slot until it clicks into place as shown in the illustration.  then appears on the toolbar.



Figure 19: Inserting the SD card

2. On the home screen, navigate to “System” and press the MFK.
3. Navigate to “Copy data” and press the MFK.

Possible messages during the copying process:


“Run...”

⇒ The copying process is being performed

“SD card full” or “Finished...error”


⇒ Format the SD card and copy again. Please
Please note: *This will delete all data on the SD card.*

⇒ Use another SD card.

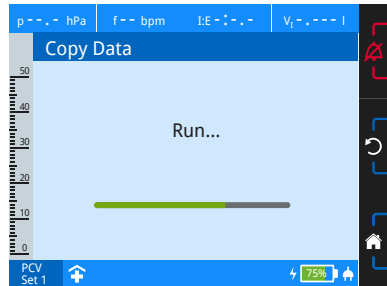

- Once successfully copied, the message “Finished...ok” will appear.

Press the MFK to close the notice window.

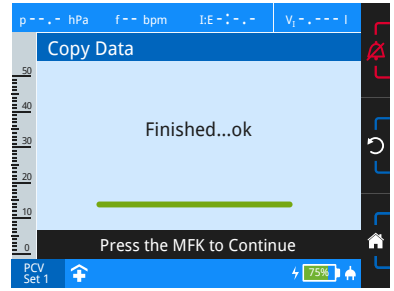

- To remove the SD card, carefully press the SD card into the SD card slot. Then, remove the SD card.



Figure 20: Removing SD card

Using oxygen

ATTENTION

Material damage due to oxygen accumulating in the device!

Oxygen may only be supplied during active ventilation.

Please note: *Before using oxygen, read the safety instructions from page 26.*

The supply of oxygen is possible in all ventilation modes. Please note that any changes to the ventilation parameters, e.g. pressure, I:E, frequency, will lead to a change of the FiO₂ concentration.

Connecting the oxygen source

ATTENTION


Material damage on the back-stop due to incorrect oxygen connection adapter!

If the wrong oxygen connection adapter is used, there is a risk that the back-stop in the connection is damaged.

⇒ Only use the supplied oxygen connection adaptor for this purpose.

Connect the oxygen supply as follows:

Connect the oxygen source to the device as shown in Figure 21.

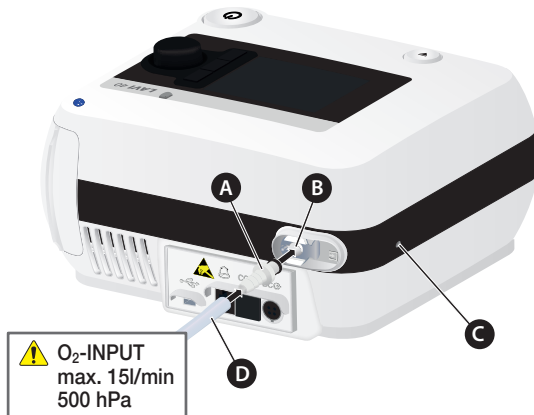

A Straight O₂ connection adapter B Oxygen connection C Oxygen output
D Tube from the oxygen source (acc. to EN ISO 5359)

Figure 21: Connecting the oxygen source

Measuring the oxygen concentration

The oxygen concentration may be inconsistent when feeding in a fixed value oxygen flow (FlowO₂). The inspirational oxygen concentration (FiO₂) can vary depending on pressure, ventilation pattern of the patient, mask or leakage. The oxygen concentration should therefore always be measured with a (FiO₂) sensor when oxygen is being supplied (see accessories on page 166). The FiO₂ sensor must be calibrated for exact results (see page 63).

Starting the supply of oxygen



Risk of injury due to incorrect or contaminated oxygen!

Oxygen that is not suited for medical purposes or is contaminated can result in serious health issues.

⇒ Use only certified and clean oxygen sources.

1. Start ventilation and wait for several respiratory cycles.
2. Start supplying the oxygen.

Stopping the supply of oxygen

1. Stop the supply of oxygen at the oxygen source.
2. Continue ventilation for a number of respiratory cycles.
3. Stop ventilation.

Calibrating the FiO₂ sensor

Calibration is performed using ambient air, assuming an oxygen content of 21%.

Automatic calibration when ventilation is off (recommended)

If ventilation is not in progress and you install the FiO₂ sensor, the FiO₂ sensor is calibrated automatically. The FiO₂ measurements are shown on the monitoring screen.


Calibration when ventilation is on

During ventilation, the FiO₂ sensor cannot be calibrated automatically, as there is oxygen in the device. First, stop the supply of oxygen and leave ventilation active for at least 30 seconds so the oxygen in the device can be broken down. Then, install the FiO₂ sensor and turn the oxygen supply back on.

Manual calibration

A manual calibration may be performed in the system settings at any time. During continuous ventilation, we recommend calibrating the FiO₂ sensor manually once a week.

To calibrate the FiO₂ sensor:

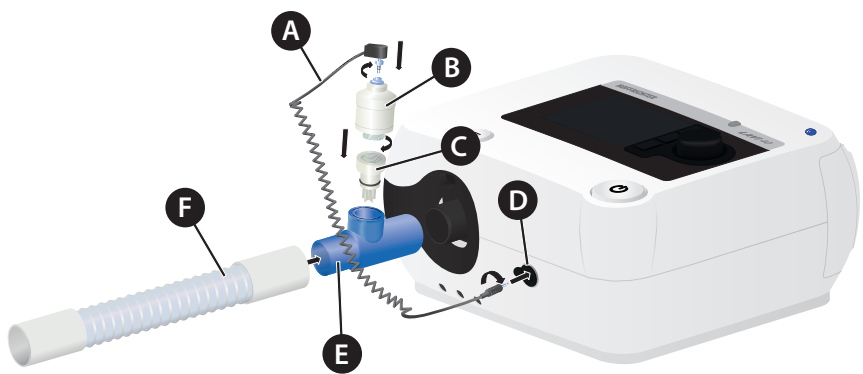


A Connecting cable B FiO₂ sensor C Gas duct housing D FiO₂ sensor cable connection
E T adapter F Leakage tube

Figure 22: Connecting the FiO₂ sensor


1. Make sure that the ventilator has been switched off.
2. Install the FiO₂ sensor according to Figure 22.
Please note: Plug and screw the straight plug of the connecting line (A) to the device and then connect the right-angled to with the FiO₂ sensor.
3. On the home screen, navigate to “System” and press the MFK.
4. Navigate to “Calibrate FiO₂ Sensor” and press the MFK. Calibration begins.

During calibration, the system shows “Run...” .

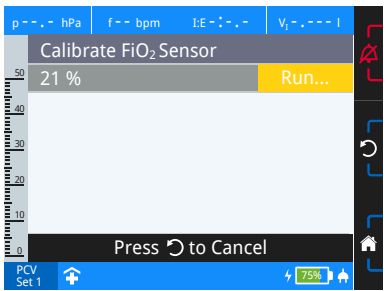


5. If the calibration was successful, "Ok" will appear after a few seconds.

To finish, press the MFK.

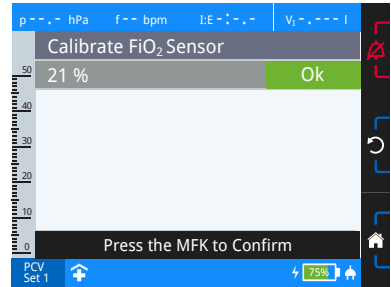

If the calibration was not successful, "Error" will appear. In the event of an error, repeat the calibration. If the calibration is still unsuccessful, replace the FiO₂ sensor.


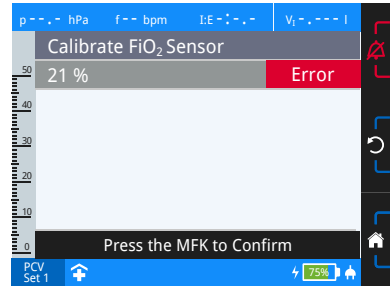


Please note:


- Depending on environmental conditions and the storage time, the sensor may take up to 15 minutes after connecting to reach signal stability again.
- FiO₂ sensors have a limited service life. The service life of the FiO₂ sensors supplied by HOFFRICHTER is approx. 1 year from the date of manufacture at an oxygen concentration of 40%. The FiO₂ sensors should not be stored for more than 6 months. You will find the date of manufacture on the FiO₂ sensor. For the longest possible sensor service life, we recommend storage at +5°C to +30°C.

Determining the oxygen concentration


The oxygen concentration (FiO_2) in the patient's respiratory system depends on the flow rate (Flow O_2) of the oxygen supply and the breathing minute ventilation (MV) of the patient. The following diagrams provide values for determining oxygen concentration at the frequencies of 15 bpm, 20 bpm and 25 bpm.¹

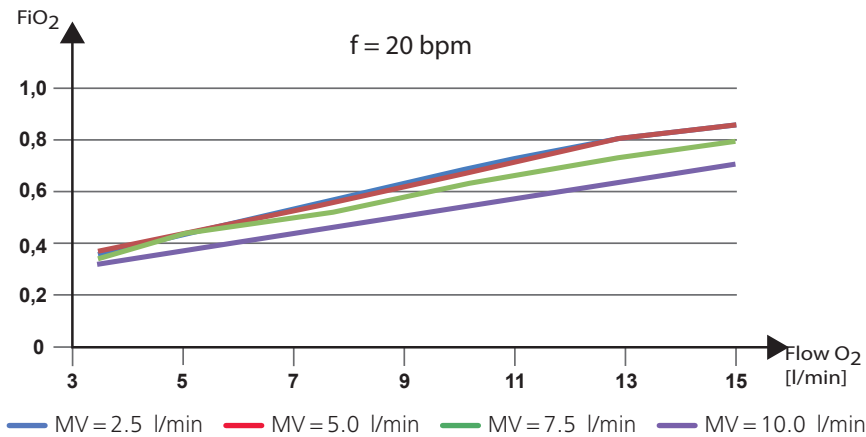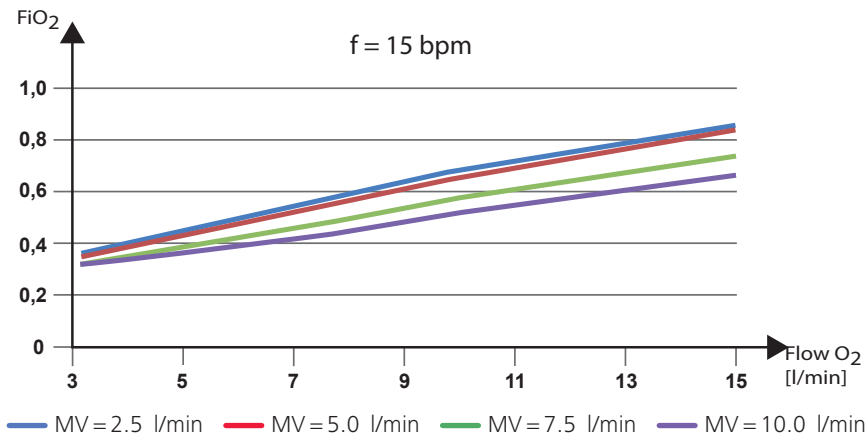

1 The measurement has been performed with PEEP = 4 mbar and using a test lung (resistance 5 mbar/l/s, tidal volume max. 1000 ml). Depending on the condition of the patient's lung, the oxygen concentration values may differ from the ones measured here. Measurement deviation can be up to max. 10 %. The measured oxygen concentration also depends on the age and condition of the FiO_2 sensor.

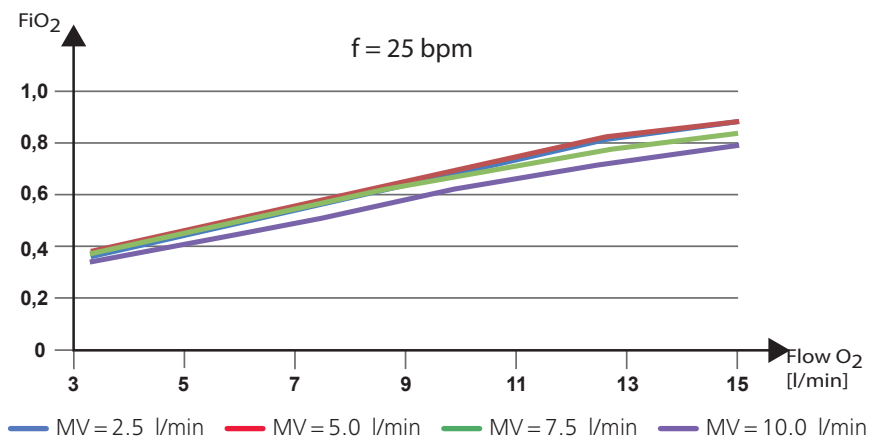

Example

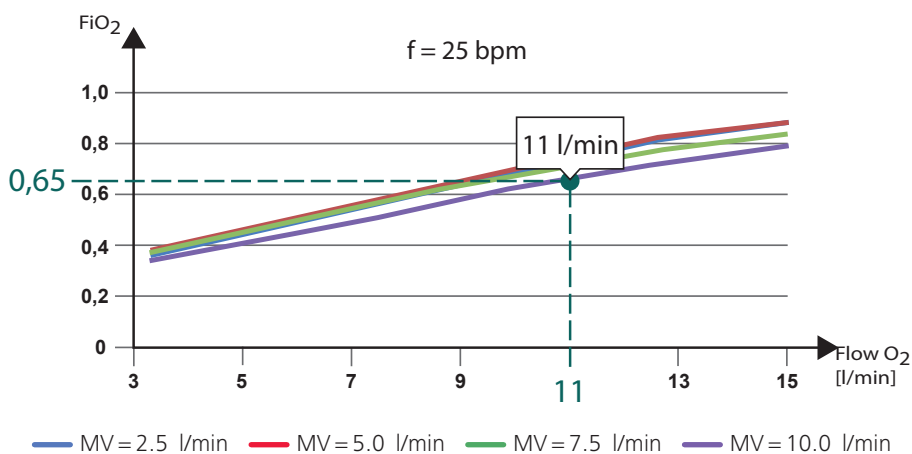

At an oxygen flow rate of 11 l/min and a minute ventilation of approximately 10 l/min, an oxygen concentration of approximately 65% can be reached.

Using the functional bag



WARNING

Risk of injury due to insufficient monitoring of the device functions!

If important device functions are not shown or available, proper operation cannot be guaranteed.

⇒ Only use the original HOFFRICHTER functional bag.

The functional bag protects the ventilator in mobile use (e.g. at the wheelchair or walker) from mechanical damage or effects of the weather. The functional bag is available as an accessory (see page 167).



Figure 23: Functional bag

When using the device in the functional bag the following instructions must be observed to ensure safe and trouble-free operation:

- Set the alarm sound to level 3.
- Make sure that all alarm messages are visible through the viewing window and that the air vents of the bag are not blocked. The air supply for the device must be guaranteed at all times.
- Use the bag in such a way that the device is protected from overheating, dust and water.
- All accessories connected, such as tube, filter, supply lines etc., must be arranged so that they cannot cause any malfunctions of the device. Accidental disconnection of the accessory parts must be avoided.
- Follow all instructions in the user's manual when using the functional bag in combination with other accessory parts.

Switching on the device

The device requires up to one minute to boot before it is ready for use.

Please note:

- *The tube circuit may be connected when the device is started up, but it must not be connected to the patient yet.*
- *If you are using oxygen therapy during ventilation, please note the section “Using oxygen” from page 62.*

Mains operation

Connect the device to the mains supply. To switch the LAVI on, press the on/off button for more than 4 seconds.

Battery operation

To switch the LAVI on in battery operation, press the on/off button for more than 4 seconds.

Switching off the device

Mains operation

Stop ventilation. To switch on the device, press the on/off button for more than four seconds.

Battery operation

Stop ventilation. The device shuts off after approximately two minutes. To switch the device off right away, press the on/off button for more than four seconds.

Please note: The device boots down when you switch it off. This process takes a few seconds. Do not disconnect the device from the mains or remove the internal battery during this time. Otherwise, data may be lost.

Start ventilation

Press the on/off button. Ventilation begins.

Stop ventilation

Press the on/off button and confirm the prompt with “Yes”.

Chapter 5

Ventilation Modes

The device offers three basic ventilation forms:

- Mandatory ventilation modes where the device controls the respiratory work for the patient completely.
- Augmented ventilation modes, where the device performs part of the respiratory work, alternating or overlapping with the patient's breathing rate.
- Spontaneous ventilation modes, where the patient does the respiratory work with the support of the device. The patient determines the inspiration time and frequency.

PCV mode

Pressure Controlled Ventilation

Main features

- Pressure controlled
- Device triggered
- Time controlled
- Fixed frequency
- No spontaneous breathing possible

In this ventilation mode, ventilation is controlled exclusively by the device. Spontaneous breathing on the patient's part is not possible. The ventilation period is based on the set frequency and a defined I:E ratio.

The inspiratory pressure (IPAP) as well as the end-expiratory pressure (PEEP) define the range of pressure for ventilating the patient. The time for pressure increase from PEEP to IPAP is set via the parameter "I flank" while the decrease is set via the parameter "E flank".

The inspiration volume is automatically adjusted to the condition of the lung (compliance and resistance).

To ensure a minimum volume with LAVI 40, you can optionally enter a target volume with the option of reaching this volume by increasing pressure (IPAP + additional pressure).




Important settings

- “High inspiration volume” alarm limit
- “Low inspiration volume” alarm limit



Tips

- Adjustable minimum volume (only LAVI 40)
- Adjustable “Additional Pressure” to reach the minimum volume (Only LAVI 40)

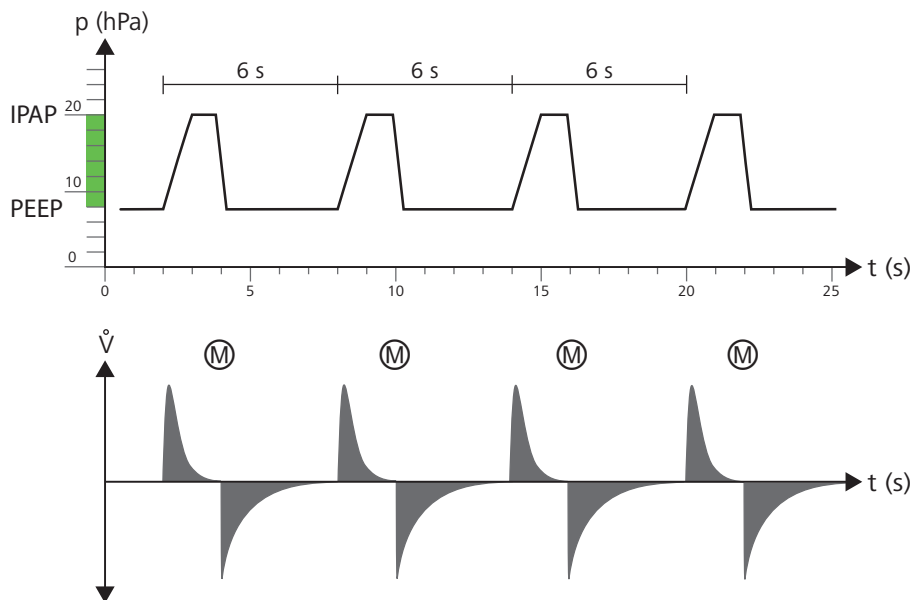

IPAP = 20 hPa PEEP = 8 hPa I:E = 1:2 (2 s : 4 s) f = 10 bpm

Ⓜ Mandatory inspiration

Figure 24: PCV mode diagram

APCV mode

Assisted Pressure Controlled Ventilation

Main features

- Pressure controlled
- Device or patient triggered
- Time controlled
- Backup frequency
- Spontaneous breathing possible

In its ventilation parameters, the pressure controlled assisted ventilation is equal to the exclusively controlled ventilation.

By setting an inspiration trigger, however, the patient can stop the expiration phase by inspiration efforts once they reach the trigger threshold and initiate the next inspiration phase. These additional breathes are only controlled by the device.

The inspiration time is defined. The patient can only shorten the expiration time by their own breathing efforts, thus increasing the set frequency.




Important settings

- “High inspiration volume” alarm limit
- “Low inspiration volume” alarm limit
- Alarm limit “High Frequency”



Tips

- Adjustable minimum volume (only LAVI 40)
- Adjustable “Additional Pressure” to reach the minimum volume (Only LAVI 40)
- Adjustable “Trigger Lock” of the inspiration trigger

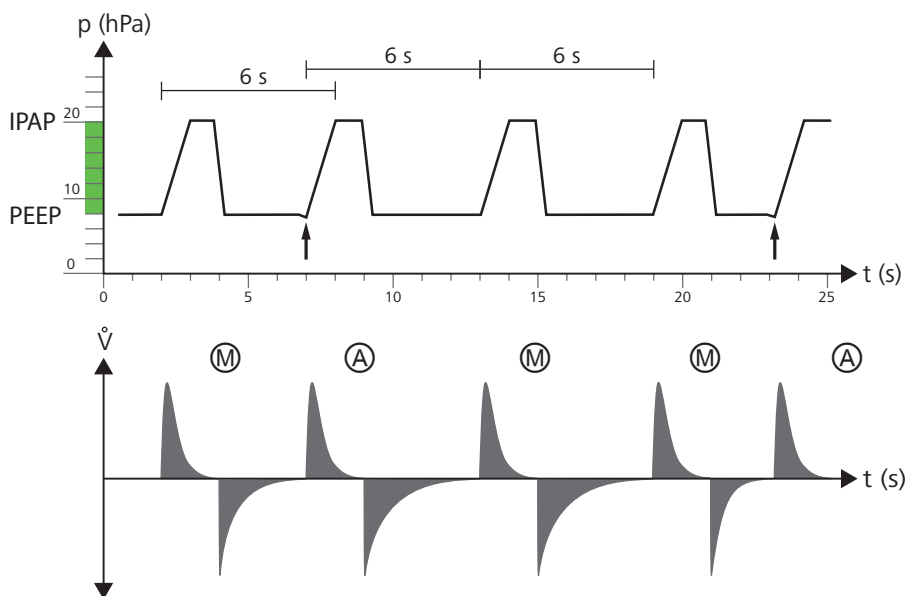

IPAP = 20 hPa PEEP = 8 hPa I:E = 1:2 (2 s : 4 s) f = 10 bpm

↑ Spontaneous respiration (M) Mandatory inspiration

(A) Pressure-controlled assisted inspiration

Figure 25: APCV mode diagram

PSV mode

Pressure Supported Ventilation

Main features

- Pressure supported
- Device or patient triggered
- Flow controlled
- Background frequency
- Spontaneous breathing possible

Pressure supported ventilation is intended to support spontaneous breathing and to initiate mechanical ventilation whenever spontaneous breathing is absent.

The inspiratory pressure (IPAP) as well as the end-expiratory pressure (PEEP) define the range of pressure for ventilating the patient. The time for pressure increase from PEEP to IPAP is set via the parameter "I flank" while the decrease is set via the parameter "E flank".

The inspiration volume is automatically adjusted to the condition of the lung (compliance and resistance). The mechanical ventilation cannot be influenced by the patient's spontaneous breathing.

The trigger thresholds of the inspiration trigger and the expiration trigger can be adjusted according to the patient's requirements.

The adjustable frequency is set as a backup frequency. As long as the patient reaches or exceeds this frequency, the device reacts with the pressure support to each spontaneous inspiration, following the patient's breathing. If the backup frequency fails to be reached, the device assumes mechanical ventilation until it registers the next spontaneous breath.

To permit respiratory pauses between the patient's breathing efforts, an apnoea time can be set to delay the start of the mechanical ventilation.

To ensure a minimum volume with LAVI 40, you can set a pressure increase (IPAP + additional pressure).




Important settings

- “High inspiration volume” alarm limit
- “Low inspiration volume” alarm limit
- Alarm limit “High Frequency”
- Alarm limit “Low Frequency”



Tips

- Adjustable minimum volume (only LAVI 40)
- Adjustable “Additional Pressure” to reach the minimum volume (Only LAVI 40)
- Adjustable “Trigger Lock” of the inspiration trigger
- Adjustable “Apnoea Time”
- Minimum (Ti Min) and maximum (Ti Max) inspiration times are adjustable

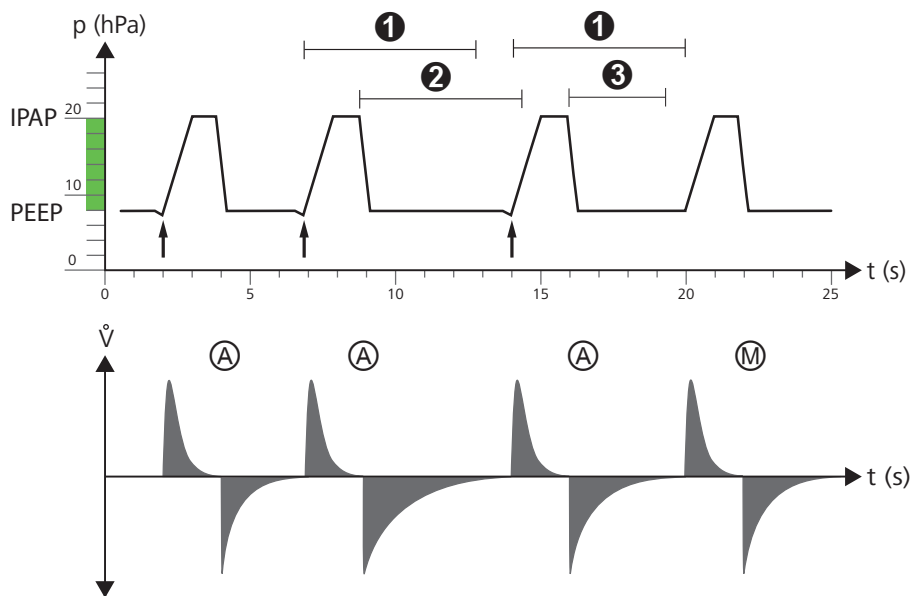

IPAP = 20 hPa PEEP = 8 hPa ↑ spontaneous respiration
 ① Background frequency ② Apnoea time > expiration time
 ③ Apnoea time < Expiration time
 (M) Mandatory inspiration
 (A) Pressure controlled supported inspiration

Figure 26: PSV mode diagram

PSV-S mode

Please note: *In PSV-S mode, the device responds only to spontaneous breathing of the patient.*

Pressure Supported Ventilation-Spontaneous

Main features

- Pressure supported
 - Patient triggered
 - Flow controlled
 - Device only responds to spontaneous breathing
 - Apnoea Alarm
-

In its ventilation parameters, the PVS-S mode is equal to the PSV mode.

The PSV-S mode ventilation parameters adjustment options are the same as those of the PSV mode. Inspiration triggers are only actuated by the patient’s spontaneous breathing, since there is no frequency setting.

The apnoea time is an alarm parameter (“Apnoea Alarm”).



Important settings

- “High inspiration volume” alarm limit
- “Low inspiration volume” alarm limit
- Alarm limit “High Frequency”
- Alarm limit “Low Frequency”
- Alarm limit “Apnoea Alarm”



Tips

- Adjustable minimum volume (only LAVI 40)
- Adjustable “Additional Pressure” to reach the minimum volume (Only LAVI 40)
- Adjustable “Trigger Lock” of the inspiration trigger
- Minimum (Ti Min) and maximum (Ti Max) inspiration times are adjustable

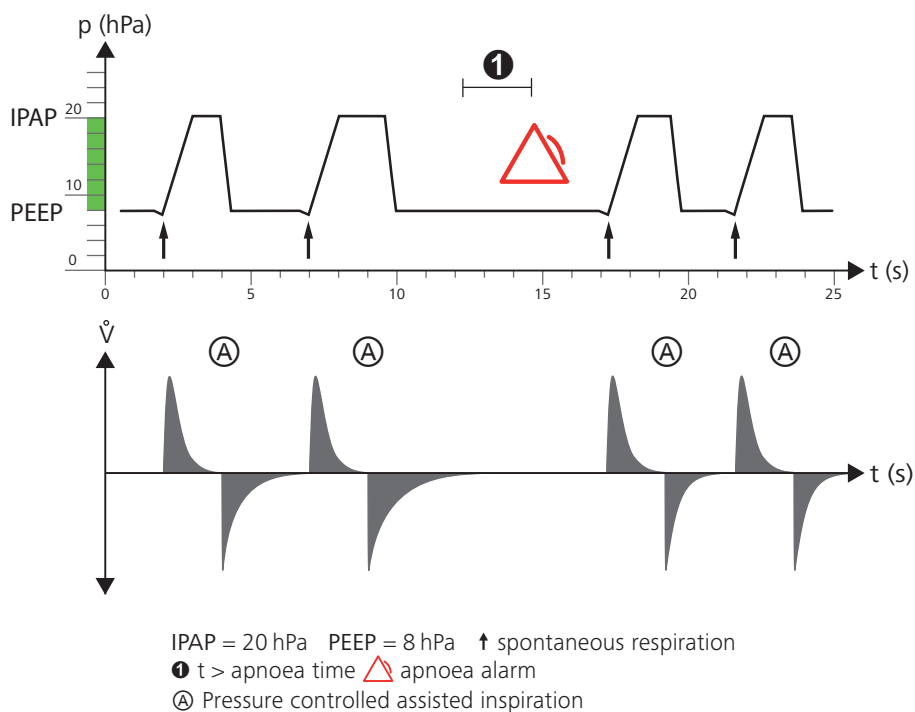


Figure 27: PSV-S mode diagram

P-SIMV mode

Please note: *P-SIMV mode is only available in LAVI 40.*

Pressure Controlled Synchronized Intermittent Mandatory Ventilation

Main features

- Pressure controlled
- Device or patient triggered
- Time controlled
- Fixed frequency
- Spontaneous breathing possible

The SIMV mode provides a combination of pressure controlled ventilation and pressure assisted spontaneous breathing.

The mechanical ventilation is based on a fixed frequency and an inspiration time.

The inspiratory pressure (IPAP) as well as the end-expiratory pressure (PEEP) define the range of pressure for ventilating the patient. The time for pressure increase from PEEP to IPAP is set via the parameter "I flank" while the decrease is set via the parameter "E flank".

The inspiration volume is automatically adjusted to the condition of the lung (compliance and resistance). The mechanical ventilation cannot be influenced by the patient's spontaneous breathing.

The patient can breathe spontaneously between mechanical ventilations when the inspiration and expiration trigger thresholds are reached. Spontaneous breathing is supported during inspiration by a predetermined and selectable pressure (PS), which is independent from the IPAP. Only the patient can predefine the length of the ventilations as well as the duration of the inspiration.

The mechanical breaths are adjusted to spontaneous breathing in terms of time. If, for example, a spontaneous inspiration occurs shortly before a SIMV period is started (within a specific expected time window = 5 seconds), the mechanical ventilation is synchronized to the breathing of the patient. Since the synchronization of the mandatory ventilation shortens controlled ventilation by Δt and the frequency would be increasing, the following breath will be extended by Δt . By contrast, patient attempts to breathe outside the expected timeframe are not supported by the machine. This ensures the set frequency is maintained.



Important settings

- “High inspiration volume” alarm limit
- “Low inspiration volume” alarm limit
- Alarm limit “High Frequency”



Tips

- Adjustable “Trigger Lock” of the inspiration trigger
- Minimum (Ti Min) and maximum (Ti Max) inspiration times are adjustable

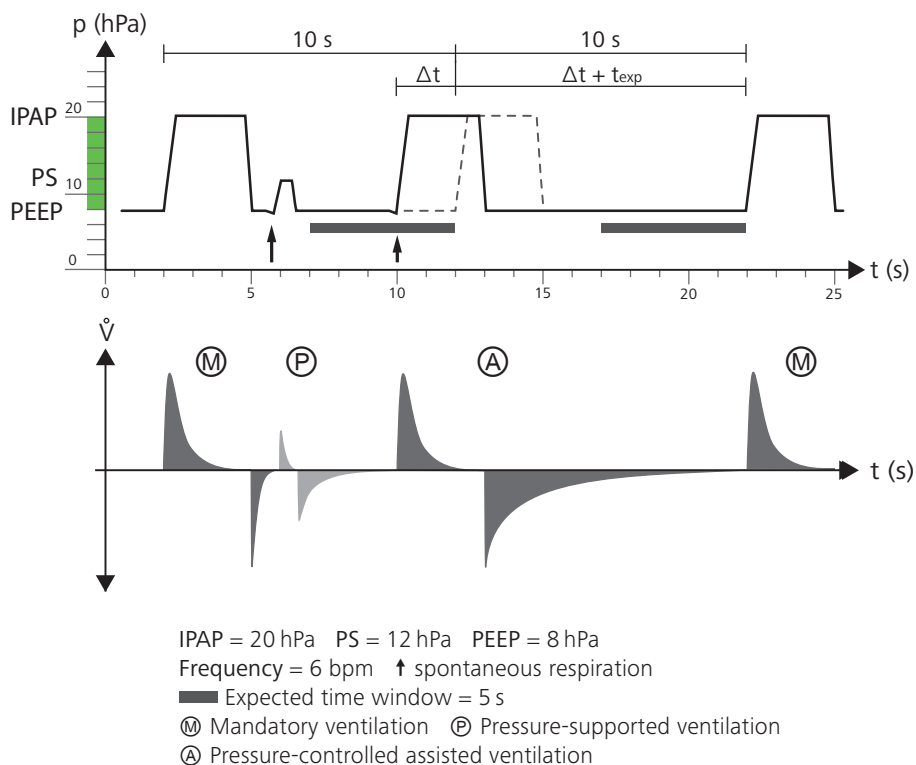


Figure 28: P-SIMV mode diagram

CPAP mode

Continuous Positive Airway Pressure

Main features

- Pressure controlled
- Independent respiration

In the CPAP mode, the device provides continuous positive pressure that helps the patient to inspire. The patient also expires against this pressure. The patient must be able to breathe spontaneously.


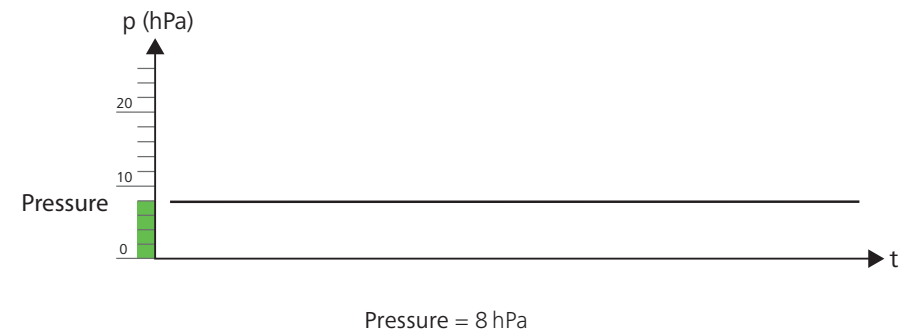


Figure 29: CPAP mode diagram


Chapter 6

Operating the Device

Menu structure

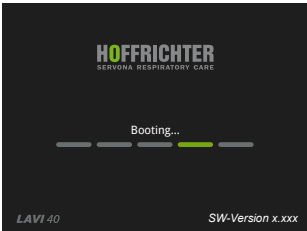


Start screen

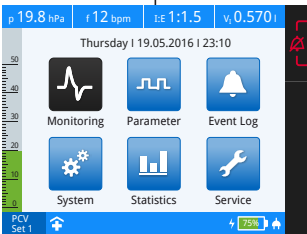

Home screen

Monitoring screen

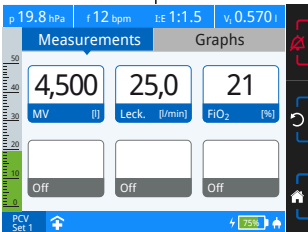


Parameter screen

	Set 1	Set 2	Set 3
Mode	PCV		
IPAP	20.0 hPa		
PEEP	5.0 hPa		
Frequency	12 bpm		
Time Inspiration	2.0 s		
I-Slope	3		
E-Slope	3		

Eventlog screen

Alarms	Alarms	Events
15.07.2016, 23:41 Set 1		
Apnoe		
09.07.2016, 22:02 Set 3		
Low Volume		
16.06.2016 21:33 Set 3		
Low Frequency		
24.05.2016 18:03 Set 3		
Leakage		

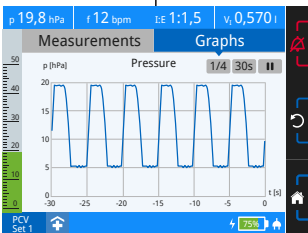

Alarms	Alarms	Events
12.07.2016, 01:53 Set 1		
Stenosis		
06.07.2016, 03:12 Set 1		
Internal Battery Empty		
06.07.2016 02:56 Set 1		
Low Internal Battery		
05.07.2016 22:24 Set 1		
Battery Operation		



Alarms	Alarms	Events
10.07.2016, 00:56 Set 1		
Therapy stopped		
10.07.2016, 00:06 Set 1		
Therapy started		
03.07.2016 22:08 Set 1	2.5 → 3.0 hPa	
High Press. Tolerance		
26.06.2016 05:12 Set 1	4.5 → 6.5 hPa	
PEEP		



Menu lock

The menu lock function can be used to protect against the accidental changing of device settings. It locks the following operating elements:

- Home key
- Escape key
- Turning and pressing the MFK
- Exception: *If a notice window is displayed that requests operation is performed via the MFK, the MFK can be used as usual.*


The following keys retain their full functionality:

- Alarm key
- Heating key
- On/off button

To activate the menu lock:

1. Hold the MFK in the home screen for longer than 1 second until the window appears at the side. Keep the MFK pressed until the progress bar has completed and  is displayed. Then release the MFK.


Please note: *If therapy is in progress, the measured values screen will automatically appear upon activating the menu lock. The home screen will be shown again once therapy is over.*

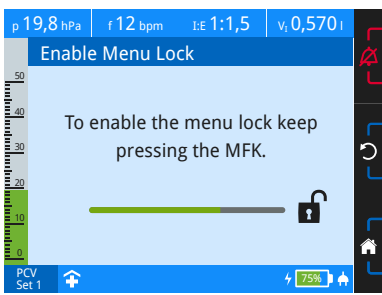

2. The symbol  will appear on the toolbar.


To deactivate the menu lock:

1. Hold the MFK for more than 1 second until the window appears at the side. Keep the MFK pressed until the progress bar has completed and  is displayed. Then release the MFK.

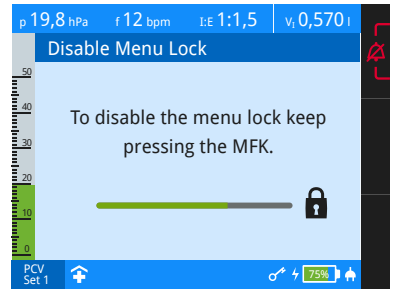

2. Select "Yes" and then press the MFK.

The symbol  will be hidden from the toolbar.

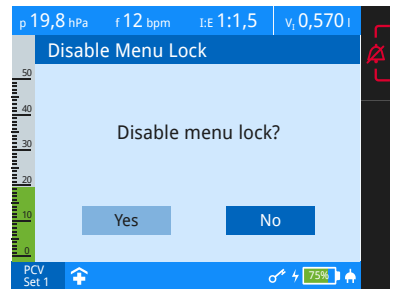

User profiles

The device can be operated in two different profiles - clinic  and home . The currently active user profile is shown on the toolbar.


In clinic mode, the user has access to all device settings. In contrast, the ventilation and alarm parameters cannot be changed in the home mode. It is possible to switch from one set to the other.

You can set the user profile on the system screen.

Home screen

You can enable individual screens from the home screen:

- Monitoring screen
 - Measurements: Monitoring measurements (numerical)
 - Graphs: Monitoring measurements (graphs)
- Parameter screen
- Selecting the ventilation set, setting the ventilation and alarm parameters
- Event log screen
 - Alarms 🧑: Displays physiologically-dependent alarms
 - Alarms 🔧: Displays technically-dependent alarms
 - Events: Displays events
- System screen
- System settings, calibrations, counters and device information
- Statistics screen
- Statistical evaluation reports
- Service screen
- System calibration and tests for service work (PIN code protected)

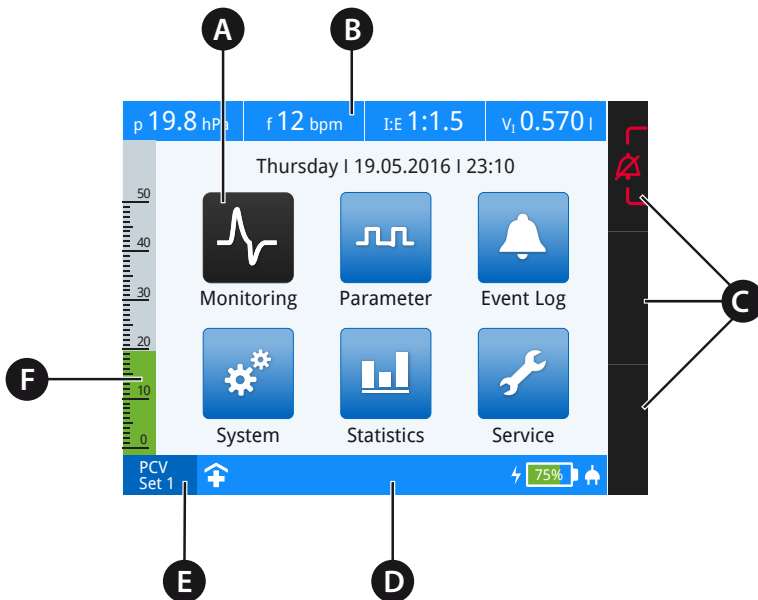

A Selected screen symbol B Measurement display C Softkeys D Toolbar
E Active ventilation mode and active set F Pressure bar

Figure 30: Home screen

To enable a screen:

1. Navigate to the desired screen by turning the MFK.




Selected screen symbol
→ black

- 2.



Symbol not selected
→ blue

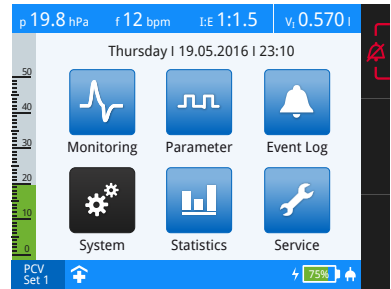

3. Press the MFK to activate the selected screen.

Change the display mode

In the factory state, the screen switches automatically 2 minutes after the last operation:

- During ventilation to the measurements/graph screen (depending on which screen was most recently active) or
- When ventilation is off, to the home screen

The time for changing the display can be modified on the system screen. If you do not want the display to change, set the parameter to "off".

Monitoring screen

The ventilation parameters are shown in real time on the monitoring screen.


The monitoring screen is divided into 2 sections:

- Measurements
- Graphs

Measurements

Depending on the configuration, the following parameters are shown under “Monitoring” > “Measurements” when ventilation is active:

- Pressure (p),
- Frequency (f)
- Inspiration to expiration ratio (I:E)
- Minute ventilation (MV),
- Volume inspiration (V_I),
- Leak rate
- FiO₂ concentration (FiO₂)

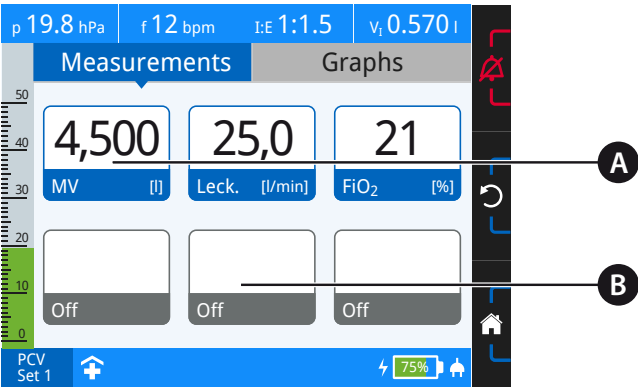


A Measurement configured B No measurement configured

Figure 31: Measurement screen, factory setting


Figure 31 depicts the default configuration as per factory settings. You can configure the measurement for each measurement box.

To configure the measured values:


1. Navigate to “Monitoring” on the home screen and press the MFK.
2. Navigate to “Measurements” and press the MFK. The first measured value box is marked.

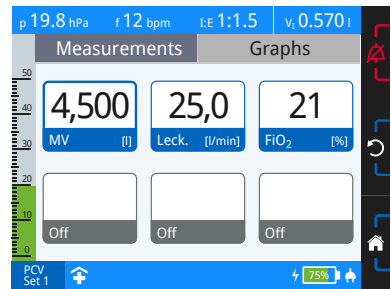

3. Navigate to the measurement box you want to configure. Press the MFK to activate the measured value box.

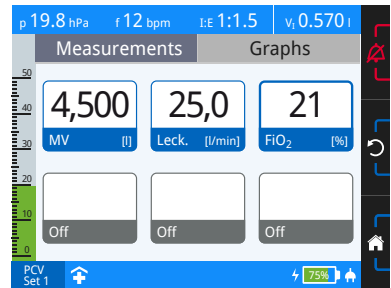


4. Change the measurement and press the MFK to confirm the setting.

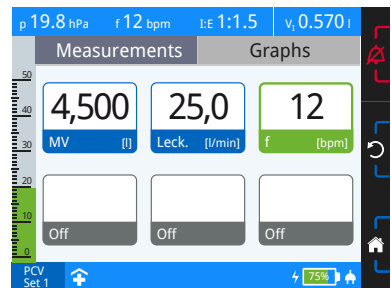

Graphs

The “Measurements” > “Graphs” section displays the following ventilation parameters when ventilation is running:

- Pressure
- Flow
- Volume
- FiO₂

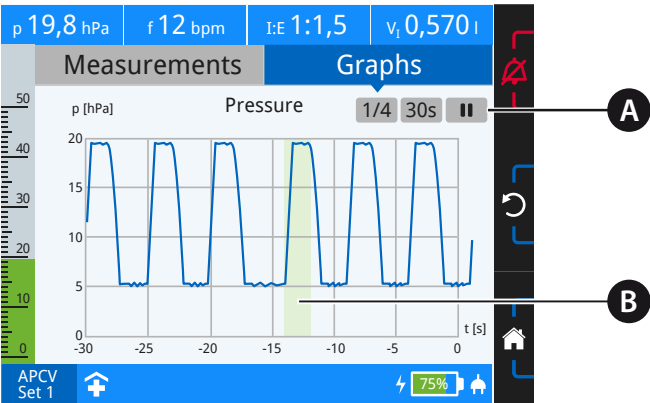

A Graph options B Inspiration triggered by patient

Figure 32: Graph screen

The graph options enable you to switch between graph views, modify scaling of the x-axis, and pause and start curve recording. You can also change the graph view and scaling of the x-axis when curve recording is stopped. Curve recording will continue in the background once you leave the graph screen.


Explanation of symbols in the graph options

Symbol	Meaning
1/4	Page selection The page selection shows how many curve views are available. You can move to the next curve view by turning the MFK.
30s	Setting the time The x-axis (t) can be scaled in 10-second steps from 10-60 seconds by turning the MFK.
⏸	Pause symbol Pressing the MFK stops curve recording.


Symbol	Meaning
	Play symbol Pressing the MFK starts curve recording.

How to change curve view:


1. Navigate to “Monitoring” on the home screen and press the MFK.
2. Navigate to “Graphs” and press the MFK.
The page selection will be marked in blue.

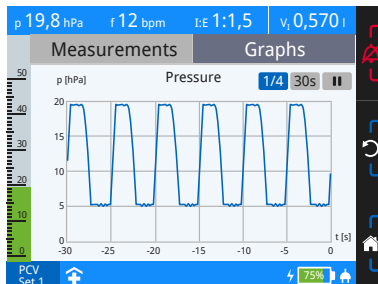

3. Press the MFK. The page selection is activated.

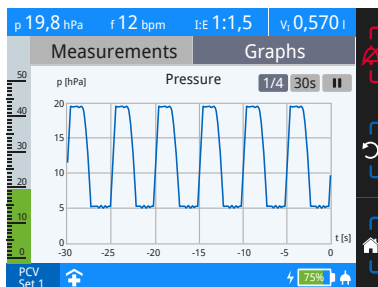


4. Turn the MFK to display further curves. Press the escape key to exit page selection.

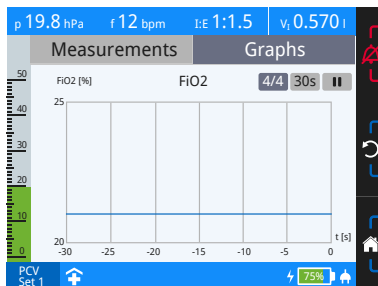


To change the scaling of the x-axis:


1. Navigate to “Monitoring” on the home screen and press the MFK.
2. Navigate to “Graphs” and press the MFK. The page selection will be marked in blue.

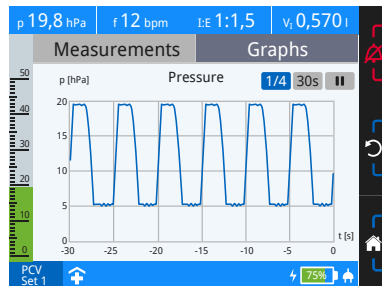


3. Turn the MFK once to the right. The time setting will be marked in blue.

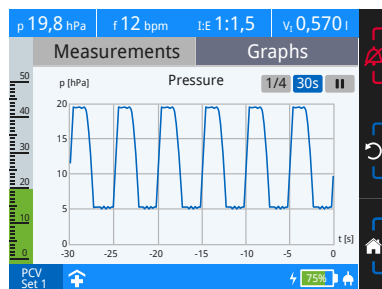

4. Press the MFK. The time setting is activated.

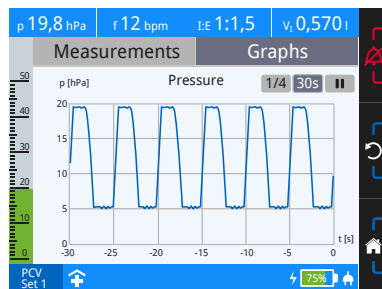


5. Turn the MFK to change the time settings. Press the MFK to confirm, or the escape key to discard.

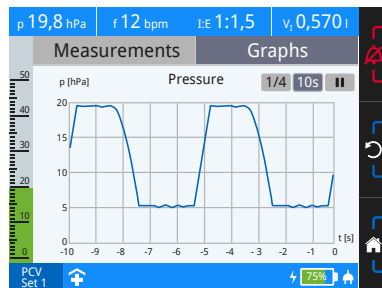


Stopping and starting curve recording:


1. Navigate to “Monitoring” on the home screen and press the MFK.
2. Navigate to “Graphs” and press the MFK. The page selection will be marked in blue.

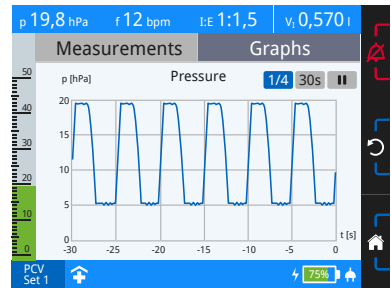

3. Turn the MFK twice to the right. The stop symbol will be marked in blue.

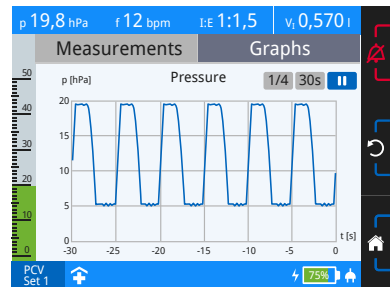

4. Press the MFK. Curve recording will be stopped and the play symbol will be displayed.
To continue curve recording, press the MFK again.

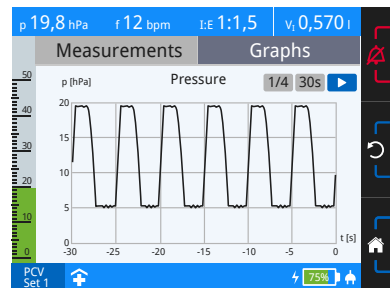


Parameter screen

You can set the ventilation and alarm parameters on the parameters screen. You can create up to 3 sets with different ventilation modes and parameters. Configure the number of sets on the system screen.

Please note: *Ventilation and alarm parameters can only be set in clinic mode.*

Active sets (green) cannot be disabled on the system to prevent accidental deletion of the configured sets.

The active set is the one that is started by pressing the on/off button. By contrast, the set being edited by the user is highlighted in blue. For example, set 1 can be edited (blue) while set 2 (green) is actively providing ventilation.

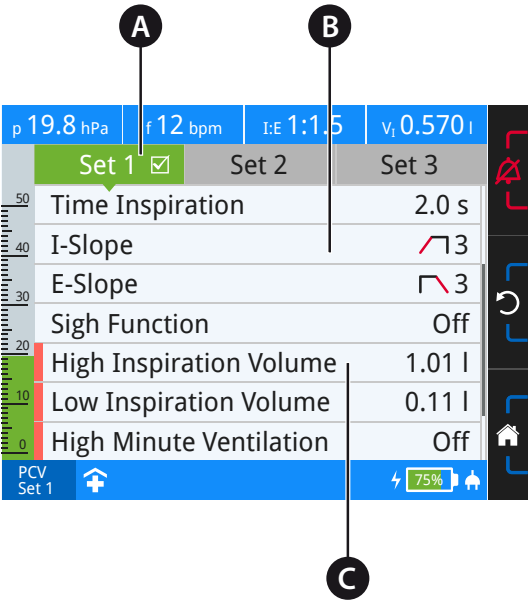

A Active set B Ventilation parameters C Alarm parameters

Figure 33: Parameter screen

Respiration parameters


Please note: *P-SIMV mode is only available in LAVI 40.*

Parameter Setting range	PCV	APCV	PSV	PSV-S	P-SIMV	CPAP	Description
Pressure 4 – 20 hPa						•	The pressure is the CPAP mode ventilation pressure, which is applied to the patient for each mechanical ventilation during inspiration and expiration.
IPAP LAVI 30: 7 – 30 hPa LAVI 40: 7 – 40 hPa	•	•	•	•	•		IPAP (= Inspiratory Positive Airway Pressure) is the ventilation pressure in the PCV-/ APCV, PSV-/PSV-S and P-SIMV mode, which is applied to the patient with each breath during inspiration. The set IPAP value is not summed up to the set PEEP, but it represents the maximum inspiratory pressure.
PS LAVI 40: 7 – 40 hPa						•	PS (= Pressure Support) is the pressure applied to the patient in the P-SIMV mode, which supports the patient in his/her own spontaneous inspiration. In the P-SIMV mode, PS is exclusively intended as pressure support of the patient's spontaneous inspiration. The set PS value is not summed up to the set PEEP, but it represents the maximum inspiratory pressure.
PEEP 4 – 20 hPa [PEEP ≤ IPAP-3 hPa]	•	•	•	•	•		The PEEP (= Positive End Expiratory Pressure) is the positive pressure which is available to the patient during expiration before starting the new inspiration. It can be spontaneous as well as mechanically controlled.

Parameter Setting range	PCV	APCV	PSV	PSV-S	P-SIMV	CPAP	Description
Frequency 4 – 50 bpm	•	•	•		•		In PCV controlled ventilation mode, the frequency is set as a defined specification by the machine. In PSV mode and in the assisted APCV mode, the set frequency is defined as the minimum frequency, which can be increased by patient spontaneous breathing. In P-SIMV mode, the frequency is defined as the frequency used to supply the mandatory breaths to the patient at the specified IPAP and during the specified inspiration time. Thus, the set frequency ensures the patient's minimum frequency. In between the mandatory breaths, the patient can increase the frequency by means of spontaneous inspiration.
Apnoea time 0 – 60 s			•				The apnoea time is intended to give the patient an additional period of time to spontaneously breath once the set frequency has lapsed. Apnoea time begins at the end of inspiration. If the apnoea time is exceeded, the device begins ventilation with the set frequency. Depending on the frequency and duration of the apnoea time, the minimum frequency may be exceeded. If the apnoea time ends before the set frequency has lapsed, the apnoea time is ignored. →Σεε αλσο page 77

Parameter Setting range	PCV	APCV	PSV	PSV-S	P-SIMV	CPAP	Description
Inspiration time 0.3 – 8 s	•	•	•		•		The inspiration time defines the duration of inspiration (in seconds). The frequency must be taken into account when setting the inspiration time. Setting up a fixed inspiration time, the I:E ratio is calculated depending on the respiratory frequency. Please note: <i>In the system settings, you can choose whether the inspiration time or the I:E ratio can be set.</i>
I:E 4.00: 1 – 1 : 4.00:	•	•	•		•		The I:E ratio is the ratio of inspiration to expiration of a breath. When setting up a fixed I:E ratio, the inspiration time depends on the set frequency. Please note: <i>In the system settings, you can choose whether the inspiration time or the I:E ratio can be set.</i>
Ti Max 1 – 10 s			•	•	•		The “Ti Max” setting limits inspiration time during spontaneous breathing, so that a late switchover to the expiration phase is prevented and a better synchronisation of patient and device can be achieved.
Ti Min off; 0.4 - 5 s			•	•	•		The “Ti Min” SETTING defines the minimum inspiration time during spontaneous breathing so that a premature switchover to the expiration phase is prevented.

Parameter Setting range		PCV	APCV	PSV	PSV-S	P-SIMV	CPAP	Description
I flank <i>Level 1 – 5</i>		•	•	•	•	•		The I flank determines the pressure increase during inspiration. When setting the I flank, the inspiration time/ I:E, the frequency and the difference between IPAP and PEEP must be considered, as unfavourable settings may prevent the IPAP from being achieved.
Flank setting	Pressure increase time ¹							
 1 (1 Pa/ms)	1.5 s ²							
 2 (2 Pa/ms)	0.75 s ²							
 3 (3 Pa/ms)	0.5 s ²							
 4 (5 Pa/ms)	0.3 s ²							
 5 (10 Pa/ms)	0.15 s ²							
1 Pressure increase time when IPAP = 20 hPa PEEP = 5 hPa 2 The values specified are reference times and depends on the pressure range set and the lung condition of the patient.								
E flank <i>Level 1 – 5</i>		•	•	•	•	•		The E flank determines the pressure decrease during expiration. When setting the E flank, the inspiration time/ I:E, the frequency and the difference between IPAP and PEEP must be considered, as unfavourable settings may prevent the PEEP from being achieved.
Flank setting	Pressure decrease time ¹							
 1 (1 Pa/ms)	1.7 s ²							
 2 (2 Pa/ms)	1.0 s ²							
 3 (3 Pa/ms)	0.6 s ²							
 4 (5 Pa/ms)	0.4 s ²							
 5 (10 Pa/ms)	0.3 s ²							
1 Pressure decrease time when IPAP = 20 hPa PEEP = 5 hPa 2 The values specified are reference times and depends on the pressure range set and the lung condition of the patient.								

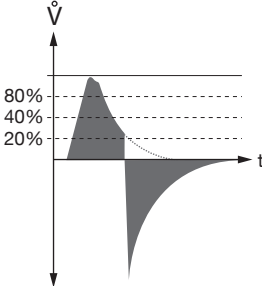Parameter Setting range			PCV	APCV	PSV	PSV-S	P-SIMV	CPAP	Description
Inspiration trigger <i>Level 1 – 10; auto</i>				•	•	•	•		The inspiration trigger is a volume trigger. It specifies the patient inspiration efforts required to obtain pressure support from the ventilator in case of spontaneous breathing. Inspiration is triggered when the inspiration volume and flow change in 30 ms are larger than the specified thresholds due to the patient's inspiration effort. Please note: <i>Please always consider the patient's clinical record when setting the trigger levels in order to avoid the risk of auto triggering.</i>
Thresholds:									
Level	Volume	Flow rate increase in 30 ms							
1	5 ml	1.2 l/min							
2	10 ml	1.5 l/min							
3	15 ml	1.8 l/min							
4	20 ml	2.1 l/min							
5	25 ml	2.4 l/min							
6	30 ml	2.7 l/min							
7	35 ml	3.0 l/min							
8	40 ml	3.3 l/min							
9	45 ml	3.6 l/min							
10	50 ml	3.9 l/min							
 1 = sensitive  10 = non-sensitive									
Expiration trigger <i>10 – 90%</i>					•	•	•		The expiration trigger is a flow trigger. The peak flow of the inspiration is measured with each breath. The expiration trigger setting defines the percentage of the peak flow at which the ventilator switches over to expiration.
 Expiration trigger threshold									

Parameter Setting range	PCV	APCV	PSV	PSV-S	P-SIMV	CPAP	Description
Trigger Lock <i>off; 0.5 - 4 s [$\leq 80\%$ of max. expiration time]</i>		•	•	•	•		The trigger lock is used specifically in the ventilation of patients with obstructive pulmonary diseases (e.g., COPD). In these patients, fluctuations often occur during the expiration phase. This leads to the device registering spontaneous respiration and the inspiration trigger is released too soon. To prevent false triggering and hyperventilation, it is possible to define a period of time (trigger lock) for the expiration phase, in which the inspiration trigger is suppressed.
Minimum Volume (Only LAVI 40) <i>off; 0.1 - 2 l [when IPAP > 37 hPa, always off]</i>	•	•	•	•			The minimum volume is the minimum tidal volume, which is considered as a volume guarantee during a pressure controlled ventilation.
Additional pressure (Only LAVI 40) <i>3 - 10 hPa [Additional pressure $\leq 40 \text{ hPa} - \text{IPAP}$]</i>	•	•	•	•			To ensure the minimum volume, the IPAP or PS is increased to the set additional pressure in gradual steps of max. 2 hPa per breath. The additional pressure value is a maximum value. To calculate the additional pressure actually required, the minimum volume is measured in relation to the actual volume and the necessary inspiration pressure is calculated from the current inspiration pressure. On the one hand, the pressure increase is limited by the set additional pressure as a maximum value, and on the other by a max. 2 hPa pressure increase per breath compared to the inspiration pressure of the previous breath.

Parameter Setting range	PCV	APCV	PSV	PSV-S	P-SIMV	CPAP	Description
Sigh function <i>on, off</i>	•	•	•	•	•		If the sigh function is activated, the IPAP is increased at each 100th breath in order to achieve approximately 150% of the measured volume of the previous breath. When in P-SIMV mode, if the 100th breath is a spontaneous pressure-supported breath, the next ventilation will be performed as a “sigh”.

Alarm parameters


Parameter <i>Setting range</i>	PCV	APCV	PSV	PSV-S	P-SIMV	CPAP	Description
Apnoea Alarm <i>off, 1 - 60 s</i>				•		•	The "Apnoea Alarm" setting establishes the time when the "Apnoea" alarm will be triggered in case of an apnoea event.
High frequency <i>off, 10 - 120 bpm</i>		•	•	•	•	•	If the measured frequency is higher than the "High Frequency" setting, the "High Frequency" alarm will be triggered.
Low frequency <i>off, 1 - 50 bpm</i>			•	•		•	If the measured frequency is lower than the "Low Frequency" setting, the "Frequency is too low" alarm will be triggered.
High Inspiration Volume <i>off, 0.2 - 2.5 l</i>	•	•	•	•	•		If the measured tidal volume is higher than the "High Inspiration Volume" setting, the "High Inspiration Volume" alarm will be triggered.
Low Inspiration Volume <i>off, 0.1 - 2 l</i>	•	•	•	•	•		If the measured tidal volume is lower than the "Low Inspiration Volume" setting, the "Low Inspiration Volume" alarm will be triggered.
High minute ventilation <i>off, 0.8 - 25 l</i>	•	•	•	•	•		If the measured minute volume be higher than the "High Minute Ventilation" setting, the "High Minute Ventilation" alarm will be triggered.
Low minute ventilation <i>off, 0.1 - 20 l</i>	•	•	•	•	•		If the measured minute volume is lower than the "Low Minute Ventilation" setting, the "Low Minute Ventilation" alarm will be triggered.
Leak rate <i>off; 20 - 200 l/min</i>	•	•	•	•	•		If the leakage measured is higher than the "Leak Rate" setting, the "Leak Rate" alarm will be triggered.
Max. FiO ₂ <i>off, 30 - 100%</i>	•	•	•	•	•		If the measured oxygen concentration is higher than the "High FiO ₂ " setting, the "High FiO ₂ " alarm will be triggered.

Parameter Setting range	PCV	APCV	PSV	PSV-S	P-SIMV	CPAP	Description
Min. FiO ₂ off, 18 – 90%	•	•	•	•	•	•	If the measured oxygen concentration is lower than the “Low FiO ₂ ” setting, the “Low FiO ₂ ” alarm will be triggered.
High pressure tolerance 1 – 10 hPa	•	•	•	•	•	•	If the measured pressure is higher than the set pressure plus the setting “High Pressure Tolerance”, then the alarm “High Pressure” will be triggered.
Low pressure tolerance 1 – 10 hPa	•	•	•	•	•	•	If the measured pressure is lower than the set pressure minus the setting “Low Pressure Tolerance” then the “Low Pressure” alarm will be triggered.


Changing the ventilator and alarm parameters

How to change a ventilation or alarm parameter:

- 1. On the home screen, navigate to “Parameters” and press the MFK.
- 2. Navigate to “Measurements” and press the MFK.

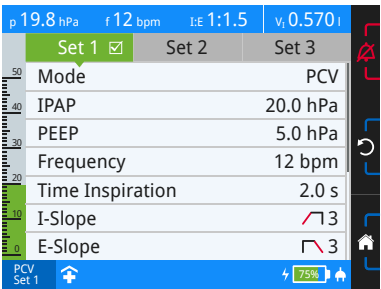


- 3. Navigate to the parameter you want to change and press the MFK.

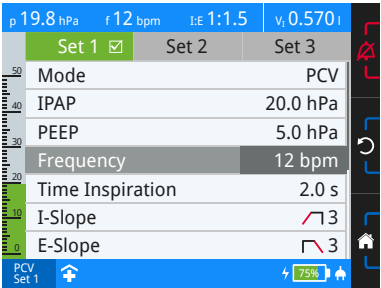

- 4. Change the setting and press the MFK to confirm.

Please note: *If the settings lead to changes in dependent parameters, these will be marked in a yellow-shaded section.*

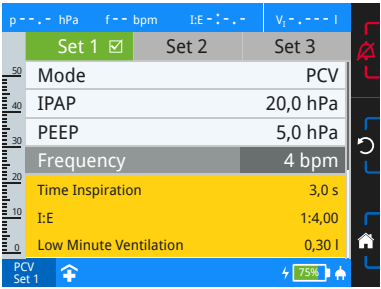

Activating another ventilation set



CAUTION Risk of injury due to activating the wrong ventilation set!
The ventilation sets can contain different ventilation and alarm parameters. Sets are not suitable for all applications.
⇒ Only activate the sets that the doctor has discussed with you.

To activate another ventilation set:


1. On the home screen, navigate to “Parameters” and press the MFK.
2. Navigate to the ventilation set you wish to activate.



3. Press the MFK twice.



4. Change the setting to “Yes” and press the MFK to confirm the setting.

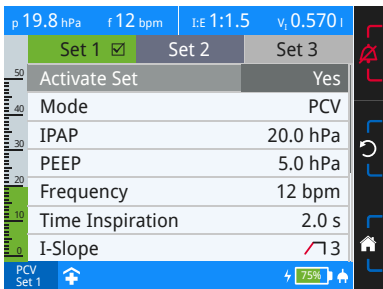

Changing a ventilation mode during ventilation

To change to ventilation mode for the active set during ventilation:

- 1. On the home screen, navigate to “Parameters” and press the MFK.
- 2. Select the active set (green) and press the MFK. The “Mode” parameter is selected.




- 3. Press the MFK.



- 4. Set the desired ventilation mode and press the MFK to confirm the settings.



- 5. Make sure the ventilation and alarm parameters are adapted to the patient. Scroll to the “Confirm Settings” at the end of the list and press the MFK.

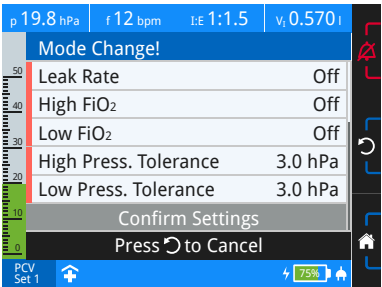

Event log screen


The event log screen is divided into three areas:

- Alarms 
- Alarms 
- Events

Alarms

Physiologically-dependant alarms whose alarm limits have been set in the parameter screen are stored under “Alarms ”.

The only technically-dependent alarms saved under “Alarms ”, by contrast, are those whose alarm conditions are implemented and built into the device. You can view the last 50 alarm events.

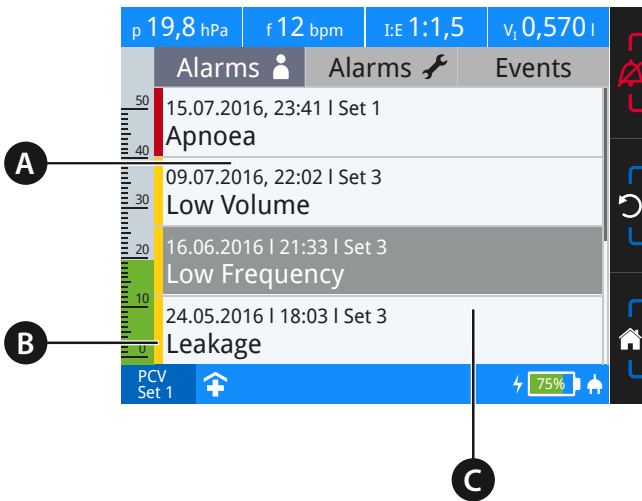


A Timestamp (date, time and set) B Alarm priority C Selected alarm

Figure 34: Event log screen (alarms)

For more information on the alarms, refer to chapter “Alarms and Messages” from page 119.

Events

Under “Events” you can see the last 50 parameter changes and events.

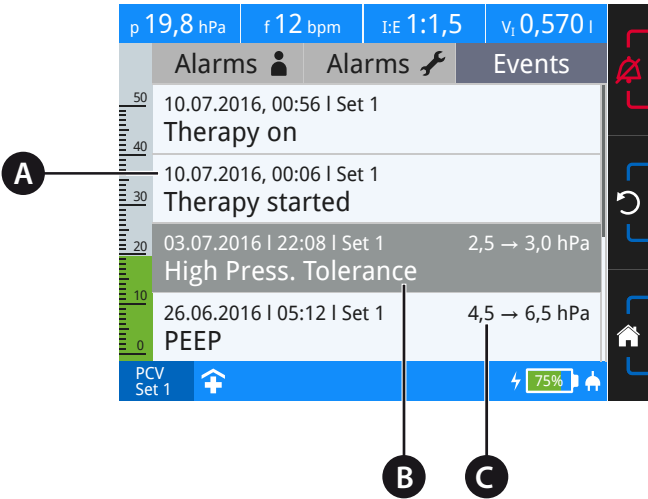

A Timestamp (date and time) B Selected event C Parameter change

Figure 35: Event log screen (events)

System screen

On the system screen, you can set the basic device settings, perform calibrations and view device information.

Please note: *Some settings cannot be changed in home mode.*

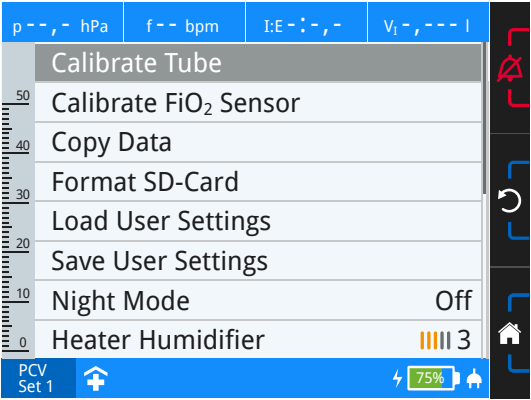


Figure 36: System screen

Parameter <i>Setting range</i>	Description	Can be changed (Home mode)
Calibrate Tube	Calibrate the connected tube circuit	✓
Calibrating the FiO ₂ sensor	Calibrating the FiO ₂ sensor	✓
Copy data	Copy data to the SD card	✓
Format the SD card	Formats the SD card. This will delete all data from the SD card.	✓
Load settings	Loads the last saved settings. <i>Please note: Settings cannot be loaded while respiration is in progress.</i>	–
Save settings	Saves all parameter- and system settings	–
Night mode <i>on, off</i>	Switch night mode on/off In night mode, the display colours are reduced and the brightness delay is lowered to 5%.	✓

Parameter <i>Setting range</i>	Description	Can be changed (Home mode)
Heater Humidifier 1 – 5	Set the heating level for the AquaTREND uni humidifier  Level 1 → lowest heating level Level 5 → highest heating level	✓
Sound volume 1 – 3	Volume of the primary sound  Level 1 → low volume Level 2 → medium volume Level 3 → high volume	–
Change the display mode <i>off, 20 s – 20 min</i>	Set when the device switches screens following the last operation action. Ventilation off: Switch to the home screen Ventilation on: Change to the measurements/graphs screen	–
Standby <i>off, 20 s – 20 min</i>	Set when the device switches to standby operation following the last operation action.	✓
Brightness display 1 – 3	Display brightness	✓
Info LED brightness 1 – 3	Info LED brightness	✓
Language	Set the device language	✓
Pressure unit <i>hPa, mbar, cmH₂O</i>	Pressure unit setting	–
Volume unit <i>l, ml</i>	Set the volume unit	–
Display time insp. <i>Seconds, I:E</i>	Set whether the inspiration time can be set in seconds or as an I:E ration on the parameter screen	–
Date	Set the date	✓
Time	Set the time	✓
Number of ventilation sets 1 – 3	Setting to establish how many ventilation sets are shown on the parameter screen	–


Parameter <i>Setting range</i>	Description	Can be changed (Home mode)
User profile <i>Clinic, home</i>	Setting the user profile  Clinic: full access to all settings  Home: restricted access to the settings (ventilation and alarm parameters cannot be changed). <i>Please note: For ventilation outside a clinic, the user mode must be set to "Home".</i>	✓
Recent Ventilation Hours	Ventilation hours since last reset	—
Ventilation total	Total hours of ventilation	—
Operating hours	Hours during which the device was connected to the mains or run in battery operation	—
Blower service in	Number of hours after which the blower must be replaced	—
SW version	Software version of the device	—
Serial number	Device serial number	—

Change the device language

To change the device language setting:

1. On the home screen, navigate to "System" and press the MFK.
2. Navigate to the 14th menu item with the MFK.

Then, press the MFK.

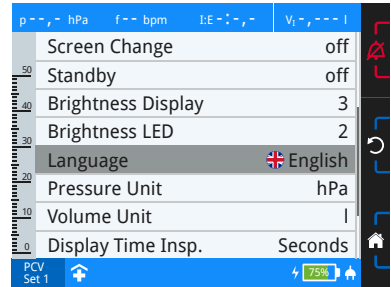

3. Change the language setting by turning the MFK, then confirm by pressing the MFK.

If your language is not set here, locate the desired language in the list below:

Language	German
Language	English
langue	Français
Lingua	Italiano
Dil	Türk


Resetting ventilation hours

LAVI features an hours counter, which you can reset to “0”. The hours counter records the period of ventilation.


To reset the ventilation hours:

Please note: *In home mode, the ventilation hours cannot be reset.*

1. On the home screen, navigate to “System” and press the MFK.
2. Navigate to “Recent Vent Hours” and press the MFK.

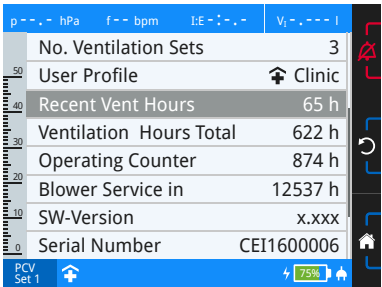

3. Select “Yes” and then press the MFK.

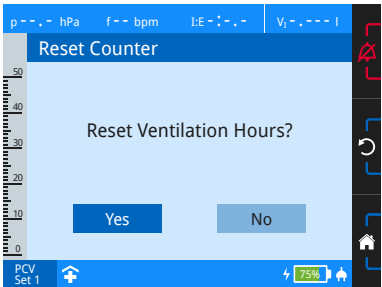


Statistics screen

The evaluation of the ventilation parameters is performed for each ventilation set based on percentiles. Percentiles are the dispersion measurement of the statistical data distribution during ventilation sessions.

Statistics are available for the following ventilation parameters:

- Minute Ventilation
- Frequency
- Tidal Volume
- I:E
- Leak Rates

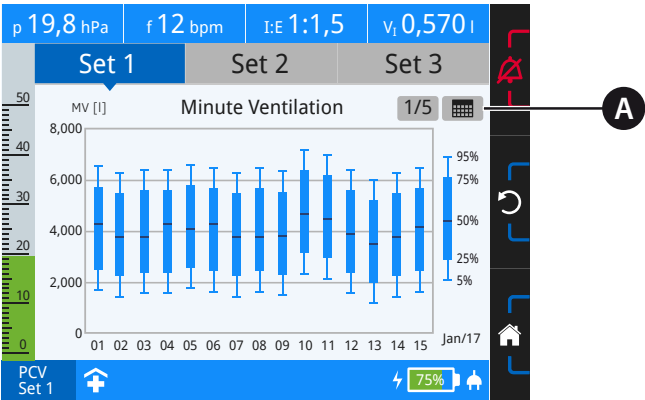

A Statistics options

Figure 37: Statistics screen


The statistics options enable you to switch between statistics views and the date of the x-axis.

Explanation of symbols in the statistics options


Symbol	Meaning
1/5	Page selection The page selection shows how many statistics views are available. You can move to the next statistic by turning the MFK.
	Set the date You can change the date (x-axis) by turning the MFK.

How to change statistics view:


- 1. Navigate to "Statistics" on the home screen and press the MFK.
- 2. Navigate to the desired set and press the MFK. The page selection will be marked in blue.

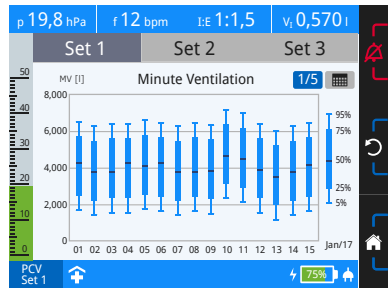

- 3. Press the MFK. The page selection is activated.

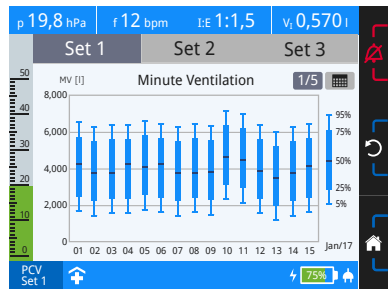


- 4. Turn the MFK to display further statistics. Press the escape key to exit page selection.

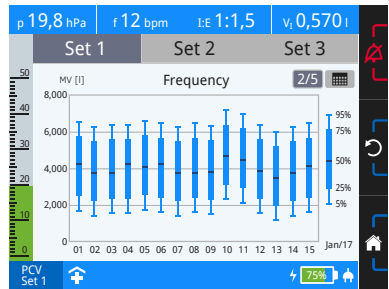


The last 15 days are always shown on the statistics screen. However, statistics are saved for the last six months. To show these statistics, proceed as follows.

To change the date of the x-axis:


1. Navigate to "Statistics" on the home screen and press the MFK.
2. Navigate to the desired set and press the MFK. The page selection will be marked in blue.

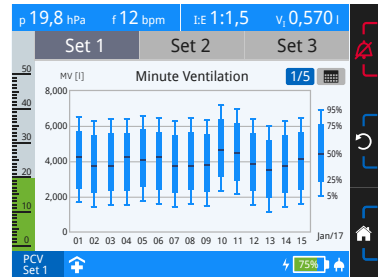


3. Turn the MFK once to the right. The date setting will be marked in blue.

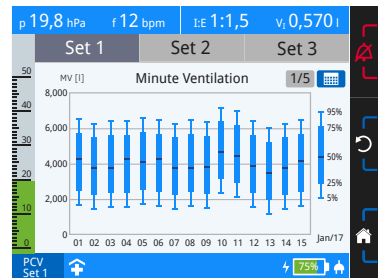


4. Press the MFK. The date setting is activated.

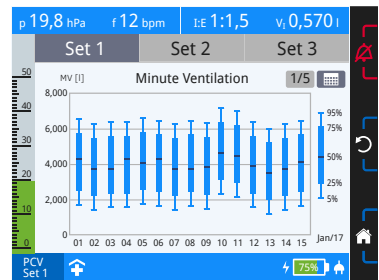

5. Turn the MFK to change the date setting. To exit date setting, press the escape key. This will not save the date setting.

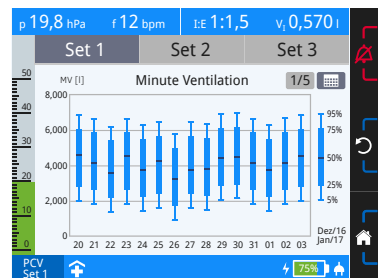

An example based on frequency

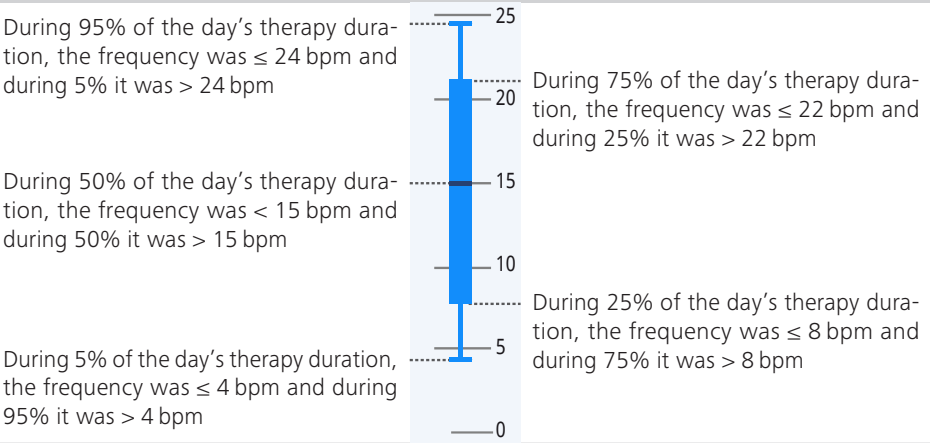

Service screen

On the service screen, authorised service companies can perform service work. Access is PIN-code protected.

Chapter 7

Alarms and Messages

General information

**CAUTION**

Risk of injury due to failure to recognise alarm!

Failure to recognise alarms can result in serious injuries for the patient.

⇒ Always operate the device so that the alarms are audible and visible by the user. Audible alarms can be forwarded using the nurse call or the alarm box.

Please note: *Alarms are switched off when the alarm event is no longer fulfilled and the priority has been lowered, but the alarm remains active until it has been confirmed with the alarm key.*

The LAVI ventilator is equipped with fixed and adjustable alarms, relating to the respective ventilation modes. There are three alarm priorities:


Alarm priority	What action is required?
HIGH	Immediate action is required. Monitor the patient and the cause of the alarm closely.
MEDIUM	Fast action is required for medium priority alarms. Correct the cause of the alarm.
LOW	User attention is requested for low-priority alarms. Low-priority alarms indicate a change at “normal” device operation. Check the cause of the alarm.

Alarm test

The user does not have to perform any alarm tests. However, the following describes an option to call the “Disconnection” alarm in order to perform a manual check of the alarm system.

To check the “Disconnection” alarm:

1. Connect the device to the mains supply.
2. Switch the device on. Do not connect the leakage tube system!
3. Start ventilation.
4. After a few seconds, the “Disconnection” alarm occurs. If this is not the case, the device must be returned for servicing.

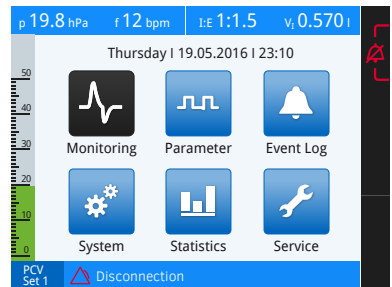

Audible alarm output (audio alarms)

Audio alarms are issued in a sequence of beeps. Alarm tones differ depending on alarm cause and priority. You will find more information from page 126.


If the alarm sound equipment is defective and emits no sound, the audible alarms will be triggered by a second alarm sound transmitter which emits only a simple audible alarm.

Please note: *Audio alarms are switched off when the alarm event is no longer fulfilled.*

Pausing audio alarms (muting audio alarms)

Audio alarms can be paused for two minutes by pressing the alarm button (mute audio alarm). If this function is enabled, the audible alarm created by new alarm events will also be suppressed as well. The info LED will visibly indicate an alarm event, even when the audible alarm has been temporarily suppressed. If the cause of the alarm is not corrected, the audible alarm will sound again after two minutes.

Please note: *Some alarms cannot be muted. You will find these in “Alarm overview” on page 126.*

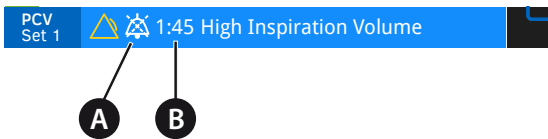

A “Audio alarm paused” symbol B “Audio alarm paused” counter

Figure 38: Alarm displays on the toolbar

The audio alarm may also be suppressed by pressing the alarm key even before an alarm event occurs, e.g. before the tube circuit is temporarily disconnected for suctioning the patient. The audio alarm can be reactivated after correcting the alarm cause, even within the two minute period, by pushing the alarm key again.

The “Audio alarm paused”  symbol indicates when the audio alarm is temporarily switched to mute. The counter tracks the time until the audio alarm will sound again.

Audio alarm volume

The audio alarm volume can be set in the system screen.

If an alarm is not confirmed by pressing the alarm key and the alarm condition remains in place, the volume levels 1 and 2 are automatically raised to a level 3 within 2 minutes.

Visible alarm output

Visible alarms are displayed as follows:

- via the info LED
- on the toolbar
- as a textbox

Alarm output via the info LED

The info LED may take on 3 different statuses, to signify the current alarm priority.

- Red, flashes rapidly (2 Hz) → high-priority alarm
- Yellow, flashes (0.5 Hz) → medium-priority alarm
- Turquoise, glows steadily → low-priority alarm


If multiple alarms are triggered simultaneously or in quick succession, the alarm with the highest priority will be displayed first. Still, all occurring alarms are listed in the event log.

Alarm output on the toolbar

Alarms are shown on the toolbar by the “Alarm active” icon and displayed with an alarm message. The icon’s colour indicates the alarm priority:

- Red symbol → high-priority alarm
- Yellow symbol → medium-priority alarm
- Turquoise symbol → low-priority alarm

If multiple alarms are triggered simultaneously or in quick succession, the alarm with the highest priority will be displayed first.

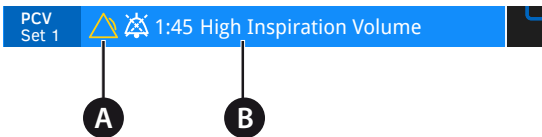

A “Alarm active” symbol B Alarm message

Figure 39: Alarm output on the toolbar


Alarm output in the textbox

120 seconds after the last performed operation the alarms will also display in a textbox as well. The textbox will disappear as soon as you press the alarm key.

The textbox colour corresponds to the highest priority alarm:

- Red textbox → high-priority alarm
- Yellow textbox → medium priority alarm
- Turquoise textbox → low-priority alarm

If multiple alarms occur at the same time the alarms are sorted and displayed in order of priority.

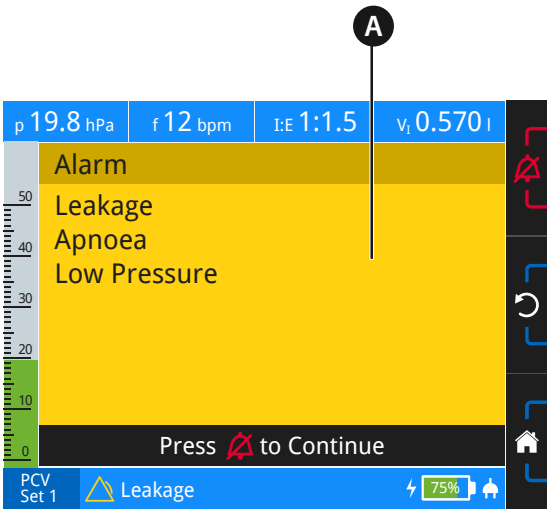

A Textbox

Figure 40: Alarm output in the textbox

Alarm log

Alarms are saved in a ring buffer. The oldest entry will be overwritten once the ring buffer is full. You can view the last 50 alarms in the event log screen under “Alarms”. All alarms saved in the ring buffer can be retrieved using the PC software “easySET”.

Alarms are permanently stored even during a complete power failure.

Forwarding alarms

Alarms can be forwarded by using a nurse call or the optionally available alarm box. This allows even better monitoring of the device at the home or clinic. The use of the remote alarm box or a nurse call is especially recommended when several ventilators are used in one room, as this allows the device generating the alarm to be easily identified. All alarms will be forwarded without delay to the nurse call or the alarm box.

Instructions on how to connect the HOFFRICHTER alarm box or nurse call can be found on page 59.



Figure 41: Alarm box

Please note: The alarm box is an optional accessory to facilitate remote output of the alarm. It does not replace monitoring of the ventilator’s primary alarm sound!

Alarm overview

Adjustable alarms

The adjustable alarms are physiologically conditional alarms. You can set the alarm limits in the parameter screen.

Alarm	Priority	Audible alarm	Info LED status	Cause	Time delay	 
Apnoea	HIGH	c a f – a f	Red flashing	Set time (“Apnoea Alarm”) has been exceeded	None	✓ -
FiO ₂ too high	MEDIUM	C b a	Yellow flashing	Measured FiO ₂ is greater than the set “Max. FiO ₂ ”	None	✓ -
FiO ₂ too low	MEDIUM	C b a	Yellow flashing	Measured FiO ₂ is less than the set “Min. FiO ₂ ”	None	✓ -
Leak rate	MEDIUM	c a f	Yellow flashing	Leakage measured greater than “Leak Rate”	3 breaths in a row	✓ -
High pressure	MEDIUM	c a f	Yellow flashing	Pressure is higher than the set “High Pressure Tolerance”	3 breaths in a row	✓ -
Low pressure	MEDIUM	c a f	Yellow flashing	Presser below set “Low Pressure Tolerance”	3 breaths in a row	✓ -
High frequency	MEDIUM	c a f	Yellow flashing	Measured frequency greater than “Max. frequency”	3 breaths in a row	✓ -
Low frequency	MEDIUM	c a f	Yellow flashing	Measured frequency less than “Min. frequency”	3 breaths in a row	✓ -
High volume	MEDIUM	c a f	Yellow flashing	Tidal volume greater than “Max. tidal volume”	3 breaths in a row	✓ -
Low volume	MEDIUM	c a f	Yellow flashing	Tidal volume is lower than the “Min. tidal volume”	3 breaths in a row	✓ -

Alarm	Priority	Audible alarm	Info LED status	Cause	Time delay		
High minute ventilation	MEDIUM	c a f	Yellow flashing	Minute ventilation greater than "Max. minute ventilation"	3 breaths in a row	✓	-
Low minute ventilation	MEDIUM	c a f	Yellow flashing	Minute ventilation is lower than the "Min. minute ventilation"	3 breaths in a row	✓	-

 = Muting
 = Can be deleted by pressing the MFK

Fixed alarms

Fixed alarms are alarms whose conditions are fixed in the device and thus cannot be changed by the user.

Alarm	Priority	Audible alarm	Info LED status	Cause	Correction			
Error flow sensor	HIGH	Ccc–Cc	Red flashing	Flow sensor defective	Device must be serviced	✓	-	-
Pressure sensor error	HIGH	Ccc–Cc	Red flashing	Pressure sensor defective	Device must be serviced	✓	-	-
Motor brake error	HIGH	Ccc–Cc	Red flashing	Motor brake defective	Device must be serviced	-	-	-
Motor controller error	HIGH	Ccc–Cc	Red flashing	e.g. over-current	Switch device on and off. If the error occurs again, return the device for servicing	✓	-	-

Alarm	Priority	Audible alarm	Info LED status	Cause	Correction			
No system flow	HIGH	Ccc–Cc	Red flashing	No leak in the ventilation system	Use an exhalation valve for the ventilation mask. Check the air outlet with a CPAP mask.	✓	-	-
Low leak rate	HIGH	c a f – a f	Red flashing	Leak rate < 50% than necessary leak rate	Check the exhalation valve for the ventilation mask. Check the air outlet with a CPAP mask.	✓	-	-
Motor speed low	HIGH	Ccc–Cc	Red flashing	Minimum motor speed exceeded	Switch device on and off. If the error occurs again, return the device for servicing	✓	-	-
Stenosis	HIGH	c a f – a f	Red flashing	No volume detected for 3 breaths	Check tube circuit and tubes for wear	✓	-	-
Ventilation off	HIGH	Ccc–Cc	Red flashing	Ventilation has been switched off due to one of the following alarms: – No system flow – Low leak rate – Motor speed low – Overcurrent – Overpressure	-	-	-	✓
Overpressure	HIGH	c a f – a f	Red flashing	Pressure switch triggered five times within 2 s	Restart device. If the error occurs again, return the device for servicing	✓	-	-
Overcurrent fuse	HIGH	c a f – a f	Red flashing	Excessive motor current	Restart device. If the error occurs again, return the device for servicing	✓	-	-

Alarm	Priority	Audible alarm	Info LED status	Cause	Correction			
Minimum volume not reached	HIGH	c a f – a f	Red flashing	Measured minimum volume less than set “Minimum Volume”	-	✓	-	-
Internal communication error	HIGH	c c c – c c	Red flashing	Communication to controller interrupted for more than 10 s	Switch device on and off. If the error occurs again, return the device for servicing	✓	-	-
Disconnection	HIGH	c a f – a f	Red flashing	Tube circuit not connected to device	Connect the tube circuit to the device	✓	-	-
Internal battery error	MEDIUM	c c c	Yellow flashing	Defective battery	Replace battery	-	✓	-
Internal Battery Empty	MEDIUM	C c c	Yellow flashing	Device operating on battery power supply, battery charge $\leq 5\%$	Battery must be charged; only 5 minutes until power supply fails completely; ventilation only possible with external power supply	-	-	-
Power failure	MEDIUM	C c c	Yellow flashing	Power supply from the mains (AC) connection has failed	Restore the power supply	-	-	✓
FiO ₂ sensor ¹ error	MEDIUM	C b a	Yellow flashing	FiO ₂ sensor has been separated or used up by device	Connect the FiO ₂ sensor to the device	-	✓	-
Low Internal Battery	LOW	e C	Glows steadily in turquoise	Device in battery operation, battery charge approx. 20% ²	Charge battery	✓	-	-

Alarm	Priority	Audible alarm	Info LED status	Cause	Correction			
Humidifier defective	LOW	e C	Glow steadily in turquoise	Humidifier defective	Humidifier must be serviced	-	-	✓

- 1 This alarm is only possible if at least one of the alarms “Max. FiO₂” or “Min. FiO₂” is activated.
- 2 Nominal battery value


 = Muting  = Inactive (can be disabled)  = Can be deleted by pressing the MFK

Fixed inactive alarms

The “Internal battery error” and “FiO₂ sensor error” alarms can be disabled by pressing the alarm button as long as the alarm condition still exists.


If one of these alarms appears on your device’s toolbar, proceed as follows:

1. First, the respective alarm is shown with the  symbol and the corresponding error text.

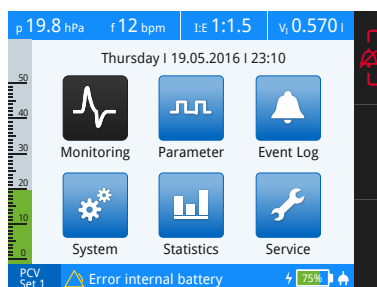

2. Press the alarm button to disable the alarm.

Please note: *After the alarm has been disabled by pressing the alarm button, it is still shown on the toolbar with the symbol  and the corresponding error message.*

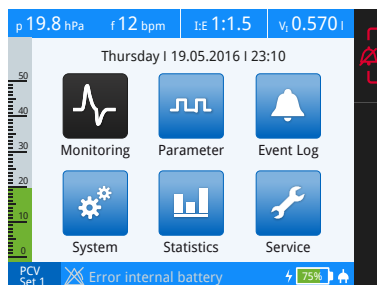


3. The alarm is only deleted after the error has been eliminated. For more information on eliminating the error, go to page 127.

Please note: *If the problem persists, please contact your service partner.*

Messages

Message output on the toolbar

Messages are shown on the toolbar. If an error occurs, the alarm is shown in place of the message, as the alarm has a higher priority.

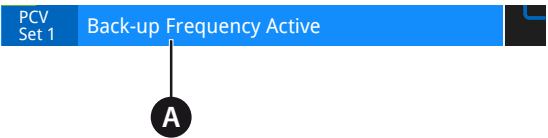

A Message

Figure 42: Messages on the toolbar

Message overview

Message	Cause	Time delay
Safety cycle active	The device is in PSV mode, the patient is not breathing spontaneously and is ventilated at the set frequency	None

Chapter 8

Cleaning and Disinfection

Important information

- Before cleaning the device, remove the power plug from the power supply.
- Cleaning the device must be performed according to this user’s manual and the applicable regulations of the hospital or nursing home.
- The device cannot be sterilised by using standard sterilisation methods. Hygienic preparation is described in the LAVI hygiene concept.
- Do not use any aggressive or abrasive cleaning agents (e. g., acetone).
- Do not immerse the device in water or solvents.
- Follow the accessory manufacturer’s instructions for cleaning and disinfection.

Overview

The following overview table describes the cleaning intervals of articles delivered by HOFFRICHTER. For articles by other manufacturers, please follow their cleaning instructions.

Component	Name	Cleaning	Disinfect	Replace
	LAVI ventilator	as needed	with every new patient	-
	Power supply unit	as needed	with every new patient	-
	Mains cable	as needed	with every new patient	-
	Leakage tube circuit	weekly	no	with every new patient in accordance with manufacturer instructions
	Mask	daily	no	with every new patient in accordance with manufacturer instructions

Component	Name	Cleaning	Disinfect	Replace
	Passive exhalation valve	weekly	with every new patient	with every new patient
	Oxygen connection adapter	as needed	no	with every new patient
	Carrying case	as needed	no	with every new patient
	Filter cassette (without filter)	as needed	with every new patient	-
	Coarse filter	weekly	no	instead of cleaning, with every new patient
	Fine filter	no	no	monthly, if severely contaminated, or with every new patient
	FiO ₂ sensor	as needed	no	in accordance with manufacturer instructions
	Bacterial filter	no	no	daily and with every new patient
	Functional bag	as needed ¹	no	with every new patient

¹ For cleaning information, please read the section “Cleaning the functional bag” on page 138.

Cleaning the device

Domestic use

For cleaning the surface of the device, use a cloth moistened with soapy water. Then wipe with a cloth moistened with clear water in order to remove any remaining of the soapy water. The device must be completely dry before commissioning.

Clinical use



WARNING

Risk of injury due to disinfectant agents!

Certain disinfectant agents can cause skin, eye or respiratory tract irritation.

⇒ Do not use disinfectant agents that cause irritation.



CAUTION

Risk of infection due to germs!

Devices afflicted with germs pose a risk of the germs spreading to other people.

⇒ Disinfect the device surface on a regular basis, or when there is any possibility of contamination.

⇒ Do not use disinfectant agents that cause irritation.

We recommend an alcohol-free surface disinfectant to disinfect the device surface, e.g. schülke mikrozyd sensitive wipes. Similar disinfectant wipes are also acceptable as well. The device must be completely dry before commissioning.

Cleaning the mask



CAUTION

Risk of injury due to damaged mask!

A heavily worn or damaged mask can result in insufficient ventilation and thus health issues.

⇒ Dispose of heavily worn or damaged masks and replace them with a new one.

Clean the mask according to the manufacturer's instructions.

Cleaning the tube circuit



CAUTION

Risk of injury due to damaged tube circuit!

A heavily worn or damaged tube circuit can result in insufficient ventilation and thus health issues.

⇒ Dispose of heavily worn or damaged tube circuits and replace them with a new one.

The supplied leakage tube circuit is intended for use on one patient only. It must not be used for other patients.

If the device is used by only one patient, the leakage tube can be cleaned. For reasons of hygiene, clean it once a week. To do so, proceed as follows:

1. Clean the tube with mild soapy water. Do not use any other agents!
2. Rinse the tube thoroughly with clear water.
3. Let the tube air-dry completely.

When using other tube circuits, the manufacturer's instructions must be observed. Tube circuits not designed for reuse must be disposed of properly.


Cleaning the exhalation valve

If the device is used by only one patient, the exhalation valve can be cleaned. For reasons of hygiene, clean the exhalation valve once a week. To do so, proceed as follows:

1. Clean the exhalation valve with mild soapy water. Do not use any other agents!
2. Rinse the exhalation valve thoroughly with clear water.
3. Let the exhalation valve air-dry completely.

When using other exhalation valves, the manufacturer's instructions must be observed. Exhalation valves not designed for reuse must be disposed of properly.

Cleaning/exchanging the air filter

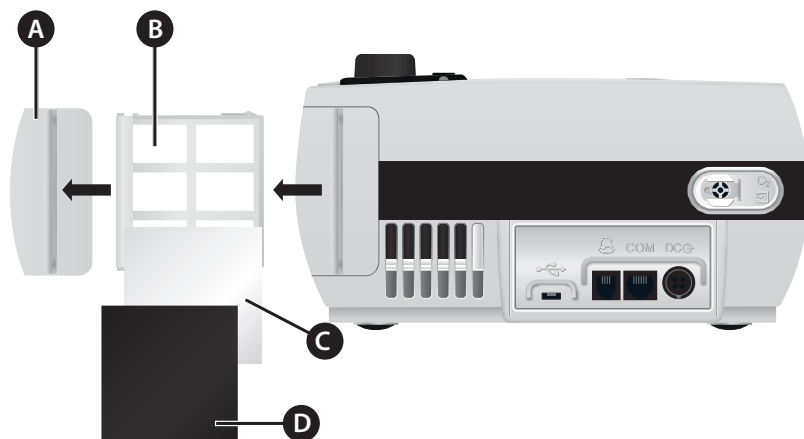

A Filter frame cover B Filter cassette C Fine filter (white) D Coarse filter (black)

Figure 43: Filter cassette structure

Cleaning the coarse filter

Clean the black course filter weekly.

1. Remove the filter cassette from the device (see Figure 43).
2. Remove the coarse filter (black) from the filter cassette.
3. Clean the filter with mild soapy water. Do not use any other agents!
4. Rinse the filter thoroughly with clear water.
5. Let the filter air-dry completely.
6. Insert the cleaned filter back into the filter cassette.
7. Slide the filter cassette into the device.

Instead of cleaning the filter, you can insert a new one or replace the entire filter cassette with a new one.

Cleaning the functional bag

Please note: *The functional bag is not suited for cleaning in a washing machine or dry-cleaning.*

To clean the functional bag, use a cloth moistened with water. A mild cleaning agent may also be used if necessary.

Replacing the fine filter

Replace the white fine filter monthly, or in case of heavy soiling.

1. Pull the filter cassette from the device.
2. Remove the coarse filter (black).
3. Remove the fine filter (white) and replace it with a new one.
4. Insert the coarse filter back into the cassette.
5. Slide the filter cassette into the device.

Replacing the filter cassette

1. Pull the filter cassette from the device.
2. Pull apart the filter cassette and the filter frame cover.
3. Reassemble the replacement cassette and the filter frame cover.
4. Slide the filter cassette into the device.

Using the device for more than one patient

If the device is intended for use by more than one patient (e.g., in operation in clinics), a suitable bacterial filter (e.g., Medisize Barr-vent S) should be used continuously to protect the device from contamination by human pathogens. The bacterial filter must be changed daily.



WARNING

Risk of infection due to germs!

- If the device has not been continuously protected by a bacterial filter, it may be contaminated with pathogenic germs. If the device is used by another patient, the patient may catch these germs.
 - ⇒ Protect the device continuously using a bacterial filter.
 - ⇒ Thoroughly clean and disinfect the device between use by different patients to ensure it is free of pathogenic germs.
 - ⇒ If in doubt, always assume that the device is contaminated, and hygienically recondition it according to the hygienic concept.
 - ⇒ If contamination with multi-resistant pathogens (e.g. MRSA) is suspected, the device must be packaged with the appropriate labelling and disinfected accordingly.
 - Accessories and consumable materials (e.g. tubes, mask etc.) can be contaminated with human pathogens. If they are used by another patient, the patient may catch these germs.
 - ⇒ Accessories and consumable materials labelled as “single use” (⊗) may only be used once for one patient.
 - ⇒ Accessories and consumable materials labelled as “single patient use” may only be used for one patient and must be cleaned in accordance with the manufacturer’s instructions.
-

Please note: *Document device preparation in the medical device book.*

Before using the device for another patient you must complete the following procedures:

Component	 What action is required?
Bacterial filter/HME filter	Replace
Coarse and fine filter	Replace
Exhalation valve	Replace or recondition according to the manufacturer's guidelines if possible
Mask	Replace or recondition according to the manufacturer's guidelines if possible
Tube circuit	Replace or recondition according to the manufacturer's guidelines if possible
Humidifier	Recondition according to the manufacturer's guidelines
Device	Disinfect the device surfaces with a disposable germicidal wipe. We recommend an alcohol-free surface disinfectant, e.g. schülke mikrocid sensitive wipes. Make sure that no liquid is drawn into the device through the openings. Let the disinfectant soak in according to the disinfectant manufacturer's instructions. Wipe any disinfectant agent residues from the device using a disposable wipe.

Chapter 9

Routine Tests and Maintenance Work

Please note: You must not perform any testing or maintenance work while the patient is still connected to the device. Provide an alternative ventilation system for the patient during that time.

The device has an expected service life of min. 5 years provided the maintenance work listed in the following overview is performed by an authorised service company and the device is used according to the user’s manual.

Overview

 When necessary?	 What action is required?	 By whom?
Before commissioning	Safety-related test (see page 145)	Operator
Weekly	Clean/replace the coarse filter (see page 138)	User
	Visually check the fine filter	User
Monthly, or before, if heavily contaminated	Replace fine filter (see page 139)	User
Every 12 months during storage	Charge battery ¹	User
Every 2 years	Maintenance 2 (refer to the service manual)	Provider/Service
	Safety-related test (see page 145)	Provider/Service
- After 15,000 h blower run time or - Every 5 years	Maintenance 5 (refer to the service manual)	Provider/Service

¹ Batteries must be charged after no more than 1 month.

Safety-related test (SRT)

To maintain and check the device functions, the device must be subjected to safety-related tests every 2 years as well as on commissioning, carried out by an authorised service technician.

The safety-related check comprises:

- A visual check for outside damage of the ventilator,
- Functional test of the alarms
- A visual inspection of the accessories according to manufacturer's instructions (tube circuit, FiO₂ sensor etc.)

All measures to be performed for the safety-related check are described in the LAVI service manual.

Battery maintenance

The LAVI internal battery is a powerful lithium-ion battery. To reach the full capacity of the battery, it is important to charge it on a regular basis.

If the "Battery maintenance" message appears, you should charge the battery to 100%, discharge it and then charge it to 100% again to recalibrate the battery indicator.

Please note: Regularly check the battery power level. Operate the device on the mains to ensure the battery is always charged.

As the number of charging cycles of lithium-ion batteries is limited, the internal battery must be replaced and disposed of after a certain period of time. This is the case when the "Replace battery" message appears.

Charging during storage

Charge the battery at least once every 12 months.

Batteries must be charged after no more than 1 month.

Chapter 10

Annex

Notes for home therapy

You should take the following into consideration before discharging a patient requiring at-home ventilation.

1. Limit the number of sets to the number that the patient should use.
2. Adjust all parameters and settings for home ventilation.
3. Finally, switch the user mode to “home”.
4. Train the patient or caretaker in the use of this device. In addition, explain which accessory equipment may be used and how to maintain and care for the device (e.g. battery care, cleaning the device).
5. For oxygen use, the patient must be instructed on how to use oxygen on site in their home environment.
6. Discuss how the patient should act, particularly if he or she is in danger or if an alarm sounds. The patient should be given a telephone number to call in case of emergency.

Device technical specifications

Main performance characteristics

The main performance characteristics of LAVI include:

- Pressure-controlled, mechanical support with a pressure precision of $\pm$ (2% of the full scale + 8% of the effective measurement),
- Alarm when the inspiratory pressure is exceeded and distribution by the remote alarm,
- Alarm when the “Min. frequency” alarm parameter is fallen below and distribution by the remote alarm,
- Alarm when the “Min. FiO₂” alarm parameter is fallen below and distribution by the remote alarm.

Please note: The main performance characteristics may be affected by electromagnetic interference (see section “Manufacturer’s declaration on electromagnetic compatibility” on page 174).

Important components

The LAVI ventilator comprises the following components:

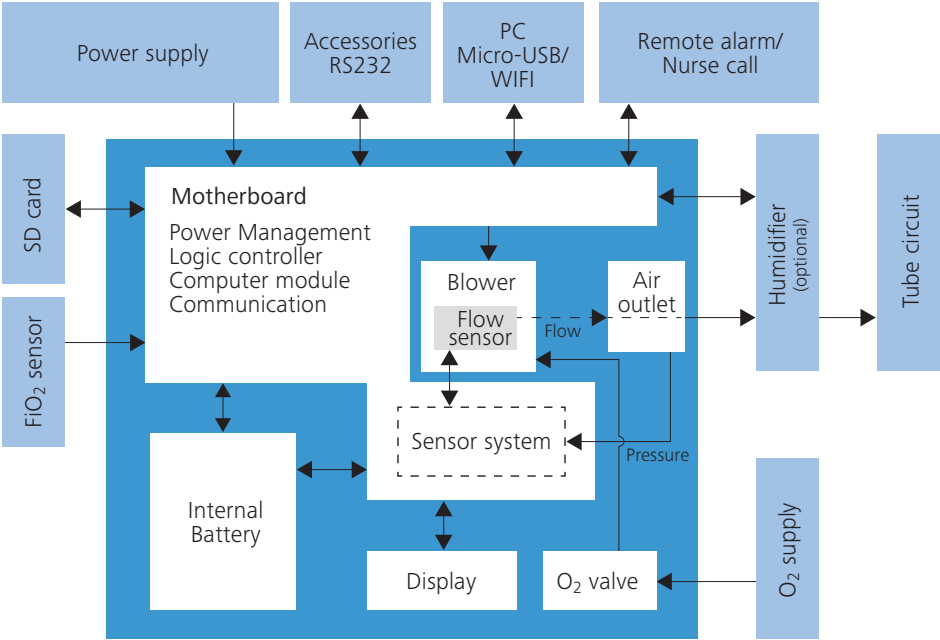


Figure 44: Block diagram for the device

Pneumatic block

The pneumatic block is the connector unit for the leakage tube circuit and consists of the following components:

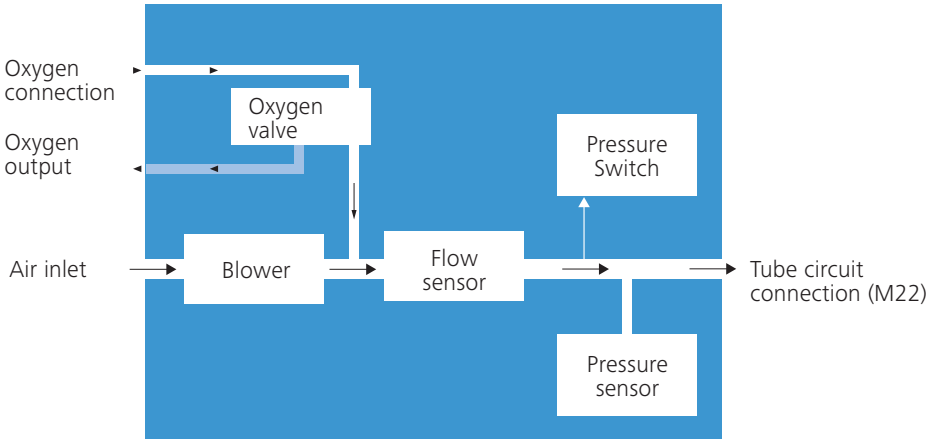


Figure 45: Pneumatic block diagram

Data management

The device has an internal memory to record data.

The following data will be saved:

- Alarms and events with date and time stamp (approx. 15000 entries)
- Statistics
- Device settings and counter
- Update files
- Initialisation files
- Measurement parameters: pressure, volume, flow, FiO_2 (approx. 50 days with 20 values recorded per second)


To be able to pass on therapy data to a doctor, you can copy them to an SD card (see page 60).

Error messages

The following error messages may arise while booting the device:

Error message	 Cause	 Solution
Flash not working	No access to the flash	Device cannot be booted and must be serviced
Default parameters are on the device	No valid parameter set available or faulty	The device uses the default parameters (factory settings) and can continue to be utilised
Calibration data corrupt	Sensor calibration data corrupt	Device can continue to be used, perform flow sensor calibration
Event log file faulty	Reading of event data failed	Device can continue to be used but must be serviced
Primary Alarm Not Working	Primary alarm sound transmitter is defective	Device can continue to be used but must be serviced
Secondary alarm sound transmitter is defective	Secondary alarm sound transmitter is defective	Device can continue to be used but must be serviced
No alarm sound transmitter available	Primary and secondary alarm sound transmitter are defective	Device can continue to be used but must be serviced
RYB LED error	Info LED defective	Device can continue to be used but must be serviced
Real-time clock error	Real-time clock defective or clock battery low/depleted	Device can continue to be used but must be serviced <i>Please note: Date recording contains incorrect date and time information.</i>
Card reader error	Card reader defective	Device can continue to be used but must be serviced
Boot error + error number	Device booting failed due to a critical error	Device cannot be booted and must be serviced
Blower servicing required	Maximum blower run time reached	Device can continue to be used but must be serviced

After booting, any errors in effect will be displayed on the toolbar via the error symbol. To show the error list detailing existing errors, press the error key on the home screen.

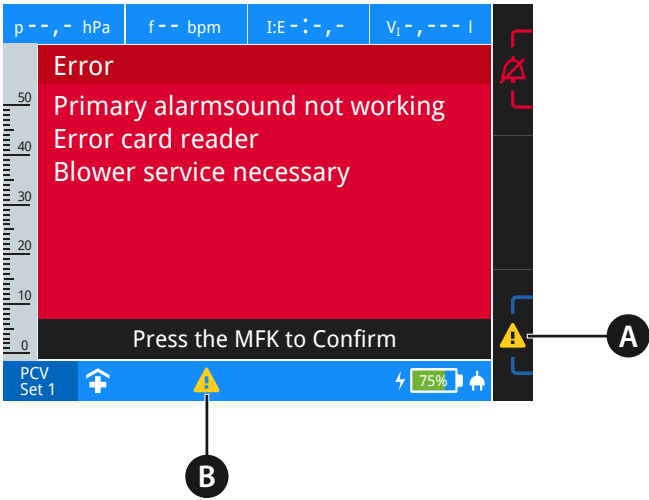

A Error button B Error symbol

Figure 46: Error list

Technical data

Supply voltage	
Mains operation	100 to 240 V AC (-20%, +10%), 50 to 60 Hz
DC operation	24 V DC / 3.75 A
Internal battery operation	Lithium-ion battery, 15 V (nominal voltage) / 3.2 Ah
Max. power consumption	75 W
Electrical protection class	Protection class II

Specifications and performance characteristics	
Dimensions (W x D x H)	215 x 203 x 115 mm
Weight	2.2 kg
Max. stable limit pressure	60 hPa
Min. stable limit pressure	0 hPa
Max. inspiratory working pressure ¹	LAVI 30: 30 hPa LAVI 40: 40 hPa
Min. inspiratory working pressure ¹	4 hPa
Max. flow at 40 hPa	With humidifier ² : 174 l/min Without humidifier: 188 l/min
Max. flow at 30 hPa	With humidifier ² : 195 l/min Without humidifier: 210 l/min
Max. flow at 4 hPa	With humidifier ² : 250 l/min Without humidifier: 267 l/min

Operating conditions	
Temperature range	+5°C to +40°C
Relative humidity	10% to 95%, non-condensing
Air pressure range ³	600 hPa to 1100 hPa (approx. 4000 to -400 m)

¹ Ensured by pressure control.

² AquaTREND uni

³ The air flow rate decreases with increasing altitude.

Storage and transport conditions		
Temperature range		
With battery	- 20°C to + 50°C	
Without battery	- 25°C to + 70°C	
Relative humidity	5% to 90%, non-condensing	
Air pressure range ¹	250 hPa to 1100 hPa (approx. 10000 to -400 m)	
Storage conditions	Store in a dry, vibration-free place, in an upright position; store device and accessories in their original packaging	
Sound pressure range of audible alarm signal (at 1 m distance)		
Minimum value	50 dBA, level 1	
Average value	55 dBA, level 2	
Maximum value	65 dBA, level 3	
Resistance the first time an error occurs		
Pressure at the patient connection opening	at 30 l/min	at 60 l/min
	0.97 hPa	4.49 hPa
Measuring conditions	Device with leakage tube circuit (art. no. 00007875) without additional accessories	
Specifications for WLAN module PAN9020U ²		
WLAN standard	IEEE 802.11/IEEE802.11b, IEEE 802.11g, IEEE 802.11n	
RF frequency range	2400 MHz	
Bandwidth	20, 40 MHz	
Max. data rate	65 Mps	
Transmission power	+ 8 dBm	
Max. reception power	-98 dBm	

¹ The air flow rate decreases with increasing altitude.

² Currently not in use

Technical requirements for accessories	
Oxygen inlet	
Connection type	Quick-connect coupling
Pressure	≤ 500 hPa
Flow	≤ 15 l/min
Bacterial filter	
Connections	22/15 mm cone (according to EN1281-1)
Resistance	< 2.3 hPa at 60 l/min
Compressible volume	< 66 ml
Internal volume	< 200 ml

The manufacturer reserves the right to make technical changes without notice.

Measured values

Parameter	Display area	Display increments	Measurement	Accuracy
Pressure ¹	0.0 – 60.0 hPa	0.1 hPa	0.0 – 60 hPa	1.5% of the measured value ³
Pressure bar ¹	0 – 50 hPa	1 Pha	0.0 – 60 hPa	1.5% of the measured value ³
Inspiration volume ²	0.000 – 9.999 l 10.00 – 99.99 l	0.001 l 0.01 l	calculated from flow measurements	0.03 l or 20% of the measured value ³
Oxygen ¹	0 – 100%	1%	0 – 100%	5%
Frequency ¹	0 – 99 bpm	1 bpm	calculated from period duration of inspiration + expiration in 0.001 s	1 bpm
I:E ¹	9.9:1 – 1:9.9 99:1 – 1:99	0.1 1	calculated from period duration of inspiration + expiration in 0.001 s	0.2
MV ²	0.000 – 9.999 l 10.00 – 99.99 l	0.1 l 0.001 l 0.01 l	calculated from flow measurements	0.03 l or 20% of the measured value ³
Leakage ¹	0 – 999 l/min	1 l/min	calculated from flow measurements	n/a

All flow and volume values were determined at 25°C (77°F) and 1030 hPa.

1

Process for determining the display value: Form average with rounding

2

Process for determining the display value: Rounding

3

The higher value applies.

Setting ranges and control precision

Respiration parameters

Parameter	Setting range	Setting increments	Accuracy
Pressure	4 – 20 hPa	0.5 hPa	1.0 hPa or 5% of the measured value ¹
IPAP	LAVI 30: 7 – 30 hPa LAVI 40: 7 – 40 hPa	0.5 hPa	1.0 hPa or 5% of the measured value ¹
PS	LAVI 30: 7 – 30 hPa LAVI 40: 7 – 40 hPa	0.5 hPa	1.0 hPa or 5% of the measured value ¹
PEEP	4 – 20 hPa <i>[PEEP ≤ IPAP-3 hPa]</i>	0.5 hPa	1.0 hPa or 5% of the measured value ¹
Frequency	4 – 50 bpm	1 bpm	1 bpm
Apnoea time	0 – 60 s	1 s	1 s
Inspiration time	0.3 – 8 s	0.1 s	0.1 s
Ti Max	1 – 10 s	0.1 s	0.1 s
Ti Min	off; 0.4 - 5 s	0.1 s	0.1 s
I:E	4.00:1 – 1:4.00	0.05	0.05
I flank	Level 1 – 5	1 level	-
E flank	Level 1 – 5	1 level	-
Inspiration trigger	Level 1 – 10; auto	1 level	-
Expiration trigger	10 – 90%	1%	1%
Trigger Lock	off, 0.5 - 4 s <i>[80 % of the max. expiration time]</i>	0.1 s	0.1 s
Minimum Volume (Only LAVI 40)	off; 0.1 – 2 l <i>[when IPAP > 37 hPa, always off]</i>	0.01 l	0.03 l or 20% of the measured value ¹
Additional pressure (Only LAVI 40)	3 – 10 hPa <i>[Additional pressure ≤ 40 hPa – IPAP]</i>	0.5 hPa	1.0 hPa or 5% of the measured value ¹
Sigh function	on, off	-	-

¹ The higher value applies.

Alarm parameters

Parameter	Setting range	Setting increments	Accuracy
Apnoea Alarm	off, 1 - 60 s	1 s	1 s
High frequency	off, 10 – 120 bpm	1 bpm	1 bpm
Low frequency	off, 1 – 50 bpm	1 bpm	1 bpm
High Inspiration Volume	off, 0.2 – 2.5 l	0.01 l	0.03 l or 20% of the measured value ¹
Low Inspiration Volume	off, 0.1 – 2 l	0.01 l	0.03 l or 20% of the measured value ¹
High minute ventilation	off, 0.8 – 25 l	0.01 l	0.03 l or 20% of the measured value ¹
Low minute ventilation	off, 0.1 – 20 l	0.01 l	0.03 l or 20% of the measured value ¹
Leak rate	on, off	-	-
High FiO2	off, 30 – 100%	1%	1%
Min. FiO ₂	off, 18 – 90%	1%	1%
High pressure tolerance	1 – 10 hPa	0.5 hPa	1.0 hPa or 5% of the measured value ¹
Low pressure tolerance	1 – 10 hPa	0.5 hPa	1.0 hPa or 5% of the measured value ¹

¹ The higher value applies.

Factory settings

Respiration parameters

Parameter	Mode		
	PCV (set 1)	PCV (set 2)	CPAP (set 3)
Pressure	-	-	4 hPa
IPAP	20 hPa	20 hPa	-
PEEP	5 hPa	5 hPa	-
Frequency	12 bpm	12 bpm	-
Inspiration time	2 s	3 s	-
Ti Max	-	4 s	-
Ti Min	-	0.4 s	-
Apnoea time	-	10 s	-
I flank	3	3	-
E flank	3	3	-
Inspiration trigger	-	3	-
Expiration trigger	-	30%	-
Trigger Lock	-	Off	-
Minimum volume (only LAVI 40)	Off	Off	-
Additional pressure (Only LAVI 40)	3 hPa	3 hPa	-
Sigh function	Off	Off	-

Alarm parameters

Parameter	Mode		
	PCV (set 1)	PCV (set 2)	CPAP (set 3)
Apnoea Alarm	-	-	10 s
High frequency	-	30 bpm	30 bpm
Low frequency	-	Off	Off
High Inspiration Volume	1 l	1 l	1 l
Low Inspiration Volume	0.2 l	0.2 l	0.2 l
High minute ventilation	Off	Off	Off
Low minute ventilation	Off	Off	Off
Leak rate	Off	Off	Off
Max. FiO ₂	Off	Off	Off
Min. FiO ₂	Off	Off	Off
High pressure tolerance	3 hPa	3 hPa	3 hPa
Low pressure tolerance	3 hPa	3 hPa	3 hPa

System parameters

Parameter	Factory setting	Setting range	Setting increments
Night mode	Off	on, off	-
Humidifier heating	3	1 – 5	1
Change the display mode	2 min	off, 20 s – 20 min	20 s – 100 s: 20 s 2 min – 20 min: 1 min
Standby	5 min	off, 20 s – 20 min	20 s – 100 s: 20 s 2 min – 20 min: 1 min
Sound volume	3	1 – 3	1
Brightness display	3	1 – 3	1
Info LED brightness	2	1 – 3	1
Language	English	German, English and others	-
Pressure unit	hPa	hPa, mbar, cmH ₂ O	-
Volume unit	l	l, ml	-
Display time insp.	Seconds	Seconds, I:E	-
Number of ventilation sets	3	1 – 3	1
User profile	Clinic	Clinic, home	-

Compliance with standards

The following standards are observed:



- **DIN EN 60601-1:2013-12**
Medical electrical equipment - Part 1: General requirements for basic safety and essential performance (IEC 60601-1:2005 + Cor. :2006 + Cor. :2007 + A1:2012); German version EN 60601-1:2006 + Cor. :2010 + A1:2013
- **DIN EN 60601-1-2:2016-05**
Medical electrical equipment - Part 1-2: General requirements for basic safety and essential performance - Collateral Standard: Electromagnetic disturbances - Requirements and tests (IEC 60601-1-2:2014); German version EN 60601-1-2:2015
- **IEC 60601-1-8:2006+A1:2012**
Medical electrical equipment - Part 1-8: General requirements for basic safety and essential performance - Collateral Standard: General requirements, tests and guidance for alarm systems in medical electrical equipment and medical electrical systems
- **DIN EN 60601-1-11:2016-04**
Medical electrical equipment - Part 1-11: General requirements for basic safety and essential performance - Collateral Standard: Requirements for medical electrical equipment and medical electrical systems used in the home healthcare environment (IEC 60601-1-11:2015); German version EN 60601-1-11:2015
- **DIN EN ISO 10651-6:2011-06**
Lung ventilators for medical use - Particular requirements for basic safety and essential performance - Part 6: Home-care ventilatory support devices (ISO 10651-6:2004); German version EN ISO 10651-6:2009

Replacement parts and accessories

Please note: *Make sure to follow all general safety guidelines when using replacement parts and accessories from page 22.*

For ordering of replacement parts and accessories, please contact a HOFFRICHTER service partner.


Replacement parts

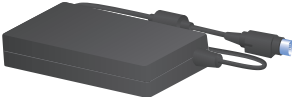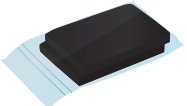Name Article number	Figure
Leakage tube (L= 180 cm, Ø 22 mm) 00007875	
Straight O ₂ connection adapter 41000104	
Mains cable 31100015	
Switched-mode power supply FSP090-RACM1 00013041	
Power supply unit holder 42101364	
Filter cassette (open) with filters 00002038	
Coarse filter, 1 pack (2 ea) 00014950	
Fine filter, 1 pack (5 ea) 00014951	

Name	Figure
Article number	
LAVI user's manual for physicians and medical professionals 50000681	
LAVI user's manual for patients 50000682	
LAVI brief instruction 50000683	
Carrying case 00002437	

Accessories

Name	Figure				
Article number					
Mask	Size XS	Size S	Size M	Size L	Size XL
With exhalation valve (vented)					
Standard nose mask	00003497	00003440	00003434	00003435	
Standard mouth/nose mask		00003441	00003436	00003437	
Cirri Mini Comfort nose mask		00003551	00003452	00003498	
Cirri Comfort nose mask		00003486	00003487	00003488	
Cirri Comfort mouth/nose mask		00003483	00003484	00003485	
Nasal Pillow 4in1			00003499		
Without exhalation valve (non-vented)					
Standard NIPPV mouth/nose mask		00003461	00003442	00003438	0003462
Standard NIPPV mouth/nose mask, autoclavable				00003439	
Cirri Comfort mouth/nose mask NIPPV		00003489	00003490	00003491	
ComfortTube system consisting of: Heated leakage tube, power supply unit, power switch 00003479 					
Humidifier AquaTREND uni 00002073 					
Passive exhalation valve 00014157 					

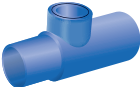Name Article number	Figure
Bacterial filter 00004932	
O ₂ connection adapter angled 41000087	
FiO ₂ measurement set consisting of: FiO ₂ sensor, T adapter, FiO ₂ sensor adapter, FiO ₂ sensor connecting line with screw connector 00004944	
FiO ₂ sensor 23000018	
T adapter 23000019	
FiO ₂ sensor adapter 23000020	
FiO ₂ sensor connecting line with screw connector 00014116	
Remote alarm box, complete including accessories 00014122	

Name Article number	Figure
Remote alarm box without accessories 00004834	
Cable for remote alarm box 00014115	
Cable for nurse call 00014117	
Battery RRC2054 00013079	
LAVI functional bag 00013082	

Internal battery

Delivery condition of replacement batteries

On delivery, the battery is in ship mode (LED status display disabled, no voltage detected at the connector). Starting the charging process activates the battery.

Safety warnings

- Do not open or take apart the battery.
- Do not subject the battery to heat or flames. Avoid storing the battery in direct sunlight.
- Do not short-circuit the battery.
- Do not store the unpacked battery in containers that could cause short-circuits with loose metallic objects.
- Do not remove the battery from its original package until you use it.
- Avoid mechanical impact.
- If liquids leak from the battery, avoid contact with your skin or eyes. Should you come into contact with it, rinse the affected area with sufficient water and consult a doctor.
- Pay attention to the plus (+) and minus (-) signs when using the battery in order to ensure proper operation.
- Do not use batteries of different manufacturers, capacities, sizes or types within one device.
- Keep the battery out of the reach of children.
- Keep the battery clean and dry.
- Only use the battery for its intended application.
- Do not store the battery in a charged state for more than 1 month.
- Do not store the battery for more than 1 year without charging it.
- The battery must be recycled or disposed of properly.

Battery description



Figure 47: Lithium-ion battery RRC2054

- 1 LED status indicator
The battery power level is shown here.
- 2 Push-button
Press to show the battery power level.

LED status indicator

The battery shows the power level when you press the push button. Each LED segment represents 25%. The LEDs light up for 4 seconds after you press the push button. If the battery voltage is too low or the battery is not operational (permanent error), the LEDs will not illuminate.

Battery power level	LED status indicator	Comment
< 10%		Flashes
10%-25%		Lights up for 4 seconds
26%-50%		Lights up for 4 seconds
51%-75%		Lights up for 4 seconds
76%-100%		Lights up for 4 seconds

Technical data

The manufacturer reserves the right to make technical changes without notice.

General	
Type	RRC2054
Voltage	15 V
Capacity	3.20 Ah
Maximum charging current	2.24 A
Maximum charging voltage	17.4 V
Max. discharge current	4 A
Dimensions (W x D x H)	77.6 x 85.4 x 23 mm
Weight	240 g

Operating temperature	
for charging	0°C to + 40°C
for discharging	- 10°C to + 55°C

Storage temperature	
< 1 year	- 20°C to + 20°C
< 3 months	- 20°C to + 45°C
< 1 month	- 20°C to + 60°C

Please note: Store the battery at temperatures under 20°C and at a low humidity. Avoid dust and corrosive gas atmospheres. Store the battery at a power level between 50 – 70%.

Symbols on the battery

Symbol	Meaning
	Observe the warning and safety instructions in the user's manual.
	Follow the operating instructions
	This symbol indicates conformity with applicable EC directives.
	Symbol indicating compliance with regulations of Australia and New Zealand
	UL recognised for Canada and the US
	Taiwan recycling symbol
	UN transport test
 Li-ion	Recycling symbol
	Do not dispose of the device with the household waste. Please contact the relevant customer services department to find out how to properly dispose of the device.
	For Canada and the US: Call 1-800-822-8837 for recycling information.
	China RoHS
	Symbol indicating compliance with regulations of Japan

Declaration of Conformity

CE: The battery meets the EC directives valid at the time of manufacture.

FCC: This product has adhered to limit values defined in section 15 of FCC regulations. Operation is subject to the following two conditions: (1) This device must not cause any harmful interference, and (2) this device must accept any received interference, including interferences that can cause undesired behaviour.

Exchanging the battery

ATTENTION

Material damage due to incorrect battery!


Using an incorrect battery can damage the device, or the power supply to the device cannot be guaranteed.

⇒ Only use battery RRC2054 approved by HOFFRICHTER.

Please note:

- *No patient may be connected to the ventilator when replacing the battery. Provide an alternative ventilation system for the patient during that time.*
- *After installation, charge the battery by running the device using the mains power.*
- *Replacement batteries can be ordered as accessories (see page 167).*

To exchange the internal battery:

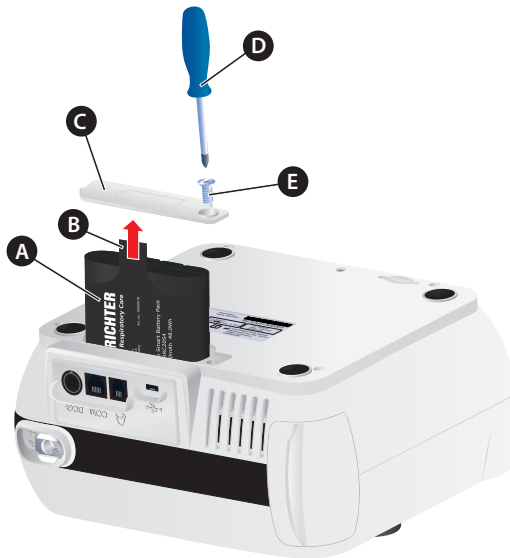

A Battery B Strap C Battery compartment cover D Phillips screwdriver PZ2
E Countersunk bolt

Figure 48: Exchanging the battery

1. Disconnect all connected components from the device, such as the mains cable, humidifier, tube circuit etc.
2. Turn the device over and carefully place the top side on a soft surface.

-
3. Unscrew the screw on the battery compartment cover using a Phillips screwdriver PZ2.

ATTENTION: Other screwdrivers could damage the screw!

4. Take off the battery compartment cover.
 5. Remove the battery out of the strap.
 6. Slide the new battery into the compartment.
 7. Put the battery compartment cover back on and tighten the screw.
-

Manufacturer’s declaration on electromagnetic compatibility

LAVI complies with standard DIN EN 60601-1-2:2016-05 and is intended for use in the electromagnetic environment described in the following. Any deviating ambient conditions can affect its main performance characteristics (e.g. pressure accuracy or alarms) or cause the device to fail.

Guidance and manufacturer’s declaration - electromagnetic emissions

The LAVI ventilator is intended for operation in environments as described below. The user of the LAVI ventilator must ensure it is operated in such an environment.

Emitted interference	Compliance	Electromagnetic environment - guidance
HF emissions according to CISPR 11	Group 1	The LAVI ventilator only uses HF energy for its internal functions. Its HF emissions are thus very low and interference with nearby electronic devices is unlikely.
HF emissions according to CISPR 11	Class B	The LAVI ventilator is suited for use in all facilities including those connected directly to a public supply network that also supplies buildings used for residential purposes.
Emission of harmonic oscillations as per IEC 61000-3-2	Class A	
Emission of voltage fluctuations/flickers as per IEC 61000-3-3	In compliance	

Guidance and manufacturer's declaration - electromagnetic immunity

The LAVI ventilator is intended for operation in environments as described below. The user of the LAVI ventilator must ensure it is operated in such an environment.

Immunity tests	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Electrostatic discharge (ESD) as per IEC 61000-4-2	± 8 kV contact discharge ± 2, 4, 8, 15 kV air discharge	± 8 kV contact discharge ± 2, 4, 8, 15 kV air discharge ± 2, 4, 8, 15 kV air discharge display	Floors should be made of wood or concrete, or be lined with ceramic tiles. If the floor is lined with a synthetic material, the relative humidity must be at least 30%.
Emitted HF disturbance as per IEC 61000-4-3	10 V/m 80 MHz – 2.7 GHz 80% AM at 1 kHz	10 V/m	Portable and mobile HF communication devices (including wires and antenna) should be used at least 0.3 m (suggested safety distance) to the LAVI ventilator.
	27 V/m 385 MHz PM: 18 Hz	27 V/m	
	28 V/m 450 MHz FM ± 5 Hz: 1kHz sine	28 V/m	
	9 V/m 710, 745, 780 MHz PM:217Hz	9 V/m	
	28 V/m 810, 870, 930 MHz PM: 18 Hz	28 V/m	
	28 V/m 1720, 1845, 1970 MHz PM: 217 Hz	28 V/m	
	28 V/m 2450 MHz PM:217Hz	28 V/m	
	9 V/m 5240, 5500, 5785 PM:217Hz	9 V/m	

Guidance and manufacturer's declaration - electromagnetic immunity

Immunity tests	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Rapid transient disturbances/bursts as per IEC 61000-4-4	± 2 kV for power cables ± 1 kV for input and output cables	± 2 kV for power cables ± 1 kV for input and output cables	The quality of the supply voltage should correspond to one of the typical business or hospital environments.
Surges/ as per IEC 61000-4-5	± 1 kV voltage outer conductor-outer conductor ± 2 kV voltage outer conductor-earth	± 1 kV voltage outer conductor-outer conductor ± 2 kV voltage outer conductor-earth	The quality of the supply voltage should correspond to one of the typical business or hospital environments.
Conducted HF disturbances as per IEC 61000-4-6	$3 V_{\text{Effective value}}$ 150 kHz – 80 MHz $6 V_{\text{Effective value}}$ in ISM and amateur radio bands between 150 kHz and 80 MHz	$3 V$ $6 V$	Portable and mobile radio devices should be used spaced at least at the suggested safety distance to the LAVI ventilator, including wires, which is calculated using the equation for transmission frequency. Suggested safety distance: 0.3 m
Magnetic fields at supply frequency (50/60 Hz) as per IEC 61000-4-8	30 A/m	30 A/m	Magnetic fields at supply frequency should correspond to typical values for business and hospital environments.

Guidance and manufacturer's declaration - electromagnetic immunity

Immunity tests	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Voltage dips, temporary interruptions and voltage fluctuations in the voltage supply lines as per IEC 61000-4-11	>95 % dip of U_T for 1/2 period	>95 % dip of U_T for 1/2 period	The quality of the supply voltage should correspond to one of the typical business or hospital environments. If the user of the LAVI ventilator requires continued functioning even in the event of an interruption to the power supply, we recommend supplying the LAVI ventilator using an uninterruptible power supply (UPS) or a battery.
	>95 % dip of U_T for 1 period	>60 % dip of U_T for 5 periods	
	>30 % dip of U_T for 25 (50 Hz) periods / 30 (60 Hz) periods	>30 % dip of U_T for 25 periods	
	>95 % dip of U_T for 250 (50 Hz) periods / 300 (60 Hz) periods	>95 % dip of U_T for 5s or 6s	

Disposal

Proper disposal saves natural resources and prevents harmful substances being released into the environment.

Device



The device must not be disposed of with the household waste. Please contact the relevant customer services department to find out how to dispose of the device etc. properly.

Batteries



Replaced batteries must be disposed in accordance with the respective local laws. Please contact the relevant customer services department to find out how to dispose of the device etc. properly.

Batteries must not be disposed of when discharged. If a battery is not fully charged, there is a risk of a short circuit. Short circuits can be prevented by insulating the contacts with adhesive strips.

Packaging



The packaging is taken back by the distributor but it can alternatively be disposed of separately with normal household waste.

FiO₂ sensor



The FiO₂ sensor must not be disposed of with the household waste. Please contact the relevant customer services department to find out how to dispose of the device etc. properly.

Disclaimer

HOFFRICHTER GmbH accepts no liability for consequences in terms of safety, reliability and performance of the product if:

- Interventions, modifications, extensions, calibration, repairs and maintenance are carried out by persons not authorised by us,
- Other manufacturers' accessories and spare parts are used that have not been approved by us for use on the product,
- The product is used for purposes other than stipulated in the user's manual or
- The hygiene and cleaning instructions stipulated in the user's manual have not been complied with.

Statutory guarantee rights remain unaffected by this disclaimer.

