Annex 3

Manual

Analysis of cyanobacterial toxins

Workshop

Analysis of cyanobacterial toxins

23.09. – 25.09.2008

National Institute of Chemical Physics and Biophysics,
Tallinn (Estonia)

Twinning Light Project EE06-IB-TWP-ESC-01 "Reducing
the health risks from algal toxins in drinking and bathing
waters"

Compendium compiled by:

Jussi Meriluoto	Åbo Akademi University, Turku (Finland)
Lisa Spoo	Åbo Akademi University, Turku (Finland)
Pia Vesterkvist	Åbo Akademi University, Turku (Finland)

Project manager and responsible for translation into Estonian:

Jüri Ruut	Health Protection Inspectorate, Tallinn (Estonia)
-----------	---

Project leader in Germany:

Wido Schmidt	Water Technology Center (TZW) Karlsruhe, branch Dresden (Germany)
--------------	--

Expert:

Risto Tanner	National Institute of Chemical Physics and Biophysics, Tallinn (Estonia)
--------------	---

Küllli Kangur (Estonian University of Life Sciences, Tartu, Estonia), Geoffrey A. Codd (University of Dundee, UK) and James S. Metcalf (University of Dundee, UK) are gratefully acknowledged for their contributions.

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ELISA in microcystin and cylindrospermopsin detection – principles and advise

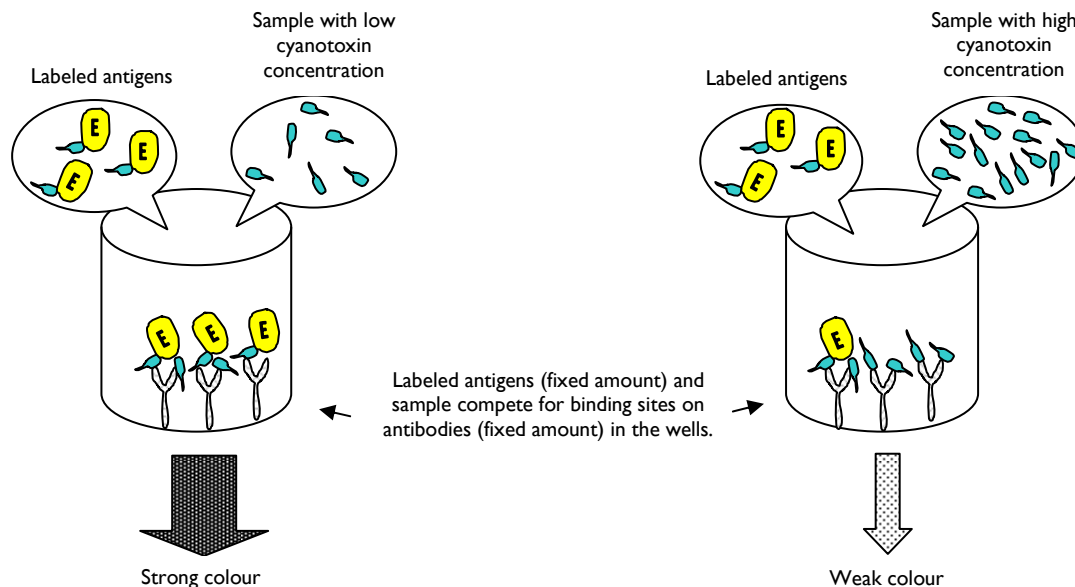
Pia Vesterkvist, Åbo Akademi University, Turku, Finland

Since cyanotoxins are threatening human and animal health, very tiny amounts in the environment should be detected. A method that do not usually need any sample pre-handling, e.g. concentration steps, are immunoassays, where specific antibodies identify the toxins, the antigens, and gives rise to a detectable signal from where quantitative results can be made. Immunoassays, like ELISA, allow water samples to be used directly and toxin concentrations of 0.05-0.1 µg microcystin or cylindrospermopsin per liter water can be detected.

Antibodies used in ELISA kits today are polyclonal or monoclonal. For research and diagnostic use, and therapeutic purposes, monoclonal antibodies are usually preferred. However, their production is much more laborious compared to polyclonal antibodies. Monoclonal antibodies originate from a single clone and are specific for one epitope, antigen-binding site, only. For example, Fischer et al., 2001 developed monoclonal antibodies specific to Adda, a unique amino acid in microcystins. Polyclonal antibodies comprise of a heterogenous mixture of different antibodies, all with their own epitope specific site, but this might increase the risk to get false positive results. Molecules related to antibodies, for example parts of an antibody, can also be used in immunoassays.

The specificity of the antibody-antigen reaction depends largely on cross reactivity, i.e. antigens with slightly different chemical structure, including toxin conjugates, will also bind to the antibodies. Since microcystins constitutes of a large family of similar structures, cross reactions can be considered as a benefit. Some microcystins can be recognized with high affinity and others with low affinity so the overall *toxicity* of the sample can usually not be correctly evaluated with ELISA. Neither can determination of individual microcystin variants be done with existing ELISA kits. That is the reason why the amount of microcystin is given as microcystin equivalents, and more specific microcystin-LR equivalents.

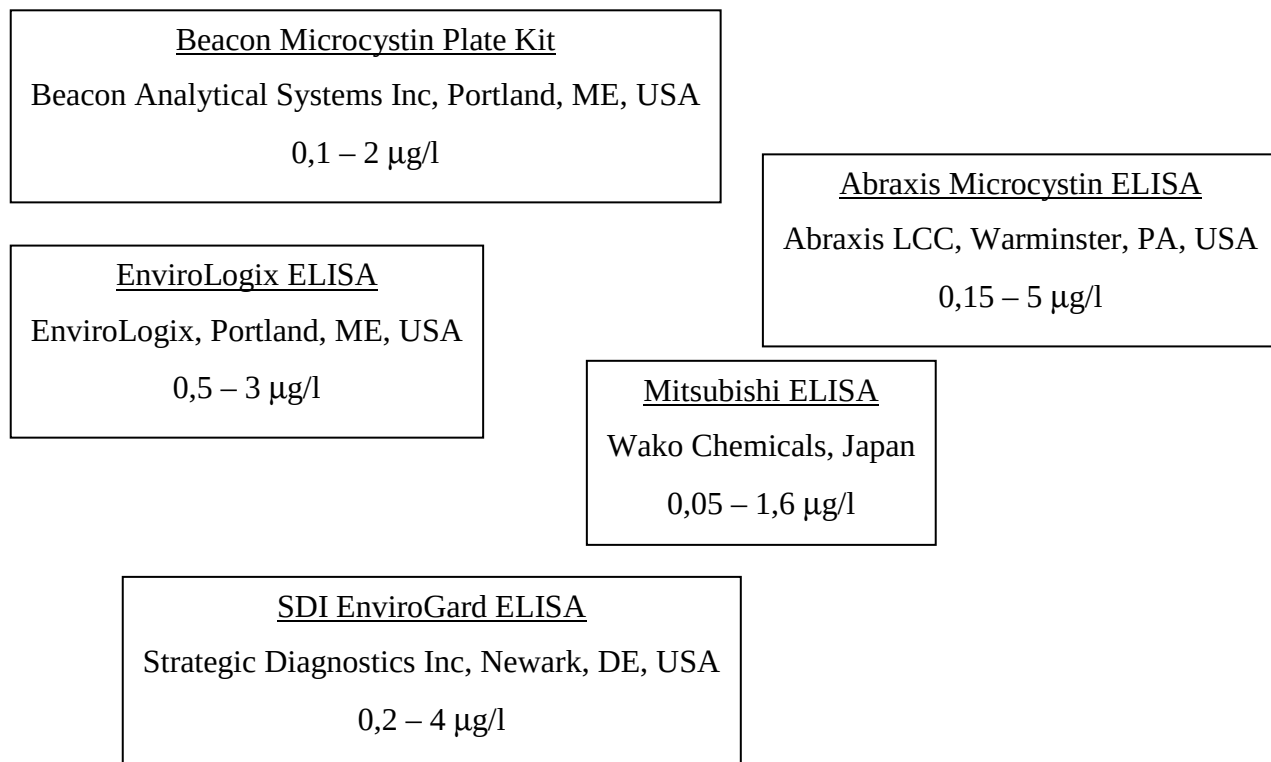
Detection of the signal in ELISA is accomplished with an enzyme (Horse Radish Peroxidase, HRP) that together with a substrate, gives rise to a coloured product. Depending on the type of assay either the antigen or the antibody is labeled with the enzyme. Four of the five commercially available microcystin ELISA kits are direct competitive immunoassays. In direct competitive immunoassays antibodies are immobilized in the microtiter wells (see figure). The sample and the labeled antigen are added to microtiter wells. These compete with each other for antibody binding sites. Washing removes non-specific compounds after which a substrate is added that reacts with the enzyme-labeled antigen providing a detectable product. The coloured product is therefore inversely proportional to the concentration of microcystin in the sample. A sample containing large amounts of microcystins or cylindrospermopsin out-compete the labeled antigen, and give rise to a less intense color since a major part of the labeled antigens are washed away. To put it short: less colour means more toxin.

Quantitative determination of microcystins and cylindrospermopsin in immunoassays is made according to a standard curve of known concentrations. The standard curve has to be done every time samples are examined for toxins since the reaction between antibodies and antigens is critical. Therefore each step when performing the assay has to be done with care, e.g. pipeting, mixing, incubation and washing.

The great advantage of ELISA is, nevertheless, that the whole assay is easy to perform, no special knowledge or training in immunology is required. It is highly sensitive and particularly suitable for screening purposes. However, the user should be aware of the methods susceptibilities. Methanol, often used for extraction of cyanobacterial cells and a common solvent for microcystins, can at concentrations higher than 5% interfere and give rise to false positive microcystin concentrations. High salinity is also troublesome and gives rise to weak false positive concentrations. Matrix effects are quite usual and give false positive results. Some plastics, like polypropen, absorbs microcystins (not cylindrospermopsin) dissolved in water and the detected concentration is lower than it should be. Certain kits are also prone to give a lower microcystin concentration with increased pH, no data exist regarding cylindrospermopsin kits and high pH.

The five commercially available microcystin ELISA kits are:

Assay procedure for Beacon Microcystin Plate Kit – in short

1. Add 50 µl of enzyme conjugate solution into each well (if possible, use a multi-channel pipette)
2. Add 50 µl of standard solutions, control and samples into the wells as duplicates
3. Add 50 µl antibody solution into each well (multi-channel pipette)
4. Cover the wells, mix the contents, incubate for 30 min
5. Wash 5 times, 250 µl/well
6. Add 100 µl substrate solution (multi-channel pipette)
7. Cover the wells, incubate 30 min
8. Add 100 µl stop solution (multi-channel pipette)
9. Read the absorbance at 450 nm within 15 min

**The following texts have been translated into Estonian language from
the publication**

**TOXIC : Cyanobacterial Monitoring and
Cyanotoxin Analysis**

ISBN 951-765-259-3

(c) Åbo Akademi University Press, 2005

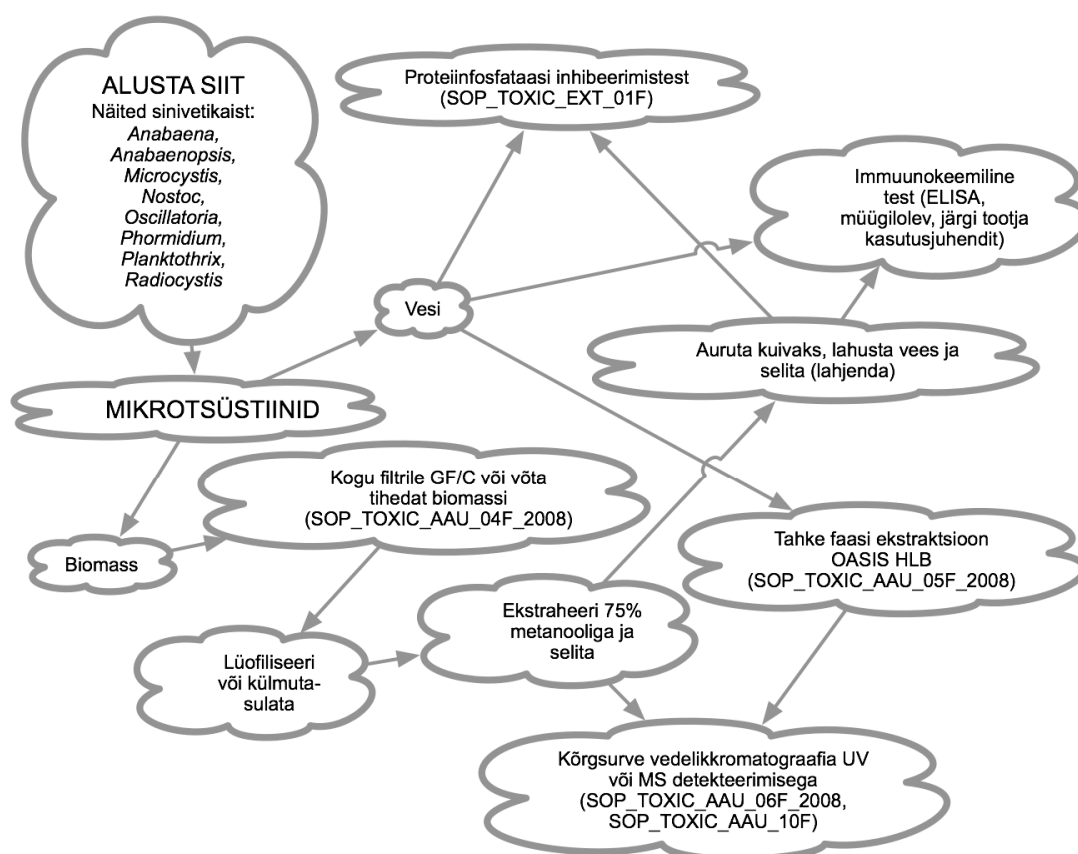
with the kind permission of the publisher.

When the original English and the translated Estonian versions are found to be conflicting, the English version must be regarded as valid.

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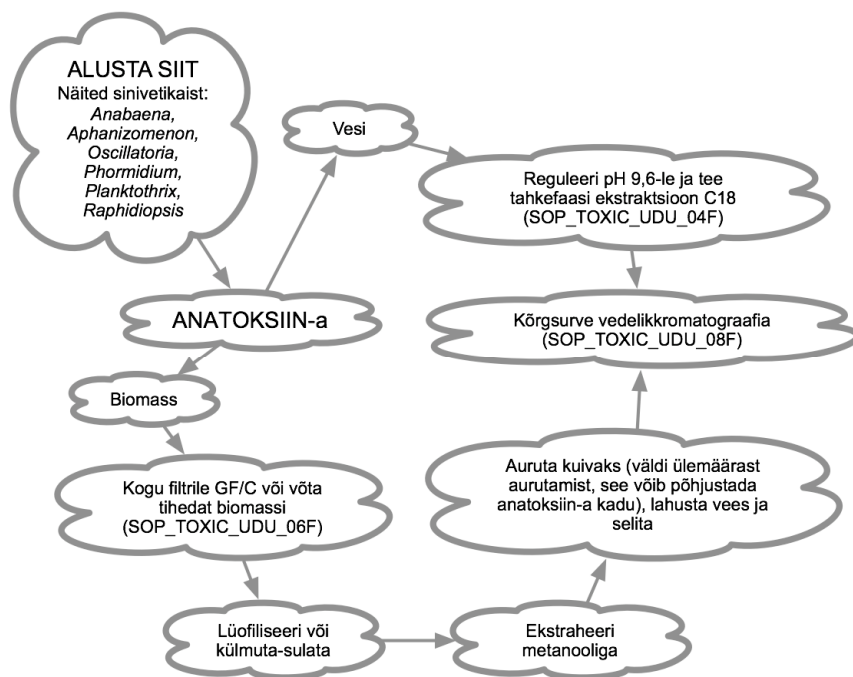
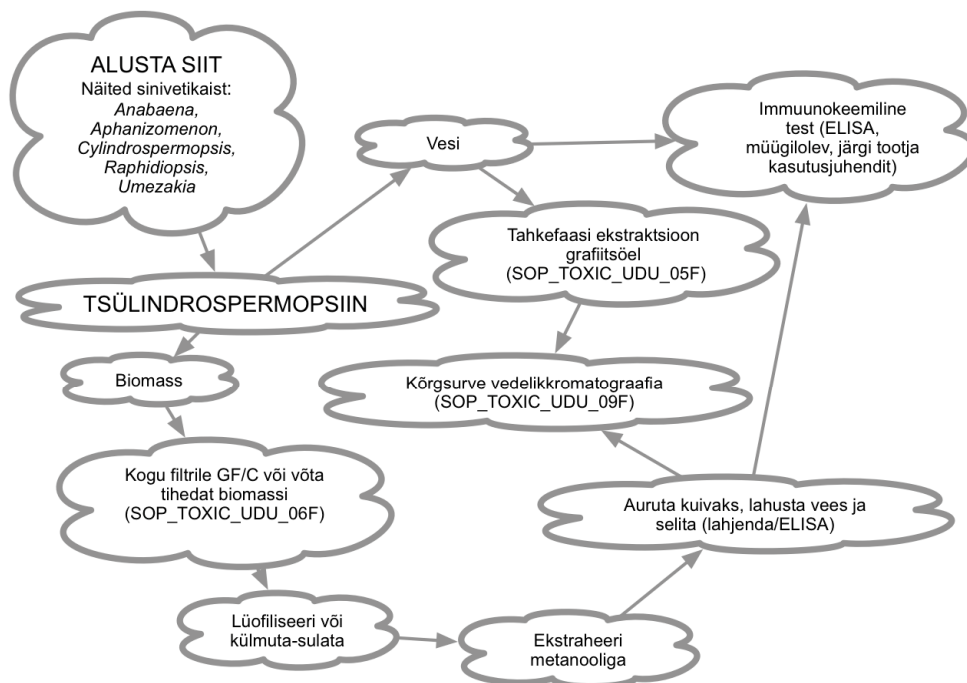
7 Selection of analytical methodology

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This chapter provides schematic workflows for the analyses of microcystins, anatoxin-a and cylindrospermopsin. High-performance liquid chromatography (HPLC) targets individual analytes and it requires access to authentic toxin standards for analyte identification. Protein phosphatase inhibition assay (PPIA) and enzyme-linked immunosorbent assay (ELISA) are group-specific, highly sensitive screening methods. In addition to microcystins and nodularins, PPIA also detects other protein phosphatase inhibitors.

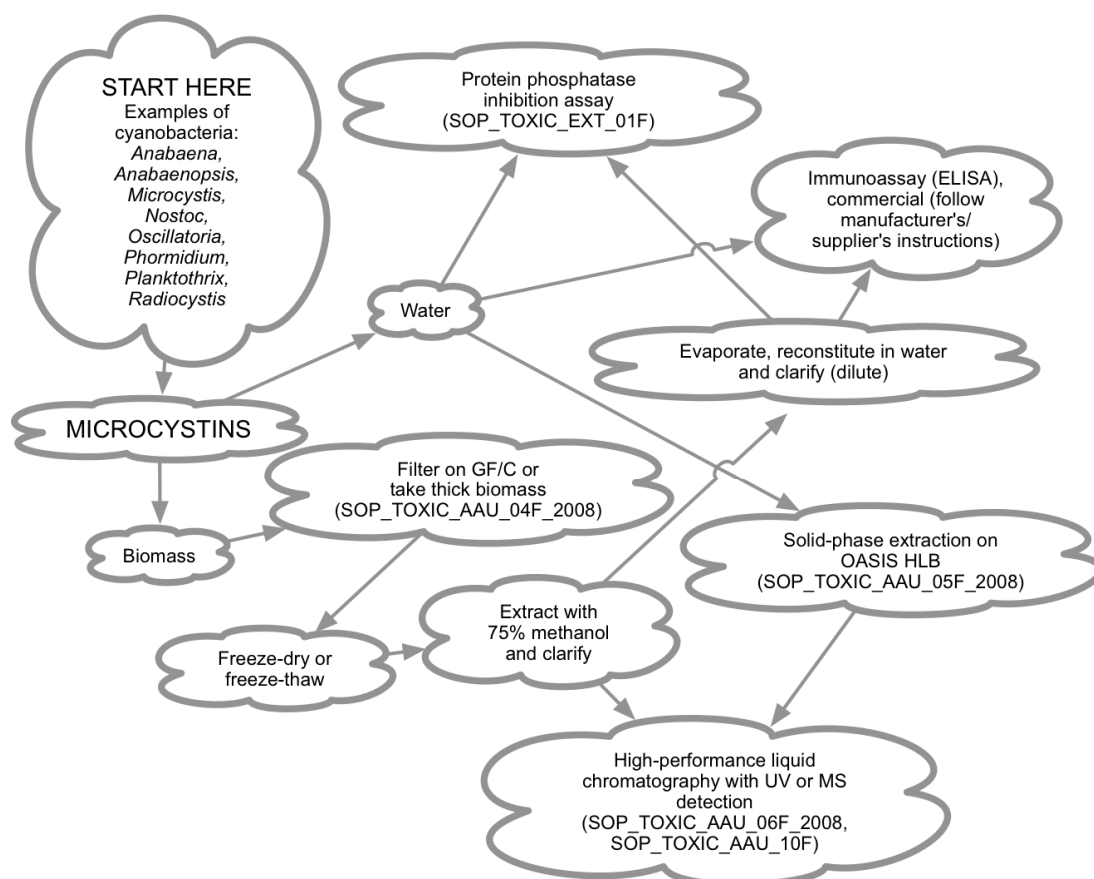
Figure 7.1: Workflow for the analysis of microcystins.

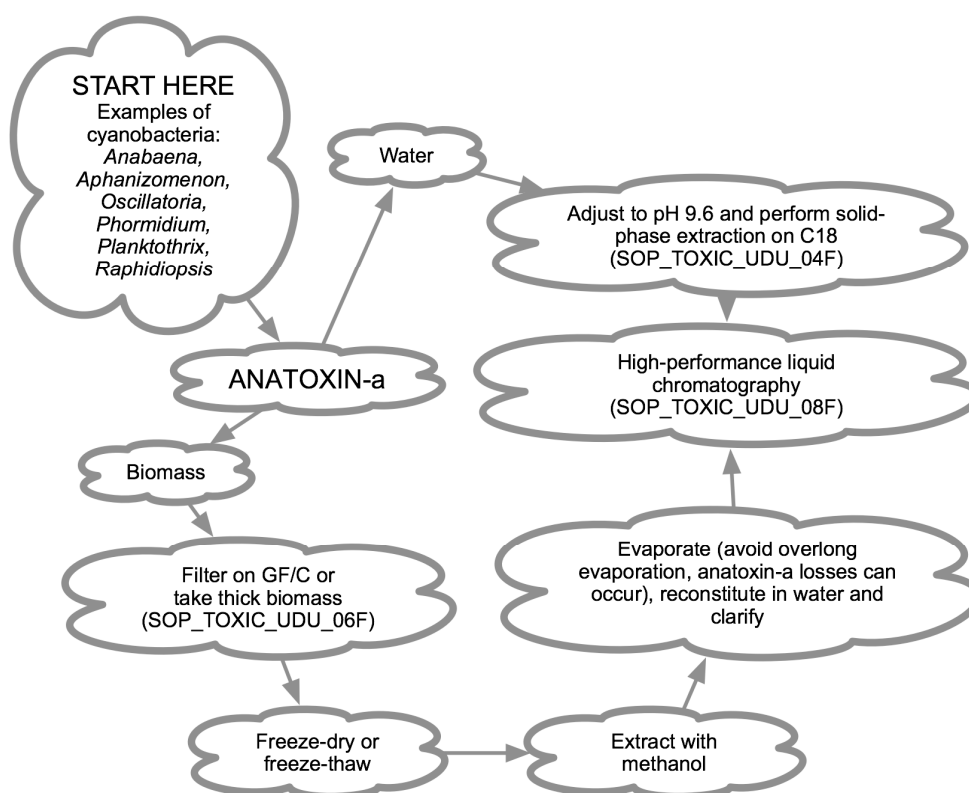
Figure 7.2: Workflow for the analysis of anatoxin-a.

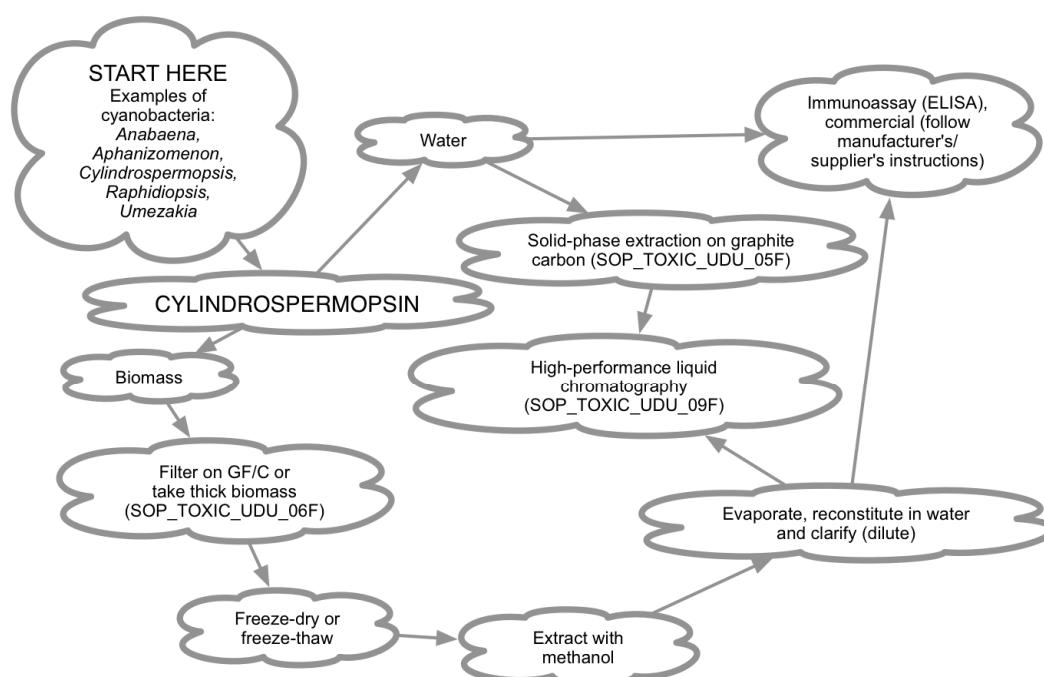
Figure 7.3: Workflow for the analysis of cylindrospermopsin.

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8.5 Kirjandusviited

1. Codd, G. A., Chorus, I., Burch, M. : Design of monitoring programmes. In Chorus, I., Bartram, J. (Eds.), Toxic Cyanobacteria in Water: A Guide to Their Public Health Consequences, Monitoring, and Management, pp. 313-328, E & FN Spon, London (1999).
2. Codd, G. A., Lindsay, J., Young, F. M., Morrison, L. F., Metcalf, J. S.: Harmful cyanobacteria: from mass mortalities to management measures. In Huisman, J., Matthijs, H. C. P., Visser, P. M. (Eds.), Harmful Cyanobacteria, pp. 1-23, Springer, Dordrecht (2005).
3. Burch, M. D.: Blue-green algal blooms - an alert levels framework for water supply contingency plans. In Steffensen D. A., Nicholson, B. C. (Eds.), Toxic Cyanobacteria. Current Status of Research and Management, pp. 133-139, Australian Centre for Water Quality Research, Adelaide (1994).
4. Bartram, J., Burch, M., Falconer, I. R., Jones, G. and Kuiper-Goodman, T.: Situation assessment, planning and management. In Chorus, I., Bartram, J. (Eds.), Toxic Cyanobacteria in Water: A Guide to Their Public Health Consequences, Monitoring, and Management, pp. 179-209, E & FN Spon, London (1999).
5. <http://www.cyanobacteria-platform.com/cbase.html>
6. Shaw, G. R., Sukenik, A., Livine, A., Chiswell, R. K., Smith, M. J., Seawright, A. A., Norris, R. L., Eaglesham, G. K., Moore, M. R. : Blooms of the cylindrospermopsin containing cyanobacterium, *Aphanizomenon ovalisporum* (Forti) in newly constructed lakes, Queensland, Australia. Env. Toxicol. **14**, 167-177 (1999).
7. Chiswell, R. K., Shaw, G. R., Eaglesham, G. K., Smith, M. J., Norris, R. L., Seawright, A. A., Moore, M. R.: Stability of cylindrospermopsin, the toxin from the cyanobacterium *Cylindrospermopsis raciborskii*. Effect of pH, temperature, and sunlight on decomposition. Env. Toxicol. **14**, 155-165 (1999).
8. Namikoshi, M., Murakami, T., Fujiwara, T., Nagai, H., Niki, T., Haigaya, E., Watanabe, M. F., Oda, T.: Biosynthesis and transformation of homanatoxin-a in the cyanobacterium *Raphidiopsis*

- mediterranea* Skuja and structures of three new analogues. *Chem. Res. Toxicol* **17**, 1692-1696 (2004).
9. *National Rivers Authority*: Toxic blue-green algae. Water Quality Series No. 2. National Rivers Authority, London, 125 pp. (1990).
 10. *New South Wales Blue-Green Algal Task Force*: Final report of the New South Wales Task Force. Department of Water Resources, Parramatta, Australia, 159 pp. (1992).
 11. *Pilotto, L. S., Douglas, R. M., Burch, M. D., Cameron, S., Beers, M., Rouch, G. R., Robinson, P., Kirk, M., Cowie, C. T., Hariman, S., Moore, C., Attewell, R. G.*: Health effects of recreational exposure to cyanobacteria (blue-green algae) during recreational water-related activities. *Aust. N. Z. J. Publ. Health* **21**, 562-566 (1997).
 12. *Kuiper-Goodman, T., Falconer, I., Fitzgerald, J.*: Human health aspects. In *Chorus, I., Bartram, J.* (Eds.), *Toxic Cyanobacteria in Water: A Guide to Their Public Health Consequences, Monitoring, and Management*, pp. 113-153, E & FN Spon, London (1999).
 13. *Pilotto, L. S., Hobson, P., Burch, M. D., Ranmuthugula, G., Attewell, R. G., Davies, R. D.*: Acute skin irritant effects of cyanobacteria (blue-green algae) in healthy volunteers. *Aust. N. Z. J. Public Health* **28**, 220-224 (2004).
 14. *Carmichael, W. W.* The cyanotoxins. In *Callow, J. A.* (Ed.), *Advances in Botanical Research*. Vol. 27, pp. 211-256, Academic Press, London and New York (1997).
 15. *Codd, G. A., Morrison, L. F., Metcalf, J. S.*: Cyanobacterial toxins: risk management for health protection. *Toxicol. Appl. Pharmacol.* **203**, 264-272 (2005).
 16. *Codd, G. A., Bell, S. G., Kaya, K., Ward, C. J., Beattie, K. A., Metcalf, J. S.*: Cyanobacterial toxins, exposure routes and human health. *Eur. J. Phycol.* **34**, 405-415 (1999).
 17. *Sipiä, V. O., Karlsson, K. M., Meriluoto, J. A. O., Kankaanpää, H. T.*: Eiders (*Somateria mollissima*) obtain nodularin, a cyanobacterial hepatotoxin, in Baltic Sea food web. *Env. Toxicol. Chem.* **23**, 1256-1260 (2004).
 18. *World Health Organization*: Guidelines for drinking water quality. Third Edition. Volume 2. Recommendations. World Health Organization, Geneva, 515 pp. (2004).
 19. *Humpage, A. R., Falconer, I. R.*: Oral toxicity of the cyanobacterial toxin cylindrospermopsin in male Swiss albino mice: determination of no observed adverse effect level for deriving a drinking water guideline value. *Env. Toxicol.* **18**, 94-103 (2003).

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SOP: Extraction of microcystins in biomass filtered on glass-fibre filters or in freeze-dried cyanobacterial biomass

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Date: 19 August 2008

2 Introduction

Water samples containing cyanobacterial cells can be filtered on glass-fibre filters (GF/C, diameter 25-47 mm) for monitoring of intracellular (biomass-bound) toxins. After freeze-thawing or freeze-drying the filters can be extracted by ultrasonication in 75% methanol. This SOP can also be applied to nodularin-R.

2 Experimental

2.1 Materials

- (a) Glass-fibre filters, equivalent to Whatman (Maidstone, UK) GF/C, diameter 25-47 mm
- (b) Methanol HPLC grade
- (c) HPLC grade water purified to 18.2 MΩ cm
- (d) Argon or nitrogen, >99.99%
- (e) Borosilicate glass tubes, 12 mm x 75 mm
- (f) Parafilm
- (g) Microcentrifuge tubes, polyethylene or polypropylene, narrow models preferred
- (h) Borosilicate glass chromatographic vials: e.g. 1.5 ml clear glass with writing surface. For small sample volumes 0.3 ml polypropylene vials (only for samples in 75% methanol) or borosilicate glass inserts can be used.
- (i) HPLC grade filters showing no adsorption of microcystins from 75% methanol, equivalent to GHP Acrodisc filters, 13 mm in diameter and 0.2 µm pore size, product number 4554, Pall Life

Sciences (Ann Arbor, MI, USA), and 2-ml all-plastic (no rubber piston) single use syringes (alternative sample clarification procedure to centrifugation)

2.2 Special equipment

- (a) Bath ultrasonicator, comparable to Branson (Danbury, CT, USA) 2510E-MT or Bandelin (Berlin, Germany) Sonorex RK 156
- (b) Probe ultrasonicator, comparable to Sonifier II W-250 from Branson/Emerson Technologies (Dietzenbach, Germany) or Bandelin (Berlin, Germany) Sonopuls HD 2070
- (c) Freeze-drying equipment (optional)
- (d) Microcentrifuge
- (e) Heater block with evaporation unit

3.3 Solutions

- (a) Mix 75 volumes of methanol and 25 volumes of water to make 75% methanol.

2.4 Extraction procedure

- (a) Place the air-dried, frozen filter samples (GF/C, diameter 25 mm) in 12 mm x 75 mm borosilicate glass tubes or in 1.5 ml borosilicate glass vials, and freeze-thaw them twice, or, to obtain maximum recovery, freeze-dry the filters. The extraction procedure described below is also suitable for the extraction of 5-8 mg samples of freeze-dried cyanobacteria.
- (b) Add 1.2 ml of 75% methanol, mix and extract in the bath ultrasonicator for 15 min. Filters with a diameter of 47 mm are extracted with 2 ml of 75% methanol in the glass tubes. The tubes should be covered with parafilm during the incubation. If vials are used, they should be fitted with caps during the bath ultrasonication.
- (c) Ultrasonicate the samples further, one at a time, with the ultrasonic disruptor equipped with a microtip probe, for 1 min. Wash the probe with 75% methanol in between samples. In the case of the Branson Sonifier II W-250, the unit is operated in a pulsed mode with a 30 % duty cycle and with an output control of 2. The Bandelin HD 2070 equipped with the MS 73 microtip is operated with a 30 % pulsed cycle and with an amplitude setting of 30 %. Avoid extensive fibre loss from the filter. Higher output effects necessitate the use of cooling around the sample with an ice bath.

- (d) Centrifuge aliquots of the extracts at $10,000 \times g$ for 10 min.
- (e) Transfer 500 μ l of the supernatants to 1.5 ml borosilicate glass vials or clean 12 mm x 75 mm borosilicate glass tubes and evaporate to dryness in a heating block at 40 °C under argon. If desired, take further aliquots.
- (f) Reconstitute the dried extracts in 100 μ l of 75% methanol, mix well, transfer to 0.3 ml polypropylene vials, centrifuge at $10,000 \times g$ for 10 min or filter through the HPLC grade filter to clarify the sample if necessary, and analyse with high-performance liquid chromatography (HPLC) according to the procedure SOP_TOXIC_AAU_06F_2008.

N.B. The efficiency of ultrasonication varies considerably between different ultrasonicator models. Therefore, the above settings and times must be adjusted according to the apparatus in use.

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Fastner, J., Flieger, I., Neumann, U.: Optimized extraction of microcystins from field samples - a comparison of different solvents and procedures. *Water Res.* **32**, 3177-3181 (1998).

Spoof, L., Vesterkvist, P., Lindholm, T., Meriluoto, J. Screening for cyanobacterial hepatotoxins, microcystins and nodularin, in environmental water samples by reversed-phase liquid chromatography-electrospray ionisation mass spectrometry. *J. Chromatogr. A* **1020**, 105-119 (2003).

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SOP: Solid phase extraction of microcystins in water samples

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Date: 19 August 2008

1 Introduction

Traditionally solid phase extraction of microcystins in natural or processed waters has been performed on C₁₈ silica cartridges. New polymeric materials, in particular the OASIS HLB consisting of a macroporous copolymer [poly(vinylbenzene-co-N-vinylpyrrolidone)], can also be useful for the purpose. Chlorinated waters are usually treated with sodium thiosulphate before concentration in order to stop the further degradation of microcystins. Use borosilicate glassware throughout the procedure whenever possible because microcystins in water solutions are easily adsorbed onto plastics commonly used in laboratories. Microcystin-LR in a solution of above 25% methanol was found not to interact with plastics to any great extent [Hyenstrand et al., 2001]. This SOP can also be applied to nodularin-R.

2 Experimental

2.1 Materials

Use analytical reagent grade reagents if not indicated otherwise.

- (a) Methanol HPLC grade
- (b) HPLC grade water purified to 18.2 MΩ cm
- (c) Glass-fibre filters, equivalent to Whatman (Maidstone, UK) GF/C, diameter 25-47 mm
- (d) Sodium thiosulphate
- (e) NH₄OH
- (f) Trifluoroacetic acid (TFA), HPLC or protein sequence analysis grade. TFA should be stored under argon in a desiccator.
- (g) Reversed-phase solid phase extraction columns suitable for microcystin work, equivalent to OASIS HLB solid phase extraction columns, size 200 mg sorbent, part number WAT106202, from

Waters (Milford, MA, USA). Waters sample volumes <<500 ml can be concentrated on smaller columns.

- (h) Argon or nitrogen, >99.99%
- (i) Borosilicate test tubes or vials, >4 ml capacity
- (j) HPLC grade filters showing no adsorption of microcystins from 75% methanol, equivalent to GHP Acrodisc filters, 13 mm in diameter and 0.2 μ m pore size, product number 4554, Pall Life Sciences (Ann Arbor, MI, USA), and 2 ml all-plastic (no rubber piston) single use syringes (alternative sample clarification procedure to centrifugation)
- (k) Microcentrifuge tubes, polyethylene or polypropylene, narrow models preferred
- (l) Borosilicate glass chromatographic vials: e.g. 1.5 ml clear glass with writing surface. For small sample volumes 0.3 ml polypropylene vials (only for samples in 75% methanol) or borosilicate glass inserts can be used.

2.2 Special equipment

- (a) Vacuum manifold, preferably transparent, equipped with stopcocks (stopcocks are optional with the OASIS HLB columns), vacuum source and vacuum control
- (b) Large volume extraction (LVE) kit for unattended loading of large sample volumes, made of PTFE tubing and adaptors for SPE column connection
- (c) pH meter
- (d) Filtration unit for 500 ml volume
- (e) Heating block, operated at 40 °C, with evaporation unit
- (f) Microcentrifuge (alternative to filtration)

2.3 Solutions

- (a) Sodium thiosulphate, 1 g in 100 ml HPLC grade water
- (b) TFA, 1% solution in HPLC grade water
- (c) NH_4OH , 2 g in 100 ml HPLC grade water
- (d) 20% methanol (20 volumes of methanol and 80 volumes of water)

- (e) 100% methanol
- (f) 75% methanol (75 volumes of methanol and 25 volumes of water)

2.4 Procedure

- (a) Filter the water sample through the glass-fibre filter if it contains particulates or cyanobacterial cells.
- (b) Transfer 500 ml of water sample into a borosilicate glass flask or beaker.
- (c) Measure the pH of the water sample, and, if necessary, adjust to pH 5-8 with dilute NH_4OH or TFA.
- (d) In case of chlorinated water: add 500 μl of sodium thiosulphate solution (1 g / 100 ml) to 500 ml water. Shake vigorously and allow to stand for 5 minutes.
- (e) Add 5 ml of methanol and mix thoroughly.
- (f) Condition the solid phase extraction cartridge, OASIS HLB, 200 mg, with 5 ml of methanol followed by 10 ml of water.
- (g) Apply the sample at a flow rate not exceeding 10 ml min^{-1} (visible drops). Regulate the flow with vacuum pressure.
- (h) Wash the cartridge with 4 ml of 20% methanol.
- (i) Dry the cartridge by drawing air through it for 2 min.
- (j) Elute microcystins with 4 ml of 100% methanol e.g. in 12 mm x 75 mm borosilicate test tube or in borosilicate glass vials. It is advisable that the elution solvent is loaded in two fractions: draw 2 ml solvent into the SPE column, soak the sorbent for 3 min, then load a further 2 ml solvent into the column and collect all the solvent (max. flow rate 4 ml min^{-1}).
- (k) Evaporate the methanol eluate at 40°C using argon or nitrogen.
- (l) Resuspend the residue in 500 μl of 75% methanol.
- (m) Centrifuge 10 min $10,000 \times g$, or, alternatively, filter, and transfer to HPLC vials.
- (n) Analyse supernatant/filtrate on HPLC according to SOP_TOXIC_AAU_06F_2008.

3 References

Hyenstrand, P., Metcalf, J.S., Beattie, K.A., Codd, G.A.: Losses of the cyanobacterial toxin microcystin-LR from aqueous solution by adsorption during laboratory manipulations. *Toxicon* **39**, 589-594 (2001).

Lawton, L.A., Edwards, C., Codd, G.A.: Extraction and high-performance liquid chromatographic method for the determination of microcystins in raw and treated waters. *Analyst (London)* **119**, 1525-1530 (1994).

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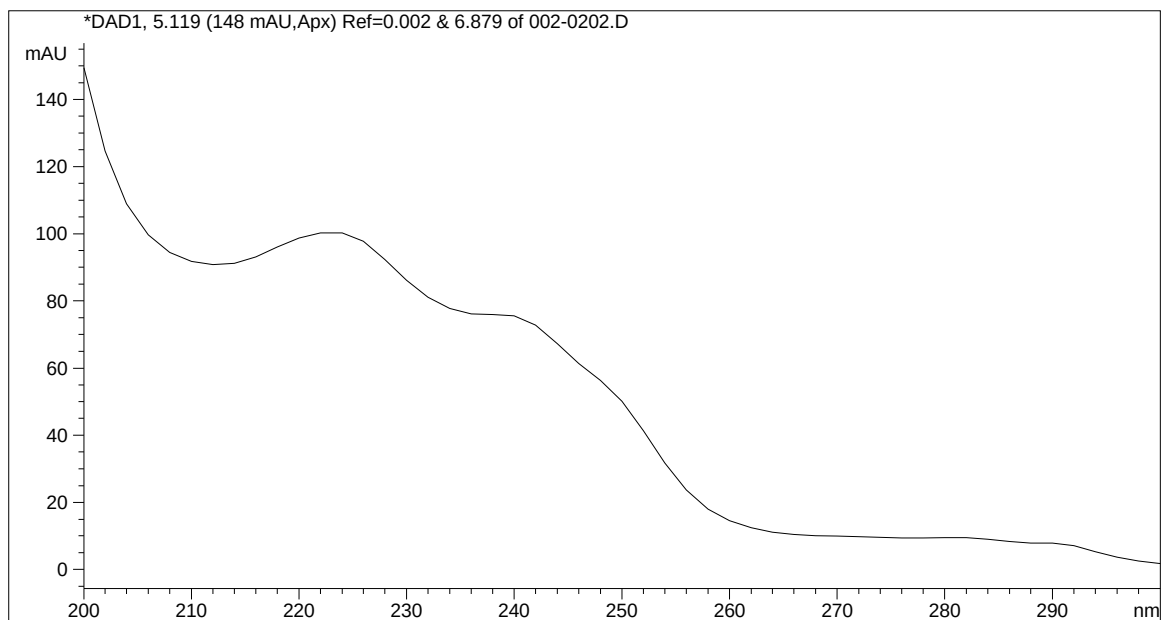
3 Kirjandusviited

Lawton, L.A., Edwards, C., Codd, G.A.: Extraction and high-performance liquid chromatographic method for the determination of microcystins in raw and treated waters. *Analyst* (London) **119**, 1525-1530 (1994).

Meriluoto, J.: Chromatography of microcystins. *Anal. Chim. Acta* **352**, 277-298 (1997).

Meriluoto, J., Lawton, L., Harada, K.-i.: Isolation and detection of microcystins and nodularins, cyanobacterial peptide hepatotoxins. In: *Holst, O.* (Ed.), *Bacterial Toxins: Methods and Protocols*, pp. 65-87, Humana Press, Totowa, NJ (2000).

Meriluoto, J.: Toxins of freshwater cyanobacteria (blue-green algae). In: *Bogusz, M.J.* (Ed.), *Forensic Science, Handbook of Analytical Separations*, Vol. 2, pp. 359-390, Elsevier, Amsterdam (2000).

Spoof, L., Karlsson, K., Meriluoto, J.: High-performance liquid chromatographic separation of microcystins and nodularin, cyanobacterial peptide toxins, on C₁₈ and amide C₁₆ sorbents. *J. Chromatogr. A* **909**, 225-236 (2001).

Spoof, L., Meriluoto, J.: Rapid separation of microcystins and nodularin using a monolithic silica C₁₈ column. *J. Chromatogr. A* **947**, 237-245 (2002).

Spoof, L., Vesterkvist, P., Lindholm, T., Meriluoto, J.: Screening for cyanobacterial hepatotoxins, microcystins and nodularin, in environmental water samples by reversed-phase liquid chromatography-electrospray ionisation mass spectrometry. *J. Chromatogr. A* **1020**, 105-119 (2003).

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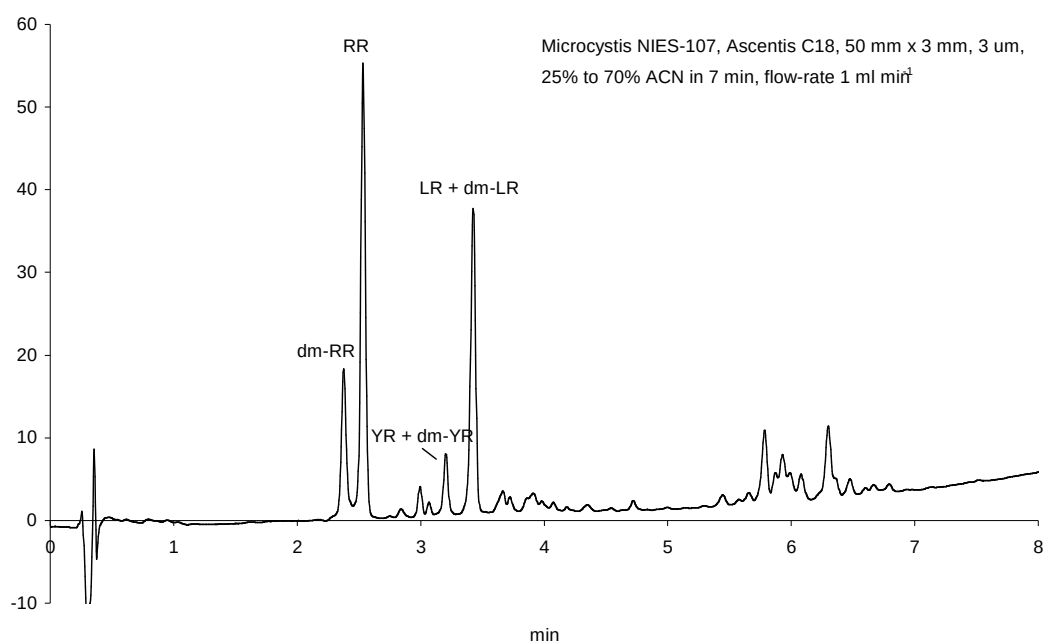
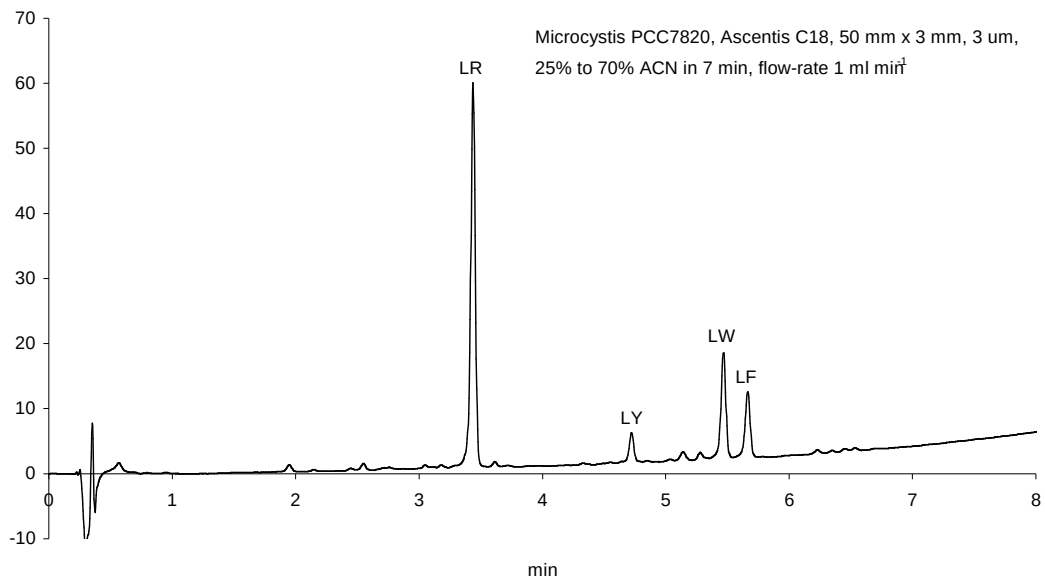
SOP: Analysis of microcystins by high-performance liquid chromatography with photodiode-array detection

Document identifier: SOP_TOXIC_AAU_06F_2008

(modified from an earlier version printed in TOXIC : Cyanobacterial Monitoring and Cyanotoxin Analysis, ISBN 951-765-259-3, (c) Åbo Akademi University Press, 2005)

Prepared by: Jussi Meriluoto and Lisa Spoof, AAU

Date: 19 August 2008

1 Introduction

Reversed-phase HPLC on C₁₈ phases is a common choice for separating smaller peptides, and the mobile phases for peptides often consist of acetonitrile gradients in the presence of perfluorinated alkyl carboxylic acids, usually trifluoroacetic acid (TFA). Microcystins make no exception, they chromatograph perfectly under these conditions. Neutral ammonium acetate - acetonitrile based eluents are commonly used in preparative toxin separations.

The rapidly growing number of known microcystins, now exceeding 80, has made it impossible to separate and quantitate all microcystins on a single chromatographic run. However, using retention time and spectrum match, the major microcystins, such as microcystin-LR, -RR, -YR with different degrees of demethylation, are fairly easy to identify by HPLC combined with photodiode- array detection. The main chromophores of microcystins, the conjugated diene in the Adda residue and the α,β -unsaturated carbonyl group in *N*-methyldehydroalanine, absorb strongly at 238 nm. UV spectra of microcystins can be divided into three main categories: normal microcystin spectra with a local absorbance maximum at 238-240 nm, the spectra of tyrosine(Y)-containing microcystins with a flatter absorbance maximum at 230-240 nm, and the spectra of tryptophan(W)-containing microcystins which have an absorption maximum at 222-223 nm and a shoulder at 238-240 nm (examples of spectra in Figs. 6 and 7).

Cyanobacterial blooms can produce very complex microcystin profiles. In order to elucidate the toxin profile and thus indirectly estimate the total toxicity of a complex sample, the chemical methods used for toxin analyses should be able to separate and quantify individual microcystins which have different toxicities. This calls for separation methods of high resolution and good selectivity, i.e. the use of long columns and long run times, or more specific detection methods such as mass spectrometry. On the other hand, many samples in e.g. water treatment trials are rather trivial chromatographically and therefore less resolution may be enough if only a few analytes are to be separated.

This SOP can also be applied to nodularin-R.

2 Experimental

2.1 Materials

- (a) Acetonitrile gradient HPLC grade
- (b) HPLC grade water purified to 18.2 MΩ cm
- (c) Trifluoroacetic acid (TFA), HPLC or protein sequence analysis grade. TFA should be stored under argon in a desiccator.
- (d) C₁₈ endcapped HPLC column (Purospher STAR RP-18 endcapped, 3 μm particles, LiChroCART 55 x 4 mm I.D., and Purospher STAR RP-18 endcapped, 5 μm particles, LiChroCART 250 x 4 mm I.D., from Merck, Darmstadt, Germany, are mentioned in this SOP as possible alternatives. Examples of chromatograms obtained with an Ascentis C18, 3μm particles, 50 x 3 mm I.D from Supelco, Bellefonte, PA, USA have been added for comparison.) A compatible guard column is also required.
- (e) Borosilicate glass chromatographic vials: e.g. 1.5 ml clear glass with writing surface. For small sample volumes 0.3 ml polypropylene vials (only for samples in 75% methanol) or borosilicate glass inserts can be used.
- (f) Microcystin standards

2.2 Special equipment

- (a) High-performance liquid chromatograph equipped with high-pressure or low-pressure gradient pump, autosampler, column oven and photodiode-array (PDA) detector
- (b) Chromatography analysis software
- (c) A positive displacement pipette with a teflon-coated piston working in a glass capillary, capable of accurately dispensing 200 μl of the strong acid TFA

2.5 HPLC mobile phase

- (a) HPLC mobile phase component A: HPLC water + 0.05% TFA.
 - i. Add 800 μl of TFA in 1.6 l of water.
 - ii. Replace every week.
- (b) HPLC mobile phase component B: acetonitrile + 0.05% TFA
 - i. Take ca 450 μl of TFA in a borosilicate glass vial. Add 400 μl of TFA in 800 ml of ACN. Discard the rest of the TFA.

- ii. Replace every week.

Please observe the following:

Work in fume hood with TFA or acetonitrile.

Dispense the TFA under the surface of acetonitrile.

Do not contaminate the original TFA bottle with acetonitrile.

Store undiluted TFA under argon.

Disposal of acetonitrile and TFA should conform to local regulations.

2.4 Chromatography

2.4.1 General procedure

- (a) The HPLC system should be set up as described in the manufacturers instructions including degassing, priming and changing columns.
- (b) Always use a guard column. Change the guard column if the back-pressure rises or peak forms deteriorate.
- (c) Set column oven at 40 °C.
- (d) Change the HPLC gradually up to starting conditions and allow to condition.
- (e) Chromatograph the samples and standards as per the recommended HPLC gradients (see below), use 10 µl injections.
- (f) Analyse the chromatogram. Compare retention times and spectra to standards.
- (g) Calculate the microcystin concentration according to the standard curve procedure described in TOXIC_SOP_AAU_03F.

Procedure for quantitation of microcystins other than microcystin-LR in HPLC-UV work: all peaks with microcystin spectra are quantified as microcystin-LR equivalents, i.e. with the same "area to ng/injection" coefficient as for microcystin-LR. If the molecular weight of the other microcystin is known, a correction for molecular weight difference can be applied in critical cases.

2.4.2 HPLC of microcystins on a short column, suitable for samples of medium complexity; column Merck Purospher STAR RP-18 endcapped, 3 µm particles, LiChroCART 55 x 4 mm I.D.

Do not use unnecessarily long HPLC columns with easy samples. The use of a short column saves time and solvents.

Table 1: Suggested gradient programme for Merck Purospher STAR RP-18 endcapped, 3 μ m particles, LiChroCART 55 x 4 mm I.D., linear gradient at a flow rate of 1 ml min⁻¹. Injection cycle about 11 minutes.

Time (min)	% A	% B
0.00	75	25
5.00	30	70
6.00	30	70
6.10	75	25
9.00	STOP	

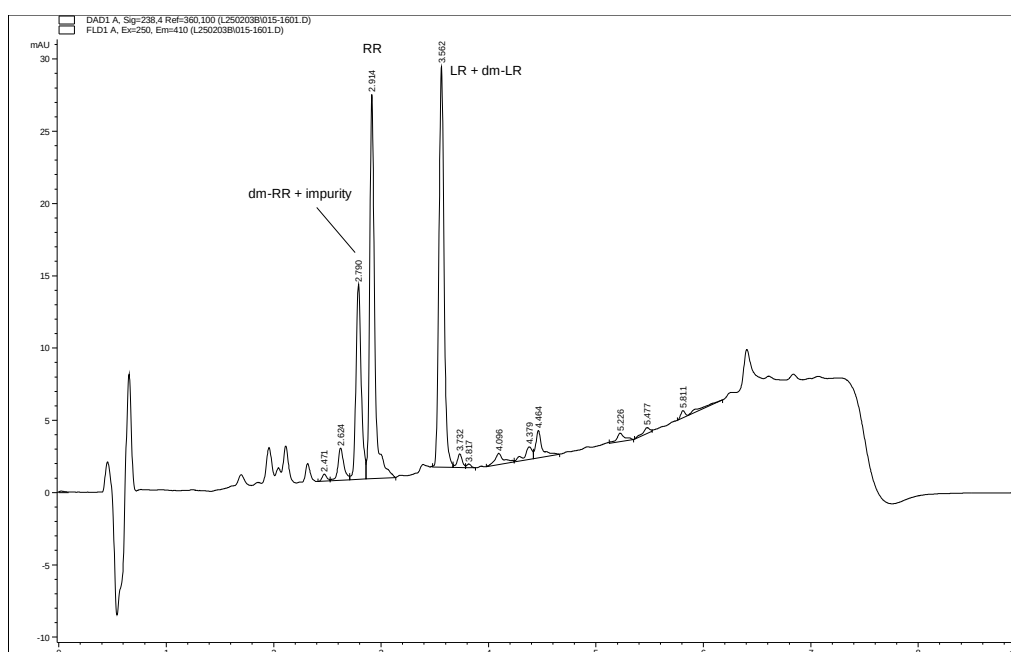
Figure 1: Trace of *Anabaena* extract. Column Merck Purospher STAR RP-18 endcapped, 3 μ m particles, LiChroCART 55 x 4 mm I.D., detection at 238 nm. Other parameters as in Table 1. Retention times: microcystin-RR 2.91 min, microcystin-LR 3.56 min.

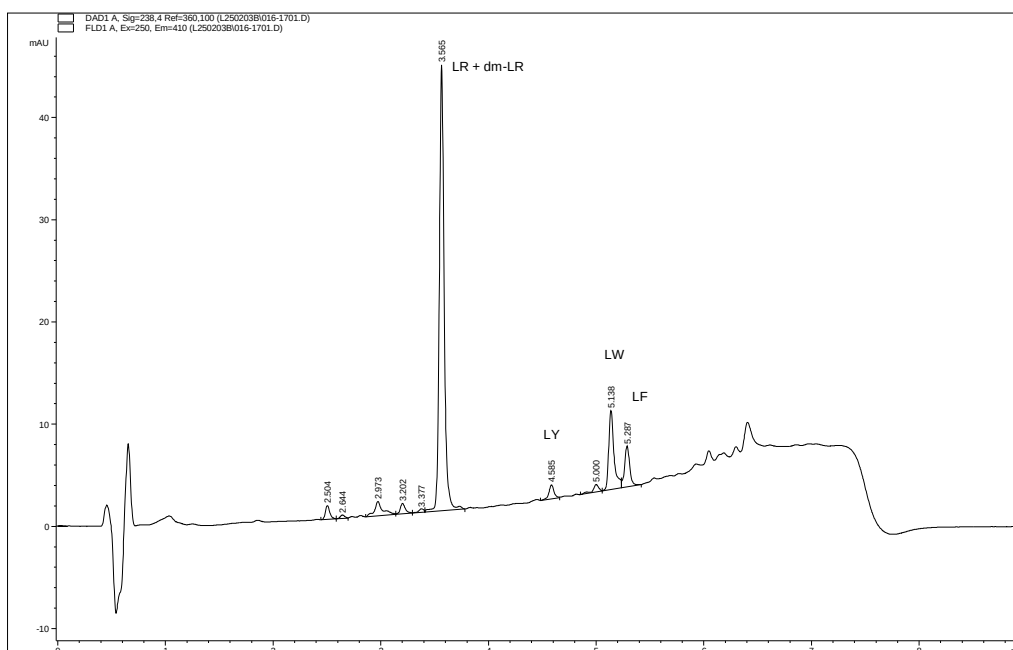
Figure 2: Trace of *Microcystis* extract. Column Merck Purospher STAR RP-18 endcapped, 3 μ m particles, LiChroCART 55 x 4 mm I.D., detection at 238 nm. Other parameters as in Table 1. Retention times: microcystin-LR 3.56 min, microcystin-LY 4.58 min, microcystin-LW 5.14 min, microcystin-LF 5.29 min.

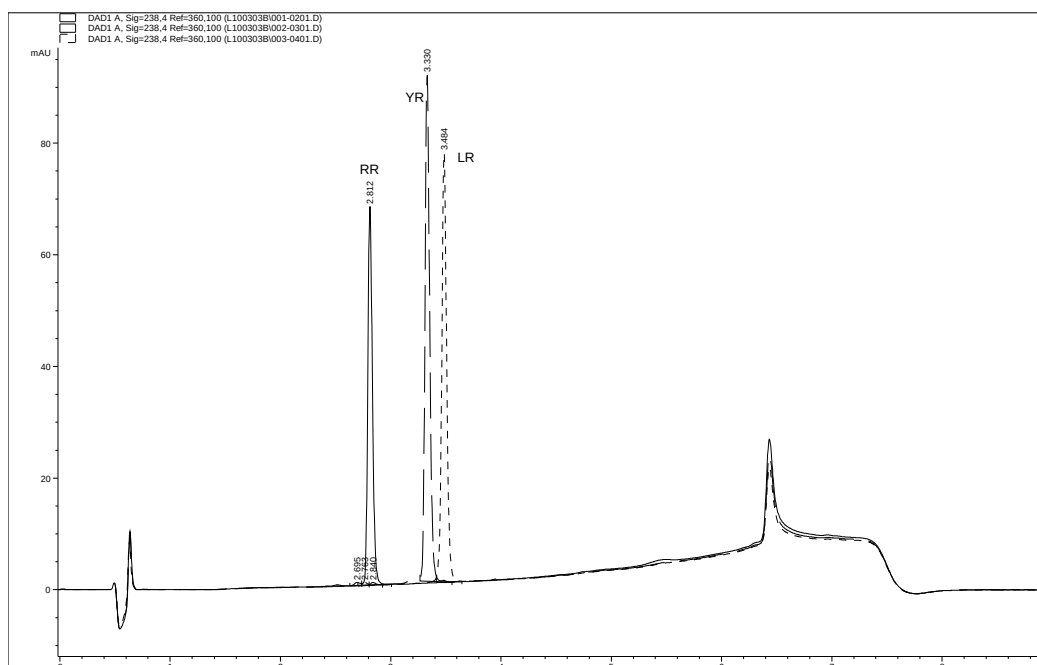
Figure 3: Overlaid traces of commercial microcystin-RR, -YR and -LR samples. Column Merck Purospher STAR RP-18 endcapped, 3 μ m particles, LiChroCART 55 x 4 mm I.D., detection at 238 nm. Other parameters as in Table 1. Retention times: microcystin-RR 2.81 min, microcystin-YR 3.33 min, microcystin-LR 3.48 min. This chromatogram

was run with an older column and the retention times have been shortened, cf. Figs. 1 and 2. The column should be replaced when the retention times have been shortened by more than 3%. Selectivity changes are possible even before this limit.

2.4.3 HPLC of microcystins on a long column, suitable for complex samples; column Merck Purospher STAR RP-18 endcapped, 5 μ m particles, LiChroCART 250 x 4 mm I.D.

Complex field or strain samples may necessitate the use of a long column.

Table 2: One possible gradient programme for Merck Purospher STAR RP-18 endcapped, 5 μ m particles, LiChroCART 250 x 4 mm I.D., linear gradient at a flow rate of 0.75 ml min⁻¹. Injection cycle about 62 minutes.

Time (min)	% A	% B
0.00	70	30
10.00	65	35
40.00	30	70
42.00	0	100
44.00	0	100
46.00	70	30
60.00	STOP	

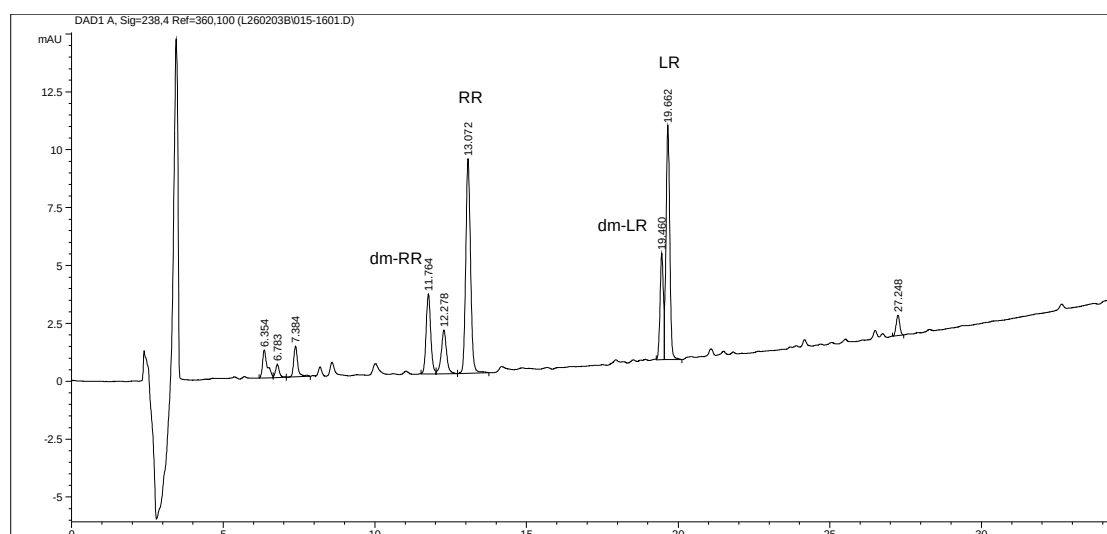
Figure 4: Trace of *Anabaena* extract. Column Merck Purospher STAR RP-18 endcapped, 5 μ m particles, LiChroCART 250 x 4 mm I.D., detection at 238 nm. Other parameters as in Table 2. Retention times: microcystin-RR 13.07 min, microcystin-LR 19.66 min.

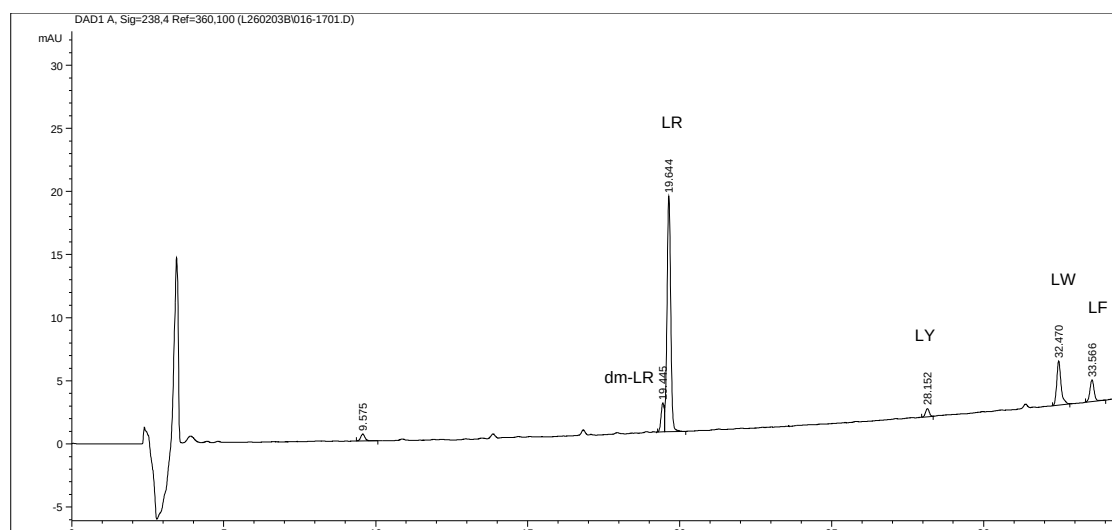
Figure 5: Trace of *Microcystis* extract. Column Merck Purospher STAR RP-18 endcapped, 5 μ m particles, LiChroCART 250 x 4 mm I.D., detection at 238 nm. Other parameters as in Table 2. Retention times: microcystin-LR 19.64 min, microcystin-LY 28.15 min, microcystin-LW 32.47 min, microcystin-LF 33.57 min.

2.5 Spectra

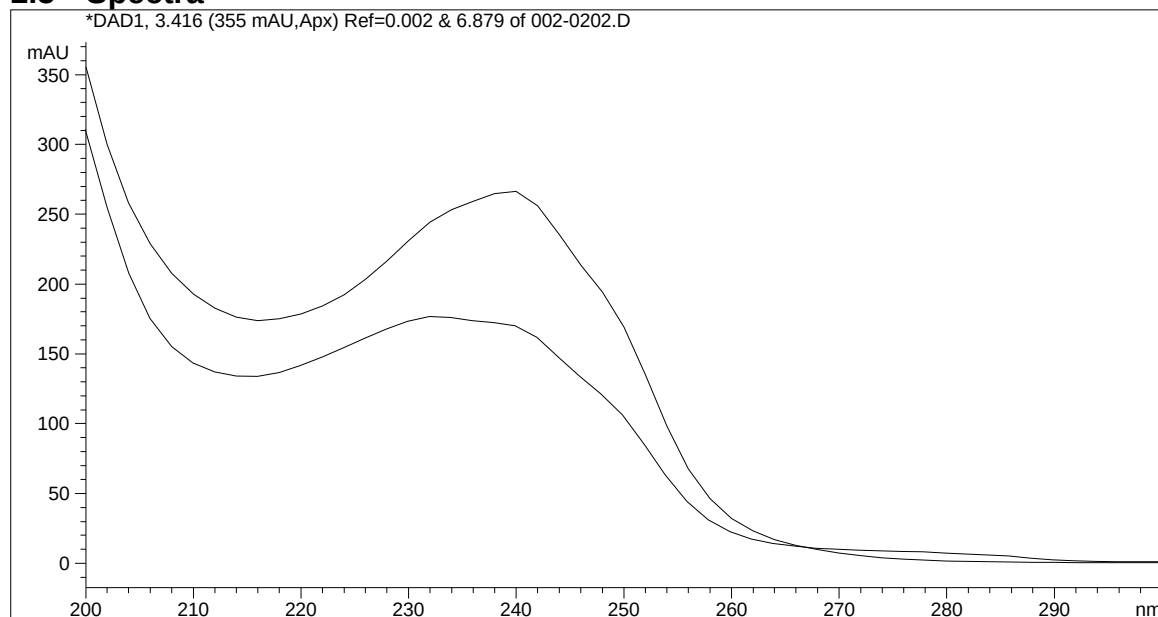
Figure 6: UV spectra of microcystin-LR (upper) and microcystin-YR (lower), chromatographic parameters according to section 2.4.2.

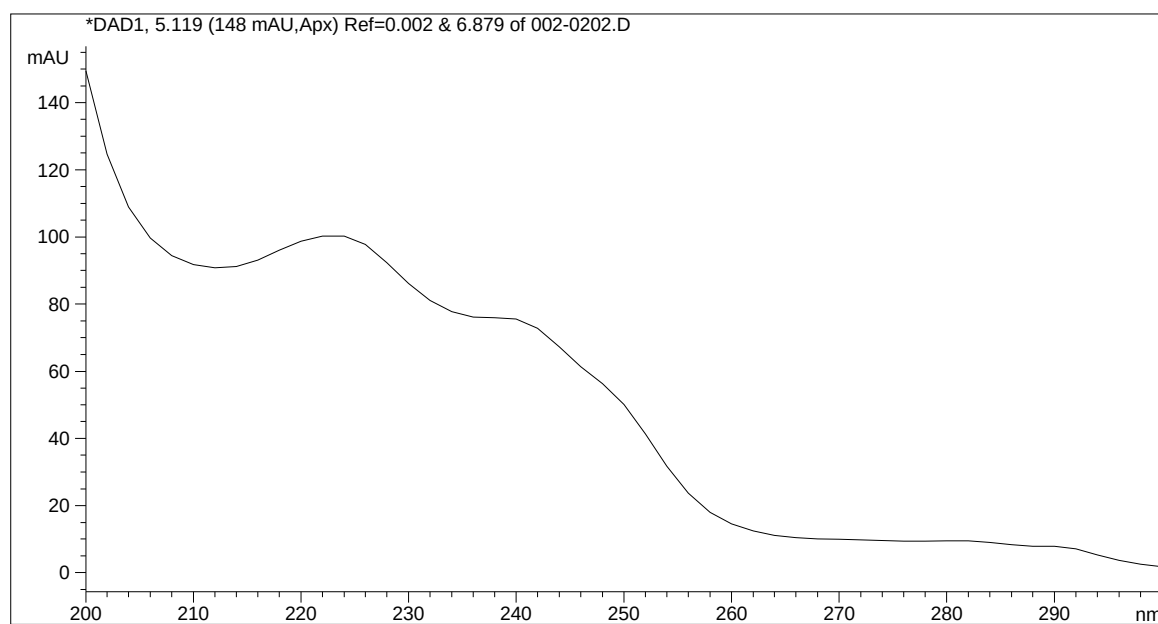
Figure 7: UV spectrum of microcystin-LW, chromatographic parameters according to section 2.4.2.

3 References

- Lawton, L.A., Edwards, C., Codd, G.A.:* Extraction and high-performance liquid chromatographic method for the determination of microcystins in raw and treated waters. *Analyst* (London) **119**, 1525-1530 (1994).
- Meriluoto, J.:* Chromatography of microcystins. *Anal. Chim. Acta* **352**, 277-298 (1997).
- Meriluoto, J., Lawton, L., Harada, K.-i.:* Isolation and detection of microcystins and nodularins, cyanobacterial peptide hepatotoxins. In: *Holst, O.* (Ed.), *Bacterial Toxins: Methods and Protocols*, pp. 65-87, Humana Press, Totowa, NJ (2000).
- Meriluoto, J.:* Toxins of freshwater cyanobacteria (blue-green algae). In: *Bogusz, M.J.* (Ed.), *Forensic Science, Handbook of Analytical Separations*, Vol. 2, pp. 359-390, Elsevier, Amsterdam (2000).
- Spoof, L., Karlsson, K., Meriluoto, J.:* High-performance liquid chromatographic separation of microcystins and nodularin, cyanobacterial peptide toxins, on C₁₈ and amide C₁₆ sorbents. *J. Chromatogr. A* **909**, 225-236 (2001).
- Spoof, L., Meriluoto, J.:* Rapid separation of microcystins and nodularin using a monolithic silica C₁₈ column. *J. Chromatogr. A* **947**, 237-245 (2002).
- Spoof, L., Vesterkvist, P., Lindholm, T., Meriluoto, J.:* Screening for cyanobacterial hepatotoxins, microcystins and nodularin, in environmental water samples by reversed-phase liquid chromatography-electrospray ionisation mass spectrometry. *J. Chromatogr. A* **1020**, 105-119 (2003).

N. B. The reader is also advised to follow up the development of the ISO standard for microcystin analysis.

ISO 20179:2005. Water quality -- Determination of microcystins -- Method using solid phase extraction (SPE) and high performance liquid chromatography (HPLC) with ultraviolet (UV) detection.

(THE ISO 20179:2005 DOCUMENT DESCRIBES A STANDARDISED METHOD FOR THE DETERMINATION OF MICROCYSTINS WHICH IS DIFFERENT FROM THE PROTOCOL DESCRIBED IN THE CURRENT SOP.)

Material added in 2008

Table 3: Suggested gradient programme for Supelco Ascentis C18, 3 μ m particles, 50 x 3 mm I.D., linear gradient at a flow rate of 1 ml min⁻¹. Injection cycle about 17 minutes.

Time (min)	% A	% B
0.00	75	25
7.00	30	70
10.00	30	70
10.10	75	25
15.00	STOP	

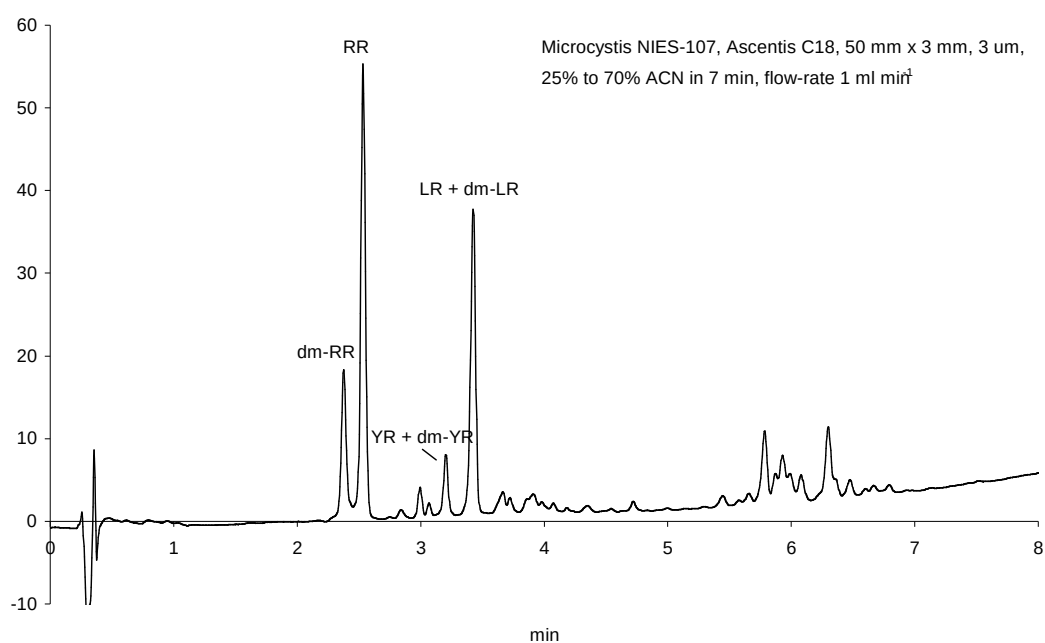
Figure 8: Trace of *Microcystis* NIES-107 extract. Column Supelco Ascentis C18, 3 μ m particles, 50 x 3 mm I.D., detection at 238 nm. Other parameters as in Table 3. Retention times: microcystin-RR 2.53 min, microcystin-YR co-eluting with demethylated microcystin-YR 3.20 min, microcystin-LR co-eluting with demethylated microcystin-LR 3.43 min.

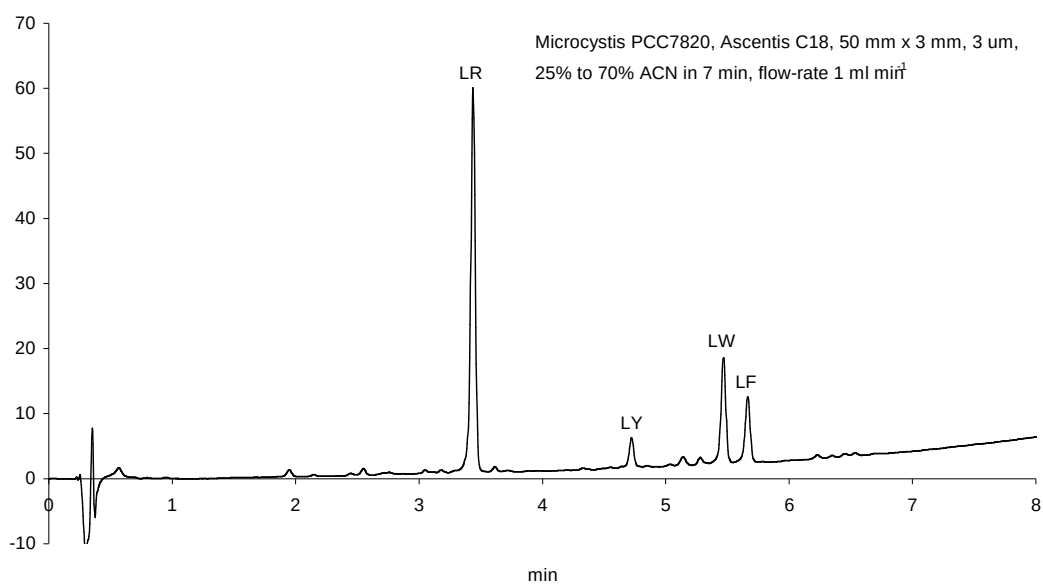
Figure 9: Trace of *Microcystis* PCC7820 extract. Column Supelco Ascentis C18, 3 μ m particles, 50 x 3 mm I.D., detection at 238 nm. Other parameters as in Table 3. Retention times: microcystin-LR 3.44 min, microcystin-LY 4.72 min, microcystin-LW 5.47 min, microcystin-LF 5.66 min.

SOP: Extraction of anatoxin-a or cylindrospermopsin from cyanobacterial biomass filtered on glass-fibre filters

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3 Kirjandusviited

Edwards, C., Beattie, K. A., Scrimgeour, C. M. Codd, G. A.: Identification of anatoxin-a in benthic cyanobacteria (blue-green algae) and in associated dog poisonings at Loch Insh, Scotland. *Toxicon* **30**,1165-1175 (1992).

Metcalf, J. S., Beattie, K. A., Saker, M. L., Codd, G. A.: Effects of organic solvents on the high performance liquid chromatographic analysis of the cyanobacterial toxin cylindrospermopsin and its recovery

from environmental eutrophic waters by solid phase extraction.
FEMS Microbiol. Lett. **216**,159-164 (2002).

SOP: Analysis of environmental cyanobacterial samples by ELISA for microcystins

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Metcalf, J. S., Codd, G. A.: Analysis of cyanobacterial toxins by immunological methods. *Chem. Res. Toxicol.* **16**, 103-112 (2003).

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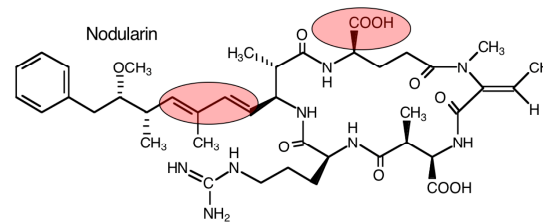
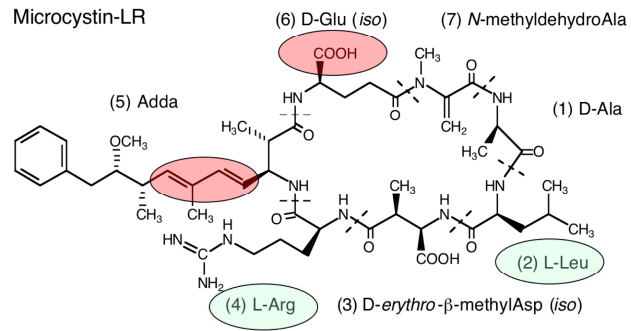
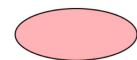
Hepatotoxic microcystins and nodularins

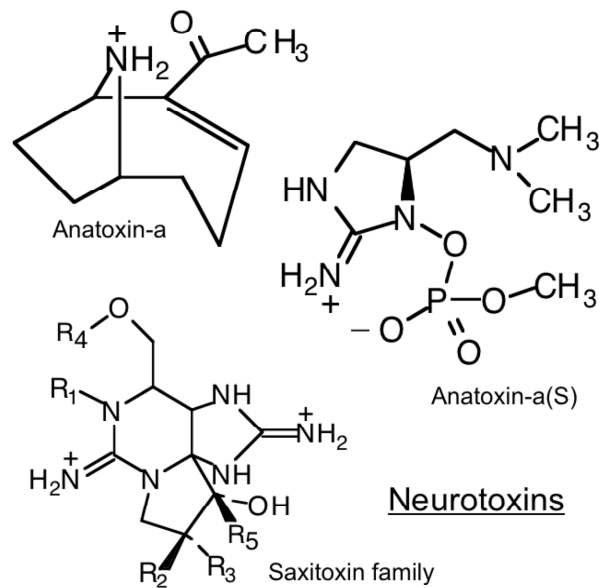
Microcystin-LR

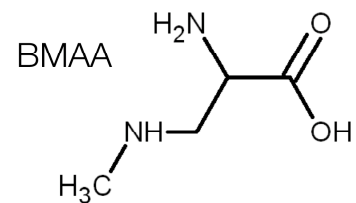
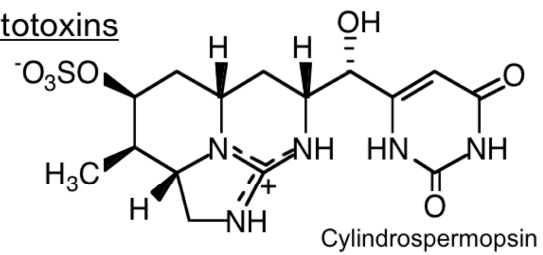
 + cyclic form: critical for toxicity

Over 80 known microcystins, ca 10 known nodularins. Hepatotoxins, protein phosphatase 1 and 2A inhibitors, tumour promoters

Carcinogens? International Agency for Research on Cancer 2006: Microcystin-LR is possibly carcinogenic to humans (Group 2B).

Cytotoxins

Acute neurotoxins: anatoxin-a, anatoxin-a(S), saxitoxin family; chronic neurotoxin: BMAA
Cytotoxins: cylindrospermopsin family

Overview of common freshwater and brackish-water cyanobacterial toxins

Jussi Meriluoto, Åbo Akademi University, Turku, Finland

Toxin	Most common producer genera (not exclusive, great variation in species/strains)	Structure	Toxic action and target organs	Toxicity	Public health concerns	Guidelines	Bioaccumulation	Typical screening and analytical methods
Microcystins (MC)	<i>Microcystis</i> , <i>Planktothrix</i> , <i>Anabaena</i> , <i>Nostoc</i>	Cyclic heptapeptide with unusual amino acids (Adda), >80 variants	Liver toxins, inhibit protein phosphatases 1 and 2A, tumour promoters, at least MC-LR possible human carcinogen	LD ₅₀ for most MCs ca 50-300 µg kg ⁻¹ (mouse, i.p.), oral toxicity weaker	Longterm exposure via drinking water (DW), recreation exposure, common in eutrophic freshwaters, traditional drinking water treatment is inadequate, chemically stable	WHO provisional guideline 1 µg per litre for microcystin-LR in drinking water (the guideline is often interpreted as total microcystin), anti-terrorism legislation	Bioaccumulate in e.g. shellfish and fish liver, and may be present in plants irrigated by contaminated water	ELISA, HPLC-DAD, LC-MS(-MS), protein phosphatase inhibition assay Standard HPLC-based method: ISO 20179:2005
Nodularins (Nod)	<i>Nodularia</i>	Cyclic pentapeptide with unusual amino acids (Adda), ca 10 variants	Liver toxins, inhibit protein phosphatases 1 and 2A, tumour promoters, at least Nod-R is carcinogenic	LD ₅₀ for most Nods 50-300 µg kg ⁻¹ (mouse, i.p.), oral toxicity weaker	Common in the Baltic Sea, recreational exposure, chemically stable	Although not officially regulated, it may be necessary to limit the intake of contaminated seafood	Bioaccumulate in e.g. shellfish, fish liver and eider (bird) liver	ELISA, HPLC-DAD, LC-MS(-MS), protein phosphatase inhibition assay
Cylindrospermopsins (Cyn)	<i>Cylindrospermopsis</i> , <i>Aphanizomenon</i> , <i>Anabaena</i>	Cyclic guanidine alkaloid with an attached uracil ring	Cytotoxins, damage many internal organs (liver, kidney, thymus, heart) inhibit protein synthesis, DNA damage	LD ₅₀ for cylindrospermopsin 2.1 mg kg ⁻¹ (mouse, i.p., 24 h)	Longterm exposure via DW, recreational exposure, present in extracellular form, chemically stable	Guideline 1 µg per litre for cylindrospermopsin in drinking water advocated by Falconer et al.	Some evidence	ELISA, HPLC-DAD, LC-MS-MS

Toxin	Most common producer genera (not exclusive, great variation)	Structure	Toxic action and target organs	Toxicity	Public health concerns	Guidelines	Bioaccumulation	Typical analytical methods
Anatoxin-a, homoana-toxin-a	<i>Anabaena</i> , <i>Aphanizomenon</i> , <i>Planktothrix</i> , <i>Oscillatoria</i>	Secondary amine, low molecular weight alkaloid	Neurotoxic, mimic the effect of acetylcholine in nerve synapses but is not hydrolysed	LD ₅₀ 200-250 µg kg ⁻¹ (mouse, i.p.)	Recreational exposure (but degrades rapidly to non-toxic in natural waters)	Anti-terrorism legislation (DW 3 µg l ⁻¹ has been suggested)		HPLC-DAD, LC-MS-MS, GC methods
Anatoxin-a(S)	<i>Anabaena</i>	Organo-phosphate	Neurotoxic, anticholin-esterase	LD ₅₀ 20 µg kg ⁻¹ (mouse, i.p.)	Recreational exposure			LC-MS-MS
Saxitoxins (known as paralytic shellfish poisons, PSPs, in the marine context)	<i>Anabaena</i> , <i>Aphanizomenon</i> , <i>Cylindro-spermopsis</i> Marine dinoflagellates	Carbamate alkaloids with different sulphation degrees. Some decarbamoyl variants. >20 variants, but not all are present in cyano-bacteria	Neurotoxic, block sodium channels in nerve axons	LD ₅₀ 10 µg kg ⁻¹ (mouse, i.p.) for saxitoxin	Recreational exposure, accumulation in marine shellfish can cause lethal poisonings	Brazil/DW: 3 µg l ⁻¹ Regulated in some seafood (EU: 80 µg / 100 g mussel meat) Anti-terrorism legislation Chemical Weapons Convention	Well-known to accumulate in marine shellfish	Complicated HPLC with derivatisation and fluorescence detection, LC-MS-MS
Beta-N-methyl-amino-L-alanine (BMAA)	All cyanobacteria???	Amino acid	Neurotoxic, associated with chronic degenerative nerve diseases (ALS, Alzheimer)	Chronic toxicity (debated)	Long-term exposure via DW		Known to bioaccumulate in a tropical food web (<i>Nostoc</i> symbiont, cycad seeds, bats, ...)	HPLC with fluorescence detector after derivatisation, LC-MS-MS
Lipopoly-saccharides (LPS, endotoxins)	All cyanobacteria	Lipid-sugar structures in cell walls of Gram neg. bacteria	Irritants, fever, gastrointestinal problems, allergenic responses	Regarded as non-lethal	Exposure to LPS may potentiate the effects of other toxins			<i>Limulus</i> amoebocyte lysate assay (LAL)

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Methods for cyanobacterial toxin analysis

Lisa Spoof, Åbo Akademi University, Turku, Finland

Method	Target toxins	Specificity	Typical sensitivity	Cost	Need of training	Suitable for	Measures	Main downsides	Other comments
Mouse bioassay	"All" acutely acting toxins	Low	LD ₅₀ MC, Nod 25-300 µg kg ⁻¹ LD ₅₀ Cyn 2.1 mg kg ⁻¹ LD ₅₀ Antx-a 200 µg kg ⁻¹ LD ₅₀ Stx 10 µg kg ⁻¹ (i.p.)	Medium	Medium	Screening	Mammalian toxicity through intraperitoneal route (i.p.)	Quick-acting neurotoxins may mask the presence of hepatotoxins, unethical, needs license	Was used frequently in the 1980s but now less frequently for cyanotoxins as better methods have become available
Invertebrate assays	"All" acutely acting toxins	Low	µg ml ⁻¹ levels (<i>Artemia</i>)	Low	Low	Screening	Invertebrate toxicity	Uncertainty about the cause of toxic agent and mechanism	E. g. <i>Artemia</i> bioassay, Thamnotox kit
Protein phosphatase inhibition assay, PPIA	MCs, Nods	High	0.2-1 µg l ⁻¹ for the colorimetric assay, direct analysis of natural waters possible	Low	Medium	Screening	Toxicity directed against a specific protein phosphatase (1 or 2A)	Other compounds inhibiting PP activity may be present in the samples, limited lifetime of reagents	Colorimetric, fluorometric and radiometric assays
Enzyme-linked immunosorbent assay, ELISA	MCs, Nods, Stxs and Cyn	High	0.1 µg l ⁻¹ , direct analysis of natural waters possible	Medium	Low	Screening/ Analysis	Presence of total toxin from a specific class with quantitation	Different cross-reactivities of MCs affect the quantitation	Most popular screening method. Commercial kits in two formats: tubes and 96-well plates
High-performance liquid chromatography with UV (or FL) detection	"All"	High	ng levels per injection	Medium	Medium	Analysis	Presence of individual toxins with quantitation, UV spectral based detection of unknown MCs may be possible	Needs authentic standards. Small toxin peaks can be overlooked due to coeluting impurities	MCs have an ISO standard, number 20179:2005. Antx-a and especially Stxs are typically derivatised to enable fluorescence detection.

Method	Target toxins	Specificity	Typical sensitivity	Cost	Need of training	Suitable for	Measures	Main downsides	Other comments
Liquid chromatography - mass spectrometry, LC-MS or LC-MS-MS	"All"	Highest	Strongly dependent on the instrument used. Direct analysis of natural waters may be possible	Expensive	High	Analysis	Presence of individual toxins with quantitation	Needs authentic standards. Ion suppression or enhancement by co-eluting matrix substances affect quantitation	Electrospray ionisation, ESI, typical. MS-MS fragmentation patterns can make the detection of unknown toxins possible
Matrix-assisted laser desorption/ionisation time-of-flight MS, MALDI-TOF MS	"All"	High	Strongly dependent on the instrument used	Expensive	Medium	Analysis	Qualitative or semi-quantitative, presence of individual toxins	Off-line technique, limited quantitation possibilities	Post-source decay for identification of unknown metabolites based on fragmentation patterns
Thin-layer chromatography TLC	"All"	Medium	In optimal conditions 10 ng /spot MCs	Low	Low	Analysis	Mainly qualitative	Needs authentic standards, complex samples difficult to analyse	
Capillary electrophoresis, CE	"All"	High	Comparable to HPLC	Expensive	High	Analysis	Presence of individual toxins with quantitation	Technical obstacles related to sensitivity and robustness	Optical and MS detection possible
3-methoxy-2-methyl-4-phenylbutyric acid, MMPB, quantitation	MCs, Nods	High	Medium	Medium	High	Screening/ Analysis	Detects total MCs and Nods. No toxin standards needed (only MMPB standard)	Time-consuming. Non-toxic variants of Adda or free Adda lead to overestimation	Oxidative cleavage of Adda followed by chromatographic detection

Antx-a = anatoxin-a, Cyn = cylindrospermopsin, MC = microcystin, Nod = nodularin, Stx = saxitoxin

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- Chorus, I., Bartram, J. (Eds.): Toxic Cyanobacteria in Water: A Guide to Their Public Health Consequences, Monitoring, and Management. E & FN Spon, London (1999).
- Fastner, J., Codd, G. A., Metcalf, J. S., Wölk, P., Wiedner, C., Utkilen, H.: An international intercomparison exercise for the determination of purified microcystin-LR and microcystins in cyanobacterial field material. *Anal. Bioanal. Chem.* **374**, 437-444 (2002).
- Hawkins, P. R., Novic, S., Cox, P., Neilan, B. A., Burns, B. P., Shaw, G., Wickramasinghe, W., Peerapornpisal, Y., Ruangyuttikarn, W., Itayama, T., Saitou, T., Mizuochi, M., Inamori, Y.: A review of analytical methods for assessing the public health risk from microcystin in the aquatic environment. *J. Water SRT - Aqua* **54**, 509-518 (2005).
- Marsálek, B., Bláha, L.: Comparison of 17 biotests for detection of cyanobacterial toxicity. *Environ. Toxicol.* **19**, 310-317 (2004).
- Meriluoto, J.: Toxins of freshwater cyanobacteria (blue-green algae). In Bogusz, M. J. (Ed.), *Handbook of Analytical Separations*, vol. 2: Forensic Science, pp. 359-390. Elsevier, Amsterdam (2000).
- Meriluoto, J., Codd, G. A. (Eds.): *TOXIC: Cyanobacterial Monitoring and Cyanotoxin Analysis*. Åbo Akademi University Press, Turku (2005).
- Metcalf, J. S., Codd, G. A.: Analysis of cyanobacterial toxins by immunological methods. *Chem. Res. Toxicol.* **16**, 103-112 (2003).
- Osswald, J., Rellán, S., Gago, A., Vasconcelos, V.: Toxicology and detection methods of the alkaloid neurotoxin produced by cyanobacteria, anatoxin-a. *Environ. Int.* **33**, 1070-1089 (2007).
- Rapala, J., Erkoma, K., Kukkonen, J., Sivonen, K., Lahti, K.: Detection of microcystins with protein phosphatase inhibition assay, high-performance liquid chromatography-UV detection and enzyme-linked immunosorbent assay. Comparison of methods. *Anal. Chim. Acta* **466**, 213-231 (2002).
- van Apeldoorn, M. E., van Egmond, H. P., Speijers, G. J., Bakker, G. J.: Toxins of cyanobacteria. *Mol. Nutr. Food Res.* **51**, 7-60 (2007).
- Spoof, L., Vesterkvist, P., Lindholm, T., Meriluoto, J.: Screening for cyanobacterial hepatotoxins, microcystins and nodularin in environmental water samples by reversed-phase liquid chromatography-electrospray ionisation mass spectrometry. *J. Chromatogr. A* **1020**, 105-119 (2003).
- Törökne, A., Asztalos, M., Bánkiné, M., Bickel, H., Borbély, G., Carmeli, S., Codd, G.A., Fastner, J., Huang, Q., Humpage, A., Metcalf, J., Rábai, E., Sukenik, A., Surányi, G., Vasas, G., Weiszfeiler, V.: Interlaboratory comparison trial on cylindrospermopsin measurement. *Anal. Biochem.* **332**, 280-284 (2004).

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Reference to the ISO (International Organization for Standardization, <http://www.iso.org>) method for determination of microcystins

- ISO 20179:2005. Water quality -- Determination of microcystins -- Method using solid phase extraction (SPE) and high performance liquid chromatography (HPLC) with ultraviolet (UV) detection.