INSTRUCTION

- Children may not be treated with this device.
- Keep the device away from children and youths below the age of 16. Children may swallow small parts and suffocate. Children may also injure themselves when using device.
- Do not apply stimulation over your neck or mouth because this could cause severe muscle spasms resulting in closure of your airway, difficulty in breathing, or adverse effects on heart rhythm or blood pressure.
- Do not apply stimulation across your chest because the introduction of electrical current into the chest may cause rhythm disturbances to your heart, which could be lethal.
- Do not apply stimulation over, or in proximity to, cancerous lesions.
- Do not apply stimulation in the presence of electronic monitoring equipment (e.g., cardiac monitors, ECG alarms), which may not operate properly when electrical stimulation device is in use.
- Do not use the stimulator in the vicinity of shortwave or microwave therapy equipment, since this may affect the output power of the stimulator.
- Never use near the heart. Gel pad should never be placed anywhere on the front of the thorax (marked by ribs and breastbone), but above all not on the two large pectoral muscles. Here it can increase the risk of ventricular fibrillation and lead to cardiac arrest.
- Keep gel pad electrodes separate during treatment, electrode in contact with other could result in improper stimulation or skin burns.

Waveform:	Biphasic square wave.
P.R (Frequency):	80Hz
Burst Frequency:	2Hz

Waveform:	Biphasic square wave.
P.R (Frequency):	2Hz/ 80Hz
P.W (Pulse width):	250us/ 200us
Output order	2Hz/250us first output 3 sec. and then output 3 sec.
Output Voltage:	80Hz/200us.
Output Current:	0 to 70mA(500 Ω Load), adjustable


Waveform:	Biphasic square wave.
Wave characteristic:	All program of EMS mode output wave characteristic
P.R. (Frequency):	are same.
For program E2: 20Hz	0 to 35Vpp(500 Ω Load), adjustable
P.W. (Pulse width):	250uS
Ramp up time	2sec.
Ramp down time	2sec.
Contract time (on time)	5sec.
Relax time (off time)	10sec.
Output Voltage:	0 to 35Vpp (500 Ω Load), adjustable
Current:	0 to 70mA(500 Ω Load), adjustable

P.W. (pulse width)	200, 250us
P.R. (Frequency)	2, 80Hz (Hz=vibration per second)
Burst frequency	2Hz
Output characteristics	Constant Voltage (CV)
P1	The pulse rate and pulse width output will be constant based on the design value.
P3	2Hz/250us first output 3 sec. and then output 3 sec. 80Hz/200us.

Problem	Possible cause	Solution/Remark
The unit has no response, or no LED lights when presses the buttons.	<ol style="list-style-type: none"> 1. If the battery has used up or wrong polarity. 2. If there's a loose connection in the battery interface. 3. If the specification of battery is incorrect. 	<ol style="list-style-type: none"> 1. Replace the new battery or install the battery with correct polarity. 2. Check and clean the battery contact points. 3. Replace the battery with the correct specification.
The stimulation is weak even the intensity is high	<ol style="list-style-type: none"> 1. If the pad is not fully on the skin or not well connected to the buckles. 2. The pad is worn 3. If the gel pad has been used too long time, and lost with adhesive. 	<ol style="list-style-type: none"> 1. Reconnect the gel pad or make sure the gel pad is placed flat against your skin. 2. Renew the gel pad. 3. Replace the new gel pad.
Halt automatically in the treatment	<ol style="list-style-type: none"> 1. If the gel pad loosen from the body 2. If the buckles are not connected well. 3. If the battery is used up. 	<ol style="list-style-type: none"> 1. Check and stick the gel pad on the well. 2. Check and connect the cable. 3. Replace the battery.

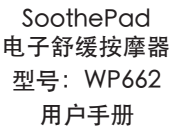	Attention: see Instructions for use!
	Applied part of type BF
	Disposal in accordance with Directive 2002/96/EC (WEEE)
	Complies with the European Medical Device Directive (93/42/EEC). Notified body is SGS United Kingdom Ltd.
	Manufacturer information: The manufacturer Famidoc Technology Co., Ltd. Address: No. 212 Yilong Road, Hexi Industrial Zone, Jinxiu Changan Town, Dongguan 523853 Guangdong Province, P.R.China
	Serial number
 	Name: Shanghai International Holding Corp. GmbH (Europe) Add: Eifflasse 80, 20537 Hamburg, Germany Tel: +49-40-2513175 Fax: +49-40-255726 E-mail: shholding@hotmail.com

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