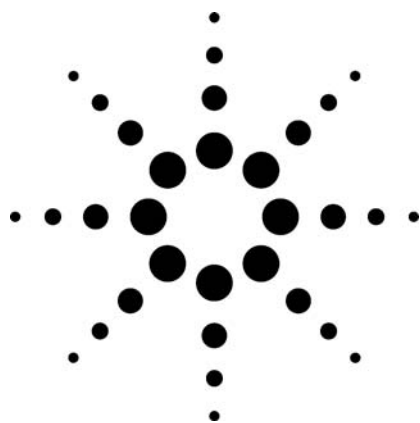




Process Monitoring of Bisphenol-A in Industrial Feedstock using High Throughput HPLC

Application

Process Control

Authors

John W. Henderson, Jr. and William J. Long
Agilent Technologies
2850 Centerville Road
Wilmington, DE 19808-1610
USA

Stephen Fetsko
Badger Licensing, LLC
1 Main Street
Cambridge, MA 02142
USA

Abstract

The chromatographic method used to monitor the Bisphenol-A manufacturing process was improved using Agilent RRHT Eclipse XDB-C18 columns. These columns use 1.8- μm particles versus conventional 3.5- μm or 5- μm particles. The improved method allowed seven times faster analyses, improved resolution, and higher sensitivity.

Introduction

Bisphenol-A (Figure 1) is a highly versatile material used to manufacture many modern products. It is also known as 4,4'-Isopropylidenediphenol, 4,4'-(1-Methylethylidene) bisphenol, or simply BPA.

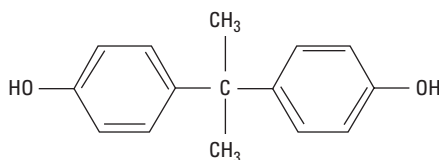


Figure 1. Bisphenol A

Every year, 2.8 million tons of BPA are produced. BPA is a building block for polycarbonate plastic and epoxy resins. Polycarbonate plastic is prized for its scratch resistance, optical clarity, and heat and electrical resistance. Because of these attributes, it is used for eyewear, CD/DVD disks, electronics, and food and drink containers. Epoxy resins are used for protective coatings because of their combination of inertness, chemical resistance, adhesion, and formability. For example, metal food cans are lined to protect taste. Epoxy resins are also used as a component in dental sealants and as a component in dental composites providing an alternative to mercury amalgam in veneers and fillings. Other uses include fungicides, polymer antioxidants, and components in automobiles and appliances.

BPA is produced through an acid-catalyzed condensation reaction of phenol with acetone. During condensation, a number of phenol-based byproducts are also formed. HPLC is used to determine the composition of many of the process streams in a commercial BPA plant.

Here we describe the use of new HPLC column technology for the possible improvement to one of the HPLC methods used in a commercial BPA facility.



Agilent Technologies

Method Optimization and Scalability

The existing HPLC method was proven and robust; however, it was complicated. We sought a similar chromatogram, based on the original method, but using simpler method parameters. Because of the challenge of changing many chromatographic parameters, essentially redeveloping the method, we chose a 4.6×50 mm, $1.8\text{-}\mu\text{m}$ Eclipse XDB-C18 column for experiments to reduce the time required. Smaller particles packed in shorter columns increase the speed of analysis and still provide enough efficiency to maintain resolution equivalent to longer columns packed with larger particles. After several trials, we developed a method that produced a chromatogram similar to the original. The short analysis time is a major advantage of Rapid Resolution High Throughput (RRHT) technology. Whereas a handful of experimental runs would take an entire work day using a typical analytical-sized column (50 min/run), the series of runs took about an hour (7.5 min/run), using an RRHT column.

We incrementally scaled up to a 4.6×250 mm column. Figure 2 shows an overlay of the sample analyzed by three 4.6-mm id columns of different lengths and particle sizes. Injection volume was also changed proportionally to length. The smaller ZORBAX particles speed up the analysis while maintaining resolution. In fact, resolution increased when using the RRHT columns despite their shorter length.

One reason this method can be easily scaled (up or down) is the uniform spherical Eclipse XDB-C18 packing. It has a proprietary engineered particle size distribution, based on ZORBAX silica with a controlled surface area and pore size. The robust proprietary packing material and proven column manufacturing techniques consistently yield reproducible columns with similar chromatographic performance, independent of the column dimensions.

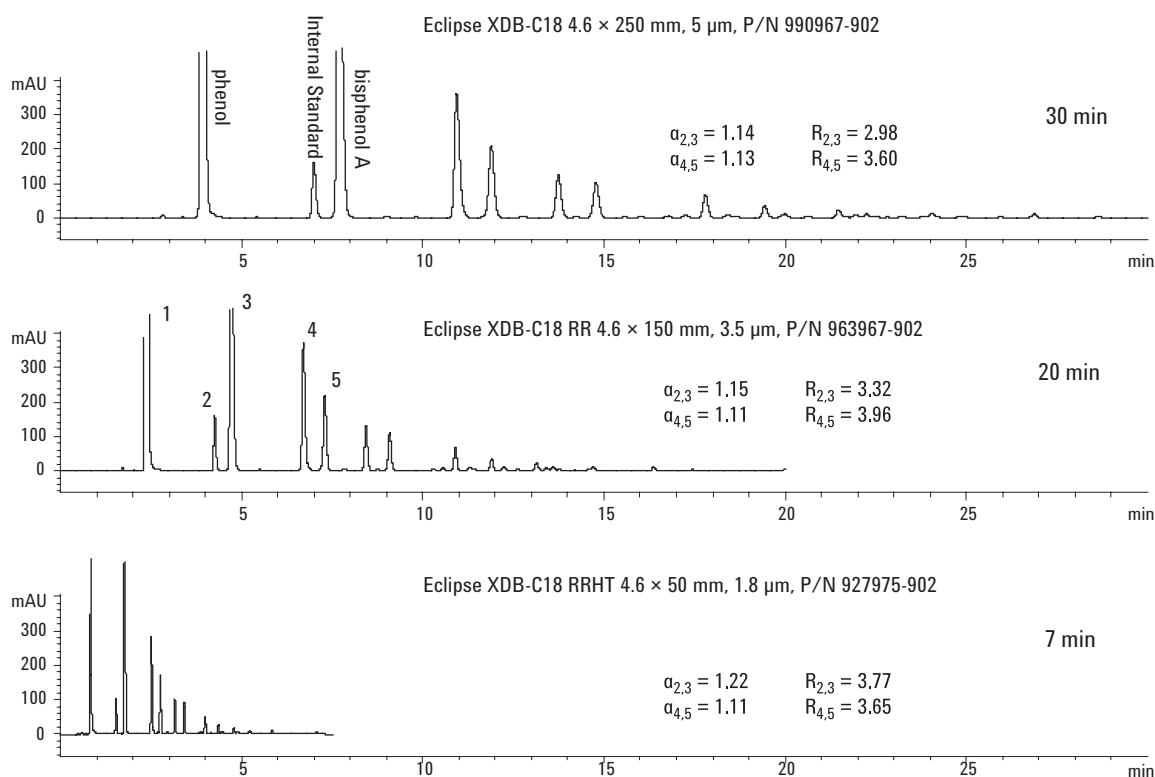


Figure 2. RRHT column configuration increased both speed and resolution.

Particle size does influence resolution. The influence can be noticed when comparing columns of identical dimensions, packed with three different particle sizes. Figure 3 shows the shortened Bisphenol-A analysis using different particle-sized Eclipse XDB-C18 columns. Resolution (R_s) is related to selectivity (α), efficiency (N) and retention (k'):

$$R_s = (1/4)(\alpha-1) \sqrt{N} [k'/(1+k')]$$

Factors affecting the selectivity term (stationary phase, mobile phase) and retention term (mobile phase, temperature) are constant for the three

chromatograms. The efficiency term is influenced by column length, linear velocity of the mobile phase (both constant), and particle size (varied in Figure 3). N increases as particle size decreases. In Figure 3 the selectivity factors (α) and retention remain about the same, but resolution actually increases. The increase in resolution due to the decrease in particle size highlights the advantage of using smaller particles. The similar selectivity and retention highlight the suitability of ZORBAX Eclipse XDB-C18 columns for scaling methods, especially to more rapid, high-throughput methods.

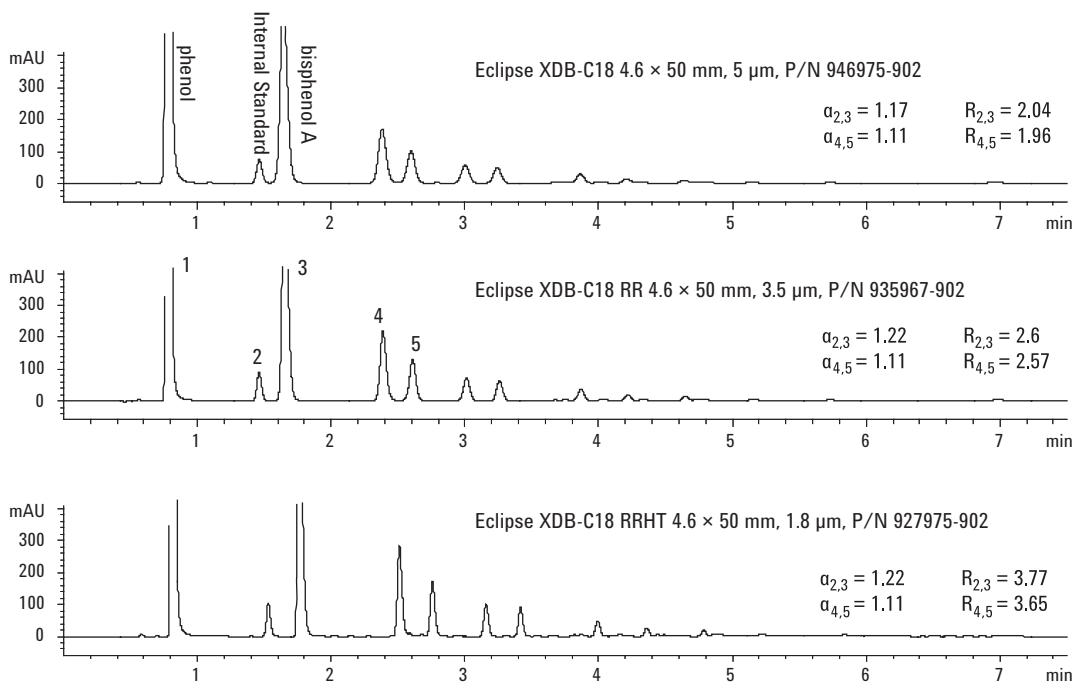


Figure 3. Effect of particle size on resolution and selectivity.

Comparing the Existing Method to the RRHT Method

Figure 4 compares the original BPA separation to the RRHT separation. The top chromatogram is an example of the analysis using the original commercial method, and the bottom is an example of the process sample analyzed with the RRHT method. The method developed with the new column technology clearly increases productivity.

Analysis time is reduced at least six-fold; solvent consumption is reduced about 12.5 times, from 100 mL/analysis to only 7.5 mL/analysis. Interestingly, the peak shape of Bisphenol-A is more symmetrical using Eclipse XDB-C18 as compared to the current C18 column used in the original analysis. The more Gaussian peak shape eluted by the Eclipse XDB-C18 column is important for accurate quantification. Other method improvements such as a simplified gradient and a binary mobile phase are listed in Table 1.

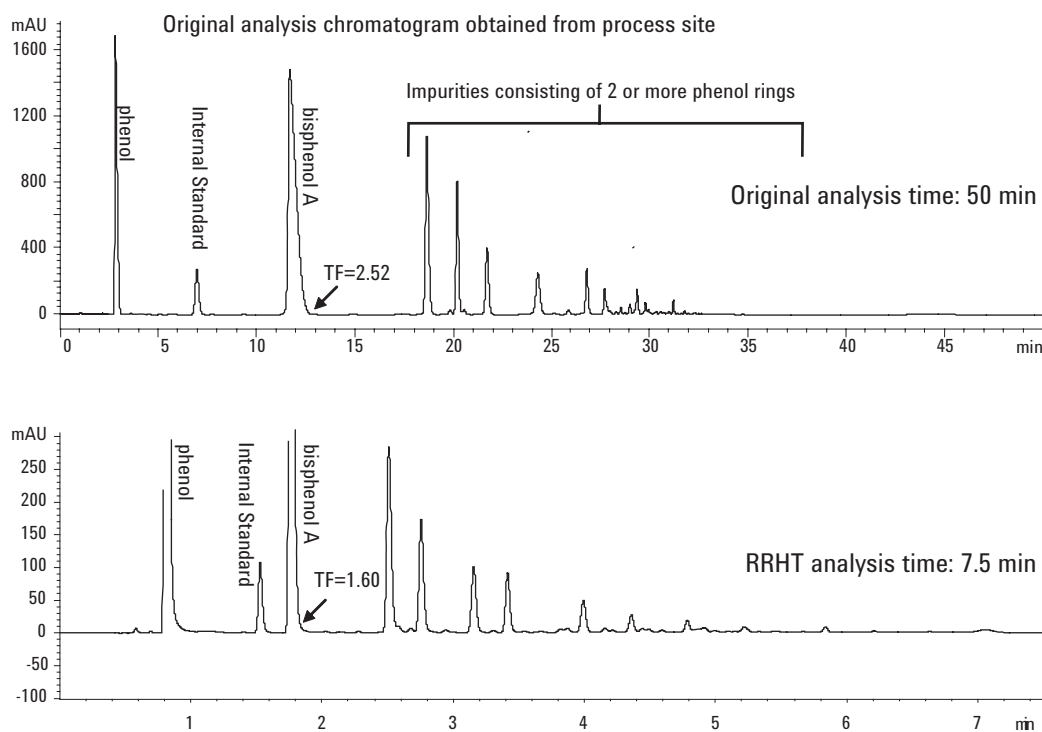


Figure 4. Comparison of methods; original to RRHT.

Table 1. Current and Improved Method Parameters

Original		RRHT	
• Column: Supelco LC-18, 4.6 × 250 mm, 5 µm • Mobile phase: A: 0.025% H₃PO₄, B: ACN, C: MeOH • Flow: 2 mL/min • Temperature: 35 °C • Sample size: 20 µL • Gradient: segmented, has isocratic holds		• Column: ZORBAX XDB-C18, 4.6 × 50 mm, 1.8 µm • Mobile phase: A: 0.1% formic acid, B: ACN: MeOH (200:800) • Flow: 1 mL/min • Temperature: 25 °C • Sample size: 2 µL • Gradient: linear, no isocratic holds	
Time	% A:B:C	Time	% B
0	65:25:10	0	60
13	65:25:10	6	95
18	50:40:10	6.01	60
23	50:40:10	8	60
27	30:50:20		
32	0:70:30		
35	0:70:30		
36	0:60:40		
40	0:50:50		
43	0:20:80		
48	65:25:10		

Conclusion

Converting an existing method to a high-throughput method is one way to improve lab productivity. Using RRHT columns initially for method development also improves productivity. Eclipse XDB-C18 RRHT columns are a good choice for converting existing C18 methods into high-throughput methods. Smaller particles packed into shorter columns provide comparable resolution to larger particles packed into longer columns in a fraction of the time. RRHT columns are advantageous for gradient method development because gradient re-equilibration is time-consuming and often overlooked in the total analysis time. Methods developed on Agilent RRHT columns can be scaled easily because of the highly uniform particles, bonded phase chemistry, and column manufacturing techniques. An existing method developed on a “traditional analytical-sized” column was easily converted to a high throughput method using an Eclipse XDB-C18 RRHT column. The method was incrementally scaled up to an analytical-sized column, and it performed with predictable results

on various column dimensions and particle sizes. The predictability of the results supports Eclipse XDB-C18 RRHT columns’ ability to easily improve applications and transfer them into high-throughput and high-resolution applications.

For More Information

For more information on our products and services, visit our Web site at www.agilent.com/chem.

Agilent shall not be liable for errors contained herein or for incidental or consequential damages in connection with the furnishing, performance, or use of this material.

Information, descriptions, and specifications in this publication are subject to change without notice.

© Agilent Technologies, Inc. 2006

Printed in the USA
July 26, 2006
5989-5231EN



Agilent Technologies