
Operating Instructions

Household
Electric Massager
EV293

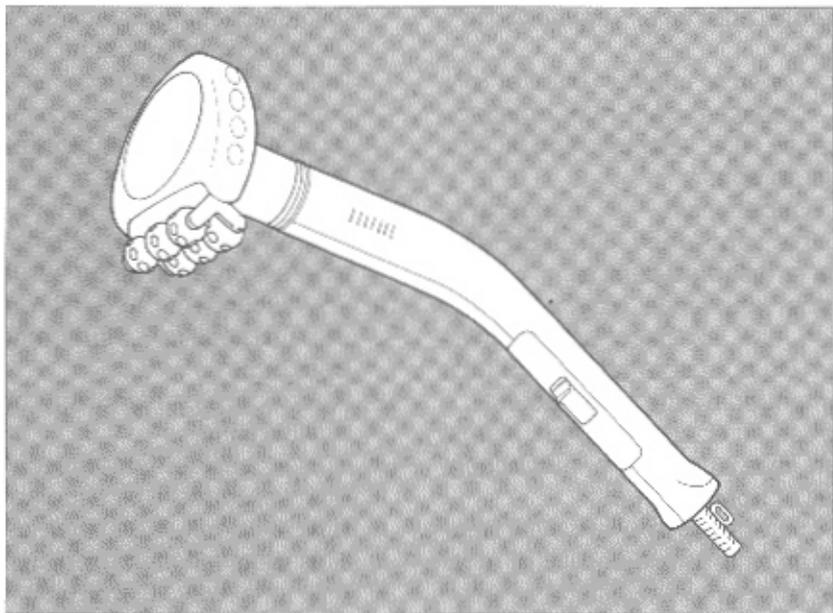Panasonic

Before operating this unit, please read these instructions completely.

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING THE MASSAGER.


DANGER – To reduce the risk of electric shock:

1. Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
2. Do not reach for an appliance that has fallen into water. Unplug immediately.
3. Do not use while bathing or in a shower.
4. Do not place or store the appliance where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.

WARNING To reduce the risk of burns, fire, electric shock, or injury to persons:

1. An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before adding or removing accessories.
2. Do not operate under blanket or pillow. Excessive heat may develop and cause fire, electric shock, or injury to persons.
3. Close supervision is necessary when this appliance is used by invalids or disabled persons.
4. Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
5. Do not use on or near children. Keep children away from the massager whether or not it is in use.
6. Never operate this appliance if it has a damaged cord or plug, if it is not working properly, or if it has been dropped or damaged, or dropped into water. Call your nearest authorized servicenter for repair.
7. Do not carry this appliance by power cord or use cord as a handle.
8. Keep the cord away from hot surfaces.
9. Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair and the like. Never operate on a soft surface such as a bed or couch where the air openings may be blocked.
10. Do not use outdoors.
11. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
12. To disconnect, turn all controls to the off position, then remove plug from outlet.
13. Never drop or insert any object into any opening.
14. To reduce the risk of electric shock, this appliance has a polarized plug (one blade is wider than the other). This plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician to install the proper outlet. Do not alter the plug in any way.

SAVE THESE INSTRUCTIONS

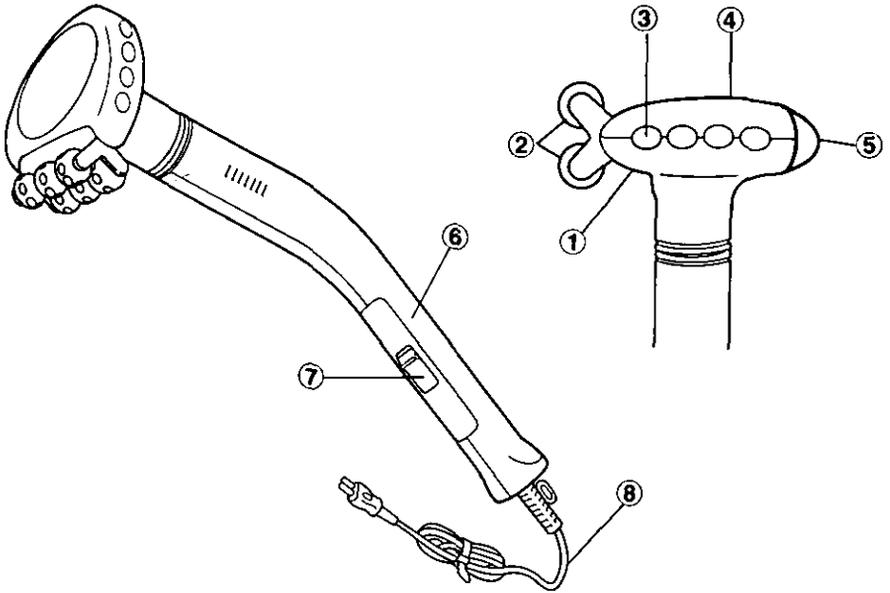Parts Identification


- | | |
|-------------------------|--------------------------|
| ① Ventilation openings | ⑤ Point massage |
| ② Double roller massage | ⑥ Grip |
| ③ Raised node massage | ⑦ 2 speed control switch |
| ④ Rounded top surface | ⑧ Power cord |

How to use

1. Make sure the switch is in the off position and hold the main body firmly before the unit is plugged into the power source. Plug the unit in and turn on the switch.
2. Begin your massage on low speed. Place the roller surface against the skin, working around the desired area. We recommend use of high speed only after becoming accustomed to lesser intensity.
3. Next, using the point massage surface, go over the stiffed or previously treated area. Set the speed controller to a comfortable level. (Do not press the massage surface too hard or bent it excessively.)
4. Finally, using the roller massage surface, gently massage for the finishing effect with the speed controller set to low.
5. Turn the switch off and unplug the unit.


(For a pleasant massage)

1. With the roller surface
move over a wide massage area
move the head around in big, wide circles
2. With the pointed surface
push down over the area of
soreness or stiffness
(for about five to seven seconds
per area)
for a stronger vibrating effect, hold
the appliance with both hands
3. With the roller surface
finish by moving the head around
in big circles uniformly all around the
area

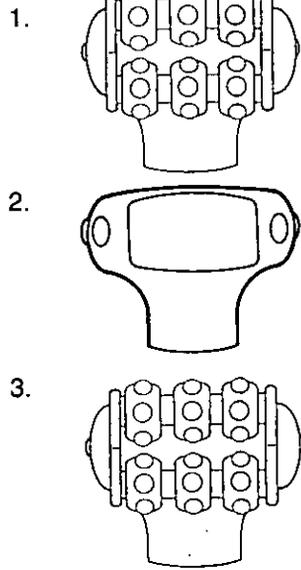Examples of How to Use

1. (Back)
With the roller surface, leisurely massage all around the area.
(Waist)
With the pointed surface, leisurely massage all around the area.
Press down on areas of muscle stiffness.
2. (Neck)
With the pointed surface, lightly press down along the neck muscles.
(Shoulder)
With the pointed surface, push down with both hands.
With the roller surface, move along the shoulder.
(Arm)
With the roller surface, move in a flowing motion down along the muscles.
3. (Leg)
With the roller surface, move in a flowing motion down along the muscles.
(Calf)
With the roller surface, move in a flowing motion down along the muscles. In case of calf pain, do not use without consulting a physician.
(Sole)
With the pointed surface, leisurely massage around the area.

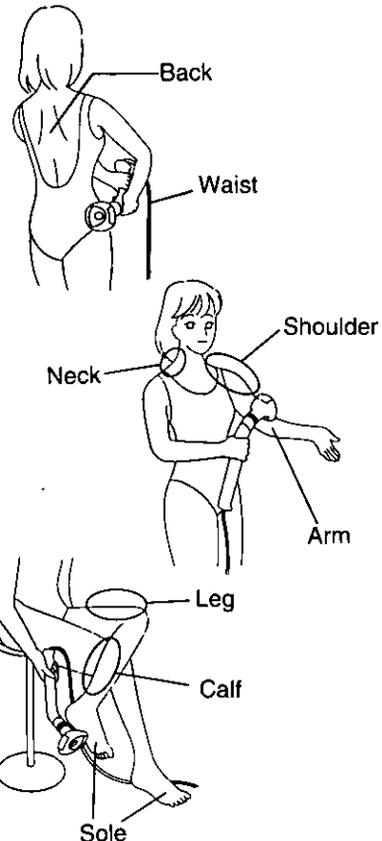Warning

1. This device should not be used on any injured, swollen, or inflamed areas, or wounds, without consulting your physician. Massage should be pleasant and comfortable. If you feel pain or discomfort, discontinue use immediately and consult a physician.
2. Avoid using on the head, abdomen, joints, or thyroid gland (just below Adam's Apple).

Caution

1. Do not hold the massager by the cord, or pull the cord sharply.
2. Do not attempt to make any repairs yourself. Do not touch the internal workings of the massager.
3. If the massager becomes soiled, clean with a cloth dampened with a mild cleanser. Do not use thinner, benzine or similar solvents. Do not immerse in liquid.
4. Avoid storage in any location where there is excessive moisture, such as the bathroom.
5. Do not bend or turn massager head, or expose it to other excessive force. Damage may result.
6. Roller can be removed for easy cleaning by pulling roller straight out. After cleaning with a damp cloth, push roller back in, making sure it firmly snaps in. Do not use massager if roller is not properly in place.

Specifications

Power source: AC120V, 60Hz

Power consumption: 10W

Speed: Approx. 6000rpm (high)

Approx. 5000rpm (low)

Dimensions: 15.2inch × 4.7inch head diameter

Weight: 19.4oz.

**MATSUSHITA ELECTRIC CORPORATION
OF AMERICA**
Home Products/Commercial Systems
Sales Group.
One Panasonic Way
Secaucus, NJ 07094

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