
Operating Instructions

Household
Electric Massager
EV328

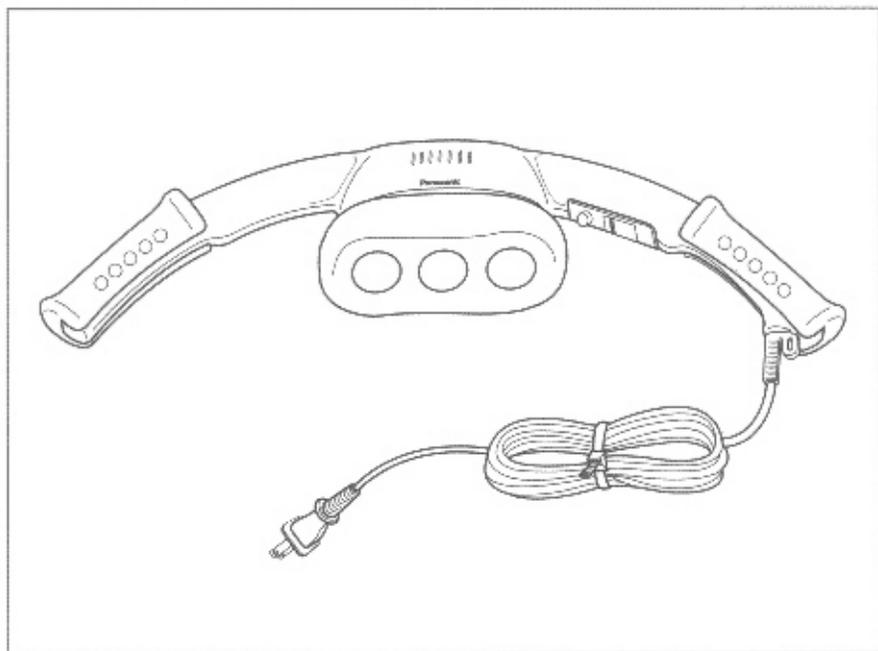Panasonic

Before operating this unit, please read these instructions completely.

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following.

Read all instructions before using the massager.


DANGER - To reduce the risk of electric shock:

1. Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
2. Do not reach for an appliance that has fallen into water. Unplug immediately.
3. Do not use while bathing or in a shower.
4. Do not place or store the appliance where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.

WARNING - To reduce the risk of burns, fire, electric shock, or injury to persons:

1. An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
2. Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
3. Close supervision is necessary when this appliance is used by invalids or disabled persons.
4. Do not use on or near children. Keep children away from the appliance whether or not it is in use.
5. Use this appliance only for its intended use as described in this manual.
6. Never operate this appliance if it has a damaged cord or plug, if it is not working properly, or if it has been dropped or damaged, or dropped into water. Call your nearest authorized service center for repair.
7. Do not carry this appliance by power cord or use cord as a handle.
8. Keep the cord away from heated surfaces.
9. Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair and the like. Never operate on a soft surface such as a bed or couch where the air openings may be blocked.
10. Do not use outdoors.
11. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
12. To disconnect, turn all controls to the off position, then remove plug from outlet.
13. Never drop or insert any object into any opening.
14. To reduce the risk of electric shock, this appliance has a polarized plug (one blade is wider than the other). This plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician to install the proper outlet. Do not change the plug in any way.
15. Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.

SAVE THESE INSTRUCTIONS

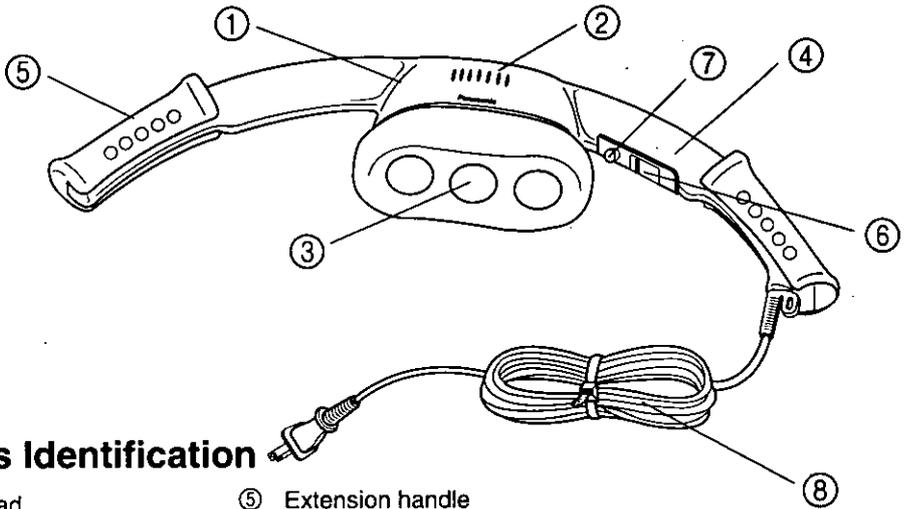Parts Identification


- | | |
|------------------------|-------------------------------|
| ① Head | ⑤ Extension handle |
| ② Ventilation openings | ⑥ Power switch (high-low-off) |
| ③ Heated surface | ⑦ Heat switch |
| ④ Arm | ⑧ Power cord |

How to use

Using the extension handle


1. Extend

Pull the extension handle until it clicks into place.

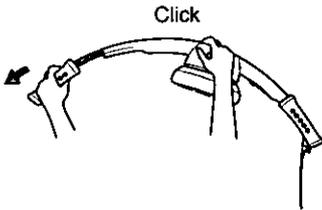2. Massage


Grip both ends.

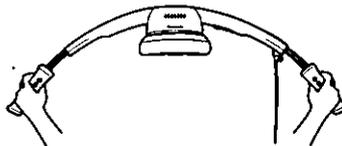3. Retract.

Caution:

Do not put the power cord and fingers between the arm and the extension handle when you retract it.

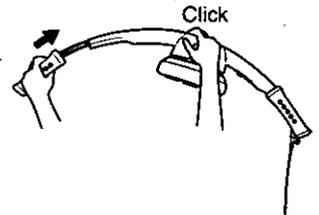Using Heat

To activate the heating mechanism, press the red heat button after turning the massager on.

Key Benefit:

The grip extends out to allow for a perfectly fitting massage in difficulty to reach areas such as the lower back, shoulders and calves.

Warning

- This device should not be used on any injured, swollen, inflamed areas, or wounds, without consulting your physician. Massage should be pleasant and comfortable. If you feel pain or discomfort, discontinue use immediately and consult a physician.

Cautions

1. If you are pregnant, ill or suffered from injury, consult your physician before using the appliance.
2. Avoid using on the head, abdomen, joints, or thyroid gland (just below Adam's Apple).
3. Do not hold the appliance by the cord or pull the cord sharply.
4. Do not attempt to make any repair by yourself. Do not touch the internal working of the appliance.
5. If the massager becomes soiled wipe with a cloth dampened with a mild cleanser. Do not use thinner, benzine or similar solvents. Do not immerse in liquid.
6. Avoid storage in any location where there is excessive moisture, such as the bathroom.
7. Do not use the appliance for more than 15 minutes at a time, and do not use it continuously for more than 3 minutes on the same parts of the body. (Prolonged usage can have reverse effects on muscles at times.)
8. When extending or retracting the extension handle, be sure to firmly grip the head of the massager.




Specifications

Power source	: AC 120V 60Hz
Power consumption	: 14W
Speed	: Approx. 6300 min ⁻¹ (high) Approx. 5000 min ⁻¹ (low)
Dimensions	: 160 X 595 X 105 mm (6.3" x 23.4" x 4.1")
Weight	: 910g (2.01lb.)

Massaging tips


(For a pleasant massage)

1. Use the massager to loosen stiff or tense muscles.
 - Move the massager in broad strokes around the stiff or tense area.
 - Keep the massager at the low setting until you become accustomed to it.
 - Avoid using the massager directly on exposed skin and, if you find it painful at any point, place a towel or other such item over whatever you already have on.
2. If you prefer a more concentrated massage on the area where you feel stiffness, use the high setting.
 - Focus on the area where you feel the stiffness or tension.
3. Finish up massaging the afflicted area by moving the massager around the area at the low setting.

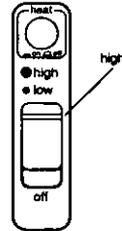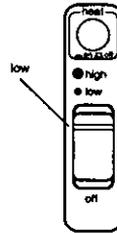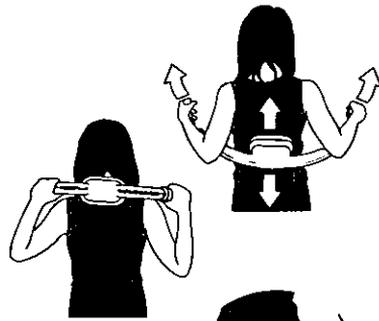Examples of How to Use

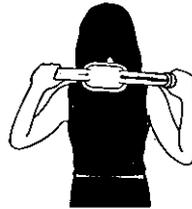1. (Shoulder)
Apply the massager steadily and evenly.
The key to achieving a satisfying massage is to place the head squarely on the afflicted.
2. (Lower Back)
Apply the massager steadily and evenly.
3. (Thigh, Calf)
Run the massager in a flowing motion down along the muscles.
(Foot)
Apply the massager steadily and evenly to the sole of the foot.

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División de Matsushita Electric Corporation of America