# IKA® RV 10 control

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	Seite		
CE-Declaration of conformity	3	Interfaces and outputs	42
Safety instructions	32	Commissioning	44
Correct use	33	Maintenance and cleaning	54
Unpacking	33	Accessories	55
Useful information	34	Error codes	57
Setting up	35	Technical data	58
Drive RV 10 control		Solvent table (excerpt)	58
Heating bath		Warranty	58
Glassware			
Hose system			

## Safety instructions

### For your protection

#### • Read the operating instructions in full before starting up and follow the safety instructions

- Keep the operating instructions in a place where they can be accessed by everyone.
- Ensure that only trained staff work with the appliance.
- Follow the safety instructions, guidelines, occupational health and safety and accident prevention regulations. **When working under a vacuum in particular!**
- Wear your personal protective equipment in accordance with the hazard category of the medium to be processed. Otherwise there is a risk of:
  - splashing liquids,
  - body parts, hair, clothing and jewellery getting caught,
  - injury as a result of glass breakage.
- **CAUTION!** Inhalation of or contact with media such as poisonous liquids, gases, spray mist, vapours, dusts or biological and microbiological materials can be hazardous to user.
- Set up the device in a spacious area on an even, stable, clean, non-slip, dry and fireproof surface.
- Ensure that there is sufficient space above the device as the glass assembly may exceed the height of the device.
- Prior to each use, always check the device, accessories and especially the glass parts for damage. Do not use damaged components.
- Ensure that the glass assembly is tension-free! Danger of cracking as a result of:
  - stress due to incorrect assembly,
  - external mechanical hazards,
  - local temperature peaks.
- Ensure that the stand does not start to move due to vibrations respectively unbalance.
- Beware of hazards due to:
  - flammable materials,
  - combustible media with a low boiling temperature,
  - glass breakage.
- **CAUTION!** Only process and heat up media that has a flash point higher than the adjusted safe temperature limit of the heating bath that has been set.  
The safe temperature limit of the heating bath must always be set to at least 25 °C lower than the fire point of the media used.
- Do **not** operate the appliance in explosive atmospheres, with hazardous substances or under water.
- Only process media that will not react dangerously to the extra energy produced through processing. This also applies to any extra energy produced in other ways, e.g. through light irradiation.
- Electrostatic discharges could pose a direct danger.
- Tasks with the device must only be performed when operation is monitored.
- Operation with excess pressure is not permitted (for cooling water pressure see "Technical Data").
- Do not cover the ventilation slots of the device in order to ensure adequate cooling of the drive.
- There may be electrostatic discharges between the medium and the drive which could pose a direct danger.
- The appliance is not suitable for manual operation.
- Safe operation is only guaranteed with the accessories described in the "Accessories" chapter.
- Refer to the operating instructions for the HB 10 control heating bath.
- Refer to the operating instructions for the accessories, e.g. vacuum pump.
- Position the positive pressure outlet of the vacuum pump under a fume hood.

- Use the HB 10.2 cover or the HB 10.1 spray guard.
- Only use the device under an all side-closed exhaust, or a comparable protective device.
- Adapt the quantity and the type of distilland to the size of the distillation equipment. The cooler must work properly. Monitor the coolant flow rate at the cooler outlet.
- The glass equipment must always be ventilated when working under normal pressure (e.g. open outlet at cooler) in order to prevent a pressure build-up.
- Please note that dangerous concentrations of gases, vapours or particulate matter can escape through the outlet at the cooler. Take appropriate action to avoid this risk, for example, downstream cold traps, gas wash bottles or an effective extraction system.
- Evacuated glass vessels must not be heated on one side; the evaporating flask must rotate during the heating phase.
- The glassware is designed for operation under a vacuum of up to 10 mbar. The equipment must be evacuated prior to heating (see chapter "Commissioning"). The equipment must only be aired again after cooling. When carrying out vacuum distillation, uncondensed vapours must be condensed out or safely dissipated. If there is a risk that the distillation residue could disintegrate in the presence of oxygen, only inert gas must be admitted for stress relief.
- **CAUTION!** Avoid peroxide formation. Organic peroxides can accumulate in distillation and exhaust residues and explode while decomposing!  
Keep liquids that tend to form organic peroxides away from light, in particular from UV rays and check them prior to distillation and exhaust for the presence of peroxides. Any existing peroxides must be eliminated. Many organic compounds are prone to the formation of peroxides e.g. dekalin, diethyl ether, dioxane, tetrahydrofuran, as well as unsaturated hydrocarbons, such as tetralin, diene, cumene and aldehydes, ketones and solutions of these substances.
- **DANGER OF BURNING!** The heating bath, tempering medium, evaporator piston and glass assembly can become hot during operation and remain so for a long time afterwards! Let the components cool off before continuing work with the device.
- **ATTENTION!** Avoid delayed boiling! Never heat the evaporating flask in the heating bath without switching on the rotary drive! Sudden foaming or exhaust gases indicate that flask content is beginning to decompose. Switch off heating immediately. Use the lifting mechanism to lift the evaporator piston out of the heating bath. Evacuate the danger zone and warn those in the surrounding area.
- **Automatic operation:** Always check evaporation processes before starting automatic operation. Automatic operation must not be used for unfamiliar evaporation processes.  
The following automatic modes are available in addition to manual mode:
  - 100% distillation,
  - Volume-/quantity-dependent distillation.
- When the device is switched off or the power supply disconnected, the internal safety lift removes the evaporator piston from the heating bath. The maximum load (evaporator piston and contents) is 1.5 kg for a permissible safety lift with no power.
- **CAUTION!** Never operate the device when the evaporator piston is rotating and the lift is raised. Always lower the evaporator piston into the heating bath first before starting the rotation drive. Otherwise hot tempering medium may be sprayed out!
- A vacuum may form inside the glassware in case of power outage. The device will vent automatically when it is switched on again.
- Set the speed of the drive so no tempering medium is sprayed out as a result of the evaporator piston rotating in the heating bath. If necessary reduce the speed.
- Do not touch rotating parts during operation.

- Imbalance may result in uncontrolled resonance behaviour of the device or assembly. Glass apparatus may be damaged or destroyed. In the event of unbalance or unusual noises, switch off the appliance immediately or reduce the speed.
- The appliance does not start up again automatically following a cut in the power supply.
- The device is only disconnected from the power supply network if the device power switch is off or the plug is pulled out.
- The socket for the mains cord must be easily accessible.

#### For protection of the equipment

- The voltage stated on the type plate must correspond to the mains voltage.
- Socket must be earthed (protective ground contact).
- Removable parts must be refitted to the appliance to prevent the infiltration of foreign objects, liquids etc.
- Protect the appliance and accessories from bumps and impacts.
- The appliance may only be opened by experts.

## Correct use

### • Use

► Together with the accessories recommended by **IKA**, the device is suitable for:

- quick and gentle distillation of liquids
- evaporation of solutions and suspensions
- crystallisation, synthesis or cleaning of fine chemicals
- drying of powder and granulate material
- recycling of solvents

Mode of operation: Tabletop device

### • Range of use

- Laboratories
- Pharmacies
- Schools
- Universities

The safety of the user cannot be guaranteed if the appliance is operated with accessories that are not supplied or recommended by the manufacturer or if the appliance is operated improperly contrary to the manufacturer's specifications.

## Unpack

### • Unpacking

- Please unpack the device carefully
- In the case of any damage a fact report must be set immediately (post, rail or forwarder)

### • Delivery scope

see table

	Drive RV 10 control	Heating bath HB 10 control	Vertical glassware RV 10.1	Vertical glassware RV 10.10 coated	Diagonal glassware RV 10.2	Diagonal glassware RV 10.20 beschichtet	Woulff bottle	Bracket	Cooler locking device, compl.	Vacuum hose (2x0,55 m)	Water discharge hose (1x1 m)	Operating instructions	Ring spanner wrench	Handhold	optional:	RV 10.4001 Vacuum valve for central vacuum system	RV 10.4002 Vacuum valve for laboratory vacuum system	RV 10.4003 Pump controller with vacuum valve for individual operation on vacuum pump
RV 10 control V	x	x	x				x	x	x	x	x	x	x	x				
RV 10 control VC	x	x		x			x	x	x	x	x	x	x	x				
RV 10 control D	x	x			x		x	x	x	x	x	x	x	x		or	or	or
RV 10 control DC	x	x				x	x	x	x	x	x	x	x	x				

## Useful information

Distillation is a thermal separating process for liquid compounds based on substance-specific, pressure-dependent boiling points through evaporation and subsequent condensation.

The boiling point temperature decreases with decreasing external pressure which means that work is usually done under reduced pressure. In this way the heating bath can be maintained at a constant temperature (e.g. 60 °C). Using the vacuum, the boiling point is set with a steam temperature of approx. 40 °C. The cooling water for the condensation cooler should not be warmer than 20 °C (60-40-20 rule).

A chemical resistant membrane pump with a vacuum controller should be used to create the vacuum. The pump is protected from solvent residue by the addition of a Woulff bottle and/or a vacuum separator.

Working with a jet pump to create a vacuum can only be recommended to a limited extent as the solvents may contaminate the environment when using these systems.

Speed, temperature, piston size and system pressure all affect the evaporator air capacity.

The optimum capacity of the flow-through condenser is approx. 60%.

This corresponds to condensation on approx. 2/3 of the cooling coil. With larger capacities there is the risk that the uncondensed solvent vapour will be extracted.

In order to prevent this, the condenser loading safety monitoring system can be activated. Before activating the safety monitoring, first carry out several distillations, noting the value given under menu "Settings", menu option "Drive", item "Maximum heat rejection rate". Now calculate the heat rejection rate using the following formula:

Heat rejection rate threshold = maximum heat rejection rate \* 1.3

Enter this value under menu "Settings", menu option "Distillation", item "Heat rejection rate threshold". The device will now stop the measurement and issue an error message if the threshold for the condenser is exceeded. If the configuration of the condenser equipment or the distillation process is modified, then it may be necessary to re-evaluate the maximum condenser loading.

The factory setting for this value is 900 W.

The device is equipped with a piston-operated safety mechanism. A vacuum may form in the glassware assembly if there is a power outage; this will be vented when the system is switched on again. If the power cuts out, the evaporator piston is automatically lifted out of the heating bath by an integrated gas spring.

**CAUTION!** The maximum load (evaporator piston and contents) is 1.5 kg for a permissible safety lift with no power for glass apparatus with diagonal or vertical condenser.

When using other types of condensers such as dry ice or intensive condensers as well as when using return distillation distributors with slip-on condensers, it may be necessary to reduce the mass of the media to compensate for the mass of this additional glassware. Thus, prior to distillation, check whether the lift goes up without power when laden with the glass and distillation material.

In addition to offering a full range of manual and semi-automatic evaporation operating modes, the RV 10 control rotary evaporator can also be used for fully automatic, quantity-controlled evaporation processes. Standard fittings for fully-automatic operation include a two-position vacuum controller, a cooling water differential temperature measurement sensor, and a cooling water flow rate meter. The device is designed for operation with a cooling water supply system (e.g. laboratory thermostat), but can also be run off a water supply line. Please refer to the Technical Data for information on cooling water pressure, temperature constancy, and flow rate.

We recommend using the RV10.5001 water regulator valve when operating directly off a water supply line.

This valve allows the cooling water flow rate to be adjusted and automatically cuts off the cooling water supply on completion of the distillation process.

The RV10.5002 water filter can be used to prevent particles of dirt passing from the water supply line into the valve mechanism.

If pressure reduction is required, we recommend fitting the RV10.5003 pressure regulator valve in the pipe immediately downstream of the pick-off point.

**Automatic volume-dependent distillation:** The apparatus must be brought up to operating temperature. This is achieved by conducting a test distillation.

Please observe the specified values during automatic distillation (see "Technical specifications").

Vacuum equipment: The RV 10 control rotary evaporator is fitted as standard with a two-position vacuum controller and a pressure measuring chamber with bleed valve.

The correct vacuum valve must be selected from the following list and installed according to the type of vacuum source used at the customer site.

Vacuum source	Code	Valve
Central vacuum system	Large-scale vacuum source for multiple consumers of different types across a site; generally a fixed-pipe installation.	RV 10.4001
Laboratory vacuum system	Laboratory-based RV vacuum pump with multiple similar consumers, e.g. 2 rotary evaporators/1 pump.	RV 10.4002
Single pump system	One rotary evaporator with one pump/desktop operation. Pump is deactivated when target pressure reached and reactivated automatically.	RV 10.4003

## Setting up

### Drive RV 10 control

**ATTENTION!** Loosen transportation lock (fig. 4a)!

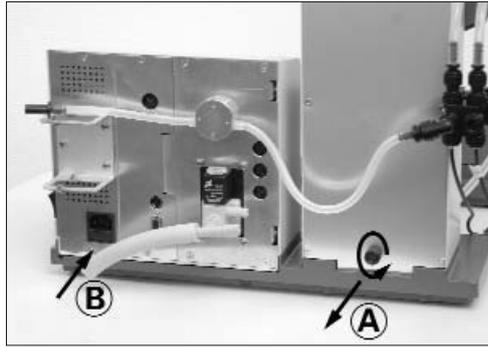
Fig. 4a

- ☞ Hold the lift with your hand on the height position and remove the thumb screw on the back of the appliance (A).
- ☞ Once the transportation lock has been removed, the lift moves slowly to its upper end position. The distance is approx. 140 mm.
- ☞ Connect the device to the power supply (B) using the supplied power cord (B).

### Adjustable base (fig. 4b)

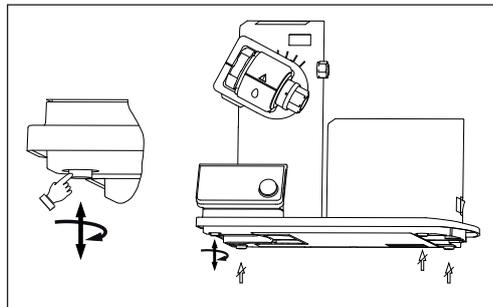
Fig. 4b

### Mounting the Woulff bottle (Fig. 5)

- ☞ Fit the hose connecting piece (C) on the left side of the lift.
- ☞ Mount the holder (D) onto the hose connecting piece (C).

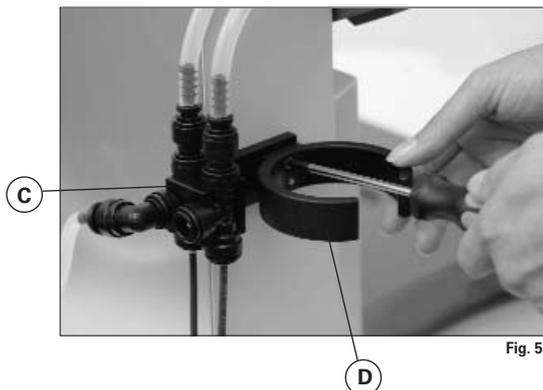
Fig. 5

- ☞ Insert the bottle and attach the supplied hose connectors to the bottle (fig. 6).

Fig. 6

- ☞ Remove the clamping device for the angle setting of the rotation drive on the right side of the lift by rotating the knurled screw counter clockwise (gently press and rotate the knurled screw (E) at the same time to move it further in or out).
- ☞ Set the drive at an angle of approx. 30° (fig. 7).

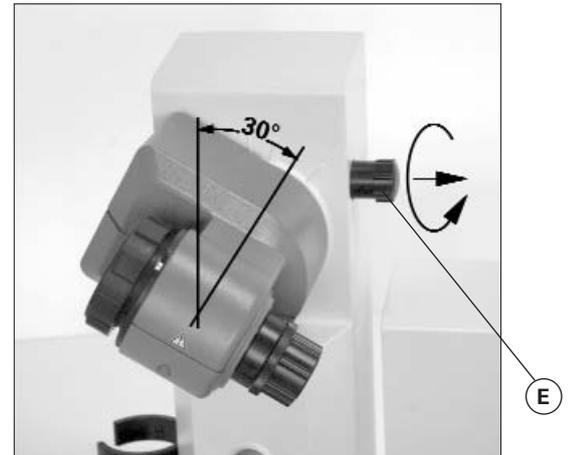
Fig. 7

- ☞ Then secure the rotation drive from being accidentally turned by tightening the knurled screw in a clockwise direction.

### Heating bath HB 10 control

**CAUTION!** Refer to the chapter entitled "Commissioning" in the heating bath instruction manual!

- ☞ Place the heating bath on the stand of the rotation drive and push it into the left position (fig. 8).

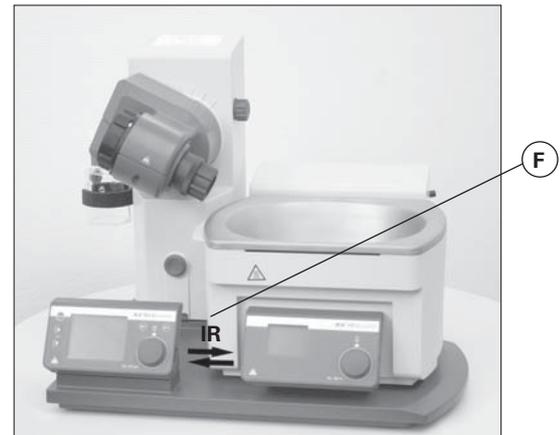
Fig. 8

*Note: Data is exchanged between the drive unit and the heating bath by means of an infrared link (F). Please note that reliable communication is only guaranteed when the infrared beam has clear line-of-sight to the detector.*

### Glassware

- ☞ Release the locking device on the steam pipe by turning it 60° counter clockwise (G), (fig. 9).

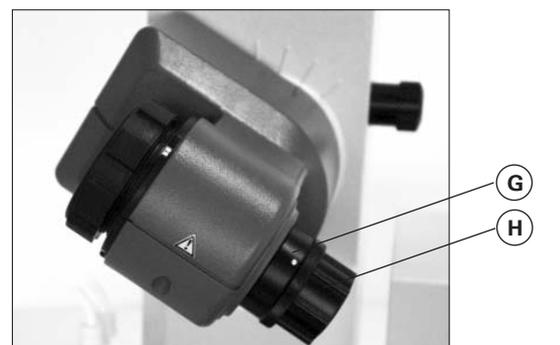
Fig. 9

- ☞ Feed the steam pipe in until it stops.
- ☞ Then lock the locking device by turning it clockwise by 60° (fig. 10).

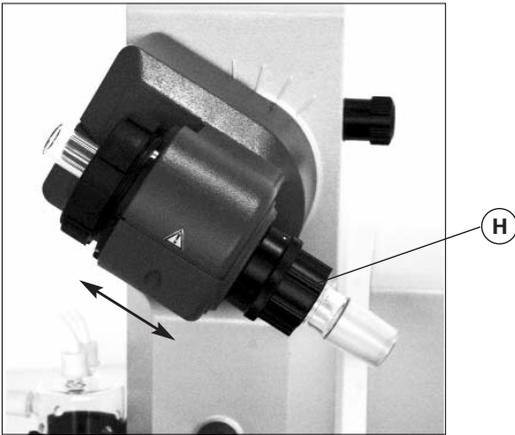
Fig. 10

- ☞ Check the axial locking device on the steam pipe.

Push-off mechanism

The plastic screw nut (H) helps loosen tight-fitting piston ground-in connections (Fig. 9).

- ☞ Hold the locking device and loosen the plastic screw nut (H).
- ☞ To do so, hold the tight-fitting evaporator piston and turn the plastic screw nut (H) until the evaporator piston neck.

*Before using the device, hand-tighten the plastic screw nut (H) against the left-hand stop. This will hold the steam pipe locking device (G) firmly in place.*

Fitting the condenser seal

- ☞ Place the RV 06.15 or RV 06.13 condenser seal in the condenser receptacle and fit the glassware to the device according to the assembly instructions (Fig. 11 a,b, c and d).

Fig. 11a

Fig. 11b

First use - Fitting the seal RV 06.15/RV 06.13

- ☞ Insert vapour tube (1).
- ☞ Insert seal RV 06.15 (2).
- ☞ Slide the union nut (3b) over the flange on the condenser (3a).
- ☞ Also slide the annular spring (3c) over the flange on the condenser (3a).
- ☞ Position the condenser (3a) on top of the seal (2).
- ☞ Tighten the union nut (3b) by hand onto the threaded flange. Run the device for 20 minutes at 120 rpm. Finally, retighten the union nut (3b) by hand.

*Note: Follow the mounting instructions for the glassware*

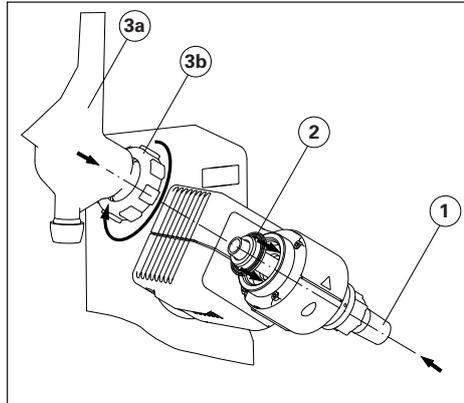
Fig. 11c

Fig. 11d

Assembling the vertical glassware cooler locking device

- ☞ Assemble the condenser locking device according to the diagram (fig. 12).
- ☞ Mount the plate (I) using the thumb screw (J).
- ☞ Put the support rod (L) on the plate (I) and attach it with the screw nut (K).
- ☞ Attach the rubber protector (M).
- ☞ Fasten the Velcro band (N) to the support rod (L).
- ☞ Secure the vertical glassware with the Velcro band (N).


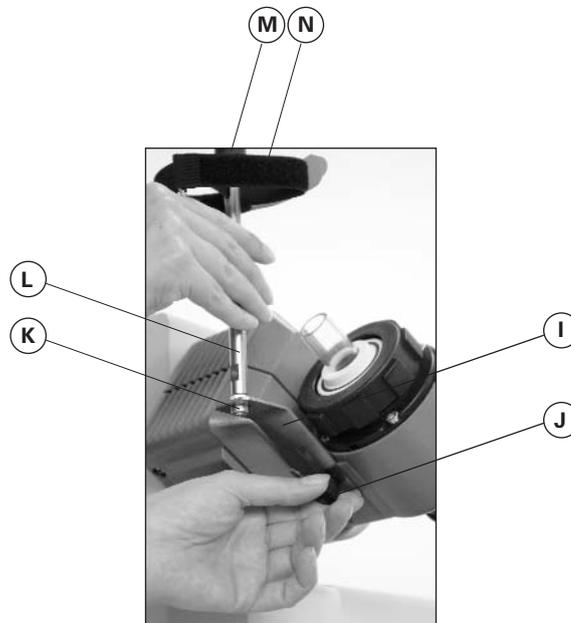
Fig. 12

## Mounting the glassware

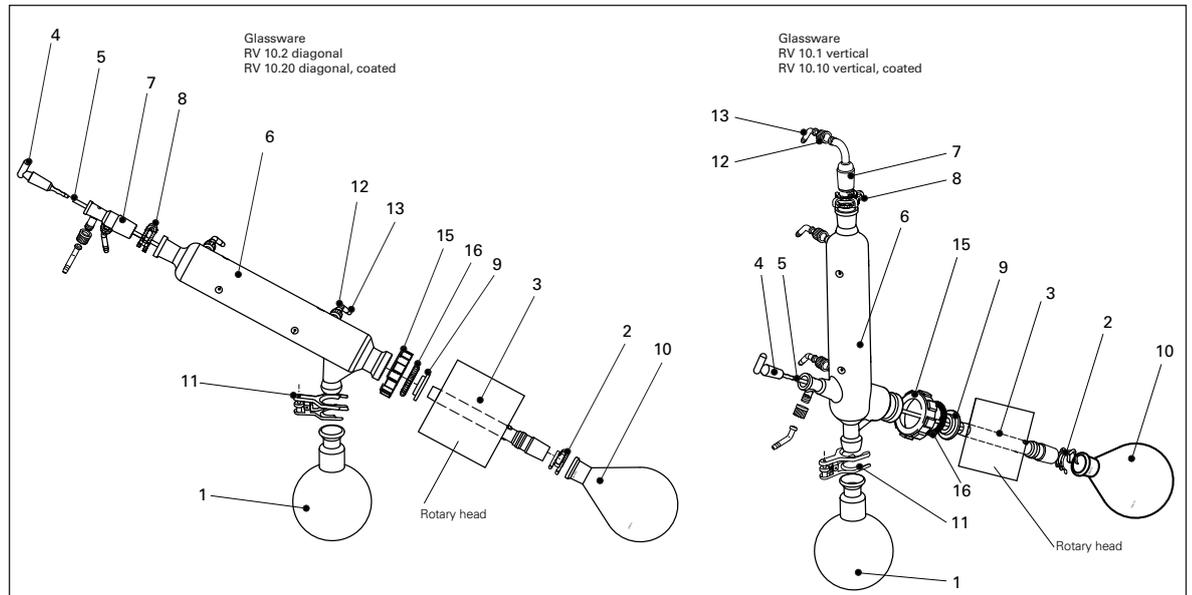Item	Designation	Quantity	
		diagonal glassware	vertical Glassware
1	Receiving flask	1	1
2	Clamp NS 29 (stainless steel)	1	1
3	Steam pipe	1	1
4	Stopcock	1	1
5	Pipe tube	1	1
6	Condenser	1 Diagonal condenser	1 Vertical condenser
7	Connection	1 Introduction sleeve	1 Vacuum connection
8	Clamp NS 29 (plastic)	1	1
9	Seal RV 06.13/ RV 06.15	1	1
10	Evaporator piston 1000 ml	1	1
11	Ball joint clamp RV 05.10	1	1
12	Screw joint cap	4	4
13	Hose connection	4	4
15	Locking nut	1	1
16	Spring ring	1	1

## Removing the condenser

- ☞ Use the ring spanner provided to loosen union nuts that are tightly fitted.
- ☞ Loosen the union nut by turning anticlockwise.

Fig.13

## Description of special condensers

- **RV 10.3 Vertical-intensive condenser with manifold**  
This vertical-intensive condenser features a double jacket design for particularly efficient condensation.  
Also available with coating (RV10.30)
- **RV 10.4 Dry ice condenser**  
Dry ice condenser for distilling low-boiling solvents.  
Cooling by dry ice, no cooling water required. Maximum condensation thanks to low temperatures.  
Also available with coating (RV10.40)  
If a dry ice condenser is used, then the flow rate control system must be switched of in menu "Distillation", menu option "Flow rate control"; see also "Commissioning" section.
- **RV 10.5 Vertical-condenser with manifold and cut-off valve for reflux distillation**  
Also available with coating (RV10.50)
- **RV 10.6 Vertical-intensive action condenser with manifold and cut-off valve for reflux distillation**  
This vertical-intensive action condenser features a double jacket design for particularly efficient reflux distillation.  
Also available with coating (RV10.60)

# Hose system

Schematic view of hose system (rear view)

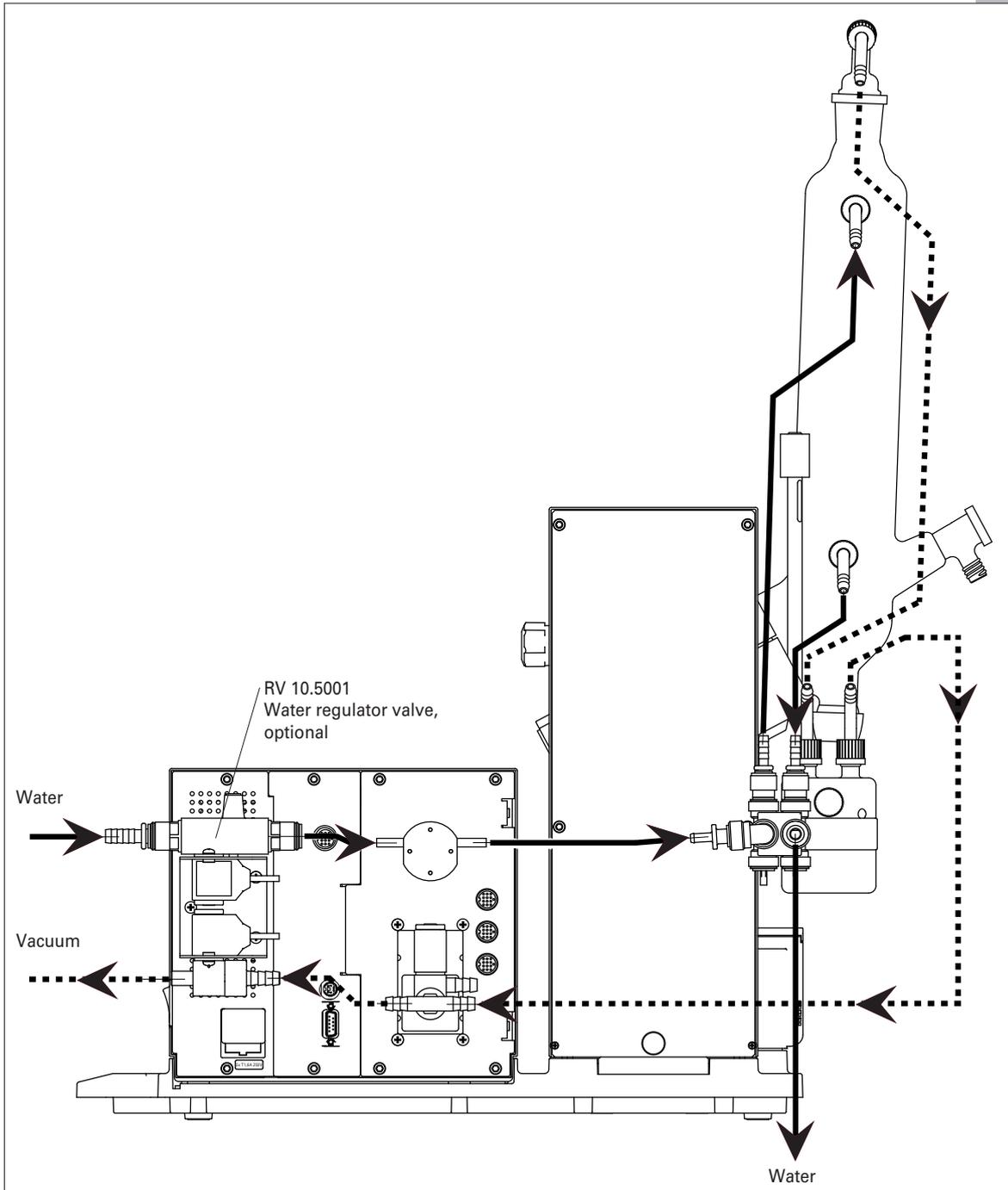
Fig. 14

Vacuum

Fit one of the optional RV10.4001, RV10.4002 or RV10.4003 vacuum valves, as appropriate for the type of vacuum system used, in the bracket provided for this purpose and connect the vacuum hose to the valve (Figure 15a, b, c).

Connection RV 10.4001


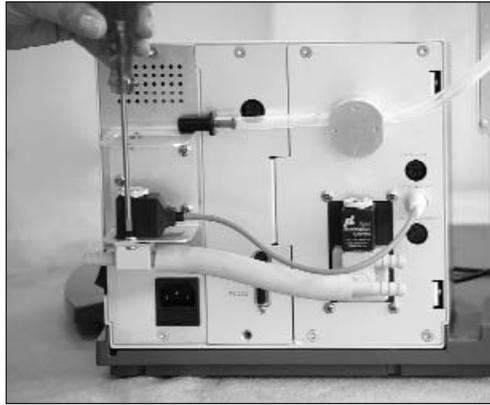
Fig. 15a

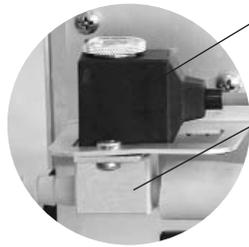Solenoid rotates freely

Direction of suction corresponds to the arrow marking. Always check for correct orientation.

Connection RV 10.4003


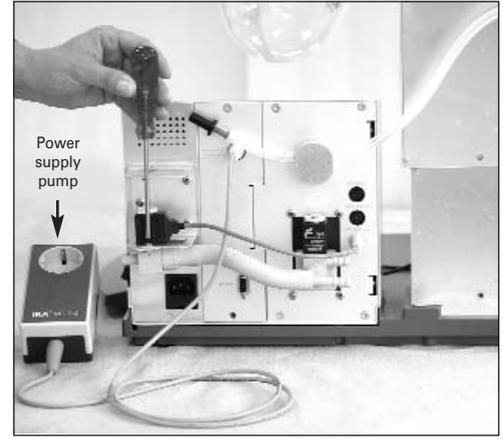
Fig. 15c

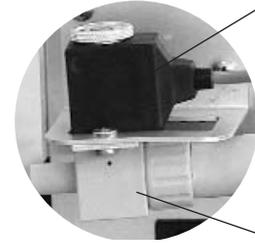Solenoid rotates freely

Direction of suction corresponds to the arrow marking. Always check for correct orientation.

Connection RV 10.4002


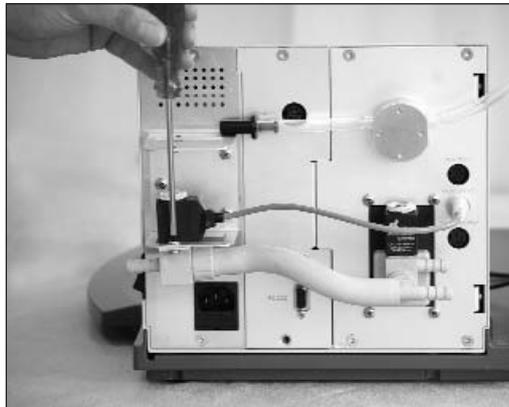
Fig. 15b

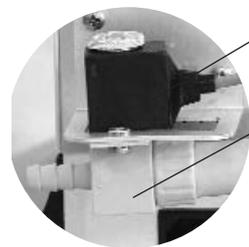Solenoid rotates freely

Direction of suction corresponds to the arrow marking. Always check for correct orientation.

Connection RV10 temperature sensor (dT)

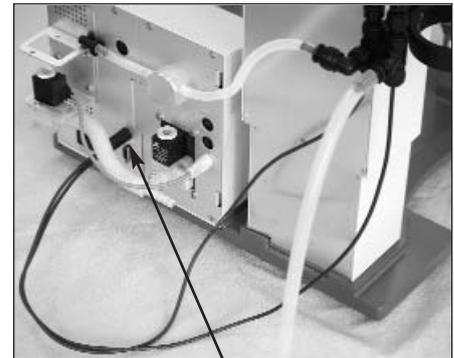
Fig. 15d

dT

Note the direction of the arrow.

## Water

- ☞ Connect the water inlet to the water supply (Figure 16a; Figure 16b with optional RV 10.5001 water regulator valve). Please refer to the technical specifications relating to the water supply. The RV 10.5001 water regulator valve is not suitable for use on the condenser unit because it would reduce the flow rate too greatly.

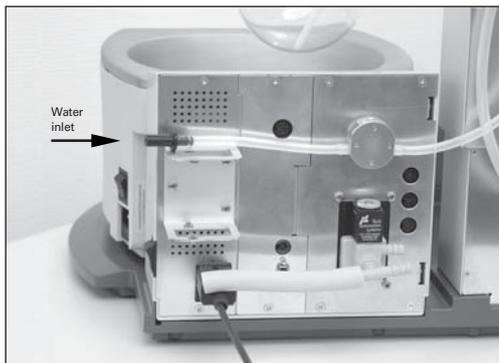
Fig. 16a

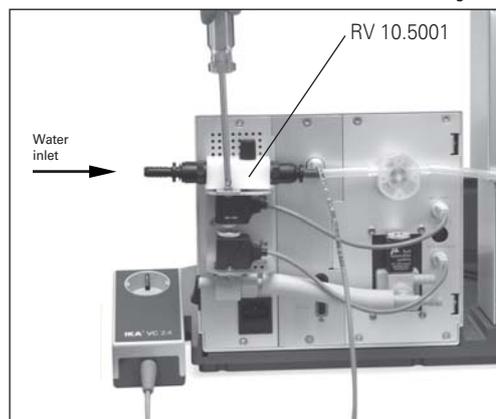
Fig. 16b

- ☞ Connect the water hoses to the glass condenser (short hose (O) = lower outlet, long hose (P) = upper inlet) and fix hoses in place (fig. 17).

Note: For special condensers see "Description of special condensers" section.

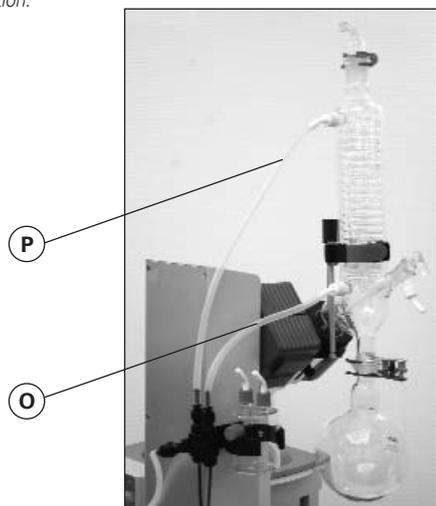
Fig. 17

- ☞ Connect the water outlet hose (supplied) by pushing it onto the nipple until the stop inside the connection socket is reached (fig. 18).

**Caution!** Ensure that the condenser inlet and condenser outlet are connected correctly.

Automatic distillation cannot be carried out if the water inlet and outlet hoses are not fitted according to the instructions and illustrations acc. fig. 14, since this would cause the temperature readings for the inlet and outlet to be measured incorrectly.

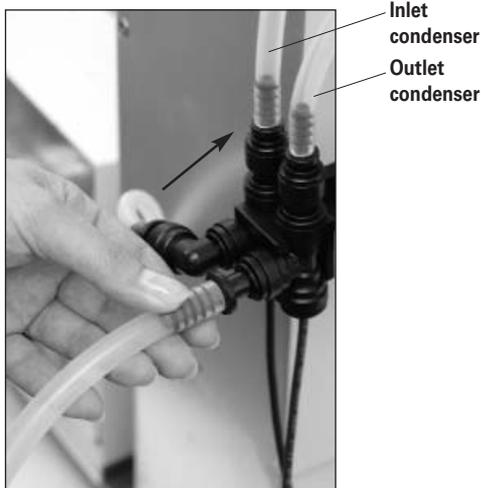
Fig. 18

- ☞ Release the hose connection using the supplied tool.

Fig. 19

Schematic view of connections (rear view)

Insert the valve connector cable (RV10.5001 or RV10.4001/2 resp. RV10.4003 Valve) resp. the network cable (RV10.4003 pump control) into the appropriate socket (Fig. 20).

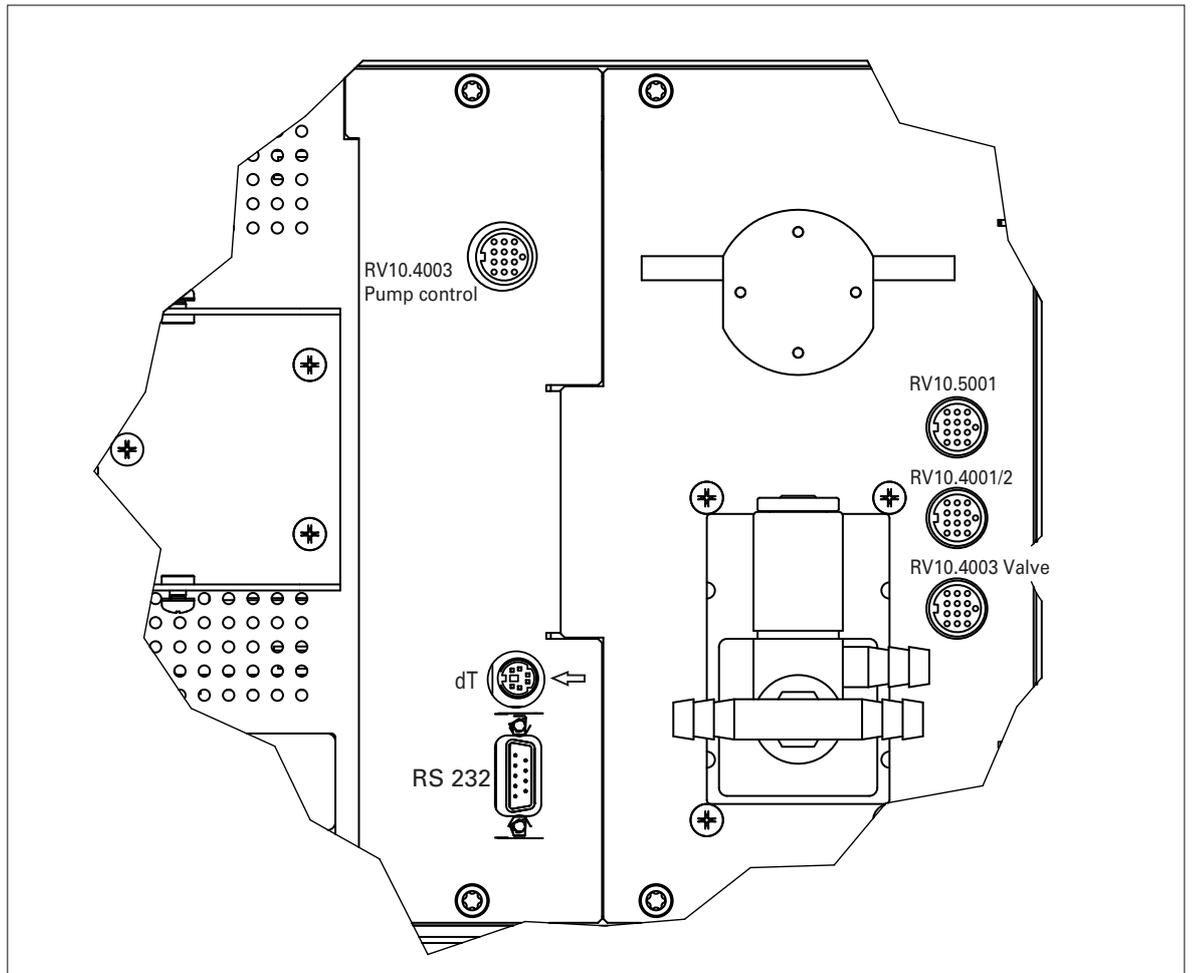
Fig. 20

Connect the pressure sensor to the Wouff bottle and the Wouff bottle to the vacuum connector on the condenser using the vacuum hoses supplied. Please note that the vacuum must always be connected to the highest connection point on the condenser (fig. 21).

Fig. 21

### HB 10 control

Data is exchanged between the heating bath and the drive unit by means of an infrared link. These are located on the left display side of the heating bath or on the right side of the drive unit. Do not place any objects between the two operating units as otherwise the data transfer may be interrupted (fig. 22)!

Fig. 22

### RV10 control

The device can be operated in "Remote" mode via an RS 232 or USB interface using labworldsoft® laboratory software.

The RS 232 interface at the back of the device, see fig. 20, is fitted with a 9-pole SUB-D jack and can be connected to a PC. The pins have serial signals. The USB interface is located on the left-hand side of the display on the drive unit and can be connected to a PC using the USB cable supplied.

*Note: Please note the system requirements as well as the operating instructions and help section included with the software.*

### USB Interface

The Universal Serial Bus (USB) is a serial bus system used to connect the RV 10 control with the PC (Fig. 23). Devices that support USB can be connected to each other whilst they are running (hot pluggable) and provide automatic recognition of the connected devices and their properties.

Use the USB interface in conjunction with labworldsoft® for operation in "Remote" mode to update the solvent library. Navigate to <http://www.ika.net/ika/lws/download/RV10.cfg> to update.

Fig. 23

### Installation

Connect the RV 10 control to the PC using the USB data cable. The RV 10 control will then transmit information to the Windows operating system to tell it which device drivers are required. Windows will then either:

- Load the driver,
- Install the driver automatically if not already installed,
- Prompt the user to perform a manual installation.

Navigate to <http://www.ika.net/ika/lws/download/stmcdc.inf>.

### Serial interface RS 232 (V24)

#### Configuration

- The functions of the interface lines between the device and the automation system are a selection from the signals specified in the EIA standard RS232 C, as per DIN 66 020 Part 1.
- For the electrical characteristics of the interface and the allocation of signal status, standard RS 232 C applies in accordance with DIN 66 259 part 1.
- Transmission procedure: asynchronous character transmission in start-stop mode.
- Type of transmission: full duplex.
- Character format: character representation in accordance with data format in DIN 66 022 for start-stop mode. 1 start bit; 7 character bits; 1 parity bit (even); 1 stop bit.
- Transmission speed: 9600 bit/s.
- Data flow control: none
- Access procedure: data transfer from the device to the computer takes place only at the computer's request.

#### Command syntax and format

The following applies to the command set:

- Commands are generally sent from the computer (Master) to the device (Slave).
- The stirrer machine sends only at the computer's request. Even fault indications cannot be sent spontaneously from the stirrer machine to the computer (automation system).
- Commands are transmitted in capital letters.
- Commands and parameters including successive parameters are separated by at least one space (Code: hex 0x20).
- Each individual command (incl. parameters and data) and each response are terminated with Blank CR Blank LF (Code: hex 0x20 hex 0x0d hex 0x20 hex 0x0A) and have a maximum length of 80 characters.
- The decimal separator in a number is a dot (Code: hex 0x2E).

The above details correspond as far as possible to the recommendations of the NAMUR working party (NAMUR recommendations for the design of electrical plug connections for analogue and digital signal transmission on individual items of laboratory control equipment, rev. 1.1).

The NAMUR commands and the additional specific **IKA** commands serve only as low level commands for communication between the rotary evaporator and the PC. With a suitable terminal or communications programme these commands can be transmitted directly to the rotary evaporator. The **IKA** software package, labworldsoft, provides a convenient tool for controlling rotary evaporator and collecting data under MS Windows, and includes graphical entry features, for motor speed ramps for example.

The following table summarises the (NAMUR) commands understood by the **IKA** control equipment.

Abbreviations used:

m	=	Numbering parameter (integer)
X	=	2 Temperature heating plate
X	=	3 Safety temperature heating plate
X	=	4 Speed
X	=	60 Interval time (1-99 seconds, 1 <= m >=99)
X	=	61 Timer (1-199 minutes, 1 <= m >=199)
X	=	62 Upper lift position (OUT_SP_62 1-> drive lift up)
X	=	63 Lower lift position (OUT_SP_62 1-> drive lift down)
X	=	66 Value vacuum controller
X	=	70 Vacuum controller hysteresis
X	=	74 Tempering medium (OUT_SP_74 0=Öl, OUT_SP_74 1=Wasser)

NAMUR Commands	Function
IN_NAME	Request designation
IN_PV_X X = 2, 3, 4, 66	Read actual value
IN_SOFTWARE	Request software Id-number, date and version
IN_SP_X X = 2, 3, 4, 60, 61, 66, 70	Read target value input
OUT_SP_X m X = 2, 3, 4, 60, 61, 62, 63, 66, 70, 74	Set target value to m
RESET	Switch to normal operation
START_X X = 2, 4, 60, 66	Switch on appliance (remote) function
STOP_X X = 2, 4, 60, 66	witch off appliance function

**PC 1.1 Cable**

This cable is required to connect the 9-pin connector to a PC. (fig. 24).

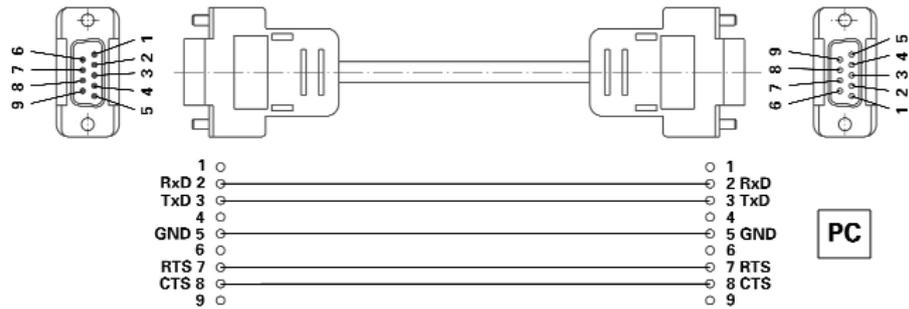
Fig. 24

# Commissioning

## Basic guide to using the menu system

### Selecting a menu

- ☞ Turn the rotary/push knob to the right/left to select the desired menu option.

*Note: The active menu item is shown with a black background on the display.*

### Editing values

- ☞ Press the rotary/push knob
- ☞ Turn the rotary/push knob to change the value (the magnitude of the change is dependent on the speed at which the knob is turned).
- ☞ Press the rotary/push knob again to complete the process.

### Exiting a menu

- ☞ Press the SET button.  
⇒ Current values are stored.  
Next-highest menu level is displayed.
- ☞ Press the ESC button.  
⇒ Current values are not stored.  
Next-highest menu level is displayed.

*Note: Some menu items do not contain values that can be stored. Select SET or ESC to exit the menu.*

*Note: Certain menu items have functionality that varies from or adds to that described above.*

### Display screen during the distillation process

Each operating mode has its own dedicated display screen. The following display characteristics are common:

- When a distillation is in progress it is not possible to leave the active display indicated.
- Press SET or ESC to exit the active display mode when no distillation is in progress. This will return the display to the main menu.

### Error messages

- ☞ Press ESC to acknowledge the fault if an error message appears on the display.
- ⇒ The error message will disappear.

*Note: In the case of serious faults, an additional display message will appear indicating that the device must be switched off and only switched back on again when the fault has been rectified.*

The unit is ready for service when the mains plug has been plugged in.

### Switching on the device

- ☞ Switch on the device on the right side of the device by the mains plug (fig.25).
- ☞ Device functions activated.

Fig.25

*Note: The heating bath must also be switched on.*

## Functional description (in as-delivered state)


### Factory settings

The values specified below correspond to the as-delivered state (display output/factory language setting is English).


- ⇒ The start-up screen will be displayed for few seconds. This shows the software and hardware versions for the device.

The main menu will appear after a few seconds

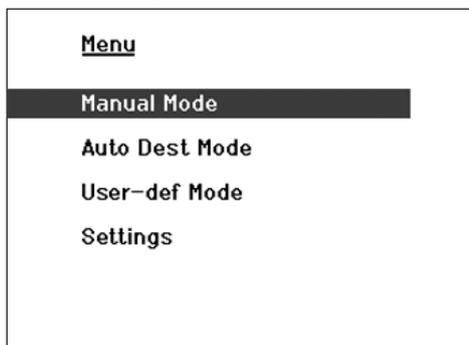- ☞ Turn the rotary/push knob to the right/left to select a menu option.
- ☞ Confirm the selection by pressing the rotary/push knob.

### "Settings" menu

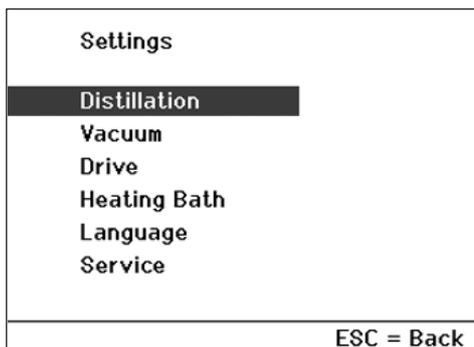- The following options are available from the "Settings" menu.
- ☞ Turn the rotary/push knob to the right/left to select a menu option.
  - ☞ Confirm the selection by pressing the rotary/push knob.

## Distillation settings

Settings Distillation	
Type of distillation	Volume
Drying	Off
Limit of cooler power	900 W
Maximum cooler power	0 W
RV10.5001 present ?	No
Flow Check	Off
Unit Amount	mL
SET = Save	ESC = Cancel

- ☞ Turn the rotary/push knob to the right/left to select a menu option.
- ☞ Press the rotary/push knob.
- ☞ Turn the rotary/push knob to change the value (the magnitude of the change is dependent on the speed at which the knob is turned).
- ☞ Press the rotary/push knob again to complete the process.

SET key: Exit menu and save changes.

ESC key: Exit menu without saving changes.

### Distillation method

Select between the automatic distillation methods "Volume" and "100%". This selection only has an effect during auto-distillation.

### Volume

The measured values of the cooling water flow rate and the cooling water temperature difference form the basis for a heat balance calculation used to determine the quantity of distillate at every stage of the distillation. The distillation process is halted when the specified quantity of distillate has been reached.

### 100%

The distillation process is stopped when the measured cooling water temperature difference falls below a threshold value, i.e. as soon as a solvent has been completely distilled off.

### Drying

The cooling water temperature difference is not monitored if this option is activated, e.g. use in processes for drying powdered media.

### Heat rejection rate threshold

A limit value can be set for the heat rejection rate that is appropriate for the condenser used; see "Worth knowing" section. The actual heat rejection rate is calculated during every distillation process. If the threshold is exceeded, the distillation will be stopped and an error message displayed.

### Maximum heat rejection rate

The maximum heat rejection rate achieved during the last distillation is displayed. This value is for information purposes only.

### RV10.5001 connected?

If no RV10.5001 is connected, then the cooling water flow rate will not be monitored.

### Flow rate control

It is also possible to deactivate the cooling water flow rate monitoring when the RV10.5001 is connected. This is necessary if, for example, a dry ice cooler is used.

### Quantity meas. unit

Select between millilitre and gram as the unit of quantity.

## Vacuum settings

Settings Vacuum	
Hysteresis	5 mbar
Scaling time	10 min
Scaling vacuum	10 mbar
Unit Vacuum	mbar
SET = Save	ESC = Cancel

- ☞ Turn the rotary/push knob to the right/left to select a menu option.
- ☞ Press the rotary/push knob.
- ☞ Turn the rotary/push knob to change the value (the magnitude of the change is dependent on the speed at which the knob is turned).
- ☞ Press the rotary/push knob again to complete the process.

SET key: Exit menu and save changes.

ESC key: Exit menu without saving changes.

### Hysteresis

The default settings for the hysteresis (as shown) are suitable for use in most situations. The (vacuum) hysteresis value describes the difference between the closing and opening pressures for the vacuum valve.

### Time scaling

The time axis of the vacuum-rotation-time graph is scaled using the value specified here. If the value = 0, then the time axis is scaled automatically.

### Vacuum scaling

The vacuum axis of the vacuum-rotation-time graph is scaled using the value specified here. If the value = 0, then the vacuum axis is scaled automatically.

### Vacuum meas. unit

Select between the units of measurement for a vacuum: mBar, Torr, and hPascal.

Settings Drive	
Right/Left interval	0 sec
Timer	0 min
Scaling time	10 min
Scaling rotation	10 rpm
SET = Save                      ESC = Cancel	

- ☞ Turn the rotary/push knob to the right/left to select a menu option.
- ☞ Press the rotary/push knob.
- ☞ Turn the rotary/push knob to change the value (the magnitude of the change is dependent on the speed at which the knob is turned).
- ☞ Press the rotary/push knob again to complete the process.

SET key: Exit menu and save changes.  
 ESC key: Exit menu without saving changes.

Right/left interval

Setting a particular value X for the right/left interval will cause the drive to change the direction of rotation every X seconds.

*Note: In interval mode the maximum speed is limited to 200 rpm.*

Timer

The timer value determines the time after which a manual distillation will be stopped

Time scaling

The time axis of the vacuum-rotation-time graph is scaled using the value specified here. If the value = 0, then the time axis is scaled automatically.

Rotation scaling

The rotation axis of the vacuum-rotation-time graph is scaled using the value specified here. If the value = 0, then the rotation axis is scaled automatically.

Settings Heating Bath	
Bath medium	Water
SET = Save                      ESC = Cancel	

- ☞ Turn the rotary/push knob to the right/left to select a menu option.
- ☞ Press the rotary/push knob.
- ☞ Turn the rotary/push knob to change the value (the magnitude of the change is dependent on the speed at which the knob is turned).
- ☞ Press the rotary/push knob again to complete the process.

SET key: Exit menu and save changes.  
 ESC key: Exit menu without saving changes.


Medium bath

Select either water or oil as the medium for the heating bath.

Language

Language	
English	
Deutsch	
Français	
Español	
中文	
SET = Save                      ESC = Cancel	

- ☞ Turn the rotary/push knob to the left/right to select the language.
- ☞ Confirm the selection by pressing SET.

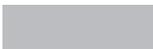Service			
<b>RV10.5001</b>	<b>Off</b>	<b>0.0L/h</b>	
<b>RV10.4001/2</b>	<b>Off</b>		
<b>RV10.4003 Valve</b>	<b>Off</b>		
<b>V4 Vacuum venting</b>	<b>Off</b>		
<b>RV10.4003 Pump</b>	<b>Off</b>		
<b>Adjustment</b>		<b>0.0 K</b>	
<b>Operation hours</b>		<b>1 h</b>	
			<b>ESC = Back</b>

- ☞ Turn the rotary/push knob to the right/left to select a menu option.
- ☞ Press the rotary/push knob.
- ☞ Turn the rotary/push knob to change the value (the magnitude of the change is dependent on the speed at which the knob is turned).
- ☞ Press the rotary/push knob again to complete the process.

ESC key: Exit menu, all changes to switch settings are reset, the calibration remains unaffected.

#### Switching valves and pump

The Service menu is used to check functionality during servicing by directly switching the valves and pump e.g. RV10.5001.

The RV10.5001 valve switches the water circuit on and off. This row also shows the current cooling water flow rate.

#### Calibrate

The "Calibrate" option is used for calibration of the temperature sensor.

The temperature difference displayed will be calibrated to zero. This action cannot be reversed.

The temperature sensor must be calibrated whenever it is replaced or changed, or if a new temperature sensor is installed.

Please contact the Service department.

The temperature sensor supplied with the system is calibrated ex-works.

#### Performing calibration

- ☞ Activate the "Calibrate" menu option by turning the rotary/push knob.
- ☞ Remove both sensors. To remove a sensor, push back the outer ring to release the locking device on the push-in connector, while at the same time pulling the temperature sensor out of the connector, see Fig. 26.

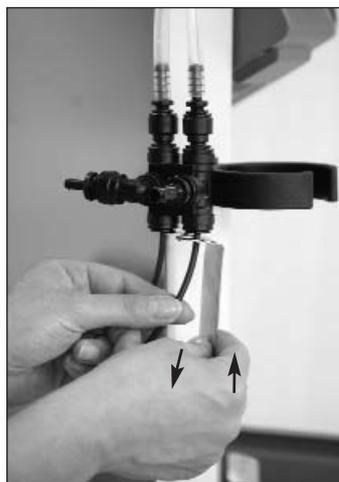
Fig.26

- ☞ To reconnect a temperature sensor, push it into the push-in connector until the limit stop is reached; an initial resistance will be felt.

- ☞ Fill water at room ambient temperature into a glass beaker (approx. 500 ml). Fully submerge both sensors into the water, see fig. 27.

Fig.27

- ☞ Wait until the temperature display in the "Service" menu, menu option "Calibrate" has stabilised, e.g. 0.2K.
- ☞ Press the rotary/push knob.
- ⇒ The temperature difference displayed will be calibrated to zero.


Calibration is only possible within the temperature range +0.5 K to -0.5 K. If the temperature is outside this range, then either a temperature sensor is faulty or connected incorrectly. Please contact the Service department.

If the characteristics of the cooling water supply deviate significantly from the specifications given in the "Technical specifications", with the result that the initial conditions for automatic distillation are not achieved, then it may be necessary to recalibrate the temperature sensors. This may, however, lead to a reduction in the accuracy with which the distillate quantity is measured in "Auto-distillation" mode.

#### Operation hours

This value serves as information and cannot be changed.

## „Manuel Mode“ menu
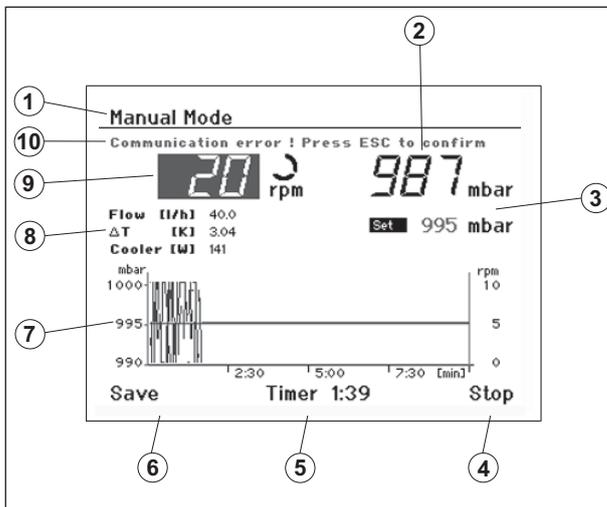
### „Manuel Mode“ display screen

1. Title
2. Vacuum display (actual value)
3. Vacuum display (target value)
4. Control field for „Start/ Stop“ distillation
5. Timer display
6. Control field for storing the distillation sequence
7. Vacuum-rotation-time graph
8. Display of flow rate, differential temperature, and condenser heat rejection rate
9. Control field for rotary speed display (target value) and speed measurement unit
10. Error message in case of fault, otherwise status indicator

### Changing the target value on the rotary speed display

- ☞ Press the rotary/push knob to change the target value.
- ☞ Turn the rotary/push knob to the right/left to change the target value. The magnitude of the change is dependent on the speed at which the knob is turned.
- ☞ Press the rotary/push knob again to save and complete the change

Note: When you select the speed > 100 rpm, smooth start is automatically activated.

### Start

- ☞ Turn the rotary/push knob to the right/left to select "Start".
- ☞ Press the rotary/push knob to start the distillation..
- ⇒ The display field will change to "Stop".


Note: If the timer has been configured, then the time will start to count backwards. If a timer has not been configured then the time elapsed since the start is shown.

### Stop

- ☞ Turn the rotary/push knob to the right/left to select "Stop".
- ☞ Press the rotary/push knob to stop the distillation.
- ⇒ The display field will change to "Start".

### Saving the distillation sequence as a procedure on completion of distillation

- ☞ Turn the rotary/push knob to the right/left to select "Save".
- ☞ Press the rotary/push knob to display the "Procedures" screen.

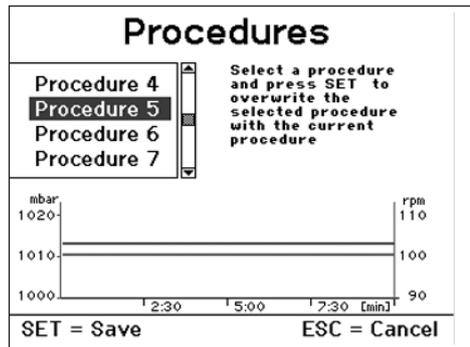
Turn the rotary/push knob to the left/right to select one of the ten procedures. The sequence will be shown graphically.

- ☞ Press SET, to overwrite the selected procedure with the current procedure (the procedure that has just been completed) and exit the menu.
- Pressing ESC will prevent the selected procedure from being overwritten.

Note: Distillation sequences saved in this manner can be repeated at a later date in the user-defined mode by selecting the corresponding procedure.

Note: If it is necessary to work in manual mode for an extended period (>4 minutes) without performing a distillation (e.g. flask not in heating bath), or if a drying process is carried out, then "Drying" mode must be activated in menu "Settings", menu option "Distillation", in order to prevent the "No temperature increase" error message being displayed.

## Auto-distillation Mode " menu

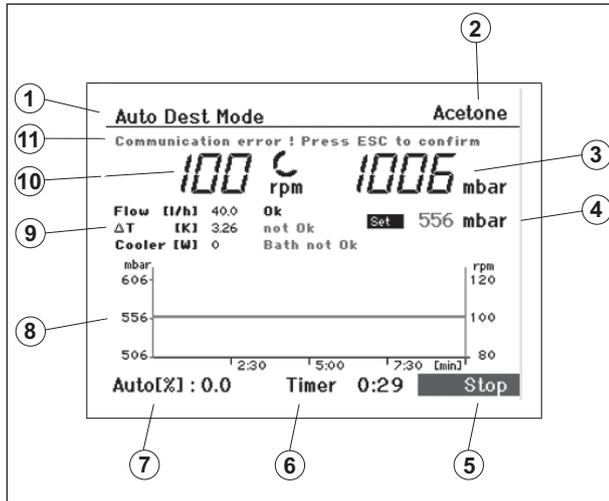
"Auto-distillation Mode" display screen

1. Title
2. Control field for solvent
3. Vacuum display (actual value)
4. Vacuum display (target value)
5. Control field for „Start/ Stop“ distillation
6. Timer display
7. Distillation progress indicator  
(% of required quantity of distillate)
8. Vacuum-rotation-time graph
9. Display of flow rate, differential temperature, and condenser heat rejection rate.
10. Rotary speed display (target value) and speed measurement unit
11. Error message in case of fault, otherwise status indicator

### Changing the selected solvent


- ☞ Turn the rotary/push knob to the right/left to select the control field "Solvent".
- ☞ Press the rotary/push knob to display a list of the various solvents available.

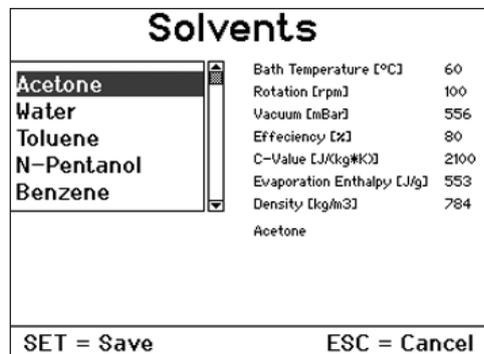
- ☞ Turn the rotary/push knob to the right/left to select a "Solvent".
- ☞ Press the SET key.
- ⇒ The solvent selected will be used for the subsequent distillation.

*Note: The display shows the parameters for each of the solvents. The parameters can only be changed for the user-definable solvents UserSolvent1 ... UserSolvent5. Please refer to the next section for instructions on how to do this.*

## User defined solvents

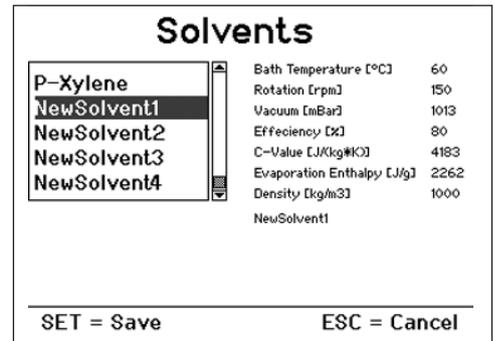- ☞ Select one of the options NewSolvent1 ... NewSolvent5 from the "Solvents" screen.
- ☞ Select one of these solvents and press the rotary/push knob.

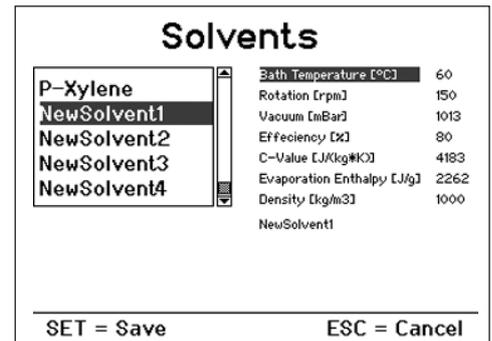All the parameters displayed for the solvent can now be modified.

- ☞ Turn the rotary/push knob to the right/left to select the parameter to be changed.
- ☞ Press the rotary/push knob.
- ☞ Turn the rotary/push knob to the right/left to change the selected parameter.
- ☞ Press the rotary/push knob.
- ☞ Press SET to save all values.
- ⇒ The name of the solvent will be changed to "UserSolvent1". This action also selects UserSolvent1 as the current solvent.

*Note: It is not possible to return to the selection list.*

## Starting automatic distillation

Press the rotary/push knob with the "Start" or "Continue" field selected. If the "100% distillation" or "Drying" options are active, then the distillation will start immediately, provided that the necessary conditions are fulfilled.

When the "Volume" distillation method is used, a further screen with additional distillation parameters will be displayed.

Acetone	
<b>Distillation :</b>	<b>Volume</b>
<b>Amount</b>	<b>500 mL</b>
<b>Distillate</b>	<b>300 mL</b>
<b>Efficiency</b>	<b>80 %</b>
<b>Start</b>	

- Turn the rotary/push knob to the right/left to select a menu option.
- Press the rotary/push knob.
- Turn the rotary/push knob to change the value (the magnitude of the change is dependent on the speed at which the knob is turned).
- Press the rotary/push knob again to complete the process.

### Distillation

it is not possible to change the distillation method anymore in this menu option.

### Amount

Enter the amount, which the flask contains.

### Distillate

Enter the amount, which you want to distill off.

### Efficiency

Balancing requires a precise thermal efficiency to be entered. The efficiency is dependent on a combination of factors that affect the distillation process, so it must initially be estimated. The first time a distillation process is carried out should therefore be used for calibration. Following completion of the distillation process, determine the quantity of distillate produced and calculate the actual heat rejection rate according to the following formula:

$$\eta^p = \frac{\eta_{th} \cdot m_{(calc)}}{m_{(gem)}}$$

$\eta_p$	Actual efficiency
$\eta_{th}$	Estimated efficiency
$m_{(gem)}$	Measured quantity of distillate
$m_{(calc)}$	Specified (calculated) distillate quantity

For subsequent distillations, the actual calculated efficiency should be entered; this will ensure that subsequent distillations under the **same distillation and ambient conditions** are completed automatically with sufficient precision.

### End of automatic distillation

Automatic distillation will stop when:

### Volume controlled

The specified quantity of distillate has been reached.


### 100% distillation

The measured cooling water temperature has fallen below a target value, indicating that the solvent is fully evaporated.

The distillation can also be stopped manually:


- Press the rotary/push knob with the "Stop" display field selected..

## "User-defined Distillation" menu
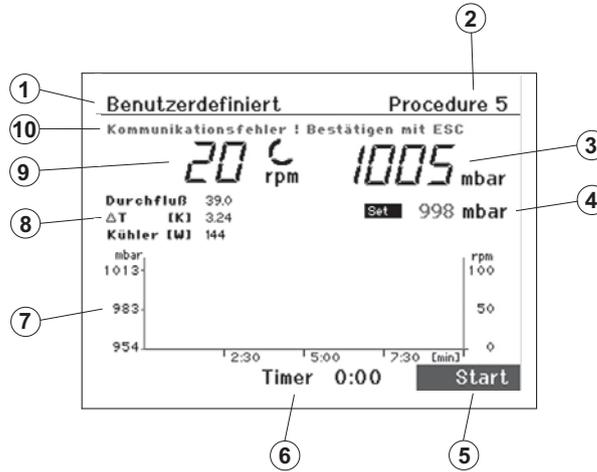
### "User-defined Distillation" display screen

1. Title
2. Control field for distillation procedure
3. Vacuum display (actual value)
4. Vacuum display (target value)
5. Control field for „Start/ Stop“ distillation
6. Timer display
7. Vacuum-rotation-time graph
8. Display of flow rate, differential temperature, and condenser heat rejection rate
9. Rotary speed display (target value) and speed measurement unit
10. Error message in case of fault, otherwise status indicator

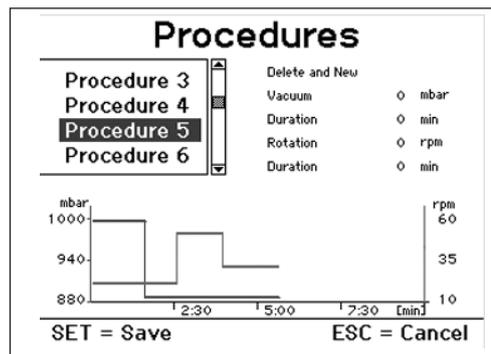- Turn the rotary/push knob to the right/left to select „Procedure“.
- Press the rotary/push knob.
- The display will show the "Procedures" screen.


- Turn the rotary/push knob to the right/left to select a procedure.
- Press the SET button.
- The procedure selected will be used for the subsequent distillation.

*Note: The vacuum-rotation-time graph is displayed for each procedure. The parameters can be changed for each procedure. Please refer to the next section for instructions on how to do this. The sequence used for a manual distillation process can also be saved as a procedure. Please refer to the "Manual Distillation" section.*

### Changing a procedure

- ☞ Select the procedure to be changed and press the rotary/push knob. The "Delete and create" display field will be activated.
- ☞ Press the rotary/push knob with the "Delete and create" field selected to delete the procedure.

New vacuum and rotation values can now be entered to define a new step in procedure.
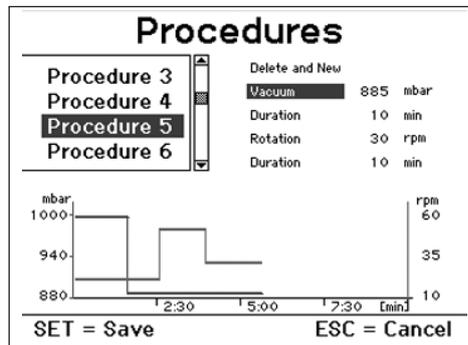
### Vacuum

- ☞ Turn the rotary/push knob to the right/left to select "Vacuum" from the on-screen options.
- ☞ Press the rotary/push knob.
- ☞ Turn the rotary/push knob to the right/left to enter a value.
- ☞ Press the rotary/push knob.

### Duration

- ☞ Turn the rotary/push knob to the right/left to select "Duration" from the on-screen options.
- ☞ Press the rotary/push knob.
- ☞ Turn the rotary/push knob to the right/left to enter a value.
- ☞ Press the rotary/push knob. This will define a new step in the procedure with the pair of values vacuum-duration and add it to the procedure. The diagram will be updated.

### Rotation

- ☞ Turn the rotary/push knob to the right/left to select "Rotation" from the on-screen options.
- ☞ Press the rotary/push knob.
- ☞ Turn the rotary/push knob to the right/left to enter a value.
- ☞ Press the rotary/push knob.

### Duration

- ☞ Turn the rotary/push knob to the right/left to select "Duration" from the on-screen options.
- ☞ Press the rotary/push knob.
- ☞ Turn the rotary/push knob to the right/left to enter a value.
- ☞ Press the rotary/push knob. This will define a new step in the procedure with the pair of values rotation-duration and add it to the procedure. The diagram will be updated.

The vacuum-duration and rotation-duration steps can be repeated a maximum of ten times.

SET key: Exit menu, store the newly created procedure and activate it as the current procedure.

ESC key: Exit menu and reset all changes.

### Stopping user-defined distillation


The user-defined distillation will stop automatically once all steps in the currently active procedure have been completed. It is also possible to stop the process manually by pressing the rotary/push knob with the "Stop" display field selected.

### Stand-by operation

- ☞ Press the "Power" button.
- ⇒ The device will enter stand-by.

- ☞ Press the "Power" button again.
- ⇒ The screen will display "Start" to indicate that the device is ready for use again.
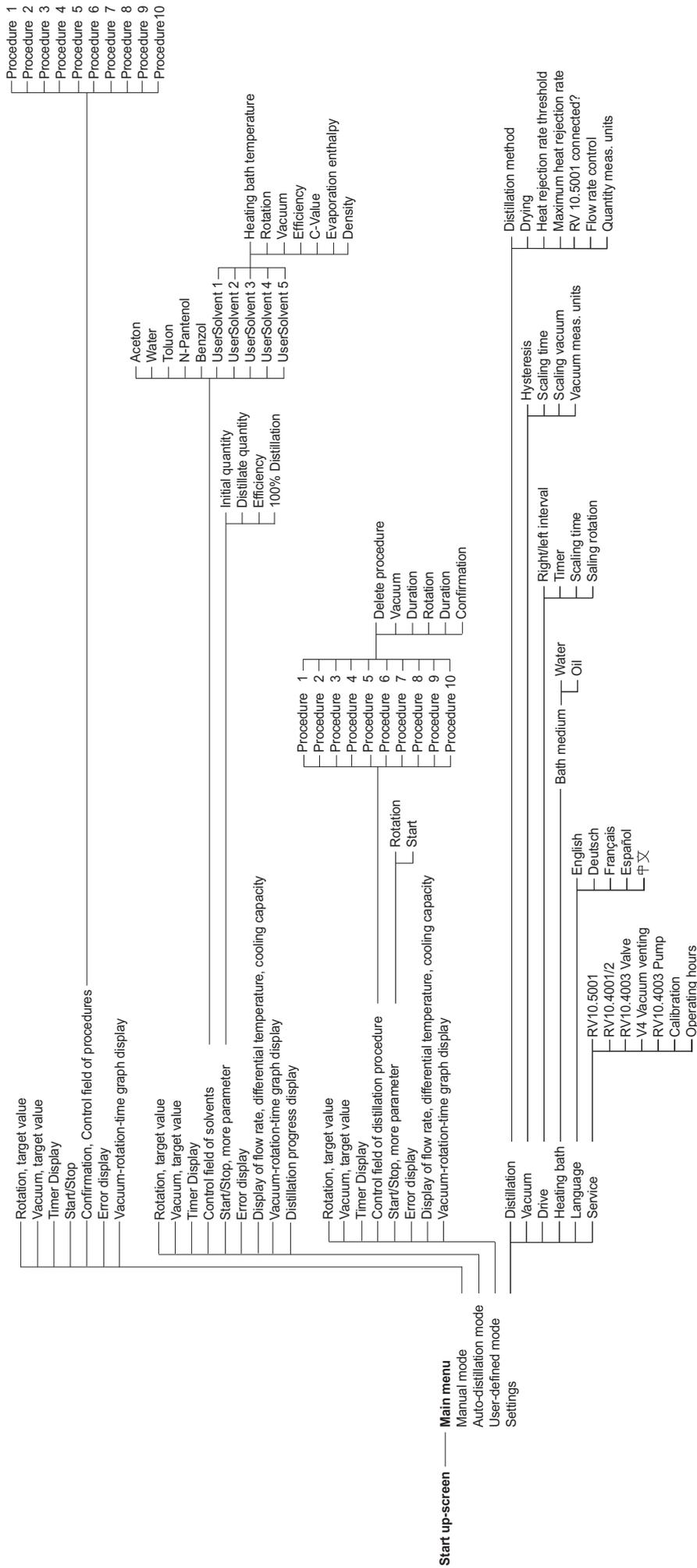
### Setting the lower end stop

**CAUTION!** Depending on the size of the piston, the setting angle of the rotation drive and the position of the heating bath and lift, the evaporator piston can be in contact with the bottom of the heating bath. CAUTION - glass may break!

Limit the lower lift position with the variable end stop.

The lift can be moved to any desired position in manual mode using the "▲" and "▼" keys. There is no automatic cut-off if a collision occurs.

- ☞ Press the "▼" key until the lift reaches the desired position.

*Note: The evaporator piston should be 2/3 immersed in the heating bath.*

- ☞ To adjust the stop element (Q), press the centre button (R) on the front of the lift (fig. 28).

- ☞ Slide the limit stop (Q) to the lowest position (Fig. 29).

- ☞ Press the "▲" key until the lift has reached the upper end stop.

*Note: The path is limited from 0 – 6 cm.*

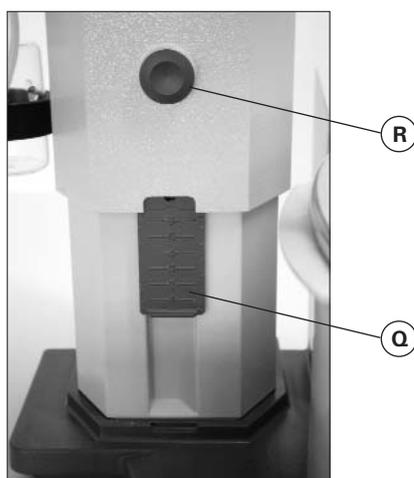
Fig.28

Fig.29

Check that the limit stop is functioning properly:

- ☞ Lower the lift by holding down the "▼" key.
- ☞ The drive stops automatically when the preset desired lower end position has been reached.
- ☞ Move the drive back to the upper position.

### Filling the evaporator piston

Prior to filling the evaporator piston, a vacuum controller is used to regulate the glass apparatus to the target pressure.

- ☞ Now fill the evaporator piston using the backfeed line, see fig. 30.
- ☞ Due to the vacuum present, the solvent is suctioned into the evaporator piston. This enables you to keep solvent loss due to suctioning to a minimum.
- ☞ You can also fill the evaporator piston manually prior to creating the vacuum. The evaporator piston should not be filled more than half its volume.

*Note: CAUTION! The maximum allowed load (evaporator piston and contents) is 3 kg.*

Fig. 30

### Setting up the heating bath

**Refer also to the operating instructions for the HB 10 control heating bath!**

- ☞ Move the lift to the bottom position and check the position of the heating bath in relation to the evaporator piston. When using larger evaporator pistons (2 or 3 litres) or depending on the angle of the rotation drive, you can move the heating bath 50 mm to the right.
- ☞ Fill the heating bath with the tempering medium until the evaporator piston is surrounded by tempering medium to 2/3 of its volume.
- ☞ Switch on the rotation drive and slowly increase the speed.

*Note: Avoid creating waves.*

- ☞ Switch on the heating bath using the main on/off switch.

*Note: Avoid stress on the glass due to different evaporator piston and heating bath temperatures when lowering the evaporator piston into the heating bath!*


*Note: If non-original accessories are used that are not supplied by IKA, then it is possible that the 50 mm travel range provided by the heating bath will not be sufficient. This applies in particular when using 3 L evaporator flasks with foam brake.*

*Use the IKA RV 10.3000 mounting plate to extend the heating bath travel range by 150 mm.*

## Maintenance and cleaning

The device is maintenance-free. It is subject only to the natural wear and tear of components and their statistical failure rate.

The seal on the glass condenser should be checked at regular intervals and replaced if necessary.

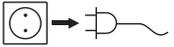For cleaning disconnect the main plug.

To clean the device use only water with a detergent that contains tensides, or use isopropylalcohol for stubborn soiling.

### Spare parts order

When ordering spare parts, please give:

- machine type,
- manufacturing number, see type plate,
- item and designation of the spare part, see [www.ika.net](http://www.ika.net), spare parts diagram and spare parts list.

### Repair

In case of repair the device has to be cleaned and free from any materials which may constitute a health hazard.

If you require servicing, return the appliance in its original packaging. Storage packaging is not sufficient. Please also use suitable transport packaging.

## Accessories

①

②

RV 10.1	NS 29/32 Vertical glassware (1)
RV 10.10	NS 29/32 Vertical glassware, coated (1)
RV 10.2	NS 29/32 Diagonal glassware (2)
RV 10.20	NS 29/32 Diagonal glassware, coated (2)
RV 10.700	NS 29/42 Vertical glassware (1)
RV 10.710	NS 29/42 Vertical glassware, coated (1)
RV 10.800	NS 29/42 Diagonal glassware (2)
RV 10.810	NS 29/42 Diagonal glassware, coated (2)
RV 10.900	NS 24/40 Vertical glassware (1)
RV 10.910	NS 24/40 Vertical glassware, coated (1)
RV 10.1000	NS 24/40 Diagonal glassware (2)
RV 10.1010	NS 24/40 Diagonal glassware, coated (2)

①

②

③

RV 10.3	Vertical-intensive condenser with manifold
RV 10.30	Vertical-intensive condenser with manifold, coated
RV 10.4	Dry ice condenser
RV 10.40	Dry ice condenser, coated
RV 10.5	Vertical-condenser with manifold and cut-off valve for reflux distillation
RV 10.50	Vertical-condenser with manifold and cut-off valve for reflux distillation, coated
RV 10.6	Vertical-intensive condenser with manifold and cut-off valve for reflux distillation
RV 10.60	Vertical-intensive condenser with manifold and cut-off valve for reflux distillation, coated

①

②

HB 10.1	Protection shield (Heating bath HB 10)
HB 10.2	Protection cover (Heating bath HB 10)


RV 10.70	NS 29/32 Steam pipe
RV 10.71	NS 24/29 Steam pipe
RV 10.72	NS 29/42 Steam pipe
RV 10.73	NS 24/40 Steam pipe

RV 10.80	NS 29/32 Evaporator piston 50 ml
RV 10.81	NS 29/32 Evaporator piston 100 ml
RV 10.82	NS 29/32 Evaporator piston 250 ml
RV 10.83	NS 29/32 Evaporator piston 500 ml
RV 10.84	NS 29/32 Evaporator piston 1000 ml
RV 10.85	NS 29/32 Evaporator piston 2000 ml
RV 10.86	NS 29/32 Evaporator piston 3000 ml
RV 10.90	NS 24/32 Evaporator piston 50 ml
RV 10.91	NS 24/32 Evaporator piston 100 ml
RV 10.92	NS 24/32 Evaporator piston 250 ml
RV 10.93	NS 24/32 Evaporator piston 500 ml
RV 10.94	NS 24/32 Evaporator piston 1000 ml
RV 10.95	NS 24/32 Evaporator piston 2000 ml
RV 10.96	NS 24/32 Evaporator piston 3000 ml
RV 10.97	NS 24/40 Evaporator piston 1000 ml
RV 10.2001	NS 29/42 Evaporator piston 50 ml
RV 10.2002	NS 29/42 Evaporator piston 100 ml
RV 10.2003	NS 29/42 Evaporator piston 250 ml
RV 10.2004	NS 29/42 Evaporator piston 500 ml
RV 10.87	NS 29/42 Evaporator piston 1000 ml
RV 10.2005	NS 29/42 Evaporator piston 2000 ml
RV 10.2006	NS 29/42 Evaporator piston 3000 ml
RV 10.2007	NS 24/40 Evaporator piston 50 ml
RV 10.2008	NS 24/40 Evaporator piston 100 ml
RV 10.2009	NS 24/40 Evaporator piston 250 ml
RV 10.2010	NS 24/40 Evaporator piston 500 ml
RV 10.2011	NS 24/40 Evaporator piston 1000 ml
RV 10.2012	NS 24/40 Evaporator piston 2000 ml
RV 10.2013	NS 24/40 Evaporator piston 3000 ml

RV 10.100	KS 35/20 Receiving flask 100 ml
RV 10.101	KS 35/20 Receiving flask 250 ml
RV 10.102	KS 35/20 Receiving flask 500 ml
RV 10.103	KS 35/20 Receiving flask 1000 ml
RV 10.104	KS 35/20 Receiving flask 2000 ml
RV 10.105	KS 35/20 Receiving flask 3000 ml

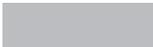RV 10.200 KS 35/20 Receiving flask, coated 100 ml  
 RV 10.201 KS 35/20 Receiving flask, coated 250 ml  
 RV 10.202 KS 35/20 Receiving flask, coated 500 ml  
 RV 10.203 KS 35/20 Receiving flask, coated 1000 ml  
 RV 10.204 KS 35/20 Receiving flask, coated 2000 ml  
 RV 10.205 KS 35/20 Receiving flask, coated 3000 ml

RV10.300 NS 29/32 Powder piston 500 ml  
 RV10.301 NS 29/32 Powder piston 1000 ml  
 RV10.302 NS 29/32 Powder piston 2000 ml  
 RV10.303 NS 24/29 Powder piston 500 ml  
 RV10.304 NS 24/29 Powder piston 1000 ml  
 RV10.305 NS 24/29 Powder piston 2000 ml

RV10.214 NS 29/32 Powder piston 500 ml  
 RV10.215 NS 29/32 Powder piston 1000 ml  
 RV10.216 NS 29/32 Powder piston 2000 ml  
 RV10.217 NS 24/40 Powder piston 500 ml  
 RV10.218 NS 24/40 Powder piston 1000 ml  
 RV10.219 NS 24/40 Powder piston 2000 ml

RV10.400 NS 29/32 Evaporator cylinder 500 ml  
 RV10.401 NS 29/32 Evaporator cylinder 1500 ml  
 RV10.402 NS 24/29 Evaporator cylinder 500 ml  
 RV10.403 NS 24/29 Evaporator cylinder 1500 ml  
 RV10.2020 NS 29/42 Evaporator cylinder 500 ml  
 RV10.2021 NS 29/42 Evaporator cylinder 1500 ml  
 RV10.2022 NS 24/40 Evaporator cylinder 500 ml  
 RV10.2023 NS 24/40 Evaporator cylinder 1500 ml

RV10.500 NS 29/32 Foam brake  
 RV10.501 NS 24/29 Foam brake  
 RV10.2024 NS 29/42 Foam brake  
 RV10.2025 NS 24/40 Foam brake

RV10.600 NS 29/32 Distillation spider with 6 sleeves  
 RV10.601 NS 29/32 Distillation spider with 12 sleeves  
 RV10.602 NS 29/32 Distillation spider with 20 sleeves  
 RV10.603 NS 24/29 Distillation spider with 6 sleeves  
 RV10.604 NS 24/29 Distillation spider with 12 sleeves  
 RV10.605 NS 24/29 Distillation spider with 20 sleeves  
 RV10.2026 NS 29/42 Distillation spider with 6 sleeves  
 RV10.2027 NS 29/42 Distillation spider with 12 sleeves  
 RV10.2028 NS 29/42 Distillation spider with 20 sleeves  
 RV10.2029 NS 24/40 Distillation spider with 6 sleeves  
 RV10.2030 NS 24/40 Distillation spider with 12 sleeves  
 RV10.2031 NS 24/40 Distillation spider with 20 sleeves

RV10.606 NS 29/32 Distillation spider with 5 flasks 50 ml  
 RV10.607 NS 29/32 Distillation spider with 5 flasks 100 ml  
 RV10.608 NS 24/29 Distillation spider with 5 flasks 50 ml  
 RV10.609 NS 24/29 Distillation spider with 5 flasks 100 ml  
 RV10.2032 NS 29/42 Distillation spider with 5 flasks 50 ml  
 RV10.2033 NS 29/42 Distillation spider with 5 flasks 100 ml  
 RV10.2034 NS 24/40 Distillation spider with 5 flasks 50 ml  
 RV10.2035 NS 24/40 Distillation spider with 5 flasks 100 ml

RV10.3000 Extension plate

RV10.4001 Solenoid valve (central vacuum system) (1)  
 RV10.4002 Solenoid valve (laboratory vacuum system) (2)  
 RV10.4003 Pump controller, including solenoid valve (3)  
 RV10.5001 Water regulator valve (4), only for use on water supply line!  
 RV10.5002 Filter (5)  
 RV10.5003 Pressure control valve (6)

①

②

③

④

⑤

⑥

①

②

RV06.13  
RV06.15Seal FKM, PTFE-coated (1)  
Seal PTFE (2)

PC 1.1

Cable

labworldsoft®

## Error codes

Any malfunctions during operation will be identified by an error message on the display.

Once a serious error message has been displayed, the lift moves to the top end position and the device can no longer be operated. The lift can be operated again.

Proceed as follows in such cases:

- ☞ Switch off device using the main switch
- ☞ Carry out corrective measures
- ☞ Restart device

Error	Effect	Cause	Correction
Communication error	Distillation stopped	IR interface obstructed Heating bath switched off	- Press ESC to confirm - Switch on heating bath - Check and clean IR interface
No rotation	Distillation stopped, no rotary drive	Rotary drive does not start or exhibits very large speed variation	- Switch device off and on at mains switch
No temperature difference	Distillation stopped	No differential temperature increase of cooling water	- Press ESC to confirm - Check flow of cooling water - Check that the cooling water hoses are fitted correctly
Temperature out of range	Distillation stopped	Differential temperature outside measurement range (> 7 K)	- Press ESC to confirm - Check that the cooling water hoses are fitted correctly
No vacuum	Distillation stopped Vacuum cannot be controlled	Vacuum target value not reached	- Press ESC to confirm - Check that the glassware and hose connections are properly sealed - Switch on the pump - Check the solenoid valve connections
No venting	Vacuum cannot be vented	No venting of condenser	- Switch off the device - Vacuum sensor faulty - Manually vent the glassware using the stop valve
Water valve not closed	Water runs continuously into cooling circuit	Water valve not closed	- Press ESC to confirm - Check the connection
Flow rate out of range	Distillation stopped	Flow rate of cooling water out of range	- Press ESC to confirm - Readjust flow rate - Switch off flow rate control (for dry ice coolers)
Condenser overload	Distillation stopped Distillate may be able to enter the vacuum pump	Condenser overload	- Press ESC to confirm - Change process parameters, e.g. - Increase pressure - Reduce bath temperature - Increase condenser flow rate within specified range
Distillation dry	Distillation stopped	The system detected that there was no more medium in the evaporating flask during the distillation.	- Press ESC to confirm - Alternatively, activate the "Drying" option in the "Settings" menu
Lift not in end position	Distillation stopped Lift motor switches off	Obstruction in path of lift	- Press ESC to confirm - Check the path of the lift. - Remove any obstructions.
Calibration error	Calibration reset to factory setting	Out of calibration range (max. +/-0.5 K)	- Press ESC to confirm - Check that the sensor is connected correctly - Carry out a new calibration procedure

If the actions described fail to resolve the fault or another error code is displayed then take one of the following steps:

- Contact the service department,
- Send the device for repair, including a short description of the fault.

## Technical Data

Operating voltage range	Vac	100 - 230 ± 10%
Rated voltage	Vac	100 - 230
Frequency	Hz	50 / 60
Power input without heating bath	W	100
Power input operation "stand by"	W	3,3
Speed	rpm	20-270
Speed display		digital
Dimensions of visible display area (W x H)	mm	70 x 52
Display		TFT-Display
Multiple languages		yes
Right and left movement/ interval operation		yes
Smooth start		yes
Lift		automatic
Lifting speed	mm/s	50
Stroke	mm	140
Setting of lower end stop		60 mm, contactless
Head angle adjustable		0° - 45°
Timer		1-199 minutes
Interval		1-60 sec
Integrated vacuum controller		yes
Interface		USB, RS 232
Ramp programming		yes
Distillation processes programmable		yes
Remote control		with accessory labworldsoft®
Cooling surface (Standard condenser RV 10.1,10.10,10.2,10.20)	cm <sup>2</sup>	1200
Minimum cooling water flow rate	l/h	30
Maximum cooling water flow rate	l/h	100
Cooling water pressure	bar	0,5
Cooling water temperature range	°C	18 - 22, constantly
Fuse		2 x T1,6A 250V 5x20
Perm. On- time	%	100
Perm. ambient temperature	°C	5-40
Perm. humidity	%	80
Protection acc. to DIN EN 60529		IP 20
Protection class		I
Overvoltage category		II
Contamination level		2
Weight (no glassware; no heating bath)	kg	18,6
Dimensions (W x D x H)	mm	500 x 440 x 430
Operation at a terrestrial altitude	m	max. 2000

Subject to technical changes!

## Solvent table (excerpt)

Solvent	Formula	Pressure for boiling point 40 °C in mbar			
Acetic acid	C <sub>2</sub> H <sub>4</sub> O <sub>2</sub>	44	Ethylacetate	C <sub>4</sub> H <sub>8</sub> O <sub>2</sub>	240
Acetone	C <sub>3</sub> H <sub>6</sub> O	556	Ethylmethylketone	C <sub>4</sub> H <sub>8</sub> O	243
Acetonitrile	C <sub>2</sub> H <sub>3</sub> N	226	Heptane	C <sub>7</sub> H <sub>16</sub>	120
N-Amylalcohol,	C <sub>5</sub> H <sub>12</sub> O	11	Hexane	C <sub>6</sub> H <sub>14</sub>	335
n-Pentanol			Isopropylalcohol	C <sub>3</sub> H <sub>8</sub> O	137
n-Butanol	C <sub>4</sub> H <sub>10</sub>	25	Isoamylalcohol,	C <sub>5</sub> H <sub>12</sub> O	14
tert. Butanol,	C <sub>4</sub> H <sub>10</sub> O	130	3-Methyl-1-Butanol		
2-Methyl-2-Propanol			Methanol	CH <sub>4</sub> O	337
Butylacetate	C <sub>6</sub> H <sub>12</sub> O <sub>2</sub>	39	Pentane	C <sub>5</sub> H <sub>12</sub>	atm.press.
Chlorobenzene	C <sub>6</sub> H <sub>5</sub> Cl	36	n-Propylalcohol	C <sub>3</sub> H <sub>8</sub> O	67
Chloroform	CHCl <sub>3</sub>	474	Pentachloroethane	C <sub>2</sub> HCl <sub>5</sub>	13
Cyclohexane	C <sub>6</sub> H <sub>12</sub>	235	1,1,2,2-Tetrachloroethane	C <sub>2</sub> H <sub>2</sub> Cl <sub>4</sub>	35
Dichloromethane,	CH <sub>2</sub> Cl <sub>2</sub>	atm.press.	1,1,1,-Trichloroethane	C <sub>2</sub> H <sub>3</sub> Cl <sub>3</sub>	300
Methylenechloride			Tetrachloroethylene	C <sub>2</sub> Cl <sub>4</sub>	53
Diethylether	C <sub>4</sub> H <sub>10</sub> O	atm.press.	Tetrachloromethane	CCl <sub>4</sub>	271
1,2,-Dichloroethylene (trans)	C <sub>2</sub> H <sub>2</sub> Cl <sub>2</sub>	751	Tetrahydrofurane (THF)	C <sub>4</sub> H <sub>8</sub> O	357
Diisopropylether	C <sub>6</sub> H <sub>14</sub> O	375	Toluene	C <sub>7</sub> H <sub>8</sub>	77
Dioxane	C <sub>4</sub> H <sub>8</sub> O <sub>2</sub>	107	Trichloroethylene	C <sub>2</sub> HCl <sub>3</sub>	183
Dimethylformamide (DMF)	C <sub>3</sub> H <sub>7</sub> NO	11	Water	H <sub>2</sub> O	72
Ethanol	C <sub>2</sub> H <sub>6</sub> O	175	Xylene	C <sub>8</sub> H <sub>10</sub>	25

## Warranty

In accordance with **IKA** warranty conditions, the warranty period is 24 months. For claims under the warranty please contact your local dealer. You may also send the machine direct to our factory, enclosing the delivery invoice and giving reasons for the claim. You will be liable for freight costs.

The warranty does not cover worn out parts, nor does it apply to faults resulting from improper use, insufficient care or maintenance not carried out in accordance with the instructions in this operating manual.

Europe  
Middle East  
Africa

---

IKA® - Werke  
GmbH & Co.KG  
Janke & Kunkel-Str. 10  
D-79219 Staufen  
Tel. +49 7633 831-0  
Fax +49 7633 831-98  
sales@ika.de

North America

---

IKA® Works, Inc.  
2635 North Chase Pkwy SE  
Wilmington NC 28405-7419  
USA  
Tel. 800 733-3037  
Tel. +1 910 452-7059  
Fax +1 910 452-7693  
usa@ika.net

China

---

IKA® Works Guangzhou  
173 - 175 Friendship Road  
Guangzhou  
Economic and Technological  
Development District  
510730 Guangzhou, China  
Tel. +86 20 8222-6771  
Fax +86 20 8222-6776  
sales@ikagz.com.cn

Japan

---

IKA® Japan K.K.  
293-1 Kobayashi-cho  
Yamato Koriyama Shi, Nara  
639-1026 Japan  
Tel. +81 743 58-4611  
Fax +81 743 58-4612  
info@ika.ne.jp

India

---

IKA® India Private Limited  
814/475, Survey No. 129/1  
Mysore Road, Kengeri  
Bangalore - 560 060  
Karnataka, India  
Tel. +91 80 26253900  
Fax +91 80 26253901  
info@ika.in

Asia  
Australia

---

IKA® Works (Asia)  
Sdn Bhd  
No. 17 & 19, Jalan PJU 3/50  
Sunway Damansara  
Technology Park  
47810 Petaling Jaya  
Selangor, Malaysia  
Tel. +60 3 7804-3322  
Fax +60 3 7804-8940  
sales@ika.com.my

Korea

---

IKA® Korea Co LTD  
1710 Anyang Trade Center  
1107 Buhung-dong,  
Dongan-gu  
Anyang City, Kyeonggi-do  
Post code: 431-817  
South Korea  
Tel. +82 31-380-6877  
Fax +82 31-380-6878  
michael@ikakorea.co.kr

Brasilia

---

IKA® Works Inc.  
Rua São Bento, 701 Sala 1  
CEP 13160-000 Centro -  
Artur Nogueira - SP  
Brasil, South America  
Tel. +55 19 3877 2399  
Fax +55 19 3877 2399  
fcabral@ika.net