

Guide to Operations

inn6va™44/44R

Large Stackable Incubator Shakers

MANUAL NO: M1282-0050

Revision L

December 5, 2007



NEW BRUNSWICK SCIENTIFIC CO., INC.

BOX 4005 • 44 TALMADGE ROAD • EDISON, NJ 08818-4005

Telephone: 1-732-287-1200 • 1-800-631-5417 Fax: 732-287-4222 • Telex: 4753012 NBSCO

Internet: http://www.nbsc.com • E-mail: bioinfo@nbsc.com

INTERNATIONAL OFFICES:

BELGIUM

New Brunswick Scientific NV-SA Stationsstraat 180/4 3110 Rotselaar België/Belgique

Tel: 32 (0)16 56 28 31 Fax: 32 (0)16 57 27 53 E-mail: sales@nbsnv-sa.be

CHINA

New Brunswick Scientific Co., Inc. A903, 904 Yin Hai Building No. 250, Cao Xi Road Shanghai 200235, P.R. China

Tel: 86 21 6484 5955 or 5966

Fax: 86 21 6484 5933 E-mail: nbschc@online.sh.cn

FRANCE

France

New Brunswick Scientific SARL 12-14, Rond Point des Champs Elysées 75008 Paris

Tel: 33 (0)1 5353 1511 Fax: 33 (0)1 5353 1557

E-mail: sales@nbssarl.fr

GERMANY

New Brunswick Scientific GmbH

In Der Au 14 D-72622 Nürtingen

Deutschland

Tel: 49 (0)7022 932490 Fax: 49 (0)7022 32486 E-mail: sales@nbsgmbh.de

THE NETHERLANDS

New Brunswick Scientific BV Kerkenbos 1101, 6546 BC Nijmegen P.O Box 6826, 6503 GH Nijmegen

Nederland

Tel: 31 (0)24 3717 600 Fax: 31 (0)24 3717 640 E-mail: sales@nbsbv.nl

UNITED KINGDOM

New Brunswick Scientific (UK) Ltd.

17 Alban Park

St. Albans, Herts. AL4 0JJ

United Kingdom

Tel: 44 (0)1727 853855 or 0800 581331

Fax: 44 (0)1727 835666 E-mail: bioinfo@nbsuk.co.uk Web: www.nbsuk.co.uk



(A) CAUTION!

This equipment *must* be operated as described in this manual. If operational guidelines are not followed, equipment damage and personal injury can occur.

Please read the entire User's Guide before attempting to use this unit.

Do not use this equipment in a hazardous atmosphere or with hazardous materials for which the equipment was not designed.

New Brunswick Scientific Co., Inc. (NBS) is not responsible for any damage to this equipment that may result from the use of an accessory not manufactured by NBS.

Copyright Notice

New Brunswick Scientific Company, Inc. Box 4005 44 Talmadge Road Edison, New Jersey 08818-4005

© Copyright 2007 New Brunswick Scientific Co., Inc.

All Rights Reserved.

Reproduction, adaptation, or translation without prior written permission from New Brunswick Scientific is prohibited.

Disclaimer Notice

New Brunswick Scientific Co., Inc. reserves the right to change information in this document without notice. Updates to information in this document reflect our commitment to continuing product development and improvement.

Manual Conventions

NOTE:

Notes contain essential information that deserves special attention.



Caution messages appear before procedures which, if caution is not observed, could result in damage to the equipment.



Warning messages alert you to specific procedures or practices which, if not followed correctly, could result in serious personal injury.

Bold

Text in boldface type emphasizes key words or phrases.



This particular *Warning* message, whether found in the manual or on the unit, means HOT SURFACE—and therefore represents a potential danger to touch.



Crush Warning messages alert you to specific procedures or practices <u>regarding heavy objects</u> which, if not followed correctly, could result in serious personal injury.

CRUSH WARNING!

Biohazard-related messages.



WARRANTY

Every Instrument manufactured by the New Brunswick Scientific Co., Inc. is warranted to be free from defects in material and workmanship.

This apparatus is warranted for three (3) years against faulty components and assembly. During the first two (2) years of this warranty period, New Brunswick Scientific Co., Inc. (NBS) will also cover labor charges associated with the repair of this instrument. Our obligation under this warranty is limited to repairing or replacing the instrument or part thereof which shall, within three (3) years after the date of shipment, prove to be defective after our examination.

This warranty does not extend to any NBS products that have been subjected to misuse, neglect, accident or improper installation or application; nor shall it extend to products that have been repaired or altered outside the NBS factory without prior authorization from New Brunswick Scientific Co., Inc.

TABLE OF CONTENTS

1	OVERVIEW	10
2	INSPECTION & UNPACKING OF EQUIPMENT	12
2.	.1 Inspection of Boxes	12
2.		
2.	.3 UNPACKING OF EQUIPMENT	12
2.	.4 Out of Box Concerns	12
2.	.5 WARRANTY REGISTRATION	13
3	PREPARING THE LOCATION	14
3.	.1 PHYSICAL LOCATION	14
3.	.2 Environment	14
3.	.3 ELECTRICAL REQUIREMENTS	14
3.	.4 SPACE REQUIREMENTS	14
4	INSTALLATION	16
4.	.1 Tools Required for Installation	16
4.	.2 Installing the Front Grille	16
4.	.3 LEVELING A SINGLE SHAKER	17
4.	.4 LEVELING THE OPTIONAL BASE	18
4.	.5 MOUNTING INNOVA 44/44R ON OPTIONAL BASE	18
4.		
	4.6.1 For Two 44/44R Shakers	
	4.6.2 For an Innova 4400/4430 on an Innova 44/44R	
4.	.,	
4.		
4.	.9 STACKING A THIRD SHAKER	23
5	FEATURES	25
5.	.1 Controls	25
5.	.2 LCD DISPLAY	26
5.	.3 Changing Screens	27
5.		
5.	.5 Alarms	28
5.		
5.		
5.		
	.9 Interior Light(s)	
	.10 Heater	
	.11 Refrigeration (44R Only)	
5.	.12 SERVICE ACCESSIBILITY	31
6	GETTING STARTED	32

	6.1	PLATFORM ASSEMBLIES	32
	6.2	INSTALLATION OF PLATFORM	32
	6.3	INSTALLING FLASK CLAMPS	34
	6.4	ELECTRICAL CONNECTIONS	35
	6.5	OPTIONAL GASSING PORT KIT	35
	6.6	OPTIONAL UV GERMICIDAL LAMP	36
	6.7	OPTIONAL PHOTOSYNTHETIC LAMPS	36
	6.8	OPTIONAL HUMIDITY MONITOR	36
	6.9	FILLING THE WATER RESERVOIR	36
	6.10	DRAINING THE WATER RESERVOIR	37
7	OPI	ERATION	39
	7.1	SAFETY PRECAUTIONS	39
	7.2	OPENING THE DOOR	40
	7.3	STARTING THE SHAKER	40
	7.4	USING THE LCD SCREENS	40
	7.4.	l Display Screen	40
	7.4.2		
	7.4	3 Setup Screen	44
	7.4.4	4 Lamps Screen	46
	7.4.	5 RS232 Screen	47
	7.4.0	6 Calibrate Screen	49
	7.4.	7 Programs Screen	49
	7.5	PROGRAMMING THE SHAKER	50
	7.5.	<i>1 Timer Only</i>	50
	7.5.2	Programmed Steps	50
	7.5	3 Creating a Program	51
	7.5.4	4 Editing a Program	55
	7.5.3	5 Running a Program	55
	7.6	MUTING THE AUDIBLE ALARM	55
	7.7	TEMPERATURE OFFSET CALIBRATION	56
	7.7.	I Calculating the Offset Value	56
	7.7.2	2 Setting the Offset	56
	7.8	USING CALSPEED	
	7.9	Power Interruption	
	7.10	SLIDE-OUT PLATFORM	58
8	MA	INTENANCE	60
	8.1	ROUTINE MAINTENANCE	60
	8.2	CLEANING EXTERNAL & INTERNAL SURFACES	60
	8.3	BIOHAZARD DECONTAMINATION	60
9	SEF	RVICE & ACCESSORIES	62
	9.1	Troubleshooting	62
	9.2	PRODUCT RETURNS	
	9.3	OPENING THE SERVICE DRAWER.	

9.4	FUSE REPLACEMENT	65
9.5	REPLACING CABINET LIGHTS	66
9.6	REPLACING OPTIONAL PHOTOSYNTHETIC LAMPS	68
9.7	Adjusting Door Tension	
9.8	REMOVING THE SERVICE DRAWER	70
9.9	REPLACING OPTIONAL GERMICIDAL LAMP	71
9.10	BELT REPLACEMENT OR ADJUSTMENT	
9.11	REPLACEMENT PARTS	
9.12	Accessories	
	2.1 Slide-Out Platforms	
	2.2 Flask Clamps for Universal Platforms	
	2.3 Replacement Clamp Hardware Kits	
	2.4 Test Tube Racks & Other Accessories	
	2.5 Optional Bases	
	2.6 Optional Stacking Kit	
	2.7 Optional Gassing Port Kit	
	2.8 Optional Photosynthetic Lamps	
	2.9 Optional UV Germicidal Lamp	
	2.10 Optional Humidity Monitor	
	2.11 Optional Handles	
9.1	2.12 Optional Remote Alarm	80
10 I	DRAWINGS & TABLES	81
10.1	SCHEMATICS	81
10.2	LIST OF DRAWINGS	
10.3	LIST OF TABLES	
11 5	SPECIFICATIONS	85
11.1	Certifications	
	APPENDIX: REMOTE PROGRAMMING	
12.1	SETTING UP WITH HYPERTERMINAL	
12.2	OVERVIEW OF COMMAND SETS	
12.3	INDEX TO COMMAND CODES	
12.4	SET COMMANDS	
12.5	PROFILE CONTROL COMMANDS	
12.6	REPORT REQUEST COMMANDS	
12.7	SET/GET DATE & TIME	
13 A	APPENDIX: PRODUCT RETURNS	95
13.1	RETURN PROCEDURE	
13.2	RETURN AUTHORIZATION AND DECONTAMINATION CERTIFICA	TF 95
		11L 73

1 Overview

The Innova 44/44R Stackable Incubator Shakers are large-capacity orbital shakers that utilize an eccentric counter-balanced drive mechanism. They provide horizontal plane rotary motion in either a 1-inch (2.54 cm) or a 2-inch (5 cm) diameter circular orbit, depending on the model. A Proportional/Integral (PI) microprocessor controller with instantaneous digital feedback controls the speed over the entire speed range.

The Innova 44R provides temperature control from 20° below ambient (as low as 4°C) to 80°C, and the Innova 44 from 5° above ambient to 80°C. Naturally, both these ranges depend on relative humidity and other ambient factors, as well as the options installed in the unit. Ambient temperature is measured at one meter from the exterior of the unit.

Erlenmeyer flasks (up to 5 liters in size) as well as a wide variety of tubes and plates can be accommodated using the New Brunswick Scientific shaker accessories described in Section 9.12. These are easily accessed on slide-out platforms.

The Innova 44/44R may be operated in the following ways:

- **Continuously:** at a set speed and temperature, until user intervention.
- **In a timed mode:** run at a set speed, time and temperature for a period of up to 99.9 hours, after which the shaker automatically shuts off.
- Via the shaker's programmable controller: run through multiple temperature and speed changes for an extended period of time.
- Via computer through an RS-232 interface.

See Section 7 for more details on these various modes of operation.

For safe operation, the Innova 44/44R shakers are designed with a safety switch that automatically stops the shaker mechanism when the door is opened.

The Innova 44/44R is equipped with visual and/or audible alarms that alert the user to the following conditions:

- The end of a timed run
- Deviations from speed setpoint
- Deviations from temperature setpoint
- Power failure
- Door open
- Unbalanced load

To accommodate customer needs, a wide variety of platforms can be used with the Innova 44/44R:

- Universal platforms are the most flexible, providing hole patterns for flask clamps, test tube racks and other accessories.
- Dedicated platforms are supplied with flask clamps attached; they are designed solely and expressly for this purpose.
- Test tube racks, microplate holders, and test tube rack holders are also available (a universal platform is needed for all test tube racks and holders).

For further information on these accessories, see Section 9.12.

2 Inspection & Unpacking of Equipment

2.1 Inspection of Boxes

After you receive your order from New Brunswick Scientific, inspect the boxes carefully for any damage that may have occurred during shipping. Report any damage immediately to the carrier and to your local NBS Customer Service Department.

2.2 Packing List Verification

Verify against your NBS packing list that you have received all of the correct materials.

2.3 Unpacking of Equipment



CRUSH WARNING!

Do not attempt to lift the Innova 44/44R by hand. Always use a lifter or other suitable equipment when raising or handling the unit.

To unpack your Innova 44/44R, you will need the following tools:

- Claw hammer
- Forklift or other lifting equipment to lift more than 570 lbs (259 kg)
- Shears to cut ³/₄-inch (19 mm) wide steel strapping
- Tool to remove 3-inch (7.6 cm) metal staples

Save all packing materials and this user's guide/operating manual.

2.4 Out of Box Concerns

If any part of your order was damaged during shipping, is missing pieces, or fails to operate properly, please fill out the *Customer Satisfaction Form 6300*, included with your equipment, and return it by fax.

2.5 Warranty Registration

Please complete and return your warranty card or register electronically at our Website: www.nbsc.com

3 Preparing the Location

3.1 Physical Location

It is essential that the instrument be situated in an area where there is sufficient space for the shaker and platform to clear walls and potential obstructions during operation. The surface on which the unit is placed must be smooth, level, and able to support the shaker under full load operating conditions.



CRUSH WARNING!

REMINDER: Do not attempt to lift the Innova 44/44R by hand. Always use a lifter or other suitable equipment when raising or handling the unit.

3.2 Environment

The shaker is designed to operate optimally in the following ambient conditions:

- 10° to 35°C
- 20 to 80% Relative Humidity (non-condensing)

3.3 Electrical Requirements

The Innova 44/44R can be equipped to run on:

- 100 Volts, 50 or 60 Hz, 1500 VA maximum
- 120 Volts, 60 Hz, 1500 VA maximum
- 230 Volts, 50 Hz, 1500 VA maximum

Check your shaker's Electrical Specification Plate (located on the back of the unit) to be sure the electrical requirements of your unit match the output of your electric supply. If they do not match, contact your NBS representative.

3.4 Space Requirements

Allow at least four inches (10 cm) around the shaker for proper ventilation and for access to Power Switch and RS-232 port accessibility on the right side. Be sure to keep the power plug and power outlet easily accessible to facilitate unplugging the unit, as needed.

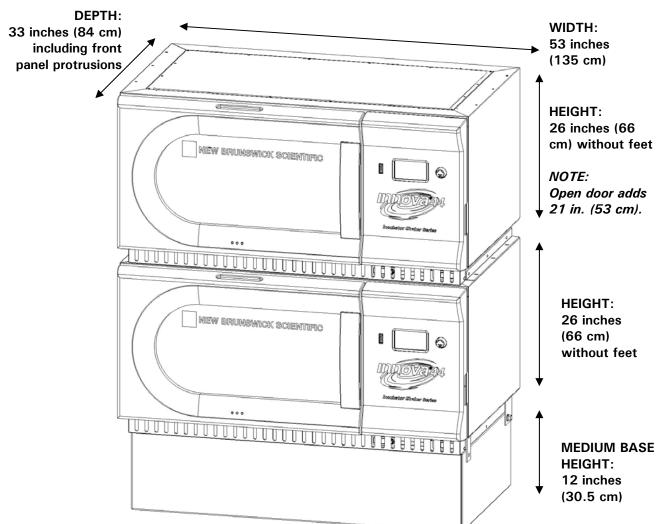


Figure 1: Innova 44/44R Front View Two Shakers Stacked on 12-Inch Base

NOTE:

Be sure to allow at least four inches (10.5 cm) around shaker for ventilation, access to power cord (rear panel), and access to power switch and RS-232 port (right side).

NOTE:

If you are stacking three shakers, use a four-inch base.

4 INSTALLATION



CRUSH WARNING!

REMINDER: Do not attempt to lift the Innova 44/44R by hand. Always use a lifter or other suitable equipment when raising or handling the unit.

NOTE:

Feet are used on the shaker ONLY when a single unit is operated without a base. At all other times (with two or three shakers stacked, or when a base is added), all feet must be removed.

4.1 Tools Required for Installation


To install the Innova 44/44R on an optional base and/or to stack Innova 44/44R units, the following tools will be needed:

- Number 2 Phillips head screwdriver
- Blade screwdriver
- 1/8-inch (3mm) Allen key
- 3/8-inch Allen key (provided)
- Level, \geq 10-inch (25.4 cm)
- Two adjustable wrenches or 7/16-inch socket wrench
- Forklift or other lifting equipment able to lift more than 570 lbs (259 kg)
- Metal leveling shims (provided in stacking kit)
- Optional: set of screw-in lifting handles, for small lifts of short duration (adding leveling shims, for example) or for minor location adjustments (see Section 9.12.11)
- Optional: stacking kit(s).

4.2 Installing the Front Grille

- 1. Remove and discard the plastic retaining clip that holds the drain hose in place for shipping.
- 2. Temporarily remove the thumbscrews (with their washers) installed at the bottom of each side panel, in the corner closest to the front of the shaker.
- 3. Hold the front grille in place (*see Figure 2 on the following page*) and reinstall the thumbscrews (and washers).

Figure 2: Installing the Front Grille

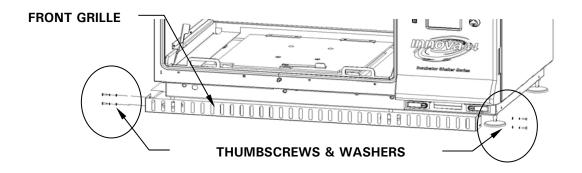

4.3 Leveling a Single Shaker

These instructions are for a single shaker without an added base.

Make sure that the shaker is placed on a level surface and that all four feet are solidly on the surface. If the shaker is not level, adjust the feet as needed to achieve leveling:

1. Immobilize the top lock nut against the unit with one wrench whenever you adjust the foot, to keep the threaded stud from falling out (see Figure 2a).

BOTTOM OF UNIT

1 7/16 inch
3.65 cm

FLATS FOR WRENCH

Figure 2a: Adjustable Foot (as shipped)

- 2. With a second wrench against the flats of the threaded stud, just above the foot. Rotate clockwise to lower the foot or counter-clockwise to raise the foot.
- 3. Place a level on the top of the unit. If necessary, make further adjustments by repeating all steps until the unit is level.
- 4. After installing a platform (*see Section 6*), fully load the shaker and do a test run at normal speed (*see Section 7*). Make additional leveling adjustments if necessary.

NOTE:

The maximum shaker lowering adjustment is $\frac{1}{2}$ inch (12.7 mm). The shaker cannot be raised above the height at which it was shipped. If more adjustment is required, you will need to add metal shims.

4.4 Leveling the Optional Base

Place the base on a sturdy (capable of bearing the weight of the combined base shakers and shaker contents), level surface, making sure that all four corners are solidly on the surface. If the base is not level, place metal shims as needed under the base until it is level.

4.5 Mounting Innova 44/44R on Optional Base

- 1. Using a forklift or lifter, raise the Innova 44/44R so that its back end is tilted toward the rear of the base.
- 2. Remove all four feet from the bottom of the shaker. They will not be needed for stacking, but you may want to keep them for future use.
- 3. With two assistants, guiding the shaker from opposite sides, lower the unit onto the base, back end first. Slowly and gently remove the forklift or lifter, lowering the front of the unit onto the base by hand.
- 4. Remove the set screws from the bottom sides of the shaker that correspond to the mounting holes for the rear corner braces (*see Figure 3 on the following page*).
- 5. Using the $\frac{1}{2}$ 13 x 1 $\frac{1}{4}$ -inch screws, lock washers and washers provided, mount the rear braces on the base. Do not fully tighten the screws yet.
- 6. Mount the rear of the Innova 44/44R to the corner braces on the base, using the $\frac{1}{2}$ 13 x 1½-inch screws, lock washers and washers provided. Do not fully tighten the screws yet.
- 7. Mount the front corner braces to the base using two sets of the $\frac{1}{2}$ 13 x 1 $\frac{1}{4}$ -inch screws, lock washers and washers provided. Do not fully tighten the screws yet.
- 8. With the remaining ½ x 13 x 1¼-inch screws, lock washers and washers, attach the front corner braces to the unit.
- 9. Tighten all screws equally.

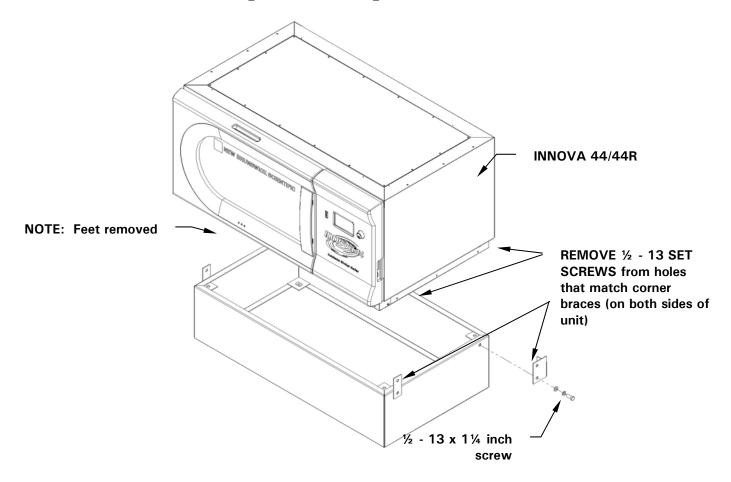


Figure 3: Mounting Shaker on Base

With the corner braces securely in place, the shaker mounted on its base will look as shown in Figure 4:

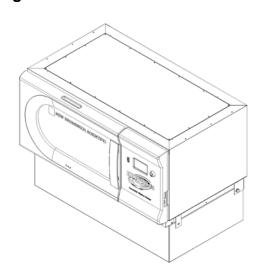


Figure 4: Shaker Installed on Base

- 10. Make sure the shaker is level; add metal shims under the base as needed to level the unit. Shims are provided with the stacking kit.
- 11. After installing a platform (*see Section 6*), fully load the shaker and do a test run at normal speed (*see Section 7*). Make additional leveling adjustments if necessary.


4.6 Installing the Stacking Kit

To stack two Innova 44/44R shakers, see Section 4.6.1. To stack an Innova 4400/4300 on an Innova 44/44R, skip to Section 4.6.2.

4.6.1 For Two 44/44R Shakers

- 1. Remove all feet from the units (putting them aside for possible future use).
- 2. If you do not plan to use an optional base, skip to Step 3. If you are using one, level the base according to Section 4.4, then mount the bottom unit on the base as explained in Section 4.5 above. If you are stacking three units, use a 4-inch base.
- 3. With reference to Figure 5a below, remove the set screws from the top of the bottom unit, on the two side edges, and remove the set screws (from the sides of the unit to be stacked) that correspond to the mounting brackets.

Figure 5a: Installing Innova 44/44R Stacking Kit

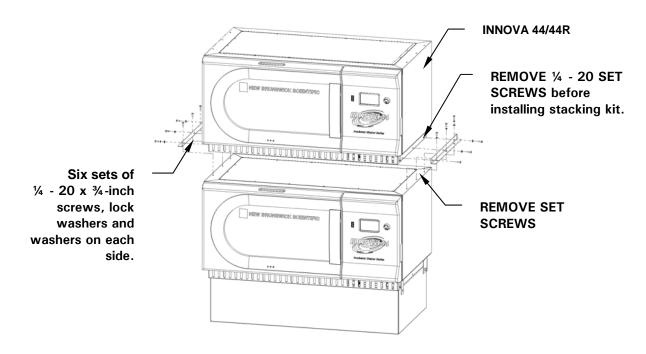


4. Using the screws and washers provided, attach the stacking kit mounting brackets to the top of the bottom unit. Secure the brackets in place as shown in Figure 5a above.

4.6.2 For an Innova 4400/4430 on an Innova 44/44R

- 1. Remove all feet from the units (putting them aside for possible future use).
- 2. If you do not plan to use an optional base, skip to Step 3. If you are using one, level the base according to Section 4.4, then mount the bottom unit on the base as explained in Section 4.5 above. If you are stacking three units, use a 4-inch base.

Figure 5b: Installing 4400/4430 to 44/44R Stacking Kit

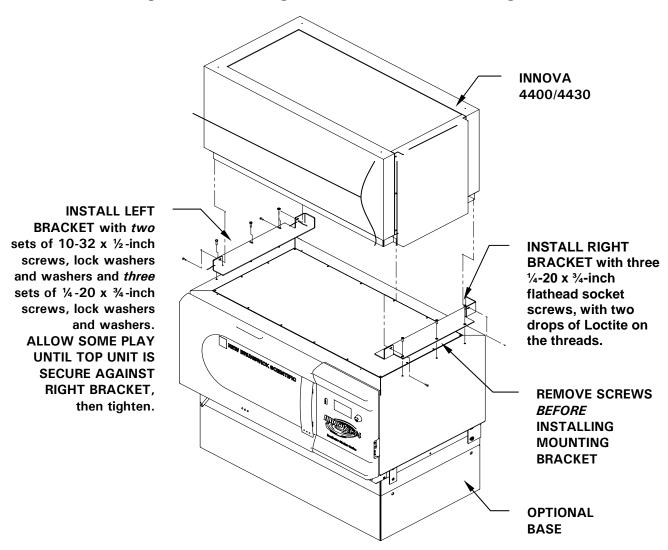

- 3. With reference to Figure 5b above, remove the screws installed at the top of the bottom unit, on the two side edges.
- 4. Use the countersunk flathead screws provided, with two drops of Loctite on their threads, to attach the righthand (facing the door) mounting bracket to the top of the unit installed on the base. Secure the bracket in place.
- 5. Use the other screws and washers provided to attach the lefthand mounting bracket to the top of the unit installed on the base. Loosely tighten these screws, to allow some play.

4.7 Stacking Two Innova 44/44R Shakers

- 1. Using a forklift or lifter, raise the Innova 44/44R to be stacked so that its back end is tilted toward the rear of the mounting brackets.
- 2. Remove the feet from the unit. They will not be needed for stacking, but you may want to keep them for future use.
- 3. With two assistants, guiding the unit from opposite sides, lower the shaker onto the mounting brackets, back end first. Slowly and gently remove the forklift or lifter, lowering the front of the unit onto the mounting brackets by hand.

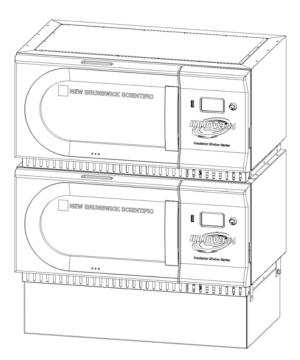


Figure 6: Stacked Shakers

- 4. As indicated in Figure 6, secure the bottom of the upper shaker to the mounting brackets (already fastened to the shaker below) with the screws, lock washers and washers provided.
- 5. Make sure the shakers are level; add metal shims under the bottom unit (or base) if needed.
- 6. After installing a platform (see Section 6), fully load the shaker and do a test run at normal speed (see Section 7). Make additional leveling adjustments if necessary.

If you are stacking a third shaker, skip to Section 4.9.

4.8 Stacking a 4400/4430 on a 44/44R

- 1. Using a forklift or lifter, raise the Innova 4400/4430 to be stacked so that its back end is tilted toward the rear of the mounting brackets.
- 2. Remove the feet from the unit. They will not be needed for stacking, but you may want to keep them for future use.
- 3. With two assistants, guiding the unit from opposite sides, lower the shaker onto the mounting brackets, back end first. Slowly and gently remove the forklift or lifter, lowering the front of the unit onto the mounting brackets by hand, and making sure it is tight against the righthand mounting bracket.
- 4. As indicated in Figure 5b, secure the bottom of the upper shaker to the mounting brackets with the screws, lock washers and washers provided. Be sure to tighten the right side first, then secure the left side.
- 5. Make sure the shakers are level; add metal shims under the base if needed.

4.9 Stacking a Third Shaker



When stacking three shakers, do not use a base taller than 4 inches (10 cm).



ldel CAUTION!

When stacking three shakers, all shaker loads must be balanced. These shakers operate best at maximum speed with a load of 34 (±3) lbs/15.5 (\pm 1.4) kg, which includes all platforms, clamps, and filled glassware.



(A) CAUTION!


For three stacked 2-inch stroke shakers, maximum operating speed must be limited to 250 rpm.

To stack a third shaker, repeat the appropriate procedures provided above to install the stacking kit, to stack the shaker, and to level the entire assembly.

5 FEATURES

5.1 Controls

Figure 7: Front Panel (detail)

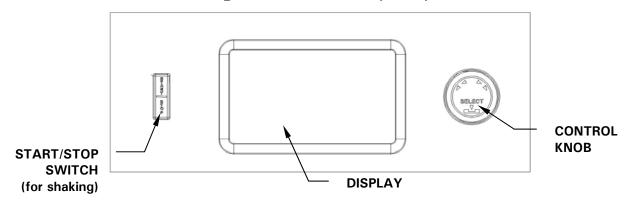


START/STOP SWITCH

This switch is used to start or stop the shaker. It will also activate the timer when a timed run is desired. If the unit is stopped and restarted, the timer automatically returns to the beginning of a run.

CONTROL KNOB

This knob is multifunctional. It is used to change screens, and to select and change operating conditions.

Figure 8: Control Panel (Right Side)

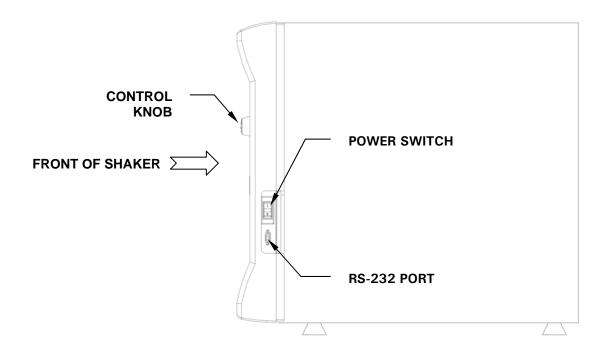

• **RS-232 PORT** See Section 5.8 for details.

• **POWER SWITCH** This rocker switch is a circuit breaker that turns power on and off to the entire Innova 44/44R.

NOTE:

In addition to the Power Switch, the Power Cord is also used to conduct power or to break the power circuit to the shaker. Whenever power to the shaker may be a hazard (during cleaning, maintenance or service work), be sure to disconnect the power cord from the electrical outlet.

5.2 LCD Display

When the unit is powered up, using the Power Switch located on the lower right side of the machine (*see Figure 8*), the initial display screen will remain for a short time while the system boots up. Then the main screen, called **DISP** for Display, appears (*see Figure 9*). This screen will indicate the same parameters that were in effect when the power was turned off.

PARAMETERS

PARAMETERS

PARAMETERS

PARAMETERS

PARAMETERS

PARAMETERS

PARAMETERS

DISP

DISP

DISP

DISP

DISP

DISP

DISP

PROGRAM RUNNING

HEATER ON

DOOR OPEN

DOOR OPEN

Figure 9: Display Screen

- For an explanation of the icons on the display screen, see Section 5.4.
- Turning the Control Knob will highlight functions and/or values that can be changed.
- Whenever you turn the Control Knob, the chamber light will go on (for details, see Section 5.9).

For more information on working in the Display Screen, see Section 7.4.1.

5.3 Changing Screens

You can change screens displayed by highlighting the screen name field in the lower left corner, pressing the Control Knob in until it clicks, rotating the Knob left or right (which also makes clicking sounds) to the desired screen and clicking the Control Knob in again. Table 1 describes the various screens:

Table 1: Screens

Screen	Meaning	Features/Modes
Name		
DISP	Display	Shows two user-selectable parameters* and actual values.
SUMM	Summary	Shows all parameters*, setpoints and actual values.
SET	Set-Up	Set day of week, set time, enable or mute alarm, lock or
	-	unlock operating parameters
LAMP	Lamps	Internal chamber light: ON (always on); OFF (always
		off); AUTO (default mode), light goes on and stays on
		when door is open, shuts off 15 seconds after door is
		closed, and goes on for 15 seconds when Control Knob is
		moved.
		¤Photosynthetic lights (GRO): ON, OFF, [★] NONE
		¤UV Light (UV) : ON, OFF, *NONE
COMM	Communication	SET: set baud rate
	(RS-232)	OFF disables RS-232
		MONITOR: PC commands shaker to read setpoints and
		actual values on a schedule determined by PC software.
		Parameters are unlocked and can be changed by program
		or manually.
		SLAVE: PC controls shaker and logs data.
		TALK: Shaker sends setpoint and actual data to PC at
		one-minute intervals.
CAL	Calibration	Allows user to enter a temperature offset.
		Allows user to calibrate speed.
PROG	Program	Allows user to set up 1- 4 programs, each with 1-15 steps.

^{*} see Table 2, Display Screen Parameters

Table 2: Display Screen Parameters

Parameter Name	Meaning
RPM	Shaking speed, in revolutions/minute
°C	Chamber temperature, in degrees Celsius
HRS	Programmed time remaining, in hours
¤ %RH	Relative Humidity, in percent.
д UV	Status of Ultraviolet germicidal lamp
¤ GRO	Status of Photosynthetic growth lamps

[□] Optional

 [□] Optional

^{*}Not installed

5.4 Display Icons

There are six icons that help identify operating conditions. Five of these icons are located at the bottom of the display, and they are visible, when applicable, in any screen you are viewing, except the Program subscreens. The sixth appears, when applicable, next to °C whenever the temperature is onscreen.

Table 3: Display Icons

Icon	Explanation
X	This icon appears when audible alarms are muted.
	This icon appears when the possibility to make manual or programmed parameter changes is disabled (locked). This is controlled by settings on the SET screen. See Section 7.4.3
鱼	This icon appears when the shaker door is open.
* *	This icon appears when the Heater is on. See Section 5.10.
	This icon appears when a user-defined Program is running. See Section 7.5.
*	This icon appears to the right of °C if the Temperature Offset feature is being used. See Section 7.7.

5.5 Alarms

If an alarm condition exists, the field in the lower right corner will alternate the Day and Time with characters indicating the nature of the alarm condition, accompanied by an audible alarm (unless muted):

Table 4: Alarms

Indication	Description
TEMP	The temperature deviates more than $\pm 1^{\circ}$ C from setpoint after achieving control temperature range. After door is opened, alarm will be disabled for 5 minutes while chamber recovers to setpoint.
RPM	The speed deviates more than ± 5 RPM from setpoint after achieving operating speed setpoint. After door is opened, alarm will be disabled for 5 minutes while chamber recovers to setpoint.
POWER	Indicates unit is powering up (both at normal power-up and after power interruption); will flash until the Control Knob is moved.
HRS	Indicates when timed run is completed.
TILT	Indicates an unbalanced shaking condition. After load is evenly distributed and/or shaker is leveled, restart by pressing START/STOP.

5.6 Glide-Up Door

The Innova 44/44R is equipped with a space-saving glide-up door. To release the door from the latch mechanism, firmly press the bottom of the large, curved push bar (located on the right side of the door), then manually move the door up to the open position.

When the door opens, the following will happen:

- Heater turns off
- Shaker stops
- Interior light goes on and, when it is in **AUTO** mode, will remain on for 15 seconds after the door is closed
- UV germicidal lamp (if so equipped) turns off

5.7 Spill Pan/Water Reservoir

The Innova 44/44R is equipped with a spill pan to protect the drive mechanism in the case of accidental spills and/or broken glassware. This pan can also be used as a water reservoir to humidify the chamber and to reduce evaporation. An optional factory-installed humidity monitor is also available.

5.8 Software Interfaces

The RS-232 port is located below the Power Switch on the right side of the control panel (see Figure 8). It can be used to interface a computer to the shaker for control of operating conditions or data logging applications (see Section 12).

The customer is responsible for securing the proper driver to interface with the RS-232.

5.9 Interior Light(s)

When the **LAMP** screen is in its default **AUTO** mode, the interior ("chamber") light is activated for 15 seconds whenever you turn the Control Knob. It will automatically shut off after 15 seconds of Control Knob inactivity.

The chamber light will also go on when the door is open.

In addition, you can set the chamber light to be continuously **ON** or **OFF** by selecting either mode in the **LAMP** screen.

There are two additional light options: interior photosynthetic growth lamps (*see Section 6.7*) and, for refrigerated units only, a UV germicidal lamp located outside the chamber but in the airflow path (*see Section 6.6*).

5.10 Heater

The chamber temperature is sensed by a 1000 ohm platinum RTD. The heater is controlled using pulse width modulation on a 2.5-second duty cycle. This cycle time is fast enough to prevent noticeable changes in air temperature due to the cycling.

Whenever the heater is on, the Heater On icon will appear in the display. The heater automatically stops running when the door is opened.

5.11 Refrigeration (44R Only)

The refrigeration system in the Innova 44R is a fixed-capacity system carefully designed to maintain the setpoint, to balance pressure within the system, and to prevent freezing on the evaporator surface.

When the shaker is powered up, there is a four-minute time delay prior to compressor start-up.

5.12 Service Accessibility

In the unlikely event that your Innova 44/44R should need service, all electronic boards, refrigeration, heating components are mounted on a pull-out drawer mechanism which is easily accessed, **by an authorized service technician**, from the front of the shaker.

6 GETTING STARTED

6.1 Platform Assemblies

The Innova 44/44R can be used with a variety of NBS platforms that will accept a wide range of clamps for flasks, test tubes, etc. A platform, which is required for operation, is a separate item, not included with the shaker assembly. Refer to Section 9.12 for details on available platforms and platform accessories.

6.2 Installation of Platform

NOTE:

In transit, there are two small plastic straps on the side of the bearing housing to secure the slide-out mechanism, and two small plastic straps that immobilize the bearing housing; all straps must be removed.

When you cut the bearing housing straps, be particularly careful not to cut the wiring that is in close proximity to the straps.

Prior to use, a subplatform and platform must be installed on the unit. The shaker is shipped with four Allen head platform screws installed in the bearing housing (see Figure 10a below, which also shows the subplatform that you must install).

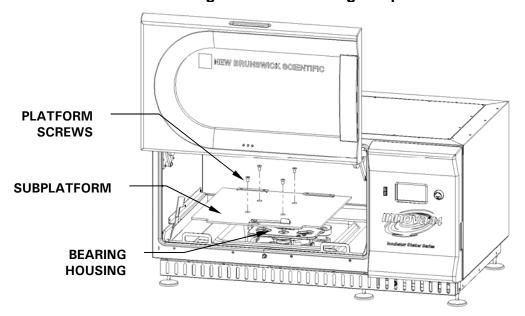


Figure 10a: Installing Subplatform

Using Figure 10a as your guide, install the subplatform:

- 1. Remove the Allen head platform screws, setting them aside.
- 2. Place the subplatform onto the bearing housing, taking care to orient it as shown in Figure 10a, with the notches and the lever toward the front of the shaker, lever facing up.
- 3. Align the subplatform holes with the tapped holes in the bearing housing, then secure the subplatform with the Allen head platform screws you previously removed.

To install the slide-out platform you have purchased:

1. With reference to Figure 10b below, place the platform on the subplatform, making sure the two slots in the back of the platform are inserted under the subplatform's blocks in the back, and the slide-out tray handles are at the front edge, facing up.

Figure 10b: Installing Slide-Out Platform

Steps 2-4 are called out in Figure 10b:

4 Rotate subplatform cam lever 180° counterclockwise (away from you), until it points to the right. This locks the platform in place. ③ Simultaneously push both side cam levers up (away from you and toward the back of the chamber). 2 Push both slide-out trav handles away from you until the platform hits its stops at the

New Brunswick Scientific Co., Inc.

back of the chamber.

6.3 Installing Flask Clamps

Flask clamps purchased for use with universal platforms (see Section 9.12.2) require installation. Clamps are installed by securing the base of the clamp to the platform with the correct type and number of screws. All clamps are shipped complete with hardware.

Clamps for 2-, 2.8- and 4-liter flasks are shipped with an additional girdle to keep the flasks in place. The girdle is an assembly of springs and sections of rubber tubing. One girdle is already in place on the clamp, the other is packed separately. To install these double girdle clamps:

- 1. Place the clamp on the platform, aligning its mounting holes with holes on the platform. Secure the clamp in place using the flat Phillips head screws provided (#S2116-3051, 10-24 x 5/16-inch). *Use Figure 11b to help you identify the proper screws, as three different types of screws are shipped with the clamps.*
- 2. With the first girdle in place, as delivered, on the upper part of the clamp body (*see Figure 11a*), insert an empty flask into the clamp.
- 3. After making sure the sections of tubing are located between the clamp legs, roll the first girdle down the legs of the clamp as far as it can go. The tubing sections will rest against the platform, and the springs will be under the clamp base.
- 4. Place the second girdle around the upper portion of clamp body (just as the first girdle was initially). Make sure that its spring sections rest against the clamp legs, while its rubber tubing sections sit against the flask, in between the clamp legs.

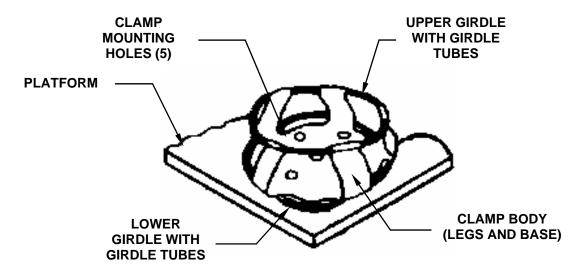


Figure 11a: Double Girdle Clamp Installation

Figure 11b: Clamp Fastener



NOTE:

The upper girdle secures the flask within the clamp, and the bottom girdle keeps the flask from spinning.

6.4 Electrical Connections

Before making electrical connections, verify that the power source voltage matches the voltage on the electrical specification plate and that the **ON/OFF** switch is in the **OFF** position. The electrical specification plate is located on the rear panel of the unit near the power connector.

Connect the power cord to the power connector, then connect the other end to a suitable, grounded receptacle. Make sure there is enough clearance to disconnect the plug whenever necessary.

POWER CONNECTOR

AIR INTAKE VENTS (do not block!)

ADJUSTABLE FOOT

6.5 Optional Gassing Port Kit

This option is factory-installed. The manifold delivers gas into the chamber via six ports. You may elect to use splitters after the manifold to increase the number of ports for your application.

6.6 Optional UV Germicidal Lamp



WARNING!

Never try to operate the UV Germicidal Lamp while the shaker door is open.

This option, a germicidal ultraviolet lamp, is placed inside the Service Drawer, outside the chamber to help reduce the risk of contamination. The lamp is identified on the display screen as **UV**.

The UV germicidal lamp is factory-installed and available on refrigerated units only. Operation instructions are provided with this option (see Section 9.12.9).

6.7 Optional Photosynthetic Lamps

This factory-installed option provides nine photosynthetic growth lamps inside the chamber. They can be turned on and off manually, by the easily set programmatble timer, or by computer. These lamps are identified on the display screen as **GRO**. The replacement bulb part number is P0300-0221.

The recommended operating temperature when using this option is 15 - 37°C. The maximum operating temperature when using photosynthetic lamps is 70°C.

6.8 Optional Humidity Monitor

This optional factory-installed accessory allows you to monitor actual humidity levels in the chamber throughout your run.

6.9 Filling the Water Reservoir

If you choose to use the spill pan as a water reservoir to reduce evaporation:

1. Temporarily remove the four thumbscrews (with their washers) that secure the front grille in place (*see Figure 13 on the following page*).

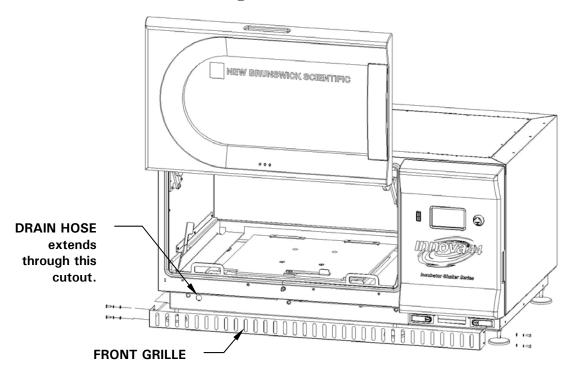


Figure 13: Front Grille

- 2. Check the drain hose to be sure it is clamped tightly.
- 3. Replace the grille and secure it with the four thumbscrews and the washers that you previously removed.

NOTE:

Never pour water directly under the subplatform. Pour water *very slowly* into the shallow area beyond the edge of the adapter plate, to protect the bearing housing.

4. Accessing the pan/reservoir from the left, the right or in front of the subplatform, **slowly** fill the reservoir with no more than 3 liters of distilled water.

6.10 Draining the Water Reservoir

To drain water from the water reservoir/spill pan:

- 1. Temporarily remove the four thumbscrews that secure the front grille (*see Figure 13*).
- 2. Disengage the drain hose, direct it to a container or drain, then unclamp it, allowing the water to gravity drain.
- 3. Reclamp the hose, tuck it back in place, and replace the grille, securing it with the thumbscrews you previously removed.

The spill pan/water reservoir drain is located in the front, on the left, under the slide-out tray.

NOTE:

The drain hose should be tightly clamped when it is not being used for draining the reservoir. Keep it clamped even if there is no water in the pan.

7 OPERATION

7.1 Safety Precautions

Before operating the shaker, verify that anyone involved with its operation has been instructed in both general safety practices for laboratories and specific safety practices for this apparatus.



The user is also responsible for following local guidelines for handling hazardous waste and biohazardous materials that may be generated from the use of this equipment.



If service should be required on a unit that is going to be returned to an NBS facility, it must be completely decontaminated and cleaned prior to its return, and a Returned Material Safety Sheet must be filled out to certify that you have complied. See Section 13, *Product Returns*.



It is the responsibility of the user to carry out appropriate decontamination procedures if hazardous material is spilled on or inside the equipment. Before using any cleaning or decontamination method other than those suggested by the manufacturer, users should check with New Brunswick Scientific that the proposed method would not damage the equipment.

This equipment is **not** "explosion-proof" and should never be used with flammable substances or used to grow organisms that produce flammable by-products.



CAUTION!

To prevent damage to the shaker and its contents, never run the shaker without a platform.

7.2 Opening the Door

Open the door by firmly pressing the bottom of the large curved push bar (located on the right side of the door) to release the latching mechanism. You can now manually move the door up to the open position, or down to the closed position (make sure it latches closed).

7.3 Starting the Shaker

To initially start the shaker, close the door and turn the power switch (located on the righthand side of the control panel) to the **ON** position. The display will come on (first showing only **New Brunswick Scientific**, then briefly displaying the model number, **44** or **44**R, and the stroke, **1 inch** or **2 inch**, and then quickly moving into the Display screen), and the audible alarm will sound. To mute it, see Section 7.6.

When the shaker begins to operate, the LCD display will track the speed as it accelerates to the last entered setpoint. The shaking action may be started or stopped by pressing the Start/Stop button on the front panel.

NOTE:

The shaker will not operate if the door is open. This is indicated by the "door open" symbol appearing in the bottom line of the display (see Figure 9 or Table 3).

7.4 Using the LCD Screens


7.4.1 Display Screen

When you turn the power on, this is the first screen to appear after the company title screen. The default display parameters are temperature (°C) and shaking speed (RPM).

You can change the displayed parameters. **To replace a parameter**:


1. Using the Control Knob, highlight the parameter that you wish to replace. For this example, we will replace **RPM** (*see Figure 14*).

Figure 14: Changing Display Parameter

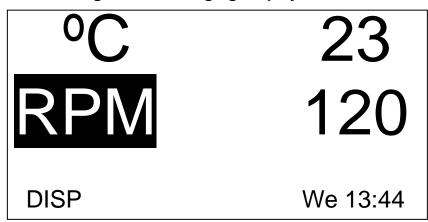

- 2. Click the Control Knob in. **RPM** will flash.
- 3. Turn the Knob until the desired parameter appears in the highlighted field. For this example, we will select **HRS**.
- 4. Click the Knob in, to set and save the parameter (see Figure 15).

Figure 15: Changed Display Parameter

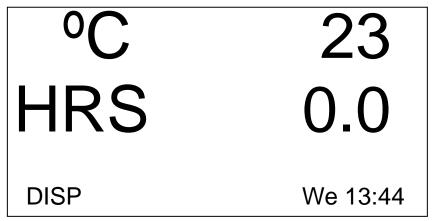

NOTE:

If you highlight an item, change it, but do not save your selection, after a few seconds the screen will revert to its previous setting.

NOTE:

"UV" (UV Germicidal lamp) and "GRO" (photosynthetic growth lamps) will appear in this screen with the word "NONE" if your shaker is not equipped with these optional features.

You can also use this screen to verify a setpoint, even though the values displayed here are actual (current) values. **To view a setpoint**:


- 1. Use the Control Knob to highlight the value (in this example, we will view the temperature setpoint, so we will highlight the current °C, which is 23).
- 2. Click the Knob in to display the current setpoint, which will flash.

At this point you can modify the setpoint or click the Knob in again to return to the normal display, which will be the actual temperature.

To modify a setpoint in this screen:

- 1. Use the Control Knob to highlight the current value (we will continue to use the temperature as our example, so we will select **23**).
- 2. Click the Knob in to display the current setpoint (in this example, **38.5**—*see Figure 16*), which will flash.

Figure 16: Changing Setpoint

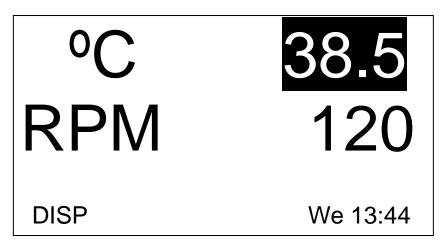

3. Turn or spin the Knob to reset the setpoint (in this example, turn the Knob left to decrease the setpoint to **37.0**).

NOTE:

If you turn the Control Knob slowly, one click left or right will change the setpoint by an increment of one tenth of a degree Celsius (0.1°C). If you spin the Knob, the value will change by larger increments.

4. Click the Knob in to set and save this new setpoint.

NOTE:

If you highlight an item, change it, but do not save your selection, after a few seconds the screen will revert to its previous setting.

5. The display will automatically return to the actual value.

To move out of this screen and into the next:

- 1. Use the Control Knob to highlight **DISP**, then click the Knob in. **DISP** begins to flash.
- 2. Turn the Knob to the right until the next screen, Summary (**SUMM**), appears. If you turn too far and enter another screen, just turn the Knob back to the left to recapture the **SUMM** screen.
- 3. Click the Knob in to select the screen and to work in it.

7.4.2 Summary Screen

In this screen (*see Figure 17*), you can see both the current **ACTUAL** readings and the **SET**points for shaking speed (**RPM**), chamber temperature (°**C**), elapsed time in a programmed run (**HRS**) and, if you are using the optional Humidity Monitor (*see Section* 9.12.10), the percentage of relative humidity (%**RH**).

PARAM ACTUAL

PARAM	ACTUAL	SET
RPM °C	100 45.1	100 45.0
HRS	0.0	0.0
%RH	50.0	N/A
SUMM		Th 16:18

NOTE:

The current day (Su, Mo, Tu, We, Th, Fr or Sa) and time always remains visible in the lower righthand corner.

The only elements you can modify in this screen are setpoints. **To change setpoints in this screen**:

- 1. Turn the Knob until the desired setpoint is highlighted, then click the Knob in. The setpoint will begin to flash.
- 2. Turn the Knob to the right to increase the number, or to the left to decrease it. One click left or right will increase the setpoint by an increment of one (one whole unit or one tenth unit, depending on the parameter). Move the Knob more rapidly (you can spin it) to change the value by larger increments.
- 3. Click the Knob in to set and save the new value.

NOTE:

If you highlight an item, change it, but do not save your selection, after a few seconds the screen will revert to its previous setting.

4. Repeat the above steps to change any or all of the other setpoints.


To move out of this screen and into the next:

- 1. Use the Control Knob to highlight **SUMM**, then click the Knob in. **SUMM** begins to flash.
- 2. Turn the Knob to the right until the next screen, Setup (SET), appears. If you turn too far and enter another screen, just turn the Knob back to the left to recapture the SET screen.
- 3. Click the Knob in to select the screen and to work in it.

7.4.3 Setup Screen

Here you can set the day of the week and the time (on a 24-hour clock). This screen also allows you to lock all of your settings from further changes, and to mute or enable the audible alarm.

Figure 18: Setup Screen

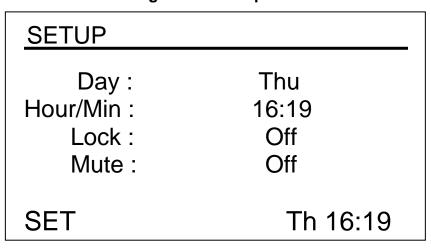

To change the day:

- 1. Turn the knob to highlight the day (**Thu** in the sample screen above), then click inward once. The day will flash.
- 2. Turn the knob left or right to select the day of choice: Sun, Mon, Tue, Wed, Thu, Fri or Sat.
- 3. Click the knob in to set and save your choice.

NOTE:

If you highlight an item, change it, but do not save your selection, after a few seconds the screen will revert to its previous setting.

To change the time (Hour/Min):

- 1. Turn the knob to highlight the time (**16:19** in the sample screen above), then click inward once. The time will flash.
- 2. Turn the knob left or right to change the time. Left moves backward, right moves forward in time. One click right or left changes by one minute; spin the knob to move more rapidly.
- 3. Click the knob once inward to set and save your choice.

To lock the settings:

- 1. Turn the knob to highlight **Lock**, then click inward once. The current status (**Off** in the sample screen above) will flash.
- 2. Turn the knob in either direction; the only other choice is **On**. Click once inward to select and save **On**, or continue turning to return to **Off**.
- 3. When you set **Lock** to **On**, a padlock icon (*see Section 5.4*) will appear at the bottom of the screen. This icon will remain on display through all main display screens until you turn the locking function off.

To mute the audible alarm:

- 1. Turn the knob to highlight **Mute**, then click inward once. The current status (**Off** in the sample screen above) will flash.
- 2. Turn the knob in either direction; the only other choice is **On**. Click once inward to select and save **On**, or continue turning to return to **Off**.
- 3. When you set **Mute** to **On**, the icon will appear at the bottom of the screen. This icon will remain on display through all screens until you turn the muting function off.

To move out of this screen and into the next:

- 1. Use the Control Knob to highlight **SET**, then click the Knob in. **SET** begins to flash.
- 2. Turn the Knob to the right until the next screen, Lamps (LAMP), appears. If you turn too far and enter another screen, just turn the Knob back to the left to recapture the LAMP screen.
- 3. Click the Knob in to select the screen and to work in it.

7.4.4 Lamps Screen

In this screen (*see Figure 19*), you can turn the chamber light (**Chamber**), the optional UV germicidal ("decontamination") lamp (**UV Decont**) and the optional photosynthetic growth lamps (**Growth**) on and off.

On means the light is always on, and **Off** means the light is always off, unless you add additional programming (see Section 7.5).

There is an additional mode for the chamber light: **Auto**. In **Auto** mode, the light will go on every time you activate the control knob or open the door. This is the default mode.

NOTE:

The Lamps screen, shown below, will always indicate Chamber light mode. If the shaker is not equipped with the optional UV germicidal lamp and/or photosynthetic growth lamps, UV Decont and/or Growth will say None.

Figure 19: Lamps Screen

Chamber: Auto
UV Decont: Off
Growth: Off

Th 16:20

To change the mode setting for any of the lamps:

- 1. Turn the Control Knob to highlight the setting for the lamp of choice, then click the knob in. The current setting will flash (in the sample screen, we will use the **Chamber** light as an example.
- 2. Turn the Knob left or right until the desired mode setting appears (**Auto** in this example).
- 3. Click the Knob in to save the new setting.

NOTE:

If you highlight an item, change it, but do not save your selection, after a few seconds the screen will revert to its previous setting.


To move out of this screen and into the next:

- 1. Use the Control Knob to highlight **LAMP**, then click the Knob in. **LAMP** begins to flash.
- 2. Turn the Knob to the right until the next screen, RS232 (COMM), appears. If you turn too far and enter another screen, just turn the Knob back to the left to recapture the COMM screen.
- 3. Click the Knob in to select the screen and to work in it.

7.4.5 RS232 Screen

This screen (*see Figure 20*) is used only if you have connected a personal computer to the RS-232 port (*see Section 5.8*). Here you can select the RS-232 port's **Mode** and the **Baud Rate** appropriate to your PC.

Figure 20: RS232 Screen

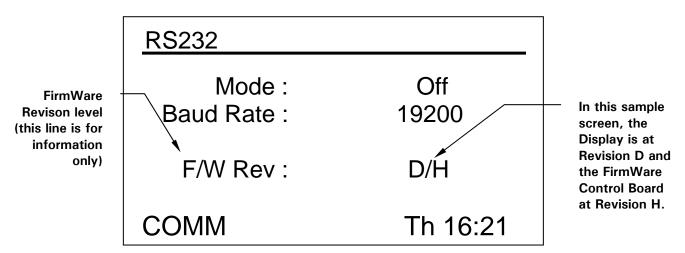

To change the Communication Mode:

- 1. Turn the Control Knob to highlight the current setting (**Off** in the sample screen above), then click the Knob in. The current setting will flash.
- 2. Turn the Control Knob left or right until the desired mode setting appears (*see Table 5: Communication Mode*):

Table 5: Communication Mode

Mode	Application	
Off	The RS-232 port is not open to communication in either direction.	
Slave	The shaker can be fully controlled from the computer.	
Talk	The shaker sends current value reports to the computer once per minute.	
Monit [Monitor]	The shaker responds only to "Report Requests".	

3. Click the control knob inward once to save the new setting.



If you highlight an item, change it, but do not save your selection, after a few seconds the screen will revert to its previous setting.

To change the Baud Rate:

- 1. Turn the control knob to highlight the current setting (**19200** in the sample screen above), then click the knob inward once. The current setting will flash.
- 2. Turn the control knob left or right until the desired setting appears: 9600, 19200 or 38400. The setting you choose should match the baud rate of your computer.

3. Click the control knob inward once to save the new setting.


To move out of this screen and into the next:

- 1. Use the Control Knob to highlight **COMM**, then click the Knob in. **COMM** begins to flash.
- 2. Turn the Knob to the right until the next screen, Calibrate (CAL), appears. If you turn too far and enter another screen, just turn the Knob back to the left to recapture the CAL screen.
- 3. Click the Knob in to select the screen and to work in it.

7.4.6 Calibrate Screen

Use this screen (see Figure 21) to create a temperature offset and to calibrate the shaking speed (for details, see Sections 7.7, Temperature Offset Calibration, and 7.8, Using Calspeed).

Figure 21: Calibrate Screen

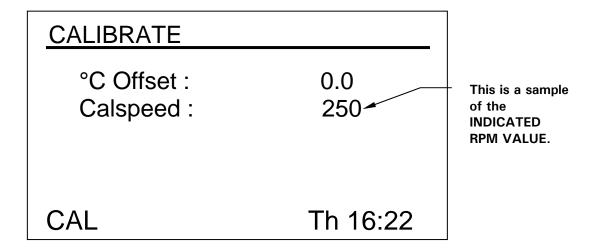

7.4.7 Programs Screen

Use this screen (*see Figure 22*) to set as many as four operating programs for the shaker. Each program can have as many as 15 steps. See Section 7.5 below for complete details.

Figure 22: Programs Screen

PROGRAMS	
Program 1 :	Off
Program 2:	Off
Program 3:	Off
Program 4 :	Off
PROG	Th 16:23

7.5 Programming the Shaker

7.5.1 Timer Only

By setting an **HRS** setpoint in the **DISP** or **SUMM** screen, the shaker may be set to automatically stop after a preset time period of 0.1 to 99.9 hours.


If the time is set to 0.0, the shaker will operate continuously until either the door is opened or the Start/Stop button is pushed.

7.5.2 Programmed Steps

The resident software for the Innova 44/44R can store up to four programs, each having as many as 15 steps. *Each step* can be programmed in one-minute increments, for total periods of one minute up to 99 hours 59 minutes each.

To enter the programming mode, use the Control Knob to select the **PROG** screen (*see Figure 22, repeated below*). At this point, you can **Run** a program, **Edit** a program, generate a **New** program or turn a program **Off**. **Off** is the default mode.

Figure 22: Programs Screen

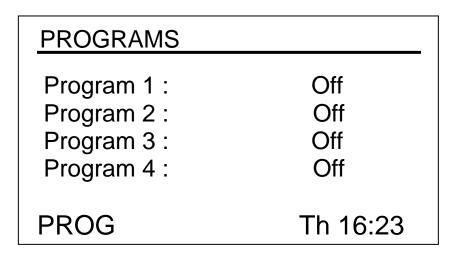


7.5.3 Creating a Program

To write a **New** program:

- 1. Use the Control Knob to highlight the mode of Program 1 (in our sample screen, Figure 22, it is **Off**), then click the Knob in. The selected field begins to flash.
- 2. Turn the Knob until the field says **New**. Click the Knob in to select this mode. The screen for Program 1 Step 1 will open (*see Figure 23*):

Figure 23: Program 1 - Step 1

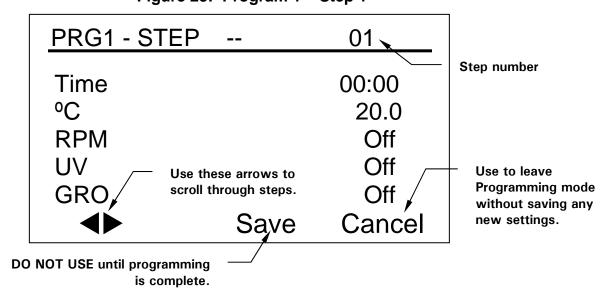

NOTE:

If your shaker is not equipped with these optional features, "UV" and "GRO" will appear with the word "Off" in this screen but will not be programmable.

- 3. Turn the Knob to highlight the Time setting (00:00 in the sample screen above), then click the Knob in. The field will flash.
- 4. Turn the Knob until the desired running **duration** for this step (from 00:01, which means one minute, to 99:59) appears, then click the Knob in to save the setting. We will set Step 1 time to eight hours for this example (*see Figure 24*).
- 5. Turn the Knob to highlight the °C temperature setting (20.0 in the sample screen above), then click the Knob in. The field will flash.
- 6. To set the temperature desired (°C from 4.0 to 80.0) for the time period you have set, turn the Knob (left to decrease, right to increase). When the desired value appears, click the Knob in to save the setting. We will set Step 1 temperature to 37.0°C (see Figure 24).
- 7. Turn the Knob to highlight the **RPM** setting (**Off** in the sample screen above), then click the Knob in. The field will flash.
- 8. Turn the Knob to select the desired shaking speed (25 to 400 RPM) for this time period, then click the Knob in to save the setting. We will set Step 1 speed to 150 RPM (*see Figure 24*).
- 9. If you do not have the optional UV lamp, skip this step. If you have the UV lamp (default setting is Off) and you wish to have it on during the time period you have set, use the Control Knob to select the field, click the Knob to select the field (it will flash), turn the Knob to select On, then click the Knob in to save the setting. In our sample Step 1 (see Figure 24), however, UV remains Off.
- 10. If you do not have the optional **GRO** lamps, skip this step. If you have the photosynthetic lamps (default setting is **Off**) and you wish to have them on during the time period you have set, use the Control Knob to select the field, click the Knob to select the field (it will flash), turn the Knob to select **On**, then click the Knob in to save the setting. In our sample Step 1 (see Figure 24), however, **GRO** remains **Off**.

NOTE:

DO NOT select "Save" yet!

11. To program Step 2 (*see sample screen in Figure 24*): use the Control Knob to highlight the arrows at the bottom left of the screen. Click the Knob to make the arrows flash, then turn the Knob to the right until Step 2 appears. Click the Knob in to work in this screen, and repeat Steps 3-10.

PRG1 - STEP 01 02 🗸 Step number Time 08:00 03:00 $^{\circ}C$ 40.0 37.0 **RPM** 150 225 UV Off Off Off Off Save Cancel

Figure 24: Program 1 – Step 2

Whatever step is shown on the right-hand side of the display is settable (in Figure 24, it is Step 2). To scroll between steps, select the arrows (bottom left), then turn the Control Knob clockwise or counterclockwise, then select the desired step.

NOTE:

The time entered for each step is <u>for that step only</u>; it is not cumulative (i.e., the elapsed time from the start of the program).

12. Continue to program as many as 15 steps in the same manner. For our sample program, we have only three steps (*see Figure 25*).

PRG1 - STEP 02 03 🗸 Step number Time 12:00 03:00 $^{\circ}C$ 40.0 4.0 **RPM** 225 UV Off Off **GRO** Off Off Cancel Save

Figure 25: Program 1 – Step 3

The three-step program we have illustrated is designed to start incubating cultures at 37°C and hold that temperature for eight hours, shaking at 150 RPM. After eight hours, Step 2 engages, increasing the temperature setpoint to 40°C for a temperature induction, holding that temperature for three hours and increasing the shaking speed to 225 RPM. After this interval, the temperature will be cooled to 4°C and held at that temperature for twelve hours; no shaking will occur during this period, as the speed is set to 0 RPM.

Please note that temperature reductions, even at temperatures above ambient (going from 40° to 30°C, for example), require refrigeration to be effective.

NOTE:

When the last step of a program is completed, the Innova 44/44R will continue operating according to the parameters of the final step programmed, until the user intervenes.

If you wish the shaker contents to experience different conditions (e.g., a different shaking speed or no shaking at all, and/or return to ambient temperature), set those parameters as the last step in your program.

To save the entire program:

- 1. Once all your steps are set, use the Control Knob to highlight **Save** at the bottom of the screen, then click the Knob in. The field will flash.
- 2. Click the Knob in again to save the program. The display will read **Process Running Saving Profile** for a few seconds, then return to the main Programs

 (**PROG**) screen.

If you wish, you can set Programs 2, 3 and 4 and save them in the same manner.


7.5.4 Editing a Program

Use the Edit function to open a program you have already created and saved, in order to modify your settings following the same procedures.

7.5.5 Running a Program

Use the Run fuction to turn a specific program on. Naturally, only one program can run at a time. When you change the mode to Run, the screen will show the Run icon:

Figure 26: Programs Screen – Running

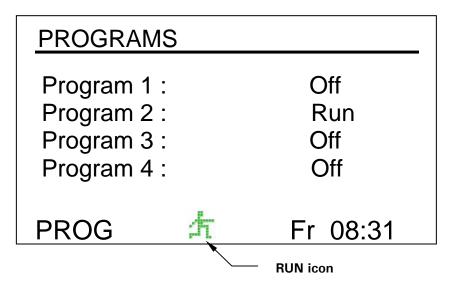

To stop a program: You can abort a run any time by changing the program's function to **Off**.

7.6 Muting the Audible Alarm

The Innova 44/44R shakers have an audible alarm that is activated under predetermined conditions (*see Table 4*). It may be muted in the following way:

- 1. Turn the Control Knob until the **SET** screen is highlighted on the display. Click the Knob in to work in this screen.
- 2. Turn the Knob to highlight the **Mute** mode (**On**), then click the Knob in. The field will flash.
- 3. Turn the Knob to change the setting to **Off**, then click the Knob to save this selection.

To reactivate the audible alarm at any time, repeat steps 1-3, reversing "off" and "on".

7.7 Temperature Offset Calibration

The temperature probe and the temperature controller are calibrated together at the factory. The temperature probe measures the temperature of the air at the probe's location, near the return vent. The controller uses the probe input to adjust air temperature, up or down, to match the temperature setpoint.

Depending on various conditions within the chamber, such as flask placement and size, the heat produced by growing organisms, heat losses due to liquid evaporation from flasks, etc., the display temperature may differ from temperatures within the flasks themselves. You can calculate the correction value for this offset and program the shaker to display a corrected temperature.

7.7.1 Calculating the Offset Value

If you wish to have the temperature display ("Indicated Temperature") match the temperature at a given point, or match the average of a series of points within the chamber ("Actual Temperature"), proceed as follows:

- 1. Let the unit equilibrate at or near the desired temperature, then record the Indicated Temperature.
- 2. Now record the Actual Temperature.
- 3. Calculate the temperature correction value using this formula:
 Actual Temperature Indicated Temperature = Temperature Offset Value.
- 4. To set the Temperature Calibration Offset, follow the procedure outlined in Section 7.7.2 below.

7.7.2 Setting the Offset


To set the temperature calibration **Offset**:

- 1. Use the Control Knob to enter the **CAL** screen (*see Figure 21, repeated below*).
- 2. Turn the Control Knob to highlight the current setting (**0.0** in the sample screen above), then click the Knob in. The current setting will flash.
- 3. Turn the Knob (left for negative settings or right for positive settings) to display the desired setting. As you turn the Knob, each click represents one tenth of a degree Celsius (0.1°C).
- 4. When you reach the desired setting, click the Control Knob in to save the new value.
- 5. When the Temperature Offset is set to any value other than 0, the icon will appear next to °C in the **DISP** and **SUMM** screens.

NOTE:

If you highlight an item, change it, but do not save your selection, after a few seconds the screen will revert to its previous setting.

Figure 21: Calibrate Screen

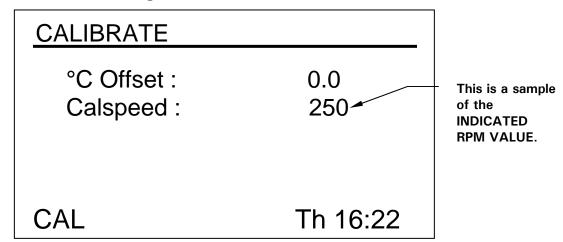

7.8 Using Calspeed

The **Calspeed** function, which is set in the **CAL** screen (*see Figure 21 above*), is used to calibrate the speed of the shaking mechanism. Calibrated at the factory, speed does not need to be recalibrated until a major operating component (e.g., drive belt) is changed.

NOTE:

Prior to engaging the Calspeed function, make sure that the platform is properly secured to the subplatform, and any flasks present are secured. Run the unit at approximately 250 rpm to ensure stable operation prior to calibration.

When the shaker is running, the **CAL** screen shows the indicated RPM value. If you wish to calibrate the speed, set the speed to a value that can be measured— a setpoint of 250 RPM works well. The use of a strobe is recommended for accuracy. If, after measuring the actual speed, you wish to adjust the indicated value:

- 1. Click the Control Knob in.
- 2. Set the new value.
- 3. Click the Knob in again to save the setting.
- 4. Turn the shaker **OFF**, wait a few seconds, then turn the shaker back **ON**.

7.9 Power Interruption

In the event of a power failure, the Innova 44/44R Shakers are equipped with an **automatic restart** function. The shaker's non-volatile memory retains all stored information.

If the shaker was in operation prior to the power interruption, the shaker will begin to operate at its last entered setpoints. The alarm/POWER display will flash, indicating that a power interruption has occurred. Turn the Control Knob in any direction to acknowledge the visual alarm. The flashing will stop.

7.10 Slide-Out Platform

The Innova 44/44R is equipped with a slide out platform mechanism as a standard feature. This allows you to easily slide the platform out of the shaker for easy access to the platform without the need for tools.

To operate the platform slide-out mechanism:

- 1. Make sure platform has stopped moving.
- 2. Rotate cam lever located on the subplatform (see Figure 27 on the following page) 180° clockwise. It should be pointing to the left.
- 3. Pull both side levers (*see Figure 27*) down. This will disengage the platform from the subplatform.
- 4. Pull **both** handles on the slide-out mechanism (*see Figure 27 on the following page*) toward you. This will move the platform to an easily accessible position.

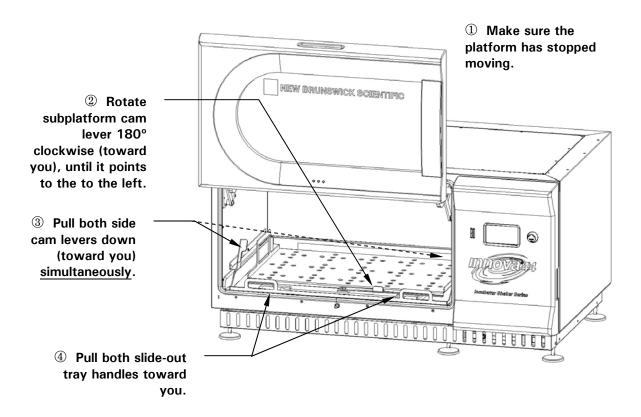


Figure 27: Disengaging Slide-Out Platform

5. To re-insert the platform, reverse this process. Make sure the two slots in the back of the platform are inserted under blocks of subplatform in the back.

NOTE:

Make sure that you have rotated the front cam lever counter-clockwise to lock the platform in position.

8 MAINTENANCE



WARNING!

When cleaning the unit, always turn off the shaker and disconnect the power cord from the power supply.

8.1 Routine Maintenance

No routine maintenance schedule is required for the Innova 44 and 44R.

To ensure that your shaker retains its attractive appearance, an occasional cleaning, using a cloth with conventional household (non-abrasive) cleaner is recommended (*see Section 8.2 below for more details*).

We also suggest that the area around the shaker be vacuumed or swept to remove dust and other debris, ensuring proper air flow in and around the shaker.

8.2 Cleaning External & Internal Surfaces

The unit may be cleaned using a damp cloth or any standard household or laboratory cleaner to wipe down its outer surfaces. Do not use abrasive or corrosive compounds to clean this instrument, as they may damage the unit and void the warranty.

If Biohazard decontamination is required, see Section 8.3 below.

8.3 Biohazard Decontamination



It is the responsibility of the user to carry out appropriate decontamination procedures if hazardous material is spilled on or inside the equipment. Before using any cleaning or decontamination method other than those suggested by the manufacturer, users should check with New Brunswick Scientific that the

proposed method would not damage the equipment.

Commercially available household bleach solutions, when diluted at a 1:10 ratio, are effective in routine decontamination of the instrument. The method for decontaminating a spill depends upon the nature of the spill.

Spills involving fresh cultures or samples known to have low concentrations of biomass should be flooded with decontamination solution and soaked for 5 min before cleanup.

Spills involving samples with high concentrations of biomass, or involving organic matter, or occurring in areas warmer than room ambient temperature, should be exposed to decontamination solution for at least one hour before cleanup.



MARNING!

Personnel involved in the cleanup of any spill should wear gloves, safety glasses, and a laboratory coat or gown during the cleanup process. Respiratory protection should be considered for spills where aerosolization is suspected.

9 Service & Accessories

9.1 Troubleshooting

If any problems occur with your shaker, do not attempt to perform any service on the unit other than specified in this manual. Unauthorized servicing may void the warranty. Please contact your local NBS Customer Service Department

In any correspondence with NBS, please refer to the model number and serial number of your unit. This information is on the electrical specification plate, located on the rear panel of the unit, above the power connector. The serial number is also located in the lower right corner panel, below the chamber door seal.

There are some problems, however, that you can investigate and correct yourself. Refer to the following Troubleshooting Guide:

Symptom(s)	Probable Cause(s) & Solution(s)	
	No power; display is not on; power cord is not plugged in and/or	
	power switch is off: plug in power cord (to working electric	
	outlet), and turn on power switch.	
	Door is open—look for Open Door icon on display: close door	
	firmly, making sure latch is engaged.	
	Door is closed but Open Door icon is on display: door magnet	
	is not adjusted correctly; call for service.	
	On/Off switch is not working: call for service.	
	Tilt Switch has been triggered. Alarm condition <i>Tilt</i> will be	
Shaker does not run.	indicated and RPM will be OFF. Check to ensure the load is	
	evenly distributed and that the shaker is level on a solid surface.	
	Restart by pressing the START/STOP button.	
	Fuse(s) burned out: check and replace as needed.	
	If you recently opened the Service Drawer to replace a fuse, it	
	may not have been seated properly: reopen and reclose	
	carefully.	
	Shaking speed has been set to Zero by program running (look	
	for Run icon in display) or by computer interface: reset shaking	
	speed.	
	Defective main board: call for service	
	Defective display controller board: call for service.	
	Jammed shaking mechanism: call for service	
	Defective motor: call for service	
	Drive belt out of alignment or worn: call for service.	

...continued...

Symptom(s)	Probable Cause(s) & Solution(s)
	If you recently opened the Service Drawer to replace a fuse, it
	may not have been seated properly: reopen and reclose
	carefully.
Shaker runs slowly	Fuse(s) burned out: check and replace as needed.
and/or no speed	Incorrect speed calibration: recalibrate shaking speed.
indication.	Defective main board: call for service.
	Defective motor: call for service.
	Drive belt out of alignment or worn: call for service.
	Shaker is running in Program mode (look for Run icon in display).
	Shaker speed has been changed by RS-232 command/
Shaker does not run at	computer interface.
set speed.	Shaker is overloaded: remove some contents & balance load.
	Defective motor: call for service.
	Drive belt out of alignment or worn: call for service.
	Load out of balance: pull out Slide-Out Platform, unload all
Operating noise	contents, reload.
	Loose component(s) in Slide-Out Platform, subplatform, and/or
	drive assembly: call for service.
	Shaker running in Program mode: Run icon in display.
	Temperature setpoint changed by RS-232/computer interface
	command.
	Heater fuse blown: replace.
	Compressor fuse blown: replace.
	Compressor over-pressure switch activated: call for service.
	Chamber fans not working: call for service
	Ambient temperature too high or too low: cool or heat the room as needed.
Incubator does not	Door is not completely closed (even though Open Door icon
reach set temperature.	may not be on display): open it and reclose it firmly.
	Sealing gasket between the Service Drawer and the Inner
	Cabinet is damaged: call for service to replace it.
	Line voltage is too low.
	Frequency on line voltage if set incorrectly: reset.
	Firmware revision does not match the Service Drawer revision:
	ask for update.
	Defective heater: call for service.
	Defective refrigeration system: call for service.
	Incorrect temperature indication (see below).
	Temperature Offset has been programmed (look for * adjacent
Incorrect temperature	to C on display).
indication.	Defective RTD assembly: call for service.
	Defective main board: call for service.

...continued...

Symptom(s)	Probable Cause(s) & Solution(s)
Incorrect relative humidity %	Defective humidity probe or broken contact: call for service.
Photosynthetic GRO lamp does not work.	Shaker is running in Program mode (look for Run icon in display), and program calls for GRO lamp to be off. GRO lamp mode has been changed by RS-232 command/computer interface. Bulb is burned out; replace (see Section 9.6). Fuse is burned out: replace (see Section 9.4) Check ballast voltage. Check cable connections. Check switch settings. GRO lamp mode has been changed by RS-232 command/computer interface.
UV germicidal lamp does not work.	Shaker is running in Program mode (look for Run icon in display), and program calls for UV lamp to be off. UV lamp is burned out; replace (see Section 9.8) Fuse is burned out: replace (see Section 9.4) Check ballast voltage. Check cable connections. Check switch settings. UV lamp mode has been changed by RS-232 command/computer interface.
Door does not stay in the fully open position.	Adjust the door tension (see Section 9.7).

9.2 Product Returns

Should you need to return your Innova 44/44R to NBS for any reason, first contact Customer Service to obtain a Returned Material Authorization (RMA) number. This number must appear on the outside of the shipping container, otherwise NBS Receiving will refuse to accept the shipment.

In addition, you must also certify that the instrument being returned has been thoroughly cleaned and decontaminated. A form for this purpose is provided in Section 13 of this manual; it can also be downloaded from our website (www.nbsc.com). A copy of this completed Return Authorization and Decontamination Certificate must be attached to the outside of the container, with a second copy packed inside with the instrument.

9.3 Opening the Service Drawer



MARNING!

Before opening the Service Drawer, always turn off the shaker and disconnect the power cord from the power supply.

The Service Drawer (see Figure 28) contains the shaker's electronic and temperature control components. It is designed to open partially (allowing access to the fuses—see Section 9.4 for Fuse Replacement instructions) or to be removed, by authorized service technicians only, from the shaker entirely, to allow easier accessibility for service.



CRUSH WARNING!

Only authorized technicians should open the Service Drawer beyond its stops, giving it proper support, or remove it entirely, which requires two people or one technician with an appropriate lifting device. The drawer weighs approximately 125 lbs/57kg.

9.4 Fuse Replacement



WARNING!

When replacing fuses, always turn off the shaker and disconnect the power cord from the power supply.

There are four fuses that can be replaced by the user. These are located inside the Service Drawer. To access this forward portion of the Service Drawer (see Figure 28):

- 1. Turn the power off and unplug the shaker.
- 2. Remove the four thumbscrews and washers that hold the front grille in place. Set them aside for reuse.
- 3. Locate and, using a screwdriver, unlock the handles that are located at the base of the Service Drawer (on the front face).
- 4. Release the handles by pushing the rectangular buttons.
- 5. Grasp the handles and pull the Service Drawer out **until it engages the stop** mechanism.



CRUSH WARNING!

Do not try to pull the drawer out beyond its stops.


- 6. Fuses are numbered 1 through 4 (*see Table 6*); access the fuse you wish to replace by using a coin or a blade screwdriver to turn and release the spring-loaded cap.
- 7. Replace the fuse with a new one of the same type and rating (see Table 6 below):

Table 6: Fuses

Fuse No.	Function	Type	Rating
F1	Heater		6.3 A
F2			3.15 A
F3			3.15 A
F4	Compressor		10 A

□ optional equipment

Figure 28: Fuse Replacement

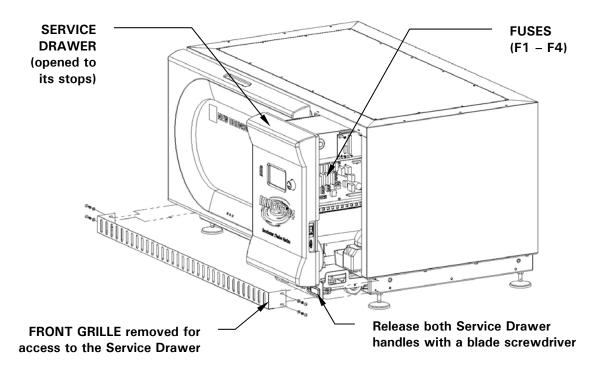

9.5 Replacing Cabinet Lights



WARNING!

Before installing or removing the quartz bulbs:

- (1) Make sure power is disconnected from the unit,
- (2) Wear plastic gloves to protect the bulb from skin oils and to protect your hands just in case the bulb should break,
- (3) Before removal, make sure the bulb has fully cooled.



Cabinet lights are connected in series. When one burns out, both lights will be off. It is therefore prudent to replace both bulbs at the same time.

Be sure to follow all the **WARNINGS** and the **NOTE** above as well as the procedure outlined below when you replace the cabinet bulbs with reference to Figure 29 below:

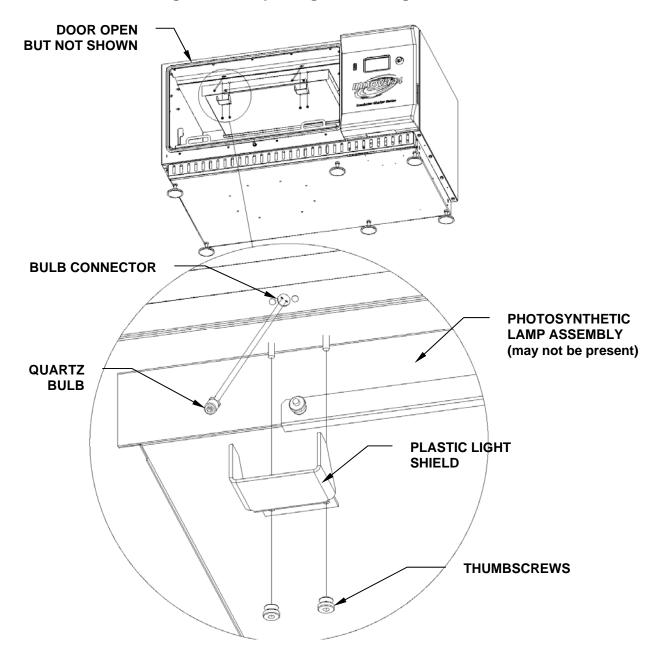


Figure 29: Replacing Cabinet Lights

- 1. Make sure power is disconnected from the shaker.
- 2. Open the door.
- 3. Unscrew the thumbscrews and remove (*with reference to Figure 29*) both plastic light shields to gain access to the bulb connectors.

- 4. Wearing plastic gloves, remove the first used bulb, then carefully insert the new bulb's prongs fully into the connector base.
- 5. Repeat step 4 for the second bulb.
- 6. Reinstall the plastic light shields with their thumbscrews.

9.6 Replacing Optional Photosynthetic Lamps



WARNING!

Before removing or installing the photosynthetic "GRO" lamp, make sure power is disconnected from the unit.

CABINET LIGHT
(see Figure 29)

PHOTOSYNTHETIC
LAMP ASSEMBLY

DIFFUSER HOLDER
THUMBSCREW

LIGHT SHIELD
THUMBSCREWS

Figure 30: Replacing Photosynthetic Lamps

To replace the photosynthetic GRO lamp bulbs, with reference to Figure 30 above:


- 1. Make sure power is disconnected from the shaker.
- 2. Open the door.
- 3. Remove the two thumbscrews from the diffuser holder.
- 4. Remove the diffuser by placing a finger underneath and rotating the diffuser in the slot at the back of the photosynthetic light assembly. All 6 lamps should now be accessible to you.

- 5. To remove the lamp for replacement, rotate the lamp 90°, taking care to support the lamp before it drops down.
- 6. To replace the lamp, orient the pins of the new lamp vertically. Simultaneously, gently press both sides of the lamp into the slots.
- 7. Rotate the new lamp 90°. The lamp should catch in the holders on both ends. Make sure both ends are secure.
- 8. Put the diffuser back in place, reinstalling the diffuser holder with its two thumbscrews.

9.7 Adjusting Door Tension

In the event the door does not stay in the fully open position, you may want to adjust the door tension. Housed in the counterbalance brackets, located on the side walls of the chamber, is an adjustment mechanism that consists of a slider controlled by the adjuster bolt at the front and a jam nut and lock washer at the rear (*see Figure 31 below*).

Figure 31: Righthand Counterbalance Bracket

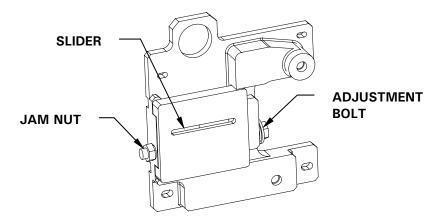

To adjust door tension, you will need two 9/16-inch open-ended wrenches and two wood blocks, each approximately 3/4 inch x 11/2 inches x 6 inches, to protect the door gasket from the upper door arms during this procedure:

- 1. Fully open the door, placing each wood block between an upper arm and the door gasket.
- 2. Using the wrenches, loosen the jam nut and lock washer from the back of both counterbalance brackets.
- 3. Turn the adjuster bolt on the front of each counterbalance bracket to force the upper door arms up. Adjust the placement of the wood blocks as needed to bring the door to the desired height without allowing the upper arms to harm the door gasket.

- 4. Adjust both sides to the same tension by comparing the position of the slider, which you can see through the slot in the black cover on each bracket.
- 5. Retighten the lock washer and jam nut on each side; tighten to lock the adjuster bolt in place.



WARNING!

The following procedures are provided <u>for your information only</u>. Do not attempt to perform these service interventions yourself unless you are an authorized service technician.

9.8 Removing the Service Drawer

When authorized service technicians need to remove the Service Drawer, they will follow this procedure:



CRUSH WARNING!

The Service Drawer is very heavy (approximately 125 lbs/57kg), so there must be at least two people to handle it, or one technician with an appropriate lifting device.

- 1. Turn the power off and unplug the shaker.
- 2. Remove the plastic fitting that positions the temperature probe inside the chamber. This is located on the right side of the inside wall.
- 3. Use a piece of tape to hold the probe flat against the inside wall.
- 4. Remove the four thumbscrews that hold the front grille in place.
- 5. Locate and, using a screwdriver, unlock the handles that are located at the base of the Service Drawer (on the front face).
- 6. Release the handles by pushing the rectangular buttons.
- 7. Grasp the handles and pull the Service Drawer out **until it engages the stop mechanism**.
- 8. Push the Service Drawer back in slightly (approximately one inch).
- 9. Push down on the Service Drawer (**not** on the plastic bezel) while pulling the drawer out again; this action will bypass the stop mechanism.
- 10. The Service Drawer is supported for *only two inches* after clearing the stops.



CRUSH WARNING!

Do not pull the Service Drawer out entirely until you are prepared to handle its weight (approximately 125 lbs/57kg).

11. Remove the Service Drawer and place it on a stable, flat surface.

When reinstalling the Service Drawer, the authorized technician will be sure to put the temperature probe back in place with its plastic fitting, keeping the probe tip 3 - 3.25inches (7.6 - 8.2 cm) from the front face of the fan cover.

9,9 Replacing Optional Germicidal Lamp



WARNING!

Before installing or removing the germicidal "UV" lamp, the authorized

- (1) Make sure power is completely disconnected from the unit,
- (2) Wear plastic gloves to protect the lamp from skin oils and to protect your hands just in case the lamp should break,
- (3) Before removal, make sure the lamp has fully cooled.

To replace the UV lamp, your service technician will follow all the **WARNINGS** above and the procedure and **WARNINGS** outlined below, with reference to Figures 32 & 33:

- 1. Make sure power is disconnected from the shaker.
- 2. Remove the Service Drawer, carefully following the instructions in Section 9.3.



CRUSH WARNING!

The Service Drawer is very heavy (approximately 125 lbs/57kg), so there must be at least two people to handle it, or one technician with an appropriate lifting device.

Do not pull the Service Drawer out entirely until you are prepared to handle its weight.

3. With reference to Figure 32 below, free the heater cable harness and the temperature probe from their tie wraps and Permagum[®].

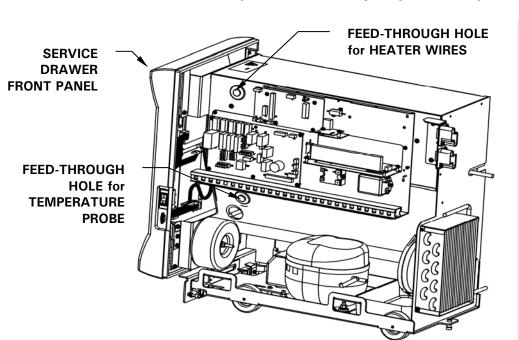


Figure 32: Service Drawer Feed-Through Holes (Germicidal Lamp Replacement)

NOTE that wires, tie wraps and Permagum® are not shown above.

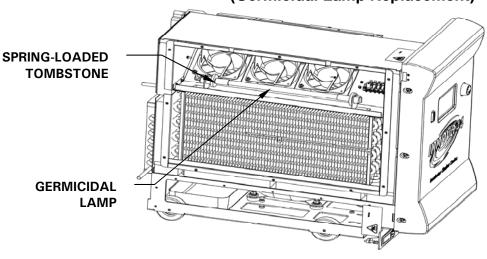


Figure 33: Service Drawer, Fan Cover Removed (Germicidal Lamp Replacement)

NOTE that the heater wire is not shown.

4. Remove the fan cover from the Service Drawer, paying careful attention to the heater cable harness and the temperature probe, both of which are attached to the fan cover. If they are freed adequately, the fan cover can be rotated and rested on top of the Service Drawer.

- 5. With reference to Figure 33, and wearing plastic gloves, remove the germicidal lamp by pressing it toward the spring-loaded tombstone. The lamp has a single pin at each end.
- 6. Still wearing plastic gloves, install a new lamp.
- 7. Pull the heater cable harness and temperature probe back through the holes in the Service Drawer (*Figure 32*).
- 8. Replace the fan cover, making sure the heater cable harness wires and temperature probe have been pulled through. Adjust the temperature probe position to ensure that the probe tip can protrude 3 3.25 inches (7.6 8.2 cm) from the fan cover face.
- 9. Secure the heater cable harness and temperature probe with tie wraps, and add Permagum[®] to fill the hole.
- 10. Reinstall the Service Drawer in the unit, following all above **WARNINGS**.
- 11. With the Service Drawer secured in the unit, use the plastic fitting to hold the temperature probe, ensuring that its tip is 3 3.25 inches (7.6 8.2 cm) from the fan cover face.

9.10 Belt Replacement or Adjustment

To gain access to the drive belt, your service technician will follow these steps with reference to Figures 34 & 35 below:

- 1. Make sure that the power is switched off and the shaker is unplugged.
- 2. Remove the platform.
- 3. Using a 7/32-inch Allen key (hex wrench), remove the four bolts that hold the subplatform in place, and remove the subplatform.
- 4. Disconnect the tachometer and motor wires by unplugging their connectors.
- 5. With a wrench, remove the nuts that hold the bearing housing to the bottom weldment.
- 6. Position the counterweight so it is closest to you, and tie it in place with tie wraps or strong cord, in order to prevent it from dropping against your fingers or causing you to lose your grip.
- 7. Remove the bearing housing with care: it weighs 47 lbs. / 21.3 kg.
- 8. Turn the bearing housing over to access the belt, motor and pulley.

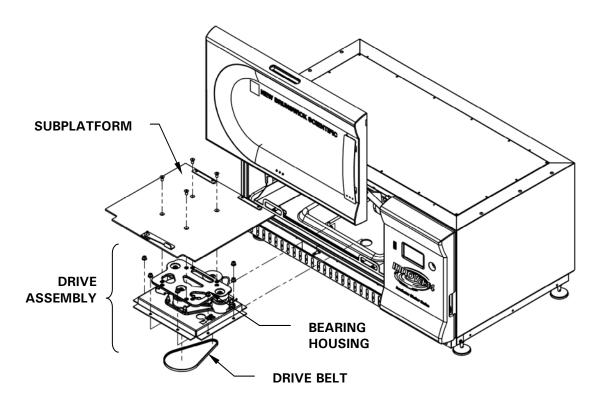


Figure 34: Drive Assembly

To replace the drive belt, follow these steps with reference to Figure 35:

- 1. Loosen the four hex nuts that hold the motor assembly in place and remove the old belt.
- 2. Guide the belt onto the two pulley grooves.
- 3. Tighten the motor nuts.
- 4. Check the belt adjustment by applying finger pressure to the belt midway between the two pulleys. The belt should deflect approximately 3/8-inch (9-10 mm).
- 5. If the deflection is not right, loosen the nuts, readjust the belt, then retighten the motor nuts.

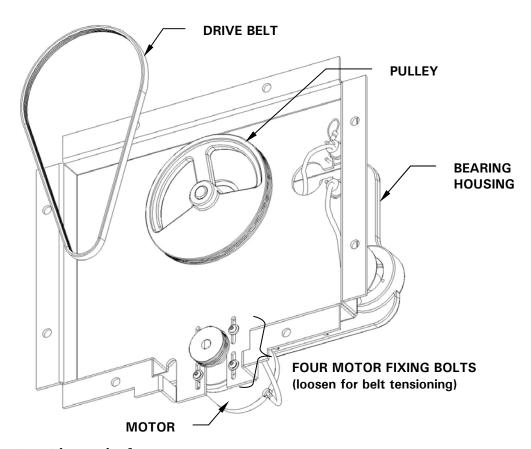


Figure 35: Belt Replacement & Adjustment

If an adjustment is required:

- 1. Loosen the four hex nuts on the motor assembly.
- 2. Move the motor assembly until the belt is tight.
- 3. Tighten the hex nuts and recheck the belt tension by exerting pressure on the belt. The belt should deflect approximately 3/8-inch (9.5 mm).

Reassembly:

- 1. In reverse order, reinstall the bearing housing. Once it is securely in place, remove the tie wraps/cord.
- 2. Secure the wiring to ensure no interference during rotation.
- 3. Reinstall the subplatform and platform.
- 4. Recalibrate speed (see Section 7.8).

9.11 Replacement Parts

In Table 7 below you will find commonly needed replacement parts for your Innova 44/44R shaker:

Table 7: Replacement Parts

Part Description	Quantity	NBS Part No.
Spare Parts Kit (includes 1 motor drive poly-V belt, and two each 10A, 6.3A & 3.15A fuses)	1	M1282-6000
Flask Clamp screws, 10-24 x 5/16-inch flat Phillips	25	S2116-3051
head (+)		
Test Tube Rack screws, 10-24 x ½-inch pan head,	10 each	M1289-0800
and washers (#10 plain and #10 lock washers)		
Motor Drive Poly-V Belt, 28 inches long	1	P0700-7070
Fuse, 10A	1	P0380-3162
Fuse, 6.3A	1	P0380-3630
Fuse, 3.15A	1	P0380-3460
Chamber Light Bulbs	2	P0300-0470
UV Germicidal Bulb (if so equipped)	1	P0300-0351
Photosynthetic Bulb (if so equipped)	1	P0300-0221
UV Germicidal Light Ballast (if so equipped)	1	P0420-5113
Photosynthetic Light Ballast (if so equipped)	1	P0420-5192

9.12 Accessories

When ordering accessories, you may be asked to provide the model number and serial number of your shaker. This information is on the electrical specification plate, located on the rear panel of the unit. The serial number is also labeled in the lower right corner of the front panel, below the chamber door seal.

9.12.1 Slide-Out Platforms

The Slide-Out Platform allows you to extend the platform in front of the shaker for easy access to all flasks. Here are the various platforms available:

Table 8: Available Slide-Out Platforms

Description	NBS Part No.
Universal Platform	M1282-9904
125 mL Erlenmeyer flask Dedicated Platform	M1282-9905
250 mL Erlenmeyer flask Dedicated Platform	M1282-9906
500 mL Erlenmeyer flask Dedicated Platform	M1282-9907
1L Erlenmeyer flask Dedicated Platform	M1282-9908
2L Erlenmeyer flask Dedicated Platform	M1282-9909
2.8L Fernbach flask Dedicated Platform	M1282-9910
4L Erlenmeyer flask Dedicated Platform	M1282-9911

Should you decide in favor of the Universal Platform, following is a list of that particular platform's flask capacity, according to flask size:

Table 9: Universal Platform Flask Capacities

Flask Type	Capacity*
50 ml Erlenmeyer Flasks	91
125 ml Erlenmeyer Flasks	39
250 ml Erlenmeyer Flasks	30
500 ml Erlenmeyer Flasks	24
1 L Erlenmeyer Flasks	14
2 L Erlenmeyer Flasks	8
2.8 L Erlenmeyer Flasks	6
4 L Erlenmeyer Flasks	6
5 L Erlenmeyer Flasks	6

^{*}Clamps for Universal Platform are sold separately.

9.12.2 Flask Clamps for Universal Platforms

The following clamps, according to flask size, are available for use with the Universal Platform:

Table 10: Flask Clamps

Clamp Size	Part Number
10 ml Erlenmeyer Flask	ACE-10S
25 ml Erlenmeyer Flask	M1190-9004
50 ml Erlenmeyer Flask	M1190-9000
125 ml Erlenmeyer Flask	M1190-9001
250 ml Erlenmeyer Flask	M1190-9002
500 ml Erlenmeyer Flask	M1190-9003
1 L Erlenmeyer Flask	ACE-1000S
2 L Erlenmeyer Flask	ACE-2000S
2.8 L Fernbach Flask	ACSB-2800S
4 L Erlenmeyer Flask	ACE-4000S
5 L Erlenmeyer Flask	ACE-5000S

9.12.3 Replacement Clamp Hardware Kits

NBS flask clamps come complete with mounting screws. Additional screws are available separately in packs of 25 (see Section 9.11).

9.12.4 Test Tube Racks & Other Accessories

Table 11: Racks & Trays

Accessory L	Description	NBS Part Number	Platform Capacity
Adjustable angle Test Tube	80 tube capacity	M1289-0100	7
Rack for tubes 8 – 11 mm	60 tube capacity	M1289-0010	9
diameter	48 tube capacity	M1289-0001	9
Adjustable angle Test Tube	60 tube capacity	M1289-0200	7
Rack for tubes 12 - 15 mm	44 tube capacity	M1289-0020	9
diameter	34 tube capacity	M1289-0002	9
Adjustable angle Test Tube	42 tube capacity	M1289-0300	7
Rack for tubes 15 –18 mm	31 tube capacity	M1289-0030	9
diameter	24 tube capacity	M1289-0003	9
Adjustable angle Test Tube	30 tube capacity	M1289-0400	7
Rack for tubes 18 – 21 mm	23 tube capacity	M1289-0040	9
diameter	18 tube capacity	M1289-0004	9
Adjustable angle Test Tube	22 tube capacity	M1289-0500	7
Rack for tubes 22 – 26 mm	16 tube capacity	M1289-0050	9
diameter	13 tube capacity	M1289-0005	9
Adjustable angle Test Tube	20 tube capacity	M1289-0600	7
Rack for tubes 26 - 30 mm	16 tube capacity	M1289-0060	9
diameter	12 tube capacity	M1289-0006	9
Microplate holder rack (stacked)	3 deep well or 9 standard	M1289-0700	16
Microplate holder rack (single layer)	5 deep well or standard	TTR-221	4
Angled Test Tube Rack Holder racks that are 4-5 in. (10-13 mr (38 mm) long.	m) wide and up to 15 in.	TTR-210	4
Angled Test Tube Rack Space to accommodate test tubes rac (13 mm) wide.		TTR-215	NA

^{*} Universal Platform Required

9.12.5 Optional Bases

Table 12: Bases

Description	NBS Part Number
Short Base, raises bottom unit 4 inches (10 cm)	M1282-0600
Medium Base, raises bottom unit 12 inches (30.5 cm)	M1282-0800
Tall Base, raises bottom unit 16 inches (40.6 cm)	M1282-0700

9.12.6 Optional Stacking Kit

Table 13: Stacking Kit

Description	NBS Part Number
Stacking hardware (to stack two Innova 44 shakers)	M1282-0500
Stacking hardware (to stack an Innova 4400 on an Innova 44)	M1282-0501

9.12.7 Optional Gassing Port Kit

Table 14: Gassing Port Kit (Factory-Installed)

Description	NBS Part Number
Gassing Port Kit	M1282-0291

9.12.8 Optional Photosynthetic Lamps

Table 15: Photosynthetic Growth Lamps (Factory-Installed)

Description	NBS Part Number
Photosynthetic Growth Lamp Assembly (9 lamps)	M1282-9922

9.12.9 Optional UV Germicidal Lamp

Table 16: Germicidal Lamp (Factory-Installed)

Description	NBS Part Number
UV Germicidal Lamp Assembly	M1282-9921

9.12.10 Optional Humidity Monitor

Table 17: Humidity Monitor (Factory-Installed)

Description	NBS Part Number
Humidity Monitor (for use with spill pan as water	M1282-5011
reservoir)	

9.12.11 Optional Handles

Table 18: Screw-In Handles

Description	NBS Part Number
Screw-In Handles (6 pieces in kit)	M1282-5042

There are three threaded holes in the front of the shaker behind the front grille, two on each side panel, and three on the rear panel. Depending on your need, you can use any of these holes for the screw-in handles (*see Figure 36 below*). Be sure to tighten the handles securely.



CRUSH WARNING!

Make sure that any persons who will use these handles (for <u>minor</u> moves or lifts <u>only</u>) take proper precautions to protect their hands, feet and backs from injury.

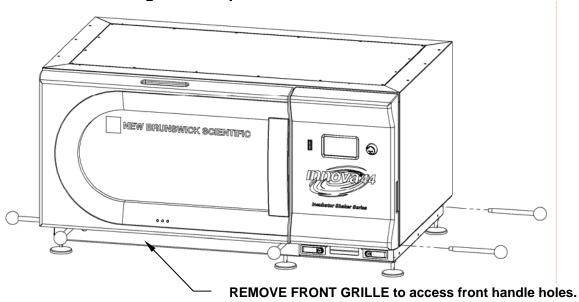


Figure 36: Optional Screw-In Handles

9.12.12 Optional Remote Alarm

The Innova 44/44R can be equipped with a factory-installed remote alarm component (NBS part number M1320-8029). When it is hooked up to your relay and receiving equipment, this device will send notification of an alarm condition to the remote location you choose.

10 Drawings & Tables

10.1 Schematics

BANK Mrzszeso

Figure 37: Control Schematics

NOTE:

For a larger version of this drawing, contact your sales or service representative.

EVAP. FAN(S) HEATER HE-1 TE A1 TEMP CONTROL HX-2PROBE EVAPORATOR FLOW CONTROL DEVICE PRESSURE SWITCH COND. FAN(S) PAH PSH 1A 1A, CONDENSER HX-1FILTER COMP. DRYER

Figure 38: Refrigeration Schematic

10.2 List of Drawings

Figure	Description	Page
1	Innova 44/44R Front View	15
2	Installing the Front Grille	17
2a	Adjustable Foot	17
3	Mounting Shaker on Base	19
4	Shaker Installed on Base	19
5a	Installing Innova 44/44R Stacking Kit	20
5b	Installing Innova 4400/4430 to 44/44R Stacking Kit	21
6	Stacked Shakers	22
7	Front Panel (Detail)	25
8	Control Panel (Right Side)	25
9	Display Screen	26
10a	Installing Subplatform	32
10b	Installing Slide-Out Platform	33
11a	Double Girdle Clamp Installation	34
11b	Clamp Fastener	34
12	Rear Panel	35
13	Front Grille	37
14	Changing Display Parameter	41
15	Changed Display Parameter	41
16	Changing Setpoint	42
17	Summary Screen	43
18	Setup Screen	45
19	Lamps Screen	47
20	RS232 Screen	48
21	Calibrate Screen	49, 57
22	Programs Screen	50, 51
23	Program 1– Step 1	51
24	Program 1– Step 2	53
25	Program 1– Step 3	54
26	Programs Screen – Running	55
27	Disengaging Slide-Out Platform	59
28	Fuse Replacement	66
29	Replacing Cabinet Lights	67
30	Replacing Photosynthetic Lamps	68
31	Righthand Counterbalance Bracket	69
32	Service Drawer Feed-Through Holes (Germicidal Lamp	72
	Replacement)	
33	Service Drawer, Fan Cover Removed (Germicidal Lamp	72
	Replacement)	
34	Drive Assembly	74
35	Belt Replacement & Adjustment	75
36	Optional Screw-In Handles	80
37	Control Schematic	81
38	Refrigeration Schematic	82

10.3 List of Tables

Table	Description	Page
1	Screens	27
2	Display Screen Parameters	27
3	Display Icons	28
4	Alarms	29
5	Communication Mode	48
6	Fuses	66
7	Replacement Parts	76
8	Available Slide-Out Platforms	77
9	Universal Platform Flask Capacities	77
10	Flask Clamps	78
11	Racks & Trays	78
12	Bases	79
13	Stacking Kit	79
14	Gassing Port Kit	79
15	Photosynthetic Growth Lamps	79
16	UV Lamp	79
17	Humidity Monitor	79
18	Screw-In Handles	80
19	COM Parameters for HyperLink	88

SPECIFICATIONS

These specifications assume a maximum load of 34 pounds (15.5 kg), including platforms, clamps, glassware and contents.

Innova 44 & 44R Incubator Shaker				
SHAKING				
Speed ¹	25-400 rpm with 1-inch (2	2.5 cm) st	roke or 25-300 r	pm with 2-inch (5.1 cm)
•	stroke; 25-250 rpm for 3	units stac	ked	
Control Accuracy	± 1 rpm			
Indication	Displayed in 1 rpm increr	nents		
Stroke/Orbit	1-inch (2.5 cm) or 2-inch			
TEMPERATURE	,	,		
Range (44)	5°C above ambient temp	erature to	80°C ²	
Range (44R)	20°C below ambient tem	perature (minimum 4°C) t	o 80°C²
Control Accuracy	± 0.1°C from 30-40 °C, ±			
Indication	Displayed in 0.1°C increr		o	
Heater	Long-life, low-watt densit		ce-type heater v	vith high temperature
	thermostat	,		
ALARMS	Visible and audible warni	ing indica	tion when speed	deviates more than 5
7.2	rpm or temperature more			
	expired. Audible alarm ca			
LCD DISPLAY	240 x 128 backlit			
RS-232	Remote control, remote monitoring, remote data logging		aaina	
SETPOINT RETENTION	All setpoints and operating status are retained in non-volatile memory			
AUTOMATIC RESTART	Automatic restart after po			
STACKING	Up to three units may be stacked; the second and third units require			
	stacking kit(s).			
DRIVE		alanced o	drive with nine pe	ermanently lubricated ball
	bearings		•	
DRIVE MOTOR	Solid-state brushless DC motor.			
SAFETY	Drive Interrupt shuts off p			
	Acceleration/deceleration	n circuit p	revents sudden s	starts and stops,
	minimizing both splashing			
	mechanical sensing swite	ch also sh	outs the motor of	fin unbalanced condition.
ELECTRICAL	■ 100V, 50Hz			
REQUIREMENTS	■ 100V, 60 Hz			00 VA per shaker
	 120 Volts, 60 Hz 		44R: 1	500 VA per shaker
	 230 Volts, 50 Hz 			
ETL REGULATORY	UL 61010A-1		CAN/CSA-C22.2 No. 1010.1	
STANDARDS	UL 61010A-2-010		CAN/CSA-C2	22.2 No. 1010.2.010
CE REGULATORY	IEC 61010-1: 2001	EN 61010-2	2-081 (2002)	EN 61000-3-2 : 1995
STANDARDS			1997 + A1 : 1998	
			2-010 (1992) + A1: 1	
DIMENSIONS ³	Single Unit		wo Units	Three Units
Width	53 inches (135 cm)		ches (135 cm)	53 inches (135 cm)
Depth (Front to Back)	33 inches (84 cm)		iches (84 cm)	33 inches (84 cm)
Height	26 inches (66 cm)	52 inc	hes (132 cm)	78 inches (198 cm)

¹ see first NOTE on following page
2 depending on ambient factors, such as relative humidity, and options installed
3 see second NOTE on following page

Innova 44 & 44R Incubator Shaker		
CHAMBER	26 inches deep X 36 inches wide X 18 inches high from top of platform	
DIMENSIONS	(66 cm deep X 91 cm wide X 46 cm high from top of platform)	
CHAMBER DOOR	Door glides out and up.	
PLATFORM	Aluminum, 18" X 30" (46 X 76 cm). Select universal or dedicated styles	
	PLEASE NOTE: Innova 44 & 44R platforms are NOT interchangeable with	
	platforms from other NBS shakers.	
GROSS WEIGHT		
44	525 lbs (238 kg)	
44R	570 lbs (259 kg)	
Short Base	36 lbs (16.4 kg)	
Medium Base	90 lbs (40.9 kg)	
Tall Base	126 lbs (57.3 kg)	

NOTE:

Use of baffled flasks will significantly reduce maximum speed for any shaker. We may be able to improve this maximum speed by using an alternative motor pulley. Contact your NBS representative for more information.


NOTE:

Add 21 in. (53 cm) to height and 6 in. (15.2 cm) to depth for door clearance.

11.1 Certifications

The Innova 44 and 44R have been tested to ETL standards, to comply with UL and CAN/CSA electrical safety standards (see "ETL Regulatory Standards" in the specifications table).

As attested in the *CE Declaration of Conformity* reproduced on the following page, the Innova 44 and 44R also conform to the appropriate CE standards (see also "CE Regulatory Standards" in the specifications table).

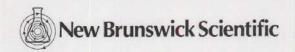

CE

CE

DECLARATION OF CONFORMITY

New Brunswick Scientific, Hereby declares that the product(s) listed below conform to the European Union directive and standards identified in this declaration.

Product(s)

Innova 44/44R

EU Directive(s)

Low Voltage (73/23/EEC/93/68/EEC) Electromagnetic Compatibility (89/336/EEC/93/68/EEC)

Standard(s)

IEC61010-1: 2001 (2nd Edition) EN61010-1: 2001 (2nd Edition) EN61010-2-010 (1992) +A1: 1996

EN61010-2-081 (2002) EN61326: 1997+A1: 1998 EN61000-3-2: 1995 EN61000-3-3: 1995

The conformity assessment procedures were performed at the following Testing Lab.:

Intertek Testing Services, 40 Commerce Way, Totowa, NJ 07512

The technical documentation relevant to the above equipment will be held at:

New Brunswick Scientific Company PO Box 4005 44 Talmadge Road Edison, New Jersey 08818-4005 U.S.A Tel. (732) 287-1200

Fax. (732) 287-4222

VP of Science & Technology

Date

12 APPENDIX: REMOTE PROGRAMMING

12.1 Setting Up with HyperTerminal

You can use the Microsoft[®] HyperTerminal program on your personal computer to communicate with your shaker:

- 1. Before you set up HyperTerminal, be sure to check the **Mode** and **Baud Rate** parameter settings (*see Section 7.4.5*).
- 2. While you are in the RS232 screen (*see Figure 20*), set the **Mode** parameter to **Slave**, and select the appropriate **Baud Rate**. The default setting, as you see in your screen, is 19200; your other choices are 9600 or 38400.
- 3. Connect the RS-232 cable between your computer and the shaker's RS-232 port (*see Figure 8*).
- 4. Start HyperTerminal, and follow the steps for Setting a New Connection. When prompted, select the COM port you will be using, then set the COM parameters (see Table 19 below):

Baud Rate	Enter the setting on your shaker (see Step 2 above)
Data Bits	8
Parity	None
Stop Bits	1
Flow Control	None

Table 19: COM Parameters for HyperTerminal

- 5. To be sure that you have established a connection with your shaker, turn the shaker **OFF**, open the HyperTerminal connection, then turn the shaker back **ON**. Once a connection has been made, the shaker will send an **OK** to HyperTerminal.
- 6. You may wish to send a trial command to the shaker (*see Sample Test below*). We suggest that you enable *Caps Lock* on your computer, because all commands must be sent as capital letters.

Sample Test:

- 1. To request the current setpoint, type **RS** at the terminal. When you type **R**, you should see the **R** on the computer display. This is an echo character, sent back by the Innova 44/44R. When you type the S, you should see the echo **S** displayed at the computer.
- 2. Use the ←Enter key to send the carriage return control character. This will move the cursor back over the first character.
- 3. **RS** is a Request command, so your shaker will send back the data corresponding to the command.

12.2 Overview of Command Sets

Using an optional RS-232 interface, the user has the capability to control various functions of the Innova 44/44R shakers by computer. A detailed command set follows this introduction.

The serial interface uses an echo character protocol that is implemented to indicate that a valid character has been received by the shaker. For every character sent by the user computer, the Innova 44/44R will return a matching character. This echo character should be read by the computer and compared to the character that was sent.

NOTE:

When you are sending a command, do not allow the pause between any two characters of your message to last longer than 10 seconds. If it takes any longer, the shaker's serial port will reset itself.

In the following example, the **Master** (your computer) will command (the **C** in **CS**) the shaker (**Slave**) to set its speed to **150 RPM**:

Master sends:	Innova 44/44R responds:
CS < <i>Space</i> > 150 < <i>CR</i> >	CS <space>150<cr></cr></space>

In the next example, the **Master** (your computer) will request that the shaker (**Slave**) report (the **R** in **RV**) the actual parameters:

Master sends:	Innova 44/44R responds:
RV <cr></cr>	First the echo: RV < <i>CR</i> >
	Then the report:
	<arg1><tab><arg2><tab><arg3><arg7><cr><lf></lf></cr></arg7></arg3></tab></arg2></tab></arg1>

NOTE:

On power-up, even if the shaker is not set in Slave mode at the time, the shaker will always send these four characters: OK < CR > < LF >. Before starting any communications, be sure to clear the computer's receive buffer.

This is the command set format:

<Command><Space><Arg1><Space><Argn...><Terminator>

Most Commands are defined as two capital characters using standard ASCII format (see the command sets provided below).

If a command has additional arguments, the ASCII space control character separates each argument. The command is terminated using the ASCII carriage return control character; an optional line feed control character may also be included but must follow the carriage return character.

When a command is sent that returns data, the characters returned are in ASCII format. Each additional argument returned will be separated by the ASCII space or tab control character; refer to the selected command for the actual format. Data returned is terminated with a carriage return line feed control character sequence.

Section 12.3 provides an Index to the 2-character command codes.

Sections 12.4-12.7 provide command sets for your use for Set Commands, Profile Control Commands, Report Request Commands and Set/Get Date & Time.

12.3 Index to Command Codes

Command Code	Meaning
CL	Command Grow Lamp
CS	Command Speed
CT	Command Temperature
CU	Command UV Lamp
PC	Clear Profile
PM	Profile Mode
PR	Read Profile Step in Memory
PS	Profile Stop/Start
PW	Write Profile Step to Memory
RI RP	Report Software Version Report Parameter List
RS	Report Setpoint Values
RV	Report Measured Values
_	
=D	Set Date and Time
?D	Get Date and Time

12.4 Set Commands

Code	Meaning	Format
CL	Command Grow	CL <space><arg1><cr></cr></arg1></space>
	Lamp	where Arg1 = Set Grow Lamp state (0 = off, 1 = on)
CS	Command Speed	CS <space><arg1><cr></cr></arg1></space>
		where Arg1 = Agitation setpoint (####)
CT	Command	CT <space><arg1><cr></cr></arg1></space>
	Temperature	where Arg1 = Temperature setpoint (###.#)
CU	Command UV	CU <space><arg1><cr></cr></arg1></space>
	Lamp	where Arg1 = Set UV Lamp state (0 = off, 1 = on)

12.5 Profile Control Commands

Code	Meaning	Format
PC	Clear Profile	PC <space><arg1><space><arg2><cr></cr></arg2></space></arg1></space>
	(clears profile	where Arg1 = Profile number (range 1 to 4) and
	step)	Arg2=Step number (range 1 to 15)
		NOTE: If Arg2 is not included, all 15 steps of the selected profile will be cleared.
PM	Profile Mode	PM <cr></cr>
	(returns status of current profile	There are no parameters to this command.
	run/hold mode)	Return data format:
		<arg1><space><arg2><space><arg3><cr></cr></arg3></space></arg2></space></arg1>
		where Arg1 =Run/hold status (0 = end, 1 = run),
		Arg2=Profile number and
		Arg3=Step number
		NOTE: If a profile has not been selected to run, a 0 will be returned for Profile number and Step number.

...continued...

Code	Meaning	Format
PR	Read	PR <space><arg1><space><arg2><cr></cr></arg2></space></arg1></space>
	Profile	where Arg1= Profile number (range 1 to 4) and
	step in	Arg2=Step number (range 1 to 15)
	memory	
		Return data format:
		<pre><arg1><space><arg2><space><arg3><space><cr><lf></lf></cr></space></arg3></space></arg2></space></arg1></pre>
		where Arg1 =Profile number (range 1 to 4),
		Arg2=Step number,
		Arg3=Temperature setpoint,
		Arg4=Agitation setpoint,
		Arg5=CO2 setpoint (future option, reads 0.0),
		Arg6=Step hours,
		Arg7=Step minutes,
		Arg8=UV lamp, and
	5 ("	Arg9=Grow lamp
PS	Profile	PS <cr></cr>
	Stop/	With no other parameters, any profile currently running will be
	Start	stopped.
	(stops or starts a	DC (Change (Augus) (Change (Augus) (CD)
	program	PS <space><arg1><space><arg2><cr> where Arg1—Profile number (range 1 to 4) and</cr></arg2></space></arg1></space>
	program profile)	where Arg1 =Profile number (range 1 to 4), and
	prome)	Arg2=Step number (range 1 to 15).
		NOTE: If Arg2 is not included, the profile will start in Step 1.
PW	Write	PW <space><arg1><space><arg2><space><arg3><space></space></arg3></space></arg2></space></arg1></space>
	Profile	<cr></cr>
	step to	where Arg1 =Profile number (range 1 to 4),
	memory	Arg2=Step number (range 1 to 15),
		Arg3=Temperature setpoint,
		Arg4=Agitation setpoint,
		Arg5=CO2 setpoint (future option); Enter 0.0,
		Arg6=Step hours (range 0 to 99),
		Arg7=Step minutes (range 0 to 59),
		Arg8=UV lamp (0 = off, 1 = on), and
		Arg9=Grow lamp (0 = off, 1 = on)

For example, the following **PW** command:

PW<Space>25.0<Space>150<Space>0.0<Space>2<Space>30<Space>0<<Space>1<CR>

will set Program 1, Step 1 to the following conditions:

- Temperature setpoint 25.0°
- Agitation setpoint 150 rpm
- CO2 setpoint 0.0%
- Step Time 2 hours 30 minutes
- UV Lamp OFF
- GRO Lamp ON

12.6 Report Request Commands

Code	Meaning	Format
RI	Report	RI <cr></cr>
	Software	
	Version	Return data format:
		< Arg1 >< <i>Tab</i> >< Arg2 >< <i>CR</i> >< <i>LF</i> >
		where Arg1=Display module software version, and
		Arg2=Control module software version
RP	Report	RP <cr></cr>
	Parameter	
	List (ASCII	Return data format:
	text header)	<arg1><tab><arg2><tab><arg3><arg6><cr><lf></lf></cr></arg6></arg3></tab></arg2></tab></arg1>
		where Arg1 ="Shaker Speed",
		Arg2="Temperature",
		Arg3="% Relative Humidity",
		Arg4="% CO2 Concentration",
		Arg5="Grow Lamp Status", and
		Arg6="UV Lamp Status"
	Report	RV <cr></cr>
RV	Measured	RS <cr></cr>
	Values	
	Б ,	Return data format:
D0	Report	<arg1><tab><arg2><tab><arg3><arg6><cr><lf></lf></cr></arg6></arg3></tab></arg2></tab></arg1>
RS	Setpoint	where Arg1 = Shaker Speed,
	Values	Arg2=Temperature,
		Arg3=% Relative Humidity,
		Arg4=% CO2 Concentration (future option, reads 0.0),
		Arg5=Grow Lamp Status, and
		Arg6=UV Lamp Status

12.7 Set/Get Date & Time

Code	Meaning	Format
=D	Set Date	=D <space><arg1><space><arg2><space><arg3><space></space></arg3></space></arg2></space></arg1></space>
	& Time	< Arg7 >< <i>CR</i> >
?D	Get Date	?D <cr></cr>
	& Time	Return data format:
		<arg1><tab><arg2><tab><arg3><tab><arg7><cr><lf></lf></cr></arg7></tab></arg3></tab></arg2></tab></arg1>
		where Arg1 =Hours (range 0-23),
		Arg2=Minutes (range 0-59),
		Arg3=Seconds (range 0-59),
		Arg4=Year (range 00-99),
		Arg5=Month (range 01-12),
		Arg6=Date (range 1-31),
		Arg7=Day (range 1-7; 1=Monday, 2=Tuesday7=Sunday)

For example, this **=D** command:

will set the real date and time on your shaker to:

1:30:00 am, 2004, July 4, Sunday.

13 APPENDIX: PRODUCT RETURNS

13.1 Return Procedure


As explained in Section 9.2, should you need to return your Innova 44/44R to NBS for any reason, first contact Customer Service to obtain a Returned Material Authorization (RMA) number. This number must appear on the outside of the shipping container, otherwise NBS Receiving will refuse to accept the shipment.

In addition, you must also certify that the instrument being returned has been thoroughly cleaned and decontaminated. A form for this purpose is provided on the following page; you can photocopy it and fill it out by hand. It can also be downloaded from our website (www.nbsc.com), if you prefer to fill it out electronically.

A copy of the completed Return Authorization and Decontamination Certificate must be attached to the outside of the container, with a second copy packed inside with the instrument.

13.2 Return Authorization and Decontamination Certificate

A sample form for you to copy and fill out is provided on the following page.

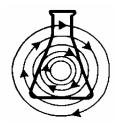

New Brunswick Scientific Return Authorization and Decontamination Certificate

Contact New Brunswick Scientific for an RMA number prior to returning any equipment, then complete this form and attach it to the outside container of the equipment being returned to our facility. In addition, please enclose a completed, duplicate copy of this form with the returned item.

Returned Material Authorization (RMA)	Number
Equipment being returned: Model Number	erSerial Number
Reason for its return:	
This equipment (check all that apply):	
New Product	Biohazards
☐ Never used	☐ Not used
	☐ Used, but decontaminated with
Hazardous Chemicals	Radioactive Materials
☐ Not used	☐ Not used
☐ Used, but decontaminated with	☐ Used, but decontaminated with
	ed above has been thoroughly cleaned and logical and radioactive contaminants and also certify protected human contact.
By:Signature	
	Print name Date:
Company:	
Address:	
Phone:	Fax: email:
1 none 1	ran tilidii
	Form 2847

14 INDEX

۸	Clamp Hardware Kits, 77
A	Cleaning, 60
Accessories, 76	Clearance, 14
Racks & Trays, 78	COM Parameters, 88
Actual Temperature, 56	COMM Screen, 47
Adjusting Door Tension, 69	Command Codes, 89
Alarms, 10, 28, 85	Command Format, 89
Arrows Icon, 28	Communication Mode, 48
ASCII Format, 90	Continuous Operation, 10
ASCII Text Headers, 93	Control Knob, 26
Audible Alarm	Control Schematic, 81
Disabling, 55	Copyright Notice, 5
Reactivating, 55	Counterbalance Bracket, 69
Audible Alarm Muted Icon, 26	Creating a Program, 51
Audible Alarms, 28	Customer Satisfaction Form 6300, 12
Automatic Restart, 58, 85	D
В	Day, 45
Dagge 79	Declaration of Conformity, 87
Bases, 78	Decontamination, 60
Baud Rate	Decontamination Certificate, 95
Setting the, 47	Decontamination Certificate, 96
Belt Replacement, 73	Dimensions
Biohazards, 39, 60	Chamber, 86
С	Overall, 85
•	Disclaimer Notice, 5
Calibrate Screen, 49	Display Icons, 28
Calibrating RPM, 57	Display Parameters, 26
Calspeed, 57	Display Screen, 40
CAUTION	Changing Parameters, 40
Symbol for, 5	Modifying Setpoints, 42
Chamber Light	Viewing Setpoints, 41
Auto, 46	Display Screens
On/Off, 46	Changing the, 27
Changing a Setpoint, 44	Door
Changing Baud Rate, 48	Opening the, 40
Changing Communication Mode, 48	Door Open Icon, 26, 28, 40
Changing Screens, 43	Door Tension
Changing the Day, 45	Adjusting the, 69
Changing the Time, 45	Double Girdle Clamps, 34

Drawings	Index of Drawings, 83
List of, 83	Index of Tables, 84
Drive, 85	Indicated Temperature, 56
Drive Belt	Inspection
Removing the, 73	Of Boxes, 12
Replacing the, 73	Installation
-	Tools for, 16
E	Installing Flask Clamps, 34
Edit Program, 55	Installing Front Grille, 16
Electrical Connections, 35	Installing Platform, 33
Electrical Requirements, 14, 85	Installing Subplatform, 32
Environment, 14	Interior Light, 30
Livironment, 14	International Offices
F	List of, 3
Firmware	L
Operation via, 10	_
Flammable Substances, 39	LAMP, 30
Flask Capacities, 77	Lamps Screen, 46
Flask Clamps, 77	LCD Display, 26, 85
Installation of, 34	Leveling
Front Grille	Base, 18
Installing the, 16	One Shaker, No Base, 17
Front View, 15	Locked Icon, 28
Fuses, 66	Locking Settings, 45
Tuses, 00	T 1 1 T 00 46
	Loudspeaker Icon, 28, 46
G	-
G Gassing Port Kit, 35, 79	М
Gassing Port Kit, 35, 79	M Manual Conventions, 5
Gassing Port Kit, 35, 79 Germicidal Lamp, 79	M Manual Conventions, 5 Missing Parts, 12
Gassing Port Kit, 35, 79	Manual Conventions, 5 Missing Parts, 12 Modifying a Setpoint, 42, 44
Gassing Port Kit, 35, 79 Germicidal Lamp, 79 Glide-Up Door, 29 Growth	Manual Conventions, 5 Missing Parts, 12 Modifying a Setpoint, 42, 44 Monitor Mode, 48
Gassing Port Kit, 35, 79 Germicidal Lamp, 79 Glide-Up Door, 29 Growth On/Off, 46	Manual Conventions, 5 Missing Parts, 12 Modifying a Setpoint, 42, 44
Gassing Port Kit, 35, 79 Germicidal Lamp, 79 Glide-Up Door, 29 Growth	Manual Conventions, 5 Missing Parts, 12 Modifying a Setpoint, 42, 44 Monitor Mode, 48
Gassing Port Kit, 35, 79 Germicidal Lamp, 79 Glide-Up Door, 29 Growth On/Off, 46	Manual Conventions, 5 Missing Parts, 12 Modifying a Setpoint, 42, 44 Monitor Mode, 48 Motor, 85 Mounting Shaker on Base, 18 Moving Screen to Screen, 43
Gassing Port Kit, 35, 79 Germicidal Lamp, 79 Glide-Up Door, 29 Growth On/Off, 46 Growth Lamps, 79 H	M Manual Conventions, 5 Missing Parts, 12 Modifying a Setpoint, 42, 44 Monitor Mode, 48 Motor, 85 Mounting Shaker on Base, 18
Gassing Port Kit, 35, 79 Germicidal Lamp, 79 Glide-Up Door, 29 Growth On/Off, 46 Growth Lamps, 79 H Handles, 80	Manual Conventions, 5 Missing Parts, 12 Modifying a Setpoint, 42, 44 Monitor Mode, 48 Motor, 85 Mounting Shaker on Base, 18 Moving Screen to Screen, 43 Mute, 46, 55 Mute Icon, 26, 28
Gassing Port Kit, 35, 79 Germicidal Lamp, 79 Glide-Up Door, 29 Growth On/Off, 46 Growth Lamps, 79 H Handles, 80 Heater, 30, 85	Manual Conventions, 5 Missing Parts, 12 Modifying a Setpoint, 42, 44 Monitor Mode, 48 Motor, 85 Mounting Shaker on Base, 18 Moving Screen to Screen, 43 Mute, 46, 55
Gassing Port Kit, 35, 79 Germicidal Lamp, 79 Glide-Up Door, 29 Growth On/Off, 46 Growth Lamps, 79 H Handles, 80 Heater, 30, 85 Heater On Icon, 26, 28	Manual Conventions, 5 Missing Parts, 12 Modifying a Setpoint, 42, 44 Monitor Mode, 48 Motor, 85 Mounting Shaker on Base, 18 Moving Screen to Screen, 43 Mute, 46, 55 Mute Icon, 26, 28 Muting the Alarm, 46
Gassing Port Kit, 35, 79 Germicidal Lamp, 79 Glide-Up Door, 29 Growth On/Off, 46 Growth Lamps, 79 H Handles, 80 Heater, 30, 85 Heater On Icon, 26, 28 HRS Alarm, 29	Manual Conventions, 5 Missing Parts, 12 Modifying a Setpoint, 42, 44 Monitor Mode, 48 Motor, 85 Mounting Shaker on Base, 18 Moving Screen to Screen, 43 Mute, 46, 55 Mute Icon, 26, 28
Gassing Port Kit, 35, 79 Germicidal Lamp, 79 Glide-Up Door, 29 Growth On/Off, 46 Growth Lamps, 79 H Handles, 80 Heater, 30, 85 Heater On Icon, 26, 28 HRS Alarm, 29 Humidity Monitor, 29, 79	Manual Conventions, 5 Missing Parts, 12 Modifying a Setpoint, 42, 44 Monitor Mode, 48 Motor, 85 Mounting Shaker on Base, 18 Moving Screen to Screen, 43 Mute, 46, 55 Mute Icon, 26, 28 Muting the Alarm, 46 N
Gassing Port Kit, 35, 79 Germicidal Lamp, 79 Glide-Up Door, 29 Growth On/Off, 46 Growth Lamps, 79 H Handles, 80 Heater, 30, 85 Heater On Icon, 26, 28 HRS Alarm, 29	Manual Conventions, 5 Missing Parts, 12 Modifying a Setpoint, 42, 44 Monitor Mode, 48 Motor, 85 Mounting Shaker on Base, 18 Moving Screen to Screen, 43 Mute, 46, 55 Mute Icon, 26, 28 Muting the Alarm, 46 N New Program, 51
Gassing Port Kit, 35, 79 Germicidal Lamp, 79 Glide-Up Door, 29 Growth On/Off, 46 Growth Lamps, 79 H Handles, 80 Heater, 30, 85 Heater On Icon, 26, 28 HRS Alarm, 29 Humidity Monitor, 29, 79	Manual Conventions, 5 Missing Parts, 12 Modifying a Setpoint, 42, 44 Monitor Mode, 48 Motor, 85 Mounting Shaker on Base, 18 Moving Screen to Screen, 43 Mute, 46, 55 Mute Icon, 26, 28 Muting the Alarm, 46 N
Gassing Port Kit, 35, 79 Germicidal Lamp, 79 Glide-Up Door, 29 Growth On/Off, 46 Growth Lamps, 79 H Handles, 80 Heater, 30, 85 Heater On Icon, 26, 28 HRS Alarm, 29 Humidity Monitor, 29, 79 HyperTerminal, 88	Manual Conventions, 5 Missing Parts, 12 Modifying a Setpoint, 42, 44 Monitor Mode, 48 Motor, 85 Mounting Shaker on Base, 18 Moving Screen to Screen, 43 Mute, 46, 55 Mute Icon, 26, 28 Muting the Alarm, 46 N New Program, 51 NOTE Symbol for, 5
Gassing Port Kit, 35, 79 Germicidal Lamp, 79 Glide-Up Door, 29 Growth On/Off, 46 Growth Lamps, 79 H Handles, 80 Heater, 30, 85 Heater On Icon, 26, 28 HRS Alarm, 29 Humidity Monitor, 29, 79 HyperTerminal, 88	Manual Conventions, 5 Missing Parts, 12 Modifying a Setpoint, 42, 44 Monitor Mode, 48 Motor, 85 Mounting Shaker on Base, 18 Moving Screen to Screen, 43 Mute, 46, 55 Mute Icon, 26, 28 Muting the Alarm, 46 N New Program, 51 NOTE

Operation Modes, 10	S
Overview, 10	Safety, 85
D.	Safety Precautions, 39
P	Safety Switch, 10
Packing List	Sample Program Screen, 53
Verification of, 12	Save Program, 54
Padlock Icon, 28, 45	Screen Parameters, 27
Parameter Locked Icon, 26	Screw-In Handles, 80
Photosynthetic Lamps, 30, 36, 79	Service Accessibility, 31
Platform	Set Commands, 91
Installing the, 33	Setpoint Retention, 85
Platform Assemblies, 32	Setting the Baud Rate, 47
Platform Screw Locations, 32	Setup Screen, 44
Platforms, 11	Changing the Day, 45
POWER Alarm, 29	Changing the Time, 45
Power Failure, 58	Locking Settings, 45
Power Interruption, 58	Muting the Alarm, 46
Problems, 62	Shaking Speed, 85
Product Returns, 64, 95	Shipping Damage, 12
Program Running Icon, 28	Slave Mode, 48
Programming, 50	Slide-Out Platform
Programs Screen, 49	Available Types, 76
-	Disengaging the, 59
R	Installing the, 33
Refrigeration, 30	Operating the, 58
Refrigeration Schematic, 82	Slide-Out Platform Accessory, 32
Remote Alarm Option, 80	Software Interface, 30
Remote Programming, 89	Spare Requirements, 14
Replacement Parts, 76	Specifications, 85
Report Request Commands, 93	SPEED Alarm, 29
Resetting the Day, 45	Spill Pan, 29
Resetting the Time, 45	Filling the, 36
Return Authorization, 95	Spills, 60
Return Authorization, 96	Stacking
Returned Material Authorization Number,	Three Shakers, 24
64, 95	Two Shakers, 22
Returned Material Safety Sheet, 64	Stacking Kit, 79
RMA, 64, 95, 96	Installation of, 20
RS-232, 30, 85, 89	Starting the Unit, 40
RS-232 Interface	Subplatform
Operation via, 10	Installing the, 32
RS232 Screen, 47	Summary Screen, 43
Run Icon, 26, 28, 55	Changing Setpoints, 44
, , ,	Symptoms, 62

T

Tables List of, 84 Talk Mode, 48 TEMP Alarm, 29 Temperature Offset, 49, 56 Setting the, 56 Temperature Offset Calibration, 56 Temperature Offset Value, 56 Temperature Range, 85 Test Tube Racks, 78 TILT Alarm, 29 Tilt Switch, 62 Time, 45 Timed Functions, 50 Timed Run, 10 Troubleshooting, 62

U

Universal Voltage Input Module, 35 Unpacking Tools for, 12, 16 UV Decont On/Off, 46 UV Lamp, 30, 79

V

Viewing a Setpoint, 41 Visible Alarms, 28

W

WARNING
Symbol for, 5
Warranty, 6, 7
Warranty Registration, 13
Water Reservoir, 29
Filling the, 36
Wavy Arrows Icon, 28
Weight, 86