+

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International certification
Excellent quality

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Sanitation license and permission No. GWX 0189 (2010)

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Production enterprises hygiene license number：GWX (2018) No. 9134th

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Made by：Shenzhen Huaide Technology Development Co., Ltd.

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Address：1st Phoenix Industrial Zone, Fuyong Town, Shenzhen, Guangdong, China

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Please keep the instruction manual carefully for future reference

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Thank you very much for purchasing this product. Please read the User Manual carefully before use.

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Warranty Card

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Special statement

- 1、若本产品日后有任何技术变更，恕不另行通知。
If there is any technical change in this product in the future, without prior notice.
- 2、产品外观、颜色若有改动，以实物为准。
If the appearance and color of the product are changed, the actual product shall prevail.

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- 1、按国家“新三包法”规定，享受三包服务,整机保修一年,主要电子零配件保修三年。
According to the national "New Three Guarantees Law", enjoy the three-pack service, the whole machine is guaranteed for one year, the main electronicParts are guaranteed for three years.
- 2、属系列情况之一的不享受三包服务，但可以付费维修（需出示发票）：
One of the series of situations does not enjoy the three-pack service,but can pay for repairs (invoices are required):
 - 1) 消费者使用、维护、保管不当造成的损坏；
Damage caused by improper use, maintenance and storage by consumers;
 - 2) 非本公司专业维修人员拆卸或操作不当造成的损坏;
Damage caused by improper disassembly or improper operation of the company's professional maintenance personnel;
 - 3) 无三包凭证及有效发票，三包凭证与维修产品型号不符合或涂改；
No three-package certificate and valid invoice, the three-pack certificate and maintenance product model do not match or be altered;
 - 4) 因不可抗力造成的损坏。
Damage caused by force majeure.
- 3、保修期的起始时间以产品发票时间为准。
The start time of the warranty period is based on the product invoice time.
- 4、超过保修期或不属于保修范围的产品，本公司维修中心竭诚为您服务。
For products that exceed the warranty period or are not covered by the warranty, the company's repair center will serve you wholeheartedly.

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International organization's assessment of chlorine dioxide

- 01 -

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Scope of application

- 02 -

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Product description

- 03 -

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Operation interface icon description

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Operation interface introduction and installation method

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The Greewin Brand series of chlorine dioxide air sterilizer is independent R&D by Shenzhen Huaide Technology Development Co., Ltd, is a new technology integrating physics, chemical engineering, biology and environmental science. It realizes 360 degree medical grade, novel and intelligent. Safe, broad-spectrum disinfection system; has obtained 5 invention patents and more than 30 utility models. Using the world's leading technology, through the synthesis, passivation, and formation of pure molecular chlorine dioxide through new materials, chlorine dioxide molecules with strong oxidizing properties quickly flow to any corner of the space, achieving 360-degree disinfection and sterilization effect without dead angle. It kills various pathogenic bacteria, pathogenic virus microorganisms, pathogens and infectious diseases in space and pipelines, and decomposes and removes indoor harmful gases such as formaldehyde, PM2.5 and TVOC, which can effectively control influenza. The spread chain of infectious diseases prevents and controls the occurrence of various chronic diseases caused by the respiratory tract; the air disinfection machine subverts the traditional disinfection mode at home and abroad, and is a revolution in traditional disinfection, opening up a dynamic disinfection to static disinfection. The revolutionary great pioneering work established a new starting point for the original disinfection to intelligent disinfection.

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International organization's assessment of chlorine dioxide:

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Chlorine dioxide is listed as an A1 safe and efficient green disinfectant by the World Health Organization (WHO) and the United States Environmental Organization (EPA) and the United Nations Food and Agriculture Organization (FAO). It is the only one in the disinfection technology. It can kill all pathogenic microorganisms and is not easy to produce drug resistance, especially for typhoid fever, hepatitis A, hepatitis B, polio and HIV. It has no carcinogenic, no teratogenic, no mutagenicity. As early as the 1970s, it was widely used in water treatment and other industries; it also has the unique function of deodorization and preservation, and has been incorporated into China's health regulations for use.

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The air sterilizer is limited to the use of green card air disinfectant chlorine dioxide, if use other disinfectants and accident happened, the company shall not be held any responsibility.

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There is still a smell of chlorine dioxide in the air after every working cycle, belongs to the normal phenomenon, after a period of time the smell will disappear automatically, has no effect on the human body.

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It is strictly prohibited to head near the outlet when a air sterilizer is turned on, please keep more than 1m.

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Prohibited for interruption the working process by force. Prohibition of non-professional to disassembly a air sterilizer.

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The ClO2 bottle should storage in accordance with the Greenwin brand ClO2 powder manual, banning children touch the air sterilizer, waste water tank and ClO2 bottle.

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The face need more than 1m from the outlet when a air sterilizer is working.

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The occurrence of suspended solids in the discharged wastewater is a normal phenomenon. The wastewater can be used for disinfection of toilet floors, sewers, etc., and avoid contact with the skin during drainage. If it is accidentally splashed on the skin or eyes, please rinse with water. If it is serious, seek medical attention immediately.

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Please do not place a air sterilizer in the position of the sunlight.

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Prohibition a air sterilizer is used in flammable and explosive substances environment.

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After the disinfectant pill is loaded into the canister, be sure to tighten it.

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Scope of application:

- 1. 医院、社康中心、药房、实验室、卫生应急中心及室内避难应急场所；
Hospitals, community health center, pharmacy, laboratory, health emergency center, and indoor refuge for emergency;
- 2. 卫生系统的应急预案措施及国家全部特殊场合的环境要求；
The environmental requirements of the department of homeland security special occasions;
- 3. 学校、幼儿园、早教中心、培训基地、室内体育馆、训练场；
Schools, kindergartens, early education center and training base, indoor stadium, training;
- 4. 家居、办公室、会议室、阅览室、新装修室内场馆；
New decorate household, office, conference room, reading room, indoor venues;
- 5. 宾馆、酒店、会所、美容院、超市；
Guesthouses, hotels, clubs, beauty salon, supermarket;
- 6. 车站、码头、地铁、口岸、机场；
The station, wharf, subway, port, airport;
- 7. 地道、防空洞、指挥中心、地下军用设施、舰艇、潜艇；
Tunnels, shelters, command center, the military facilities, ships and submarines;
- 8. 养老院、康复中心、制药厂、食品厂、无菌车间；
Nursing homes, rehabilitation centers, pharmaceutical factory, food factory, aseptic workshop;
- 9. 宠物医院、宠物爱好者家庭、育婴家庭及场所。
Pet hospital, pet lovers family, parental family and place.

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

The table below lists some of microbiologists the chlorine dioxide can kill:

微生物名称 Microbiologists	易引起的病名 Easy cause of disease	微生物名称 Microbiologists	易引起的病名 Easy cause of disease
大肠杆菌 Staphylococcus aureus	尿道炎、腹膜炎、败血症等 Urethritis, peritonitis, septicaemia, etc.	禽流感病毒 (AIV) Avian influenza virus	禽流感 Bird flu
金黄色葡萄球菌 Escherichia coli	引起局部化脓感染，也可引起肺炎、伪膜性肠炎、心包炎、红眼病等，甚至败血症、脓毒症等全身感染 Cause local purulent infection, also can cause pneumonia, pseudo membranous enteritis, pericarditis, pink eye, etc., or sepsis, sepsis systemic infection, etc	非典病毒(SARS病毒) SARS	非典型性肺炎 Atypical pneumonia
		霉菌 Mucedine	引起发霉 Typhoid fever caused by mold
		肝炎双球菌 Hepatitis diplococcus	肝炎 Hepatitis b
结核杆菌 Mycobacterium	结核病 Tuberculosis	链球菌 Streptococcus	中耳炎 Otitis media
乙型肝炎病毒 Hepatitis B	乙型肝炎 Hepatitis B	肠病毒 Enterovirus	病毒性脑膜炎 Virus meningitis
霍乱菌 Cholera bacteria	烈性肠道传染病 Strong intestinal infectious diseases	手足口病毒hand, foot and mouth virus virusvirus tuberculosis	手足口病 Hand-foot-and-mouth disease
沙氏菌 Salmonella choleraesuis fungal	胃肠炎 Gastroenteritis	H1N1流感病毒 H1N1tuberculosis	甲型H1N1流感(猪流感) Influenza a H1N1 influenza (swine flu)
真菌 Fungus	皮肤癣菌病、痢疾 Skin tinea disease, dysentery	牙斑菌 Plaque bacteria	牙龈炎、口臭 Gingivitis, bad breath,
军团杆菌 Legionella	肺炎 Pneumonia	伤寒杆菌 Typhoid bacillus	伤寒 Pneumonia
百日咳杆菌 Bordetella pertussis	百日咳 Whooping cough	流感病毒 Influenza	流感 Flu
炭疽杆菌 Bacillus anthraci	炭疽病 Anthrax	破伤风菌 Clostridium tetani	破伤风 Tetanus

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Product description

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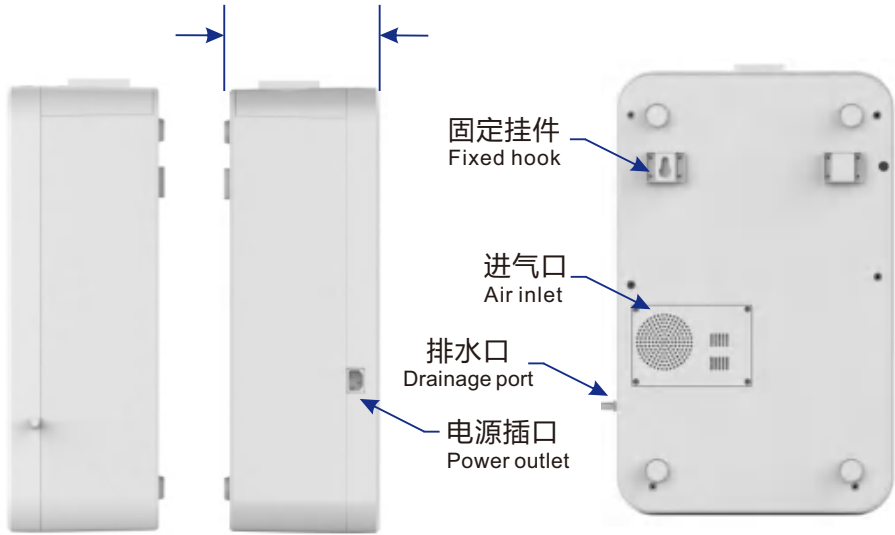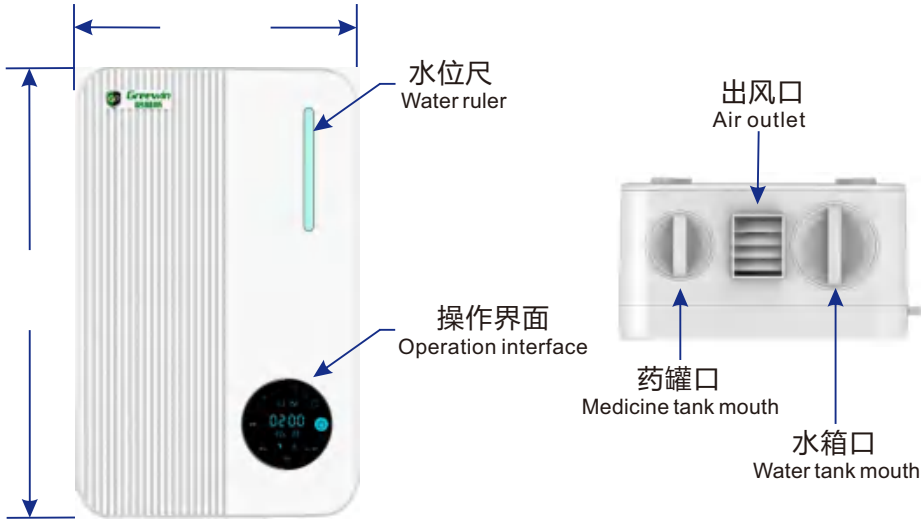Operation interface icon description

图标 Icon	名称 Name	图标含义 Icon meaning
	缺水报警 Water shortage alarm	水位达到最低水位时此图标会闪烁报警,机器会停止工作,请打开水箱盖加水到正常水位,方可正常使用。 When the water level reaches the minimum waterlevel, the icon will flash alarm and the machine will stop working. Please open the water tank coverand add water to the normal water level before it canbe used normally.
	缺丸剂报警 Lack of pills alarm	药罐内无消毒丸剂时此图标会闪烁报警,机器会停止工作,请打开药罐盖加入药丸,方可正常使用。 When there is no disinfectant pill in the tank, the icon will flash and the machine will stop working. Please open the can lid and add the pill before it can be used normally.
	故障报警 Fault alarm	机器内部出现故障时此图标会闪烁,此时请与客服联系,我们会有专业的售后人员为您服务。 This icon will flash when there is a malfunction inside the machine. Please contact customer service at this time, we will have professional after-sales personnel to serve you.
	排水提示 Drainage prompt	机器内部排水箱水满时此图标会闪烁,机器会停止工作,请按键排完水箱内的水,方可正常使用。 When the water inside the machine is full, the icon will flash and the machine will stop working. Please press the water in the water tank to use it normally.
	单次工作模式 Single working mode	此模式图标点亮时,机器处于单次工作状态,机器完成一个工作周期后停止工作,进入待机状态。 When this mode icon is lit, the machine is in a single working state, and the machine stops working after completing one working cycle and enters the standby state.
	定时模式 Timing mode	此模式图标点亮时,机器会遵照定时设置时间工作。 When this mode icon is lit, the machine will work according to the timing setting time.
	连续模式 Continuous mode	此模式图标点亮时,机器会连续工作。 When this mode icon is lit, the machine will work continuously.

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Operation interface introduction

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

Wall-mounting method: Fix the screws provided at the appropriate position on the wall, the spacing is the same as the pendant. Check whether the products are parallel, and hang the product on the pendant. Check whether it is firm and safe after the installation.

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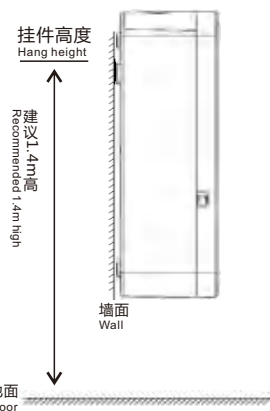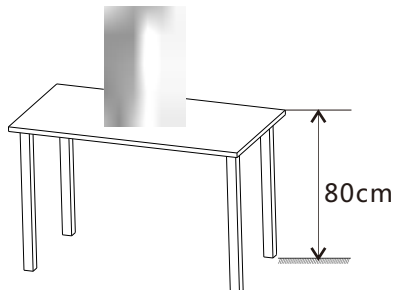Operation method

1. 打开药罐盖子,投入格林威消毒丸剂,药罐最大容量为50粒,并将药罐盖子拧紧。
Open the lid of the canister and put it into the Greenway Disinfecting Pill. The maximum capacity of the canister is 50 and tighten the lid of the canister.

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Open the lid of the tank and add water. The maximum capacity of the tank is 3800ml. Please be careful not to exceed the maximum water level and tighten the lid of the tank.

3. 接通电源(输入电压100-240V),单击"电源键"开机,显示屏点亮,进入待机界面。
Turn on the power (input voltage 100-240V), click "Power button" to turn on, the display lights up, enter the standby interface.