



PHOENIX LTC

MED CARTS

Contents:

Page 2:	Drawer Removal
Page 3:	Drawer Installation
Page 4:	Drawer Front Removal and Installation
Page 5:	Slide Removal and Installation
Page 6:	Side Shelf Removal and Installation
Page 7:	Caster Removal and Installation
Page 8:	Lock Core Removal and Installation
Page 9:	Cabinet Lock Removal
Pages 10-11:	Cabinet Lock Installation
Page 12:	Internal Narc Lid Removal and Installation
Pages 13-14:	Lock-Bar Removal and Installation
Page 15:	Locking Pin Removal and Installation
Page 16:	Narc Box Removal and Installation
	Narc Box Hinge Removal and Installation
Page 17:	Narc Box / Drawer Lock Removal
Page 18:	Panel Removal and Installation
Pages 19-21:	EMAR Arm Installation
Page 22:	Wood Top Replacement
Pages 23-32	Parts List
Page 33:	Misc. Information and Measurements

Drawer Removal

Tools Needed: None


1. Lock front casters by pushing caster levers down with the front of your foot.



2. Fully extend the drawer to be removed.



3. Stand in front of the drawer.
4. Locate the release tabs on both the slides
5. Push each release tab to the opposite side of the slide channel. (see enlarged detail below)
6. Hold the tabs in the pushed (down/up) position and carefully pull the drawer from the cart.

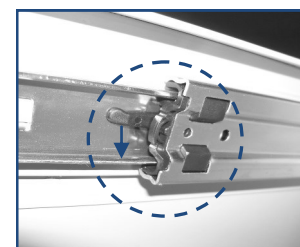

Drawer Installation

Tools Needed: None


1. Extend slides out of the cart. **IMPORTANT:** Move the ball bearings to the very front of the slide.



THE ABOVE PHOTO IS INCORRECT!
(the ball bearing section is NOT at the front)

*The photo to the right—shows:
Top slide ball bearing section is in the wrong
position for drawer installation.*

*Lower slide has ball bearing section in the
correct position for drawer installation.*

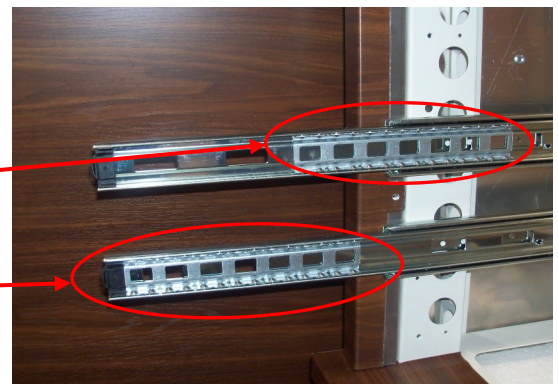

2. Insert the drawer left and right male section of the slide into the female section of the slide. Make sure that the male section feeds into the ball bearings. If this is NOT completed before attempting to push the drawer closed, the slide may be damaged.
3. Push the drawer evenly into the cart until it is in the closed position.

IMPORTANT: Do not force the drawer closed. It is normal to feel resistance from the point when the drawer is half closed until the closed position. *If you feed the drawer stopping (not resistance), remove and reinsert the drawer, confirming the ball bearings are at the very front of the slide.*



Drawer Front Removal and Installation

Tools Needed: Phillips Head Screw Driver


1. Open the drawer.
2. Remove the handle by removing the screws on the inside front of the drawer. Save/set aside both handle and screws. If you have two(2) knob type handles, they must both be removed.
3. Remove the old drawer front, by slightly lifting the top lip, and pulling forward. When the top clears the front of the drawer, it will disengage from the bump-outs on the bottom.
4. Place the new drawer front on the top edge of the drawer.
5. Slide the new drawer front under the bottom of the drawer. It will need a little pressure until the holes in the bottom of the new front engage the bump-outs on the bottom of the drawer. The last little bit, may take some firm pressure.
6. Make sure the holes for the handle are aligned in both the drawer and the new Front.
7. Insert the two(2) screws through the drawer and drawer front into the threaded holes in the handle on each side.
8. Leave the screws a little loose, double check alignment of new drawer front with other drawers. The front may be shifted a tiny amount to align with other drawer fronts.
9. Finish tightening both of the screws in the handle.



Slide Removal and Installation

Tools Needed: Phillips Screwdriver


1. Remove the Drawer. (see page 2)
2. Move the slide out until the back screw is exposed.
3. Remove the screw by turning counterclockwise.
4. Let the rear of the slide rest in the down position when the screw is removed
5. Move the ball bearings section until you can see the front screw head.
6. Remove the back screw by turning it counterclockwise.
7. Reverse the above process to install new slide, using built in pre-threaded holes.

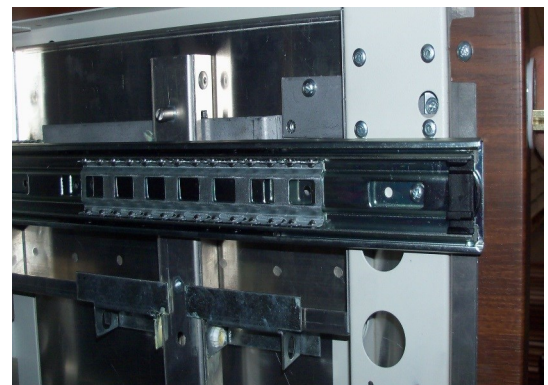

Side Shelf Removal

Tools Needed: None


1. Remove the top drawer. (see page 2)
(***Note: you may need to remove the 2nd drawer if you need more space to work*)
2. Pull out the side shelf until it stops.
3. Reach inside the cart on the right-hand side and find the two(2) buttons (see photo).
4. Depress and hold both the buttons in as you pull the side shelf out. You will feel it slip into the safety position.
5. Again, depress and hold the button in as you pull the side shelf out. **NOTE: Be Careful as you pull the shelf out, as it will fall free of the cart.**

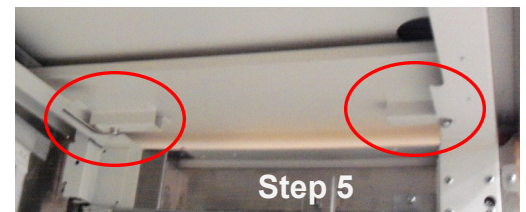

Side Shelf Intallation

Tools Needed: None

1. Carefully insert the end of the shelf into the slot on the right hand side of the cart, until the buttons hit the edge of the opening. (*see lower photo*)
2. Depress and hold both the buttons in as you push the side shelf in. You will feel it slip into the safety position.
3. Reaching inside the cart, again, depress and hold the button in as you push the side shelf in.
4. Continue pushing side shelf until it is in the closed position.

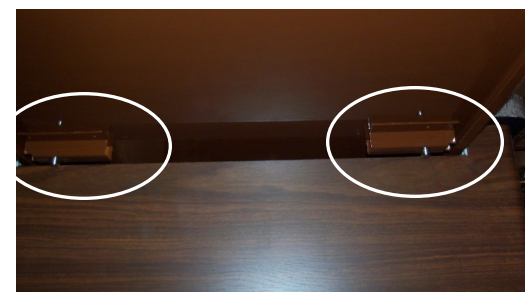

Caster Removal and Installation

Tools Needed: 3/8" or Adjustable Wrench

1. Remove all the Drawers. (see page 2)
2. With two or more people, carefully lay the cart on its back.
3. Replace the proper caster in the proper position.
 - a. Two front casters are braking. Identified by light gray levers.
 - b. Back-right caster is tracking. Identified by the dark gray lever.
 - c. Back-left caster is a swivel. It does not have a lever.
4. Using a 3/8" or adjustable wrench, turn the four bolts holding the caster plate counterclockwise until the plate can be removed.
5. To install a new caster, reverse the above process.



Lock Core Removal and Installation

Tools Needed: WHC2 Control Key

1. Insert the WHC2 Control Key into the lock.
2. Turn clockwise approximately 1/8 of a turn and pull.
3. To install a new core, reverse the above process.



Cabinet Lock Removal

Tools Needed: Phillips Screwdriver

1. Remove the top 4" Drawer (see page 2).
2. Lock the cabinet by pushing the lock in until you feel it latch.



3. Remove the two(2) small screws in the Lock Link.
4. Unlock the cart.



5. Pull the lock assembly from the cart.




Lock Link and Lock

Cabinet Lock Installation

Tools Needed: Phillips Screwdriver

1. Snap the plastic ring onto the front of the lock chamber, or make sure the existing one is secure.
2. Insert the lock assembly into the lock chamber, with the brass button on the top. Make sure that the Lock Link goes thru the slot in the back of the lock chamber. *(pictured without plastic ring)*
3. Push the brass pin down, and push the lock assembly into the locking position.

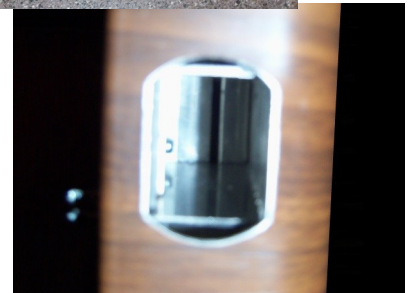


Cabinet Lock Installation continued

Tools Needed: None

4. Line up the holes in the Lock Link with the matching holes in the Lock Bar.



5. Insert two(2) Pan Head Phillips Screws 8-32 x 3/8.

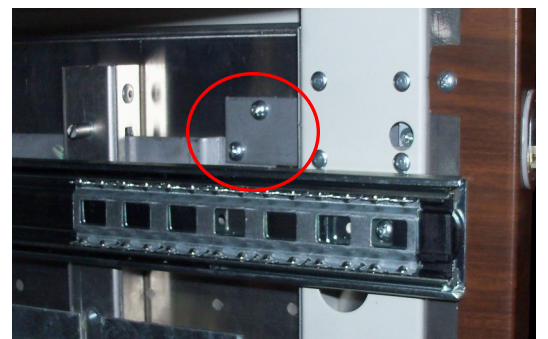

Internal Narc Lid Removal

Tools Needed: Pliers & Phillips Screwdriver

1. Remove the Drawer above the Narc Drawer(see page 2).
2. Remove the Narc Drawer (see page 2).
3. Remove four Phillips screws on the bottom of the internal narc lid and remove the lid.

NOTE: The lid may need a gentle nudge or bump to get it to move. It is a very tight installation.



Internal Narc Lid Installation

Tools Needed: Pliers & Phillips Screwdriver

Lock Bar Removal

Tools Needed: Phillips Screwdriver


1. Remove all the drawers (see page 2).
2. Remove the cabinet lock (see pages 9)
3. Remove the slides on the right-hand side, facing the cart front (page 5)



4. Remove both the Lock Bar Guides supporting the Lock Bar, by removing two screws in each one.
5. Remove both the upper Lock Bar Retainers, by removing two screws in each one.



6. Remove both the lower Lock Bar Retainers, by removing two screws in each one.

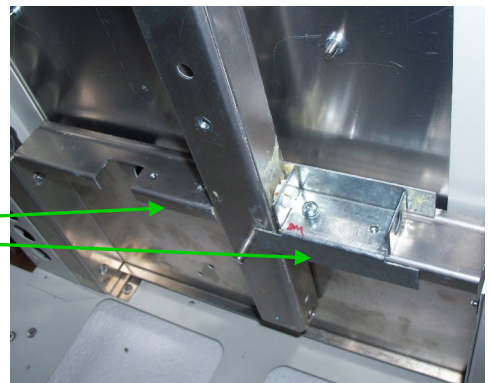


Lock Bar Removal Continued

7. Remove the center screw in the Lock Rail and Lock Bar.



8. Remove the Lock Bar and Lock Rail Assembly, by pulling the front side of the Lock Bar out past the edge of the upper Lock Plate. This will allow the Spring in the back to decompress and make it possible to remove the back part of the Lock Bar from the slot.

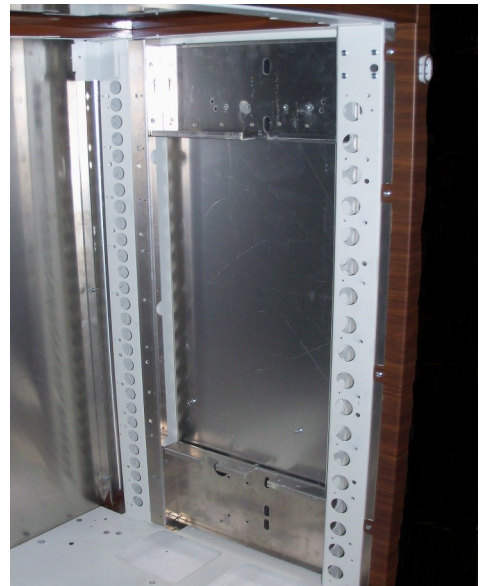

Lock Pin Installation

Tools Needed: Flat Screwdriver

1. Parts Needed:
 - a. Locking Pin, one
(refer to part list for correct part number)



2. Screw the new Lock pin into the appropriate pre-threaded hole in the lock rail.



Narc Box Removal

Tools Needed: Phillips Screwdriver


1. Open the narc box.
2. Remove the two screws in the back of the box.
3. Lift the box out of the drawer.



Narc Box Installation

Tools Needed: Phillips Screwdriver

1. Place narc box in the drawer on the right-hand side.
2. Line it up with the threaded openings.
3. Using the 1/4-20 screws provided, screw the narc box into the back of the drawer to secure it to the cart.

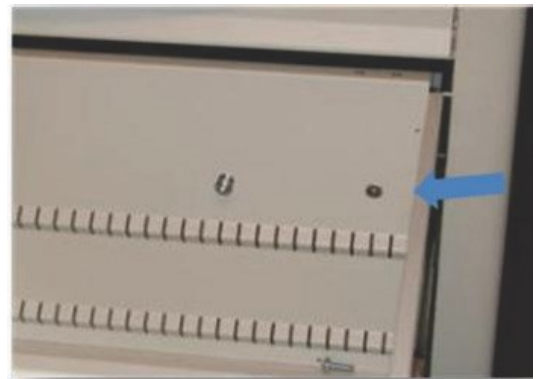

Narc Box Hinge Removal and Installation

Tools Needed: Small adjustable wrench

1. Remove the eight nuts on the hinge.
2. Remove the hinge.
3. Reverse the process to install the hinge.



Narc Box / Drawer lock Removal and Installation

Tools Needed: Phillips Screwdriver


1. Remove four screws
2. Remove the lock.
3. Reverse the process to install the lock.



Panel Removal and Installation

Tools Needed: Phillips Screwdriver


1. Remove all drawers (see page 2).
2. Remove corner bumpers (four screws each).
3. Completely remove cabinet lock (pages 9).
4. Remove the side shelf (see page 6)
5. Remove screws from panels. Go through the holes in the frame to remove screws located in the back of the panels.
6. Remove the top.
WARNING: *Be careful not to let the back fall off when the top is removed.*

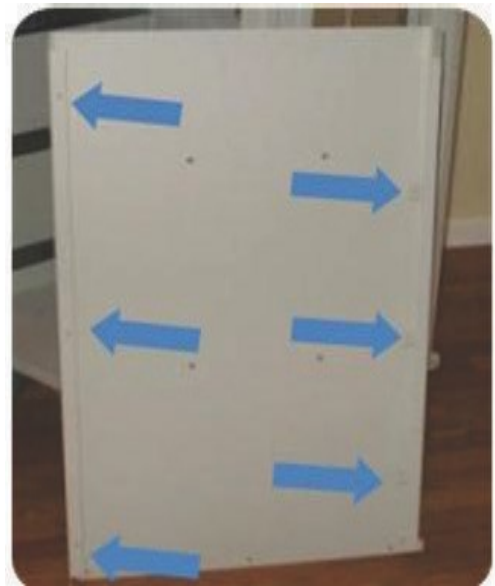

EMAR Installation Continued

Package Contents

- 1 each — Ergotron Arm Mount Kit
(used only for V3 carts, "PLXMNT-1")
- 1 each — Ergotron LX Arm
- 1 each — EMAR Installation Hardware Kit
(used only for V2 carts, "EMAR KIT")

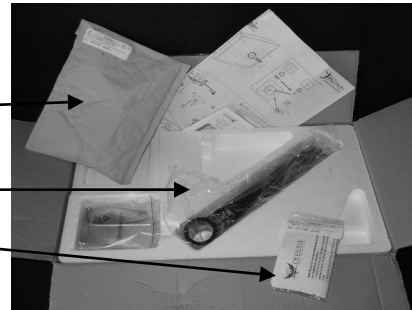


***** NOTE *****

- ⇒ If the Serial number on your cart begins with **"B"**, you have a V2 cart — Discard the manila envelope of parts labeled "PLXMNT-1". You will use the Ziploc bag labeled "EMAR KIT" and instructions included in the bag.
- ⇒ If the Serial number on your cart begins with **"P"**, you have a V3 cart — Discard the Ziploc bag of parts labeled "EMAR KIT". You will use the manila envelope labeled "PLXMNT-1" and the following instructions on this page.
- ⇒ Additional helpful hints and instructions may be found @ <http://phoenixltc.com/parts-store/refurbish-manuals>
- ⇒ Be sure to coordinate these instructions with your EMAR installation instructions. No additional bracing is needed since the cart frame is manufactured with strength to accommodate the EMAR Arm.


Additional Tools Needed: Marker, Phillips Screwdriver, Center Punch, Power Drill, 1-1/8" Hole Saw



1. Remove the top drawer.
(see Refurbish Manual for help with this. You may need to remove the additional drawer(s) if you need more space to work)
2. If your cart top was drilled and plugged at the factory, or if you ordered a new top that was drilled at the factory, you can skip to step 9.

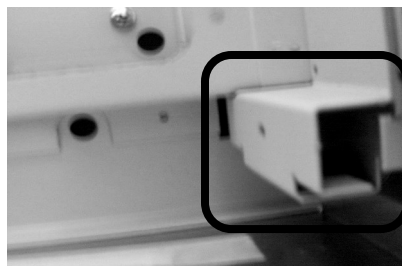


3. Remove the Side Shelf. (see Refurbish Manual for help with this.)
4. Find the square tube in the rear inside corner of the cart.
5. Using a marker, carefully mark the approximate shape of the square tube on the under side of the cart top.


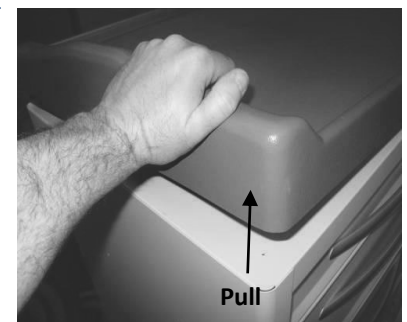


6. Remove the Cart Top
*Plastic Top — Pull the top off the cart. It will have some resistance because it is attached by Velcro.
*Wood Top — Unscrew from inside cart.

EMAR Installation Continued

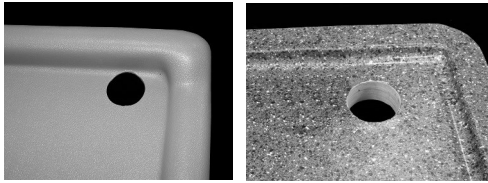


7. Using a 1 1/8" hole saw, Cut out a circular hole, by centering the square marked on the underside of the top in step # 5 in the circular hole.

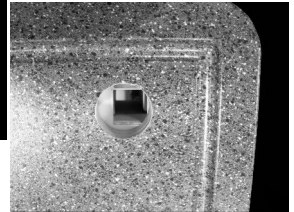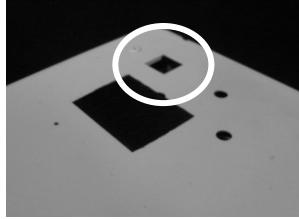


8. Replace the top on the cart. The round cut hole should be centered over the square hole/opening in the top of the cart.



9. If you have a top that is drilled/plugged at the factory, Remove the plug, and discard.

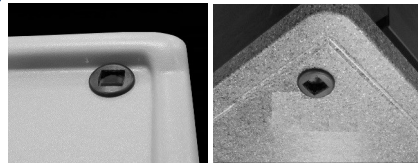




10. Insert the new plug cut with the square hole.
*Plastic top—short plug.
*Wood top—Tall plug.

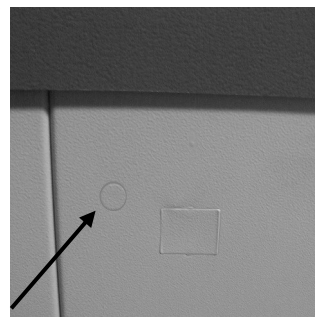

11. On the back of the cart, locate the small round knockout. This should be directly under the corner hole and plug with square cutout.




12. Knock out or Drill out the hole in the rear cart skin. (if drilling—it will be helpful to use a center punch to start the drill hole, use a power drill and a 1/4" drill bit).

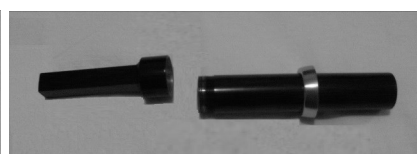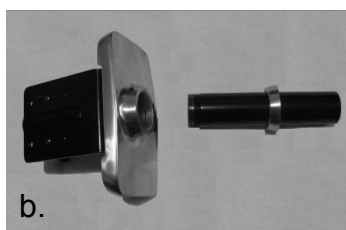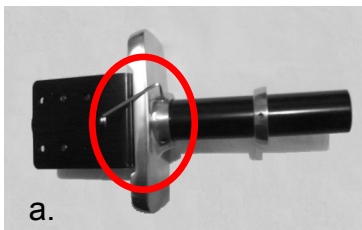


13. Remove the Post from the base of the ERGOTRON LX Desk Mount LCD Arm: (pictured from the back)

- Loosen the two(2) set screws in the back of the base with the allen wrench included.
- Unscrew the post from the square base. (Note: this may be easily accomplished by fastening it to a solid surface)
- Screw the post into the Arm Mount.

EMAR Installation Continued

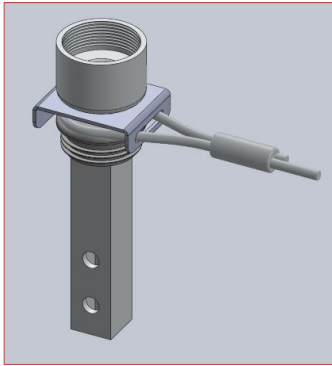

14. Above view is for reference only, of Arm Mount with Cable Lock Plate, Plastic Plug with Square hole, and Locking Cable.

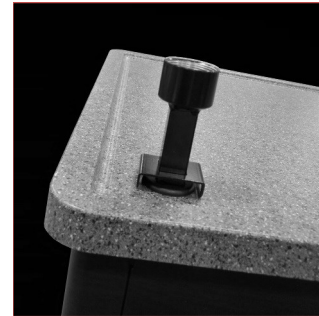

15. Insert the square tube of the Arm Mount through the square hole in the cart top and into the cart. The holes in the square tube need to face towards the rear of the cart.



16. Insert the Pan Head Phillips Machine Screw 1/4-20 X 1.5" through the Flat 1/4" washer and into the hole in the back of the cart , and into the Arm Mount.






17. Tighten the screw.

Wood Top Replacement

Tools Needed: Screwdriver

1. Remove at least top two drawers from the cart.
(see page 2)
(**Note: you may need to remove the additional drawer(s)
if you need more space to work)
2. Remove the Side Shelf
3. Find the existing holes in the underside of the cart top as pictured. There are three in the left end where the shelf is, and two in the right side.
4. Line up the holes in the bottom of the top with the rivet head on the top surface of the cart.
5. Using the screws provided, attach the wooden top to the cart, through the existing holes identified in step 3.
6. Re-install the drawers removed in step 1.

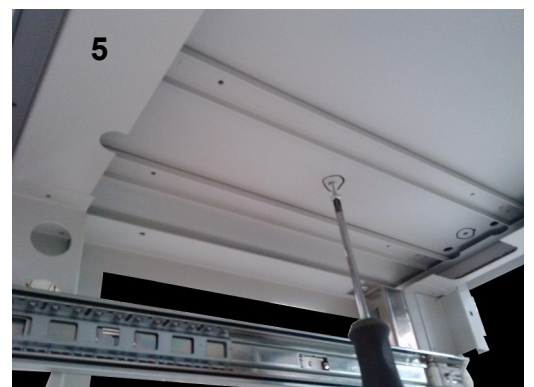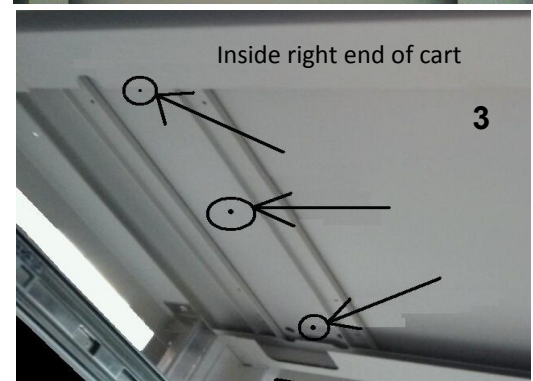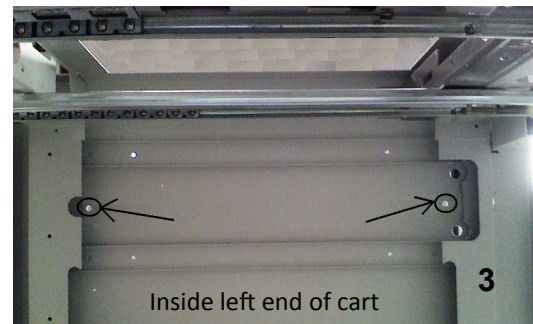

Replacement Parts

Item #	Description	Cart Version	
CP-2	Replacement Caster Package (6" Plastic casters)	Set	All
CP-1	Replacement Caster Package (5" Plastic Casters)	Set	All
M12982A7	5" caster with brake	Each	All
Y700069A7	5" tracking caster	Each	All
M12983A7	5" swivel caster	Each	All
63-1200	Double-Sided Foam Tape	Foot	All
T19557	Waste Bin Bracket (No Hardware)	Each	V1 & V2
PDRDTNT-1	Drawer Detents with Rivets (one pair)	Each	V1 & V2
PSMTOP-1	SM Gray Plastic Top	Each	All
PMEDTOP-1	MED Gray Plastic Top	Each	All
PLGTOP-1	LG Gray PlasticTop	Each	All
PXLTOP-1	XL Gray Plastic Top	Each	All
PSMWDTOP-1	SM Maui Grantie Top (predrilled with screws)	Each	All
PMEDWDTOP-1	MED Maui Grantie Top (predrilled with screws)	Each	All
PLGWDTOP-1	LG Maui Grantie Top (predrilled with screws)	Each	All
PXLWDTOP-1	XL Maui Grantie Top (predrilled with screws)	Each	All
PSLIDELS-1	V2 Complete Slide, Left (rivets & screws)	Each	V1 & V2
PSLIDERS-1	V2 Complete Slide, Right (rivets & screws included)	Each	V1 & V2
PSLIDEKIT-3	V3 Complete Auto-Close Slide, Left (rivets & screws included)	Each	V3
P-NARCDSP	Narcotics Box for Phoenix LTC Carts (DisPill Package)	Each	All
P-NARCBOX	Narcotics Box for Phoenix LTC Carts (single wide)	Each	All
T19706A3	Exchange Tray for X-Large cart	Each	All
T19706A4	Exchange Tray for Large cart	Each	All
P-WASTEBIN-CRM	Waste Bin with lid - Cream	Each	All
P-WASTEBIN-D	WASTE BIN WITH LID/BRACKET LTC - DK	Each	All
P-WASTEBIN-L	WASTE BIN WITH LID/BRACKET LTC - LT	Each	All
P-WASTEBIN-M	WASTE BIN WITH LID/BRACKET LTC -MED	Each	All
P-WASTEBIN-O	WASTE BIN WITH LID, LTC - OAK WOOD	Each	All
P-WASTEBIN-W	WASTE BIN WITH LID, LTC - WALNUT WOOD	Each	All
P-WASTEBIN-Y	WASTE BIN WITH LID, LTC - CHERRY WOOD	Each	All
M14045	WHITE-PLASTIC LOCK BUSHING, LTC CART	Each	V1 & V2
41-THB-8-10	Thumb Screws	Each	All
43-FLT-M8-10	Washer for Horizontal Best push bar	Each	V1 & V2
46-1070	E-CLIP 5/16 5133-31 (for horizontal Best push bar ONLY)	Each	V1 & V2
49-1475	Lock Bar Spring	Each	All
YM12991	Phoenix Bubble Label	Each	All
EMAR KIT	SCREWS & WASHERS FOR ERGOTRON ARM	Each	All

	Item #	Description	QTY	Version
Keys & Locks	40-2230	Best Core, 3 Keys, & Housing	Each	All
	40-2235	Best Housing (no core or keys)	Each	All
	40-2230S	Best Housing with Turned Down Pin (no core or keys)	Each	All
	CORE 2-2	Replacement Best Core with 3 Keys	Each	All
	KEY 2 CONTROL	Best Control Key for All Phoenix LTC Locks	Each	All
	KEY 2 MASTER	Best Master Key for All Phoenix LTC Locks	Each	All
	KEY 2-##	Best Central Lock Key, fill in 10-59 (Min Purchase 2)	Each	All
	CENTRAL LOCK KEY	Central Lock Key General # To Add Add'l Keys to Cart	Each	All
	PNARCLOCK	Narcotic Lock with Three Keys and four screws	Each	All
	NARC LOCK KEY	Narc Box Key General # To Add Add'l Keys to Cart	Each	All
	40-2216-1	Master Key for Narc Locks	Each	All
	PKENDALLKEY	Kendall Sharps Cabinet Key, all keyed same (min 2)	Each	All
	PSHARPKEY	Bemis Sharps Cabinet Key, all keyed same (min 2)	Each	All
Plastic Tops	PSMTOP-1	SM Gray Plastic Top	Each	All
	PMEDTOP-1	MED Gray Plastic Top	Each	All
	PLGTOP-1	LG Gray PlasticTop	Each	All
	PXLTOP-1	XL Gray Plastic Top	Each	All
Maui Granite	PSMWDTOP-1	SM Maui Grantie Top (predrilled with screws)	Each	All
	PMEDWDTOP-1	MED Maui Grantie Top (predrilled with screws)	Each	All
	PLGWDTOP-1	LG Maui Grantie Top (predrilled with screws)	Each	All
	PXLWDTOP-1	XL Maui Grantie Top (predrilled with screws)	Each	All
Slides & Related	40-1690	Slide Rubber Retention Insert	Each	V1 & V2
	PDRDTNT-1	Drawer Detents with Rivets (one pair)	Each	V1 & V2
	PSLIDELS-1	V2 Complete Slide, Left (rivets & screws)	Each	V1 & V2
	PSLIDERS-1	V2 Complete Slide, Right (rivets & screws included)	Each	V1 & V2
	PSLIDEKIT-3	V3 Complete Auto-Close Slide, Left (rivets & screws included)	Each	V3
Drawer Handles	M10191	Drawer Handle - Silver Knob (2x per drawer)	Each	All
	M10195	Drawer Handle - Silver - Small Cart (1x per drawer)	Each	All
	M10193	Drawer Handle - Silver - Medium Cart (1x per drawer)	Each	All
	M10197	Drawer Handle - Silver - Large Cart (1x per drawer)	Each	All
	M10196	Drawer Handle - Silver - Small side of Xlarge Cart (1x per drawer)	Each	All
	M10193	Drawer Handle - Silver - Wide side of Xlarge Cart (1x per drawer)	Each	All

Narc Parts	M14056-A	Narc Box Hinge - Standard, Global and Double (2)	Each	All
	M14056A3	Narc Box Hinge - Dispill	Each	All
	M14056-A	HINGE, NARC BOX W/SPRINGS, SCREW ON	Each	All
Bumpers	56-2080	Corner Bumpers, Gray (no hardware)	Each	All
	56-2082	Corner Bumpers, Brown (no hardware)	Each	All
	41-PHMS-10-75	Screws for Corner Bumpers	Each	All
Pull-Out Side-Shelves	A17096A6	PLTC, Pull Out Side Shelf (also order 46-1220)	Each	V1 & V2
	46-1220	SNAP BUTTON CLIP, .375 HD SOLID SS	Each	V1 & V2
	A17096A9-CRM	PLTC, Pull Out Side Shelf for MED, LRG, XLRG Cart (also order 46-1220, F-0216)	Each	V3
	A17096A9-DK	PLTC, Pull Out Side Shelf for MED, LRG, XLRG Cart (also order 46-1220, F-0216)	Each	V3
	A17096A9-MED	PLTC, Pull Out Side Shelf for MED, LRG, XLRG Cart (also order 46-1220, F-0216)	Each	V3
	A17096A9-WNT	PLTC, Pull Out Side Shelf for MED, LRG, XLRG Cart (also order 46-1220, F-0216)	Each	V3
	A17096A8-CRM	PLTC, Pull Out Side Shelf for SMALL Cart (also order 46-1220, F-0216)	Each	V3
	A17096A8-DK	PLTC, Pull Out Side Shelf for SMALL Cart (also order 46-1220, F-0216)	Each	V3
	A17096A8-MED	PLTC, Pull Out Side Shelf for SMALL Cart (also order 46-1220, F-0216)	Each	V3
	A17096A8-WNT	PLTC, Pull Out Side Shelf for SMALL Cart (also order 46-1220, F-0216)	Each	V3
Caster Plates	T19535A1	V2-Caster Plate (16 rivets & 2 Metal Supports need to be ordered)	Each	V1 & V2
	48-1062	Rivets for Caster Plate	Each	V1 & V2
	T19536	Metal Supports for V2 Caster Plate	Each	V1 & V2
	T19535A2	CASTER MOUNT PLATE, PXLTC CART (V3) - cart takes 4x	Each	V3
	M2163A1	HEX BOLT SERRATED 1/4-20 X 1/2 ZINC (for casters and plate)	Each	All
Odds and Ends	99-5979-3	Rivet Gun	Each	All
	63-1200	Double-Sided Foam Tape	Foot	All
	41-THB-8-10	Thumb Screws	Each	All
	YM12991	Phoenix Bubble Label	Each	All
	EMAR KIT	SCREWS & WASHERS FOR ERGOTRON ARM	Each	V1 & V2
	92-1000	Phoenix LTC Key Chain	Each	All
	28-1556	PLTC, Touch Up Paint, Cream Texture	Each	All
	WLOO-RTNSTRP	Strap Kit for Returning Carts (2 straps and 3 clips for returning a cart)	Each	All
	T19706A3	Exchange Tray for X-Large cart	Each	All
	T19706A4	Exchange Tray for Large cart	Each	All
	P-NARCBOXSTAY	Keeps lid down inside cart when lid not locked....	Each	All

Lock Parts	M14045	WHITE-PLASTIC LOCK BUSHING, LTC CART	Each	V1 & V2
	43-FLT-M8-10	Washer for Horizontal Best push bar	Each	V1 & V2
	46-1070	E-CLIP 5/16 5133-31 (for horizontal Best push bar ONLY)	Each	V1 & V2
	49-1475	Lock Bar Spring	Each	All
	46-1040	E-CLIP .25 (for lock pin on vertical lock bar ONLY)	Each	V1 & V2
	M14049	Lock Pin for Single Drawer Set (SM MED & LG Shell)	Each	V1 & V2
	M14049A6	LOCK PIN (V3 - RT OF DRWR)	Each	V1 & V2
	M14049A7	LOCK PIN (V3 - LEFT OF DRWR)	Each	V1 & V2
	M14049A8	LOCK PIN, AVK LEFT	Each	V1 & V2
	M14049A9	LOCK PIN (SPECIAL 4" DRAWER)	Each	V1 & V2
Waste Basket	T19557	Waste Bin Bracket (No Hardware)	Each	V1 & V2
	P-WASTEBIN-CRM	Waste Bin with lid - Cream	Each	All
	P-WASTEBIN-D	WASTE BIN WITH LID/BRACKET LTC - DK	Each	All
	P-WASTEBIN-L	WASTE BIN WITH LID/BRACKET LTC - LT	Each	All
	P-WASTEBIN-M	WASTE BIN WITH LID/BRACKET LTC - MED	Each	All
	P-WASTEBIN-O	WASTE BIN WITH LID, LTC - OAK WOOD	Each	All
	P-WASTEBIN-W	WASTE BIN WITH LID, LTC - WALNUT WOOD	Each	All
	P-WASTEBIN-Y	WASTE BIN WITH LID, LTC - CHERRY WOOD	Each	All
Drawer Install	44-AVK-ALS4-832-80	AVK for Drawers	Each	V1 & V2
	44-AVK-AA170	AVK FASTENER MANUAL INSTALL TOOL	Each	V1 & V2
	44-AVK-AA271-420	AVK Fastner Install Set 1/4-20	Each	V1 & V2
	44-AVK-AA271-832	AVK Fastner Install Set 8-32	Each	V1 & V2
Crème SM Drawers	TBD (V1,2 OR V3)	3" Drawer SM Shell	Each	Note Version
	TBD (V1,2 OR V3)	4" Drawer SM Shell	Each	Note Version
	TBD (V1,2 OR V3)	6" Drawer SM Shell	Each	Note Version
	TBD (V1,2 OR V3)	9" Drawer SM Shell	Each	Note Version
	TBD (V1,2 OR V3)	12" Drawer SM Shell	Each	Note Version

Crème MED Drawers	TBD (V1,2 OR V3)	3" Drawer MD Shell	Each	Note Version
	TBD (V1,2 OR V3)	4" Drawer MD Shell (XL LH SIDE)	Each	Note Version
	TBD (V1,2 OR V3)	6" Drawer MD Shell (XL LH SIDE)	Each	Note Version
	TBD (V1,2 OR V3)	9" Drawer MD Shell (XL LH SIDE)	Each	Note Version
	TBD (V1,2 OR V3)	12" Drawer MD Shell (XL LH SIDE)	Each	Note Version
Crème LG Drawers	TBD (V1,2 OR V3)	3" Drawer LG Shell	Each	Note Version
	TBD (V1,2 OR V3)	4" Drawer LG Shell	Each	Note Version
	TBD (V1,2 OR V3)	6" Drawer LG Shell	Each	Note Version
	TBD (V1,2 OR V3)	9" Drawer LG Shell	Each	Note Version
	TBD (V1,2 OR V3)	12" Drawer LG Shell	Each	Note Version
Crème XL Drawers	TBD (V1,2 OR V3)	3" Drawer RH XL Shell	Each	Note Version
	TBD (V1,2 OR V3)	4" Drawer RH XL Shell	Each	Note Version
	TBD (V1,2 OR V3)	6" Drawer RH XL Shell	Each	Note Version
	TBD (V1,2 OR V3)	9" Drawer RH XL Shell	Each	Note Version
	TBD (V1,2 OR V3)	9" Drawer RH XL Shell With Narc Lock and Internal Lid	Each	Note Version
	TBD (V1,2 OR V3)	12" Drawer RH XL Shell	Each	Note Version
Wood SM Drawers	TBD (V1,2 OR V3)	3" Drawer SM Shell - Wood	Each	Note Version
	TBD (V1,2 OR V3)	4" Drawer SM Shell - Wood	Each	Note Version
	TBD (V1,2 OR V3)	6" Drawer SM Shell - Wood	Each	Note Version
	TBD (V1,2 OR V3)	9" Drawer SM Shell - Wood	Each	Note Version
	TBD (V1,2 OR V3)	12" Drawer SM Shell - Wood	Each	Note Version
Wood MED Drawers	TBD (V1,2 OR V3)	3" Drawer MD Shell - Wood	Each	Note Version
	TBD (V1,2 OR V3)	4" Drawer MD Shell - Wood (XL LH SIDE)	Each	Note Version
	TBD (V1,2 OR V3)	6" Drawer MD Shell - Wood (XL LH SIDE)	Each	Note Version
	TBD (V1,2 OR V3)	9" Drawer MD Shell - Wood (XL LH SIDE)	Each	Note Version
	TBD (V1,2 OR V3)	12" Drawer MD Shell - Wood (XL LH SIDE)	Each	Note Version
Wood LG Drawers	TBD (V1,2 OR V3)	3" Drawer LG Shell - Wood	Each	Note Version
	TBD (V1,2 OR V3)	4" Drawer LG Shell - Wood	Each	Note Version
	TBD (V1,2 OR V3)	6" Drawer LG Shell - Wood	Each	Note Version
	TBD (V1,2 OR V3)	9" Drawer LG Shell - Wood	Each	Note Version
	TBD (V1,2 OR V3)	12" Drawer LG Shell - Wood	Each	Note Version
Wood XL Drawers	TBD (V1,2 OR V3)	3" Drawer RH XL Shell - Wood	Each	Note Version
	TBD (V1,2 OR V3)	4" Drawer RH XL Shell - Wood	Each	Note Version
	TBD (V1,2 OR V3)	6" Drawer RH XL Shell - Wood	Each	Note Version
	TBD (V1,2 OR V3)	9" Drawer RH XL Shell - Wood	Each	Note Version
	TBD (V1,2 OR V3)	9" Drawer RH XL Shell Narc Lock & Lid - Wood	Each	Note Version
	TBD (V1,2 OR V3)	12" Drawer RH XL Shell - Wood	Each	Note Version

V2 Crème Panels	T19660A2-CRM	V2 RH Panel MED & LG, Cream Texture	Each	V2
	T19660A3-CRM	V2 RH Panel XL, Cream Texture	Each	V2
	T19659-CRM	V2 LH Panel MED, LG & XL, Cream Texture	Each	V2
	T19658A1-CRM	V2 Back Panel MED Cream Texture	Each	V2
	T19658A2-CRM	V2 Back Panel LG Cream Texture	Each	V2
	T19658A3-CRM	V2 Back Panel XL Cream Texture	Each	V2
	Ti9534A5-CRM	V2 Column Cover	Each	V2
V2 Wood Panels	T19660A2-xxx	V2 RH Panel MED & LG, Wood (Specify color)	Each	V2
	T19660A3-xxx	V2 RH Panel XL, Wood (Specify Color)	Each	V2
	T19659-xxx	V2 LH Panel MED, LG & XL, Wood (Specify Color)	Each	V2
	T19658A1-xxx	V2 Back Panel MED, Wood (Specify Color)	Each	V2
	T19658A2-xxx	V2 Back Panel LG, Wood (Specify Color)	Each	V2
	T19658A3-xxx	V2 Back Panel XL, Wood (Specify Color)	Each	V2
	Ti9534A5-xxx	V2 Column Cover, Wood (Specify Color)	Each	V2
V3 Creme Panels	T19658A7-CRM	V3 Large Back Panel, Cream Texture	Each	V3
	T19659A5-CRM	V3 Small, Medium, Large, X-Large Left Side Panel, Cream Texture	Each	V3
	T19660A7-CRM	V3 Small, Medium, Large Right Side Panel, Cream Texture	Each	V3
	T19658A6-CRM	V3 Medium Back Panel, Cream Texture	Each	V3
	T19534A7-CRM	V3 X-Large Front Column Cover, Cream Texture	Each	V3
	T19658A8-CRM	V3 X-Large Back Panel, Cream Texture	Each	V3
	T19660A8-CRM	V3 X-Large Right Panel, Cream Texture	Each	V3
V3 Wood Panels	T19658A7-XXX	V3 Large Back Panel, Wood	Each	V3
	T19659A5-XXX	V3 Small, Medium, Large, X-Large Left Side Panel, Wood	Each	V3
	T19660A7-XXX	V3 Small, Medium, Large Right Side Panel, Wood	Each	V3
	T19658A6-XXX	V3 Medium Back Panel, Wood	Each	V3
	T19534A7-XXX	V3 X-Large Front Column Cover, Wood	Each	V3
	T19658A8-XXX	V3 X-Large Back Panel, Wood	Each	V3
	T19660A8-XXX	V3 X-Large Right Panel, Wood	Each	V3
	V3-LOCKFIX	Parts to fix original V3 locking issues	Each	V3

Replacement Parts Continued

Item #	Description	Cart Version	
40-1690	Slide Rubber Retention Insert	Each	V1 & V2
T19557	Trash Can Bracket (does NOT include thumb screws)	Each	V1 & V2
44-AVK-ALS4-832-80	AVK for Drawers	Each	V1 & V2
44-AVK-AA170	AVK FASTENER MANUAL INSTALL TOOL	Each	V1 & V2
44-AVK-AA271-420	AVK Fastner Install Set 1/4-20	Each	V1 & V2
44-AVK-AA271-832	AVK Fastner Install Set 8-32	Each	V1 & V2
46-1040	E-CLIP .25 (for lock pin on vertical lock bar ONLY)	Each	V1 & V2
M14049	Lock Pin for Single Drawer Set (SM MED & LG Shell)	Each	V1 & V2
M14049A6	LOCK PIN (V3 - RT OF DRWR)	Each	V1 & V2
M14049A7	LOCK PIN (V3 - LEFT OF DRWR)	Each	V1 & V2
M14049A8	LOCK PIN, AVK LEFT	Each	V1 & V2
M14049A9	LOCK PIN (SPECIAL 4" DRAWER)	Each	V1 & V2
40-2090	Drawer Catch (2 Per Drawer)	Each	V1 & V2
48-1030	Rivets for Drawer Catch (8 per Drawer)	Each	V1 & V2
M10191	Drawer Handle - Silver Knob (2x per drawer)	Each	All
M10195	Drawer Handle - Silver - Small Cart (1x per drawer)	Each	All
M10193	Drawer Handle - Silver - Medium Cart (1x per drawer)	Each	All
M10197	Drawer Handle - Silver - Large Cart (1x per drawer)	Each	All
M10196	Drawer Handle - Silver - Small side of Xlarge Cart (1x per drawer)	Each	All
M10193	Drawer Handle - Silver - Wide side of Xlarge Cart (1x per drawer)	Each	All
M14056-A	HINGE, NARC BOX W/SPRINGS, SCREW ON	Each	All
Y700209A3	DisPill Divider	Each	All
56-2080	Corner Bumpers, Gray (no hardware)	Each	All
56-2082	Corner Bumpers, Brown (no hardware)	Each	All
41-PHMS-10-75	Screws for Corner Bumpers	Each	All
92-1000	Phoenix LTC Key Chain	Each	All
46-1220	SNAP BUTTON CLIP, .375 HD SOLID SS	Each	V1 & V2
A17096A6	PLTC, Pull Out Side Shelf (also order 46-1220)	Each	V1 & V2
28-1556	PLTC, Touch Up Paint, Cream Texture	Each	All
A17096A9-CRM	PLTC, Pull Out Side Shelf for MED, LRG, XLRG Cart (also order 46-1220, F-0216)	Each	V3
A17096A9-DK	PLTC, Pull Out Side Shelf for MED, LRG, XLRG Cart (also order 46-1220, F-0216)	Each	V3
A17096A9-MED	PLTC, Pull Out Side Shelf for MED, LRG, XLRG Cart (also order 46-1220, F-0216)	Each	V3
A17096A9-WNT	PLTC, Pull Out Side Shelf for MED, LRG, XLRG Cart (also order 46-1220, F-0216)	Each	V3
A17096A8-CRM	PLTC, Pull Out Side Shelf for SMALL Cart (also order 46-1220, F-0216)	Each	V3
A17096A8-DK	PLTC, Pull Out Side Shelf for SMALL Cart (also order 46-1220, F-0216)	Each	V3
A17096A8-MED	PLTC, Pull Out Side Shelf for SMALL Cart (also order 46-1220, F-0216)	Each	V3
A17096A8-WNT	PLTC, Pull Out Side Shelf for SMALL Cart (also order 46-1220, F-0216)	Each	V3

Replacement Parts Continued

Item #	Description	Cart Version	
PLXMNT-1	ERGOTRON ARM MNT WLDMNT, LTC CART (A17072A3)	Each	V3
T19535A1	V2-Caster Plate (16 rivets & 2 Metal Supports need to be ordered)	Each	V1 & V2
48-1062	Rivets for Caster Plate	Each	V1 & V2
T19536	Metal Supports for V2 Caster Plate	Each	V1 & V2
T19535A2	CASTER MOUNT PLATE, PXLTC CART (V3) - cart takes 4x	Each	V3
M2163A1	HEX BOLT SERRATED 1/4-20 X 1/2 ZINC (for casters and plate)	Each	All
99-5979-3	Rivet Gun	Each	All
WLOO-RTNSTRP	Strap Kit for Returning Carts (2 straps and 3 clips for returning a cart)	Each	All
P-NARCBOXSTAY	Keeps lid down inside cart when lid not locked....	Each	All
	3" Drawer SM Shell	Each	
	4" Drawer SM Shell	Each	
	6" Drawer SM Shell	Each	
	9" Drawer SM Shell	Each	
	12" Drawer SM Shell	Each	
	3" Drawer MD Shell	Each	
	4" Drawer MD Shell (XL LH SIDE)	Each	
	6" Drawer MD Shell (XL LH SIDE)	Each	
	9" Drawer MD Shell (XL LH SIDE)	Each	
	12" Drawer MD Shell (XL LH SIDE)	Each	
	3" Drawer LG Shell	Each	
	4" Drawer LG Shell	Each	
	6" Drawer LG Shell	Each	
	9" Drawer LG Shell	Each	
	12" Drawer LG Shell	Each	
	3" Drawer RH XL Shell	Each	
	4" Drawer RH XL Shell	Each	
	6" Drawer RH XL Shell	Each	
	9" Drawer RH XL Shell	Each	
	9" Drawer RH XL Shell With Narc Lock and Internal Lid	Each	
	12" Drawer RH XL Shell	Each	
	3" Drawer SM Shell - Wood	Each	
	4" Drawer SM Shell - Wood	Each	
	6" Drawer SM Shell - Wood	Each	
	9" Drawer SM Shell - Wood	Each	
	12" Drawer SM Shell - Wood	Each	

Replacement Parts Continued

Finish	Item #	Description	Cart Version	
		3" Drawer MD Shell - Wood	Each	
		4" Drawer MD Shell - Wood (XL LH SIDE)	Each	
		6" Drawer MD Shell - Wood (XL LH SIDE)	Each	
		9" Drawer MD Shell - Wood (XL LH SIDE)	Each	
		12" Drawer MD Shell - Wood (XL LH SIDE)	Each	
		3" Drawer LG Shell - Wood	Each	
		4" Drawer LG Shell - Wood	Each	
		6" Drawer LG Shell - Wood	Each	
		9" Drawer LG Shell - Wood	Each	
		12" Drawer LG Shell - Wood	Each	
		3" Drawer RH XL Shell - Wood	Each	
		4" Drawer RH XL Shell - Wood	Each	
		6" Drawer RH XL Shell - Wood	Each	
		9" Drawer RH XL Shell - Wood	Each	
		9" Drawer RH XL Shell Narc Lock & Lid - Wood	Each	
		12" Drawer RH XL Shell - Wood	Each	
V2 Replacement skins				
Cream	T19660A2-CRM	V2 RH Panel MED & LG, Cream Texture	Each	V2
	T19660A3-CRM	V2 RH Panel XL, Cream Texture	Each	V2
	T19659-CRM	V2 LH Panel MED, LG & XL, Cream Texture	Each	V2
	T19658A1-CRM	V2 Back Panel MED Cream Texture	Each	V2
	T19658A2-CRM	V2 Back Panel LG Cream Texture	Each	V2
	T19658A3-CRM	V2 Back Panel XL Cream Texture	Each	V2
	Ti9534A5-CRM	V2 Column Cover	Each	V2
Wood	T19660A2-xxx	V2 RH Panel MED & LG, Wood (Specify color)	Each	V2
	T19660A3-xxx	V2 RH Panel XL, Wood (Specify Color)	Each	V2
	T19659-xxx	V2 LH Panel MED, LG & XL, Wood (Specify Color)	Each	V2
	T19658A1-xxx	V2 Back Panel MED, Wood (Specify Color)	Each	V2
	T19658A2-xxx	V2 Back Panel LG, Wood (Specify Color)	Each	V2
	T19658A3-xxx	V2 Back Panel XL, Wood (Specify Color)	Each	V2
	Ti9534A5-xxx	V2 Column Cover, Wood (Specify Color)	Each	V2

Replacement Parts Continued

Finish	Item #	Description	Cart Version	
V3 Replacement skins				
Cream	T19658A7-CRM	V3 Large Back Panel, Cream Texture	Each	V3
	T19659A5-CRM	V3 Small, Medium, Large, X-Large Left Side Panel, Cream Texture	Each	V3
	T19660A7-CRM	V3 Small, Medium, Large Right Side Panel, Cream Texture	Each	V3
	T19658A6-CRM	V3 Medium Back Panel, Cream Texture	Each	V3
	T19534A7-CRM	V3 X-Large Front Column Cover, Cream Texture	Each	V3
	T19658A8-CRM	V3 X-Large Back Panel, Cream Texture	Each	V3
	T19660A8-CRM	V3 X-Large Right Panel, Cream Texture	Each	V3
Wood	T19658A7-XXX	V3 Large Back Panel, Wood	Each	V3
	T19659A5-XXX	V3 Small, Medium, Large, X-Large Left Side Panel, Wood	Each	V3
	T19660A7-XXX	V3 Small, Medium, Large Right Side Panel, Wood	Each	V3
	T19658A6-XXX	V3 Medium Back Panel, Wood	Each	V3
	T19534A7-XXX	V3 X-Large Front Column Cover, Wood	Each	V3
	T19658A8-XXX	V3 X-Large Back Panel, Wood	Each	V3
	T19660A8-XXX	V3 X-Large Right Panel, Wood	Each	V3
	V3-LOCKFIX	Parts to fix original V3 locking issues	Each	V3

Misc. Information and Measurements

Drawer Height

Drawer	Height
3" Drawer	
4" Drawer	
6" Drawer	
9" Drawer	
9" with Narc Lock	
12" Drawer	

Wood Finish - Add

All Drawers include:

- Dividers
- Backers
- Slides
- Handles
- Lock Ramp
- Detents