Advance Table

User Manual

Product No. A-71101-EU

80028189
Version C

Enhancing outcomes for
patients and their caregivers.®

Hill-Rom

ENGLISH	7
中文简体	67
DANSK	133
NEDERLANDS.....	193
FRANÇAIS	255
DEUTSCH	318
ITALIANO	382
日本語版	445
LATVISKI	508
POLSKI	570
PORTUGUÊS	635
ROMÂNESC	698
SRPSKI	760
SLOVENSKY	823
SLOVENŠČINA.....	884
ESPAÑOL	945
SVENSKA	1008
TÜRK.....	1068

Manufactured by:

ALLEN MEDICAL SYSTEMS, INC.
100 DISCOVERY WAY
ACTON, MA 01720
800-433-5774 (NORTH AMERICA)
+1-978-266-4200 (INTERNATIONAL)

Authorized European Union Representative:

HILL-ROM SAS
B.P. 14 - Z.I. DU TALHOUE
56330 PLUVIGNER
FRANCE
TEL.: +33 (0)2 97 50 92 12

No part of this text shall be reproduced or transmitted in any form or by any means, electronic or mechanical, including photocopying, recording, or by any information or retrieval system without written permission from Allen Medical Systems, Inc. (Allen Medical).

The information in this manual is confidential and may not be disclosed to third parties without the prior written consent of Allen Medical.

This manual (D-720722-A1) was originally released and supplied in English. For a list of available translations, contact Allen Medical Technical Support.

First Edition, December 2014

First Printing, 2014

Product may be covered by one or more patents. Please consult listing at www.hill-rom.com/patents for any patent(s).

Mizuhosi® is a registered trademark of Mizuho Orthopedic Systems, Inc.

Velcro® is a registered trademark of Velcro Industries, BV.

Allen®, Bow®, C-Flex®, Easy Lock® are registered trademarks of Allen Medical Systems, Inc.

Ultra Comfort™ is a trademarks of Allen Medical Systems, Inc.

Quick View™ is a trademark of Hill-Rom Services, Inc.

The information contained in this manual is subject to change without notice. Allen Medical makes no commitment to update or keep current, the information contained in this manual.

Allen Medical reserves the right to make changes without notice in design, specifications, and models. The only warranty Allen Medical makes is the express written warranty extended on the sale or rental of its products.

NOTE:

Replace this manual (D-720722-A1) if it is damaged and/or can not be read.

To order additional copies of this manual (D-720722-A1), contact Allen Medical.

This manual was created from template D-FACT-7201.

Index

INTENDED USE	7
Model Identification	8
Document Symbols	8
Product Symbol	8
Safety warnings and cautions	14
Quick View™ list of features	17
View 1	17
View 2	18
Controls and indicators	19
Instructions for use	21
Turn the unit on and off	21
Battery charge level	22
Potential Equalization Conductor	22
Lock and Unlock the Brakes	22
Manually Release the Brakes	23
Reset the Electronic Brakes after being Manually Released	23
Adjust the length of the Allen® Advance Table	23
Install a support top	24
Remove a support top and/or bracket	27
Install the Prone Support Set (A-71300) on the rail top	28
Pre-use check	31
Adjust the height of the Allen® Advance Table	32
Put the Allen® Advance Table in a Trendelenburg position	33
Tilt the support top	34
Level the table and the support top	35
Rotate the support top 180°	36
Transfer a patient	41
Access the cervical traction assembly	42
Shut down the Allen® Advance Table	44
Clean and disinfect	45
Clean	45

Clean difficult spots	46
Steam clean	46
Disinfect	46
Troubleshooting	46
The Allen® Advance Table will not go into pendant control mode from the manual flip mode.....	46
The pendant control and manual flip indicators are flashing	47
The brakes locked and unlocked indicators are flashing	47
The Allen® Advance Table “beeps” when a pendant control is pressed	47
Service calls	47
Preventive maintenance.....	47
Maintenance checks	48
Storage and transport.....	48
Expected life	49
Discard the unit	49
Parts replacement	49
Pendant.....	50
Fuses	51
Accessories.....	51
Specifications	54

INTENDED USE

The intended users of the Allen® Advance Table are healthcare employees who have been trained to use the product, and who have the physical strength and cognitive skills to operate and control the product. Follow facility safety protocols if an intended user does not have the physical strength or cognitive skills to operate and control the product safely. The Allen® Advance Table is intended to be used in an operating room environment to support and position patients during surgical procedures. Several types of specialty tops can be attached to the Allen® Advance Table to permit surgical access and 360° radiolucency for various types of surgical operations. The patient can be secured into various operable positions from the supine, prone, or lateral position. Electronic controls on the Allen® Advance Table pendant permit the adjustment of table height, angle, tilt, and auxiliary function for the lateral top. The Allen® Advance Table contains electronic locking casters that permit the Allen® Advance Table to be moved and locked into position by use of an electronic control on the Allen® Advance Table pendant.

This manual provides the information required for normal operation of the Allen® Advance Table. Before you operate the Allen® Advance Table, be sure that you have read and understood in detail the contents of this manual. It is important that you read and strictly obey the aspects of safety contained in this manual. Any reference to a side of the Allen® Advance Table is from the patient's view lying on an installed patient support top on his or her back.

NOTE:

When we discuss the installation of the Prone Support Set, any position reference is with the patient lying face down on the patient support top.

This manual includes the different features and options for the Allen® Advance Table, some of which may not apply to your unit. To identify the model that you have, look at the model number label that is on the inside of the head-end base, near the power outlet.

Model Identification

Document Symbols

This manual contains different typefaces and symbols to make the content easier to read and understand:

- Standard text—used for regular data.
- **Boldface text**—emphasizes a word or phrase.
- **NOTE:** —sets apart special data or important instruction clarification.
- WARNING or CAUTION




- A WARNING identifies situations or actions that may have an effect on patient or user safety. To ignore a warning could cause patient fall, patient injury, or user injury.
- A CAUTION identifies special procedures or precautions that persons must obey to help prevent equipment damage.

Product Symbol

Symbol	Description
Table Frame Symbols	
	Conforms to European Medical Device Directive 93/42/EEC. ^a
	Medical Electrical Equipment Evaluated in accordance with UL60601-1, EN 60601-1, IEC 60601-1, IEC 60601-2-46, and IEC 60601-1-2.
	Type B applied part.
IPX 4	Degree of protection against the ingress of water.

Symbol	Description
	Manufacturer.
	Equipotential.
	WARNING symbol in accordance with IEC 60601-1. (This symbol is black with a yellow background. The yellow background indicates a warning).
	Follow the Operating Instructions symbol in accordance with IEC 60601-1. (This symbol is black with a white background.)
	Failure to Follow the Operating Instructions could put the patient or operator at risk, symbol in accordance with IEC 60601-1. (This symbol is white with a blue background. The blue background indicates a mandatory action.)
	Shows the maximum patient weight for the Allen® Advance Table.
	The Allen® Advance Table is not suitable for use in oxygenenriched environments or where there is a flammable anesthetics mixture with air, oxygen, or nitrous oxide.
	Battery charge level indicator.
	Power indicator.

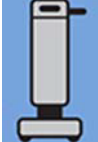Symbol	Description
	Identifies the emergency brake releases.
	Shows that when the flip handle is in the lowered position, the Allen® Advance Table is in flip mode (the top can be turned manually).
	Shows that when the flip handle is in the raised position, the Allen® Advance Table can be operated by the pendant.
	Identifies the location of the cervical traction assembly.
 Pendant symbols	
	Identifies the controls in the area as tilt controls. This is not a control.

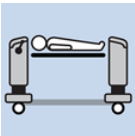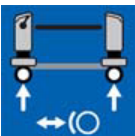Symbol	Description
	Tilt left control—turns a support top up to 25° counterclockwise.
	Tilt right control—turns a support top up to 25° clockwise.
	Identifies the controls in the area as column height controls. This is not a control.
	Up control—raises both columns at the same time.
	Down control—lowers both columns at the same time.
	Trendelenburg control—lowers the head-end column as it raises the foot-end column.
	Reverse Trendelenburg control—raises the head-end column as it lowers the foot-end column.
	Identifies the controls in the area as brake controls. This is not a control.

Symbol	Description
	Lock control—locks the brakes.
	Unlock control—when pressed along with the key control, unlocks the brakes.
	Key control—when pressed along with the unlock control, unlocks the brakes.
 Display symbols	
	Tilt indicator area—shows how far the support top is tilted to the left or right.
	Tilt level indicator—when the Tilt level indicator is green , the patient is not tilted to the left or right.
	Brake status indicators—the Locked indicator is green when the brakes are locked; the Unlocked indicator is orange when the brakes are unlocked.
	Pendant control indicator—is green when the Allen® Advance Table is in Pendant mode (the Allen® Advance Table can be operated by the pendant controls).

Symbol	Description
	Manual flip indicator—is orange when the Allen® Advance Table is in flip mode (the top can be turned manually).
Prone Support Set^b	
	Identifies the chest pad.
	Identifies the hip pads.
	Identifies the right thigh pad.
	Identifies the left thigh pad.

A. The CE mark was first applied in 2013.

b. The labels show the patient lying face down on the supports.

Safety warnings and cautions

WARNING:

Obey these safety instructions to help prevent patient fall, injury, and/or equipment damage:

- Read and understand all warnings in this manual and on the unit itself prior to use with a patient.
- This device should be operated by trained persons only.
- The procedures in this manual are only manufacturer's suggestions. The final responsibility for patient care with respect to this device remains with the attending physician.
- No modification of this equipment is allowed.
- Support tops and accessories from Allen Medical are designed to industry standard dimensions. Before you use products that are **not** from Allen Medical, you must verify with the product manufacturer that the product can be used safely with the Allen® Advance Table.
- Before you use Allen™ Patient Support Tops with a table not manufactured by Allen Medical, you must verify with the product manufacturer that the support tops can be used safely with their product.
- There are no known contraindications for use of the Allen® Advance Table.
- Except for patient transfer, always use safety straps to help secure the patient to the support top.
- Never remove a cross pin from a support top while a patient is on the support top.
- Before each use of the unit, make sure that all parts are correctly installed and that the unit operates correctly.
- If the unit has been dropped or damaged, before you use the unit, make sure facility maintenance has examined the unit completely and determined that it is safe to operate.
- To avoid the risk of electric shock, the unit must only be connected to a supply mains with protective earth.
- The potential for electrical shock exists with electrical equipment. Failure to follow facility protocols may cause death or serious injury.


- If the Allen® Advance Table is to be used near high-frequency surgical equipment, cardiac defibrillators, and cardiac defibrillator-monitors, refer to the product manufacturer's instructions to make sure the equipment will operate as intended.
- Incorrect use or handling of the power cord may cause damage to the power cord. If damage has occurred to the power cord or any of its components, immediately remove the Allen® Advance Table from service, and contact the applicable maintenance persons.
- Do not expose the unit to excessive moisture.
- Fluid spills on to the electronics can cause a hazard. If such a spill occurs, unplug the unit and remove it from service. When fluid spills occur outside that seen in normal use, immediately do as follows:
 - a. Unplug the unit from its power source.
 - b. Remove the patient from the unit.
 - c. Clean the fluid spill from the unit.
 - d. Have maintenance examine the unit completely.
 - e. Do not put the unit back into service until it is completely dry, tested, and determined to be safe to operate.
- Be aware of these pinch points on the Allen® Advance Table:
 - Between the caster and base skins
 - Between an H-bracket and the column
 - Between the head of the cross pin and the end cap during an extreme tilt and Trendelenburg, or Reverse Trendelenburg position.
 - Between the cross pin drop lock and the bowtie.
 - Between an H-bracket and the column during the removal of a cross pin.
- The unit is not intended for patient transport. However, if you must use the unit to transport a patient, make sure that the unit is in its lowest position and the patient is in a level position (not tilted).
- Before you transport the unit, make sure of these:
 - The unit is in its lowest position.
 - The power cord is disconnected from the facility power outlet and stored on the unit.
 - Accessories, supports, and any equipment that is to be transported with the unit are within the perimeter of the unit and secure on the unit.

Safety warnings and cautions

- During transport, make sure to do these:
 - Use caution so that the unit does not tip or overbalance.
 - For inclines and thresholds, approach them as you move forward or backward, rather than sideways.
 - Push the unit safely; do not make sharp corners or turn the unit at high speeds.
- Do not use tape on the Allen® Advance Table. Adhesive residue could affect the functions of the Allen® Advance Table or the support tops.
- However, you may use tape on the lateral top. Use an alcohol wipe to remove adhesive residue.


Quick View™ list of features

View 1

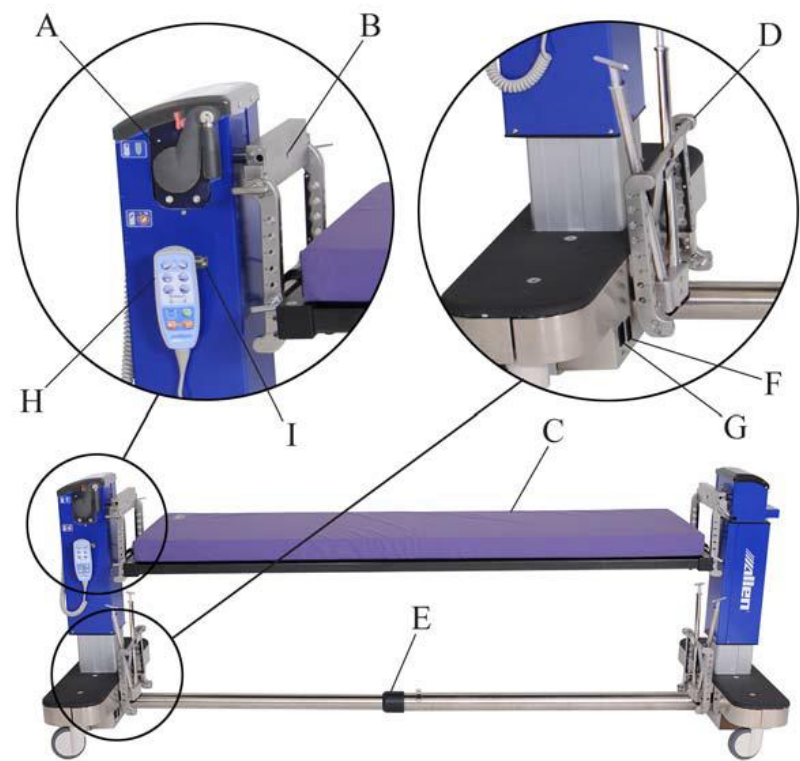Item	Feature	Page	Item	Feature	Page
A	Flip handle	36	F	On/Off switch	21
B	Cervical traction rope insertion	42	G	Power cord outlet	21
C	Patient support top with pad	24	H	Pendant	20
D	Cross pin and H-bracket storage	28	I	Pendant holder	44
E	Collapsible cross tube	23			

View 2

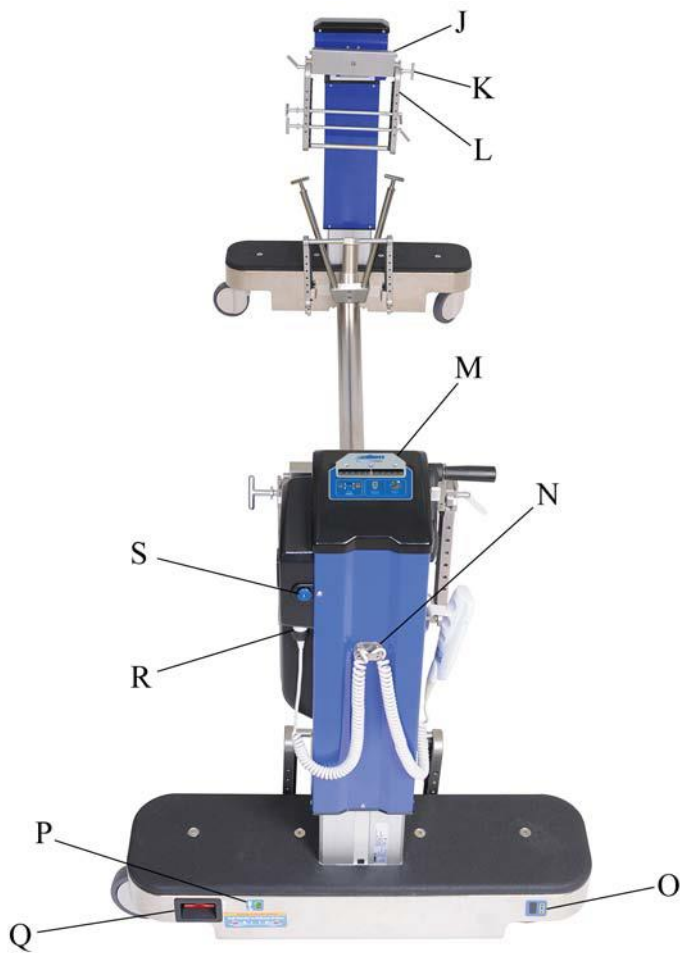Item	Feature	Page	Item	Feature	Page
J	Bowtie holder	24	O	Battery charge level	22
K	Cross pin	24	P	Power indicator	21
L	H-bracket	24	Q	Emergency brake release	23
M	Display	20	R	Pendant control connector	50
N	Pendant cord holder	44	S	Cervical traction release	42

Controls and indicators

You will find the different controls and indicators for the Allen® Advance Table on the cross tube, on both columns, and on the pendant.

You will find these controls and indicators on the head-end base:

- Model number label
- Serial number label
- Plug outlet
- On/Off switch
- Equipotential ground point

- Manual Brake release (red handle on both columns)
- Power indicator—when **on**, the Allen® Advance Table is plugged in.
- Battery charge level

- Manual flip handle—when released, the flip feature can be used.

- Cervical traction release—when pressed, the traction assembly moves out from the column.

Controls and indicators

You will find these controls on the pendant:

- Left/Right tilt—turns a support top approximately 25° clockwise or counterclockwise.
- Height Adjust—raises or lowers the table columns. The height can be adjusted from 67.3 cm to 127.0 cm (26.5" to 50"; measured from the support top to the floor).
- Trendelenburg / Reverse Trendelenburg—raises one column in relation to the opposite column.
- Electronic brake lock/unlock—locks and unlocks the brakes.

You will find these important information indicators on the display:

- Tilt level indicator (when the indicator is green, the patient is not tilted to the left or right)
- Tilt indicator and tilt indicator area
- Brakes locked and unlocked indicators
- Pendant Control indicator
- Manual flip indicator

You will find the retract/extend handle at the center of the cross tube.

Lift the retract/extend handle to adjust the overall length of the table from 175.3 cm to 276.9 cm (69" to 109"). Lower the handle to lock the cross tube in position.


Instructions for use

Turn the unit on and off

1. Make sure the flip handle is in the **up** and in **locked** position.

2. Plug the power cord into the table outlet and into an applicable power outlet.

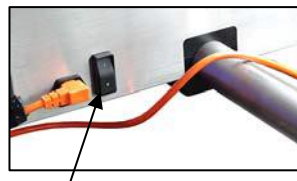On/Off switch

NOTE:

When the Table is plugged in, the Power indicator will be on.

3. Press the **On/Off** switch to the applicable position:

- 1 = On
- O = Off

NOTE:

When the On/Off switch is **on** and the table is **not** plugged into a power outlet, the table operates on battery power. Battery power is for short term use only and should be used for emergency power situations only.

Battery charge level

⚠ WARNING:

If the batteries are not charged before you use the Allen® Advance Table, the Allen® Advance Table may not function in an emergency situation. Patient fall or injury could occur.

When the Battery Charge Level indicator shows a **low** level of charge or if you hear a “beep” from the Allen® Advance Table when you try to use the pendant controls, plug the Allen® Advance Table into an applicable power outlet to charge the batteries. To keep the batteries charged, plug the Allen® Advance Table into an outlet nightly and when you store the Allen® Advance Table.

NOTE:

A “beep” from the Allen® Advance Table lets you know that the battery charge level is low or depleted.

NOTE:

When the Allen® Advance Table is **not** plugged in: the Battery Charge Level indicator will be **on** when the On/Off switch is **on**; the Battery Charge Level indicator will be **off** when the On/Off switch is **off**.

Potential Equalization Conductor

There is a plug on the Allen® Advance Table for you to connect an equipotential bonding cable. This will equalize any potential difference when the Allen® Advance Table is to be used with a patient who has intravenous, cardiac, or other medical equipment.


Lock and Unlock the Brakes


⚠ CAUTION:

Do not press the lock control while the Allen® Advance Table is being moved. Equipment damage could occur.

- **Lock**—press and hold the **lock** control until the **locked** indicator on the display is **green**.
- **Unlock**—press and hold the **key** control as you press the **unlock** control until the **unlocked** indicator on the display is **orange**.

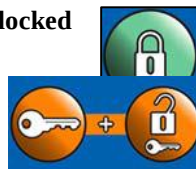When the brakes are locked, the **locked** indicator on the display is green; when the brakes are not locked, the **unlocked** indicator is **orange**.

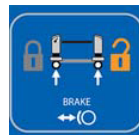During a loss of power, if the brakes are locked and it is necessary to move the Allen® Advance Table, use the emergency brake release on both columns to manually unlock the brakes.

Manually Release the Brakes

Pull the emergency brake release on **both** columns.

NOTE:

Each column has a release. You must pull both releases to unlock the brakes. The handles do not have to be pulled at the same time.


Reset the Electronic Brakes after being Manually Released

After you use the emergency brake releases to unlock the brakes, you must do these steps to lock or reset the brakes:

1. Make sure there is power to the Allen® Advance Table, and its On/Off switch is **on**.

2. Press and hold the **key** control as you press the **unlock** control. You will hear the brakes move.

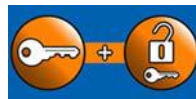3. Press and hold the **lock** control until the **locked** indicator on the Allen® Advance Table display is **green**.

4. Try to move the Allen® Advance Table to make sure the brakes are locked. If the brakes did not lock, repeat 2 and 3.

Adjust the length of the Allen® Advance Table

1. If there is a support top on the Allen® Advance Table, remove the top (see page 27).

2. **Unlock** the brakes.

3. Lift the cross tube handle to unlock the cross tube.

4. With a person at each end of the Allen® Advance Table, adjust the length as necessary:

Instructions for use

- **Retract**—push the columns toward each other.
 - **Extend**—pull the columns away from each other.
5. When the cross tube is fully extended, lower the handle to lock the cross tube into position.

NOTE:

To lock the cross tube, it must be fully extended.

Install a support top


⚠ WARNING:

Only use Allen support tops with the Allen® Advance Table, the Mizuhosi® Flat Top, Rail Top or Orthopedic Top. Patient fall, injury, or equipment damage could occur.

⚠ WARNING:

Use patient support tops that are designed for the simultaneous use of tilt and Trendelenburg or Reverse Trendelenburg only. Otherwise, patient fall, injury, or equipment damage could occur.

1. Make sure the patient support top is designed for the simultaneous use of tilt and Trendelenburg or Reverse Trendelenburg.
2. Make sure the brakes on the Allen® Advance Table are locked (the **locked** indicator on the Allen® Advance Table display is **green**).

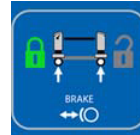3. If the lower H-brackets are not installed on each end of the Table, proceed as shown below. Otherwise, go to Step 4.


- a. Hold the H-bracket so that the open end is down, the J ends are toward you, and the bracket is angled toward you.

- b. Install the J ends of the H-bracket on the bowtie so that the H-bracket fits into the notches on the sides of the bowtie.

- c. Insert the cross pin into the H-bracket, through the bowtie, and through the other end of the H-bracket.

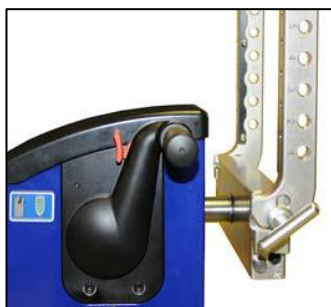⚠ WARNING:

Do not use the H-bracket without the cross pin installed correctly. To do so could cause patient fall, injury, or equipment damage.


- d. Make sure the drop lock on the end of the cross pin moves to the side to hold the H-bracket in position.

- e. Repeat these steps to install an H-bracket on the other column.

Instructions for use

- f. Lift the end (A) of the flip handle until it locks in position.
- g. Push and hold the red safety lever (B) forward as you pull the flip handle to the **released** or down position.

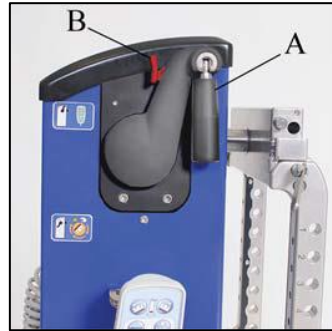- h. With the help of another person, turn the H-bracket on **each** column so that it is toward the bottom of the column.

- i. Push the flip handle up to the locked position. You do not need to press the red safety lever.

- 4. Determine which hole location on the H-bracket to install the support top.

NOTE:

The hole locations are numbered on the side of the H-brackets.

⚠ WARNING:

Failure to use correct lifting methods could cause patient fall, injury, or equipment damage.

5. With the help of another person, do as follows:
 - a. Align the ends of the support top with the pre-determined hole location. If you are to install a flat top without a pad installed, make sure the side with the Velcro® tape is up (patient side).


NOTE:

You may have to adjust the length of the Allen® Advance Table so that the columns are at the correct distance apart to fit the support top (see “Adjust the length of the Allen® Advance Table” on page 23).

⚠ WARNING:

Never remove a cross pin from a support top while a patient is on the support top. Patient fall or injury could occur.

- b. At each end of the support top, insert a cross pin into the H-bracket, through the support top connector assembly, and through the other end of the H-bracket.

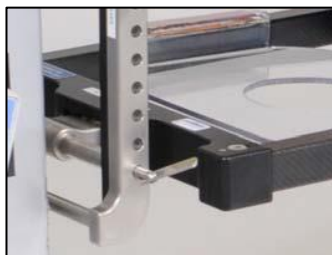**⚠ WARNING:**

Do not use the support top without the cross pins installed correctly through the H-brackets. To do so could cause patient fall, injury, or equipment damage.

- c. Make sure the drop lock on the end of the cross pin moves to the side to hold the support top in position.

Remove a support top and/or bracket**⚠ WARNING:**

Never remove a cross pin from a support top while a patient is on the support top. Patient fall or injury could occur.

1. Make sure the patient has been moved from the support top.
2. Do the applicable steps of the install procedure in reverse order.

NOTE:

To remove a cross pin, align the drop lock with the pin, and hold the lock in as you remove the cross pin from the H-bracket.

NOTE:

When the H-brackets and cross pins are not in use, they can be stored on the inside of the columns.

Install the Prone Support Set (A-71300) on the rail top

The Prone Support Set includes one chest support, six hip supports (two each of small, medium, and large), two thigh supports, one sling, and one head plate.

Sling

⚠ WARNING:

The sling is intended to provide a support for the patient's legs. We recommend that you use additional pressure reducing material between the patient's legs and the sling to help prevent patient injury.

1. Lay the closed part of the sling over the support top rails.
2. Connect the buckles.
3. Adjust the straps as necessary for the amount of slack you want in the sling.

Hip and Thigh Supports

NOTE:

The thigh supports are side specific.

1. Determine which size hip support to use, and then do as follows to install each support on the support top:

- a. Put the support in position on the top. Make sure the clamp side of the support is toward the outside of the top.
 - **Thigh support**—make sure the support will be in a position to support the patient's upper thigh area.
 - **Hip support**—make sure the support will be in a position that is in between the patient's iliac crest and the top of the pelvic bone.
- b. Close the blue clamp on to the locking posts.

- c. Slide the blue clamp to the locked position. When locked, only the green area will show on the pin.

NOTE:

Look under the rail to make sure the clamp is correctly engaged on the locking posts.

2. Install Allen® Ultra Comfort™ Covers on to the supports. Make sure that—

- The foam side of the cover is on the outside.

Instructions for use

- The covers are pulled over the blue section of the support only. The black section should not be covered.

NOTE:

Allen® Support Pads are designed to work with Allen® Ultra Comfort™ Covers only.

3. To remove a support from the Allen® Advance Table, slide the blue clamp to the unlocked position. When unlocked, only the red area will show on the pin.

Chest Support

1. Put the chest support in position on the rail top at a location dependent on the patient's gender:
 - Female—install the support so that it will be above the breasts and below the clavicle paying close attention to avoid the brachial plexus.
 - Male—install the pad so that it will be across the nipple area.

NOTE:

When the support is installed, the indent should be toward the foot end of the Allen® Advance Table.


2. Attach the Velcro® straps to hold the support in position.
3. Install an Allen® Ultra Comfort™ Cover on to the support. Make sure that—
 - The foam side of the cover is on the outside.
 - The cover completely covers the blue material of the support.

NOTE:

Allen® Support Pads are designed to work with Allen® Ultra Comfort™ Covers only.

Head Plate

1. Put the head plate in position on the rail top at a location to support the patient's head.
2. Attach the Velcro® straps to hold the head plate in position.
3. Install the foam head block on to the head plate.

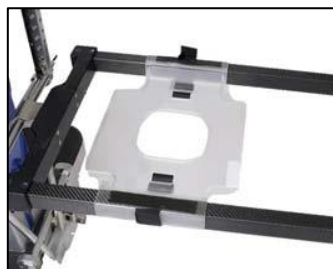NOTE:

The Allen® Head Plate is designed to work with the Allen® Foam Head Block only.

Pre-use check

Before each use of the Allen® Advance Table, make sure all necessary accessories are attached and the Allen® Advance Table is ready for use. If any of the controls or indicators do not operate correctly, do **not** use the Allen® Advance Table. Contact Allen Medical or your local distributor.

Condition	Associated Indicator(s)
The Allen® Advance Table power cord is connected to the Allen® Advance Table and to an applicable power outlet.	The power indicator is on .
The On/Off switch is on .	The battery charge level indicator is on .
The brakes are locked .	The locked indicator is green ; the unlocked indicator is off .
The Allen® Advance Table is in pendant mode.	The pendant control indicator is green ; the manual flip indicator is off .
The columns are level (at the same height).	Near the bottom on the outside of the columns, there are numbers to help identify the height level.
The support top is horizontal (not tilted).	The tilt level indicator is green ; the tilt indicator is in the center of the tilt indicator area.

Adjust the height of the Allen® Advance Table

⚠ WARNING:

Before you adjust the Allen® Advance Table or use the tilt feature, make sure that the brakes are locked, the patient is secured to the support top and is fully prepared for the movement, and the Allen® Advance Table is free of objects that could interfere with the adjustment. Otherwise, patient fall, injury, or equipment damage could occur.

1. Make sure that—

- The brakes are locked (the **locked** indicator on the Allen® Advance Table display is green).

- The patient is secured to the support top and is fully prepared for the adjustment.

NOTE:

We recommend that you use the two large security straps to secure the patient to the support top.

- The Allen® Advance Table is free of objects that could interfere with the adjustment.

⚠ WARNING:

Observe lines closely during articulations. Always use good line management techniques. Failure to do so could cause patient fall, injury, or equipment damage.

2. Press and hold the applicable raise or lower control until the Allen® Advance Table is at the correct height.

NOTE:

Near the bottom, on the outside of the columns, there are numbers to help identify the height level.

Put the Allen® Advance Table in a Trendelenburg position

⚠ WARNING:

Before you adjust the Allen® Advance Table or use the tilt feature, make sure that the brakes are locked, the patient is secured to the support top and is fully prepared for the adjustment, and the Allen® Advance Table is free of objects that could interfere with the adjustment. Otherwise, patient fall, injury, or equipment damage could occur.

1. Make sure that—

- The brakes are locked (the **locked** indicator on the Allen® Advance Table display is green).

- The patient is secured to the support top and is fully prepared for the adjustment.

NOTE:

We recommend that you use the two large security straps to secure the patient to the support top.

- The Allen® Advance Table is free of objects that could interfere with the adjustment.

⚠ WARNING:

Observe lines closely during articulations. Always use good line management techniques. Failure to do so could cause patient fall, injury, or equipment damage.

2. Press and hold the applicable **Trendelenburg** or **Reverse Trendelenburg** control until the support top is at the correct angle.

3. To return to the horizontal position, press and hold the opposite Trendelenburg control until the lowered column is at the same height as the raised column.

NOTE:

Near the bottom, on the outside of the columns, there are numbers to help identify the height level.

Tilt the support top

NOTE:


For bariatric patients, we recommend that you install three supine access tops (A-70404) to help support the patient during the tilt procedure.

⚠ WARNING:

Before you adjust the Allen® Advance Table or use the tilt feature, make sure that the brakes are locked, the patient is secured to the support top and is fully prepared for the movement, and the Allen® Advance Table is free of objects that could interfere with the adjustment. Otherwise, patient fall, injury, or equipment damage could occur.

1. Make sure that—

- The brakes are locked (the **locked** indicator on the Allen® Advance Table display is green).

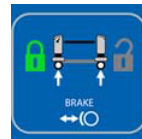- The patient is secured to the support top and is fully prepared for the adjustment.

NOTE:

We recommend that you use the two large security straps to secure the patient to the support top.

- The Allen® Advance Table is free of objects that could interfere with the adjustment.

⚠ WARNING:

Observe lines closely during articulations. Always use good line management techniques. Failure to do so could cause patient fall, injury, or equipment damage.

- #### 2. Press and hold the applicable **tilt control** until the support top is at the correct angle. When you release the control, the support top will lock at the selected angle.

The tilt indicator on the display will show how far the support top is tilted to the left or right.

- #### 3. To return the support top to the flat position, press and hold the opposite **tilt control** until the top is no longer tilted. The tilt indicator on the display will be in the center of the tilt indicator area and the tilt level indicator will be green.

Level the table and the support top

To determine if the Allen® Advance Table is level, look at the display:

- The **tilt level indicator** will be green.
- The **tilt indicator** will be in the center of the tilt indicator area.


⚠ WARNING:

Before you adjust the Allen® Advance Table or use the tilt feature, make sure that the brakes are locked, the patient is secured to the support top and is fully prepared for the movement, and the Allen® Advance Table is free of objects that could interfere with the adjustment. Otherwise, patient fall, injury, or equipment damage could occur.

1. Before you adjust the Allen® Advance Table, make sure that—

- The brakes are locked (the **locked** indicator on the Allen® Advance Table display is green).
- The patient is secured to the support top and is fully prepared for the adjustment.

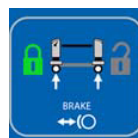NOTE:

We recommend that you use the two large security straps to secure the patient to the support top.

- The Allen® Advance Table is free of objects that could interfere with the adjustment.
2. If the Allen® Advance Table is in a Trendelenburg position, press and hold the opposite **Trendelenburg** control until the lowered column is at the same height as the raised column and the **tilt level indicator** on the display is **green**.

NOTE:

Near the bottom, on the outside of the columns, there are numbers to help identify the height level.

3. If the support top is tilted, press and hold the opposite **tilt control** until the top is no longer tilted. The **tilt indicator** on the display should be in the center of the tilt indicator area, and the **tilt level indicator** should be **green**.


Rotate the support top 180°

⚠ WARNING:

Before you adjust the Allen® Advance Table or use the flip or tilt features, make sure that the brakes are locked, the patient is secured against the prone support set, the patient is fully prepared for the adjustment; the Allen® Advance Table is free of objects that could interfere with the adjustment, and all loose or unsupported objects have been secured. Otherwise, patient fall, injury, or equipment damage could occur.

1. Make sure that—

- The brakes are locked (the **locked** indicator on the Allen® Advance Table display is **green**).
- The patient is fully prepared for the adjustment.
- The Allen® Advance Table is in the correct position for you to install a second support top:
 - The columns are at the same height.
 - The Allen® Advance Table and the installed support top are level. See “Level the table and the support top” on page 35.
 - The Allen® Advance Table is free of objects that could interfere with the adjustment.
 - All loose or unsupported objects have been secured.
- The rail top has the prone supports installed. If it does not, install them (see page 28).

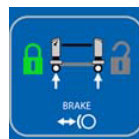2. Install the second support top as follows:

- a. At one end of the Table, put the upper H-bracket in position on the bowtie. Hold the H-bracket so that the open end is down, the J ends are toward you, and the H-bracket is angled toward you.
- b. Install the J ends of the H-bracket on the bowtie so that the H-bracket fits into the notches on the sides of the bowtie.

- c. Insert the cross pin into the H-bracket, through the bowtie, and through the other end of the H-bracket.

⚠ WARNING:

Do not use the H-bracket without the cross pin installed correctly. To do so could cause patient fall, injury, or equipment damage.

- d. Make sure the drop lock on the end of the cross pin moves to the side to hold the H-bracket in position.
- e. Repeat steps Step a through Step d on the opposite bowtie.

- f. Determine which hole location on the H-bracket to install the support top.

NOTE:

The hole locations are numbered on the side of the H-brackets.

⚠ WARNING:


Failure to use correct lifting methods could cause patient fall, injury, or equipment damage.

- g. With the help of another person, align the ends of the support top with the pre-determined hole location.

⚠ WARNING:

Never remove a cross pin from a support top while a patient is on the support top. Patient fall or injury could occur.

- h. At each end of the support top, insert a cross pin into the H-bracket, through the support top connector assembly, and through the other end of the H-bracket.

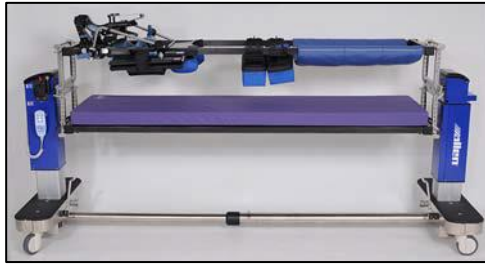⚠ WARNING:

Do not use the support top without the cross pins installed correctly through the H-brackets. To do so could cause patient fall, injury, or equipment damage.

- i. **Make sure that—**
- All cross pins are correctly installed into the H-brackets and support tops.
 - The drop lock on the end of the cross pins has moved to the side to hold the support tops and H-brackets in position.
3. Loosely install a minimum of four safety straps evenly spaced and around both support tops.

NOTE:

We recommend that you use the four Allen® Safety Straps that were included with the support top.

4. Make sure the support tops are in the correct position, and then tighten the safety straps.
5. Make sure the superior support top is progressively compressed against the patient and once again tighten the safety straps.
6. Look at the display, and make sure the **tilt level indicator** and the **brakes locked indicator** are **green**.

⚠ WARNING:

We recommend that a minimum of four caregivers are available to help flip the patient. However, the number of caregivers will vary based on the patient weight and should be determined by the facility staff for each flip procedure. Failure to do so could cause patient fall, injury, or equipment damage.

⚠ WARNING:

Observe lines closely during articulations. Always use good line management techniques. Failure to do so could cause patient fall, injury, or equipment damage.


⚠ WARNING:

Failure to observe the patient's head while the support top is being flipped could cause patient fall or injury.

7. Make sure of these, and then follow facility protocols to manually flip the support tops 180°:
 - There is a sufficient number of persons available for the patient to be flipped.
 - Persons and equipment will not be in the way of the support tops when they get flipped.
 - To watch the lines so that they don't get pinched or caught during the flip.
 - The patient's head stays in the correct position during the flip.

Instructions for use

8. Do as follows to unlock the support tops for the flip:
 - a. Lift the end (A) of the flip handle until it locks in position.
 - b. Push and hold the red **safety lever** (B) forward as you pull the flip handle (A) down to the **released** position.

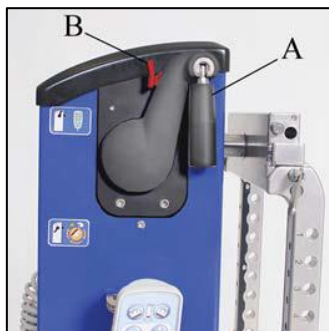NOTE:

The **manual flip indicator** on the display turns **orange** to show that the support tops are no longer locked in position and they can be flipped. The **pendant indicator** on the display turns **off** to show that you can **not** use the pendant controls while the Allen® Advance Table is in flip mode.

9. After the patient has been flipped, and the tops are level (the **tilt level indicator** on the display is **green**), push the Flip handle up to the **locked** position. You do not need to press the red **safety lever**.

NOTE:

The **flip indicator** on the display turns **off** to show that the support tops are locked in position. The **pendant indicator** on the display turns **green** to show that you can use the pendant controls.

10. Make sure the support tops are locked in position. If the tops did not lock, repeat Step 9.
11. Remove the straps from around the support tops.

⚠ WARNING:

Failure to use correct lifting methods could cause patient fall, injury, or equipment damage.

12. With the help of another person, carefully remove the upper support top.
13. Carefully remove the upper H-brackets.
14. Make sure the patient is in the correct position for the surgical procedure.


Transfer a patient

⚠ WARNING:

Before you transfer the patient, make sure the brakes on the Allen® Advance Table and the other transfer device are locked and the patient and staff are fully prepared for the movement. Otherwise, patient fall, injury, or equipment damage could occur.

1. Make sure that:

- The brakes on the Allen® Advance Table and the other transfer device are locked. (The **locked indicator** on the Allen® Advance Table display is **green**.)
- The Allen® Advance Table is in **pendant control mode**. (The **pendant control indicator** on the Allen® Advance Table display is **green**.)

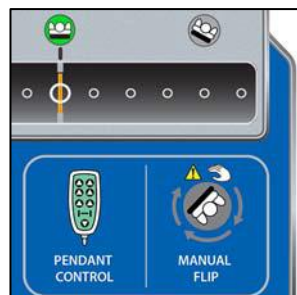- The patient and staff are fully prepared for the transfer.
2. Make sure the Allen® Advance Table support top and the surface of the transfer device are at the same height.

⚠ WARNING:

Observe lines closely during the transfer. Always use good line management techniques. Failure to do so could cause patient fall, injury, or equipment damage.

3. Observe lines closely as you transfer the patient per facility protocols.

Access the cervical traction assembly

The cervical traction release control is a push button on the upper-left side of the head-end column.

1. Push the **cervical traction release** button. The assembly moves out from the column.

2. Pull the assembly out from the column, and make sure the assembly is in the locked-out position. The cervical traction release button will return to its original position.


⚠ WARNING:

The maximum diameter of the traction rope should be 6.35 mm (0.25"). The use of a traction rope that is greater than 6.35 mm (0.25") in diameter could cause binding inside the cervical traction assembly. Patient injury could occur as a result of inaccurate traction.

3. Put the traction rope through the hole in the center of the bowtie. The rope will go through the traction assembly and out the left side of the column, near the release.
4. Put the traction rope over and around the pulley.

5. If desired, install the optional cervical traction bar on the support top-side of the head-end H-bracket as follows:
 - a. Determine the hole location for the cervical traction bar.

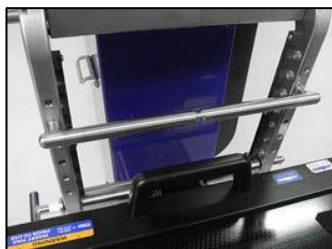NOTE:

The hole locations are numbered on the side of the H-brackets.

- b. Put the bar into the H-bracket openings. Make sure that the bar is level by verifying that each end of the bar is installed into the H-bracket at the same numbered hole location.
- c. Make sure the drop lock on each end of the bar moves to the side to hold the bar in position.
- d. Pull the traction rope from the bowtie down and underneath the bar.

NOTE:

The bar has an indentation to help keep the rope in position during cervical traction.

⚠ WARNING:

The maximum weight load for the traction assembly is 22.6 kg (50 lbs). To exceed this weight could cause patient fall, injury, or equipment damage.

⚠ WARNING:

Make sure the traction weight is in a safe position, free from possible unintended movement. Failure to do so could cause patient fall, injury, or equipment damage.

⚠ WARNING:

Do not press the cervical traction release button while traction weight is being applied. To do so could cause the assembly to retract and the traction weight to move. Patient fall, injury, or equipment damage could occur.

6. Use traction weight per facility protocols. Make sure that—
 - The traction weight is in a safe position such that it is free from possible unintended movement.
 - The traction rope stays in the indented section of the cervical traction bar.

Instructions for use

- Do not press the **cervical traction release button** while traction weight is being applied.

NOTE:

The cervical traction bar is an optional accessory that is supplied to help with the angle of pull for the cervical traction. If used, install the bar at the applicable location on the H-bracket.

- After use of the traction weight, remove the weight and the rope from the assembly.
- Press the **cervical traction release button**, and then push the assembly into the column so that it is in the locked-in position.

Shut down the Allen® Advance Table

- Level the Allen® Advance Table and the support top “Level the table and the support top” on page 35.
- Press and hold the **lower control** until the Allen® Advance Table is in its lowest position.

- Hang the pendant and its cord on their holders.

- Turn the Allen® Advance Table **off**. This removes all power from the Allen® Advance Table.
- Unplug the power cord from the facility power outlet, and store the power cord on the Allen® Advance Table.

NOTE:

We recommend that you keep the Allen® Advance Table plugged in when it is not in use so that the batteries will stay charged. Know and follow your facility's protocol.

Clean and disinfect

We recommend that you clean and disinfect the Allen® Advance Table and the support tops between patient use. Refer to your facility's cleaning and disinfection policies, and follow the guidelines below.

⚠ WARNING:

The potential for electrical shock exists with electrical equipment.
Failure to follow facility protocols may cause death or serious injury.

⚠ WARNING:

The potential for electrical shock exists with electrical equipment.
Failure to follow facility protocols may cause death or serious injury.

⚠ WARNING:

Always wear the applicable personal protective equipment when you clean the unit. Otherwise, injury could occur.

⚠ WARNING:

Do not put the unit in an equipment wash room. Patient fall, injury, or equipment damage could occur.

⚠ WARNING:

Do not expose the unit or its components to excessive moisture. Patient fall, injury, or equipment damage could occur.

⚠ CAUTION:

Do not use harsh cleaners, solvents, or detergents. Equipment damage could occur.

NOTE:

Do not use iodophors; they will stain the product.

⚠ CAUTION:

If you use a quaternary-type disinfectant, make sure to rinse it off the product immediately. Otherwise, the disinfectant could cause corrosion on painted and plated components.

Clean

⚠ CAUTION:

Failure to use cleaning and disinfectant products in accordance with the manufacturer's instructions could cause equipment damage.

We recommend that you clean the unit with detergent and warm water. Do not use excessive liquid or harsh cleansers.

Clean difficult spots

To remove difficult spots or stains, we recommend that you use standard household cleansers and a soft bristle brush. To loosen heavy, dried-on soil, you may first need to saturate the spot.

Steam clean

Do not use any steam cleaning device on the unit or put any Allen® Advance Table components in an autoclave or steam cleaning system. Excessive moisture can damage mechanisms in this unit and its components.

Disinfect

CAUTION:

Failure to use cleaning and disinfectant products in accordance with the manufacturer's instructions could cause equipment damage.

When there is visible soil and also between patient use, we recommend that you disinfect the unit with an EPA registered (US only), tuberculocidal, disinfectant such as Morning Mist Neutral Disinfectant.

Dilute the disinfectant as shown on the manufacturer's label.

Troubleshooting

WARNING:

Only Allen Medical-authorized persons should service the Allen® Advance Table. Service by unauthorized persons could cause patient fall, injury, or equipment damage.

If the troubleshooting information below does not fix the problem, contact Allen Medical Technical Support “Service calls” on page 47.

The Allen® Advance Table will not go into pendant control mode from the manual flip mode

Make sure that:

- The brakes are locked.
- The columns are at the same height.
- The Allen® Advance Table and the installed support top are level. See “Level the table and the support top” on page 35.
- The Allen® Advance Table is free of objects that could interfere with the adjustment.
- The flip handle is in the **up** and **locked** position.

If all of the conditions above are met, and the Allen® Advance Table still will not go into the manual flip mode, reset the power: press the **On/Off** switch to **off** and then to **on**.

The pendant control and manual flip indicators are flashing

Make sure that the flip handle is fully in the **up** and **locked** position or **down** and **released** position. If the indicators continue to flash, reset the power: press the **On/Off** switch to **off** and then to **on**.

The brakes locked and unlocked indicators are flashing

Press and hold either the brakes lock or unlock pendant control until the applicable brakes locked/unlocked indicator on the display stays on; the indicators are no longer flashing.

The Allen® Advance Table “beeps” when a pendant control is pressed

A “beep” from the Allen® Advance Table lets you know that the battery charge level is low or depleted. Plug the Allen® Advance Table into a power outlet to charge the batteries. To keep the batteries charged, plug the Allen® Advance Table into an outlet nightly and when you store the Allen® Advance Table.

Service calls

When you call Allen Medical about your unit, be prepared to give the serial number from the product identification label. You will find the label on the inside of the head-end base near the on/off switch.

Preventive maintenance

⚠ WARNING:

Only Allen Medical-authorized persons should service the Allen® Advance Table. Service by unauthorized persons could cause patient fall, injury, or equipment damage.

We recommend that you do annual preventive maintenance on the Allen® Advance Table. An effective way to reduce downtime and make sure the Allen® Advance Table operates correctly is to do preventive maintenance.

The suggested maintenance checks below include specific items for you to examine.

Maintenance checks

If any of these checks fail, repair or replace the part as applicable. If the repair or replacement does not correct the problem, remove the unit from service and contact Allen Medical Technical Support (refer to “Service calls” on page 47).

Overall condition
The unit is clean and has been disinfected.
The unit is not damaged: the welded assemblies are in good condition and there are no impacts or corrosion.
All attaching hardware is installed and correctly tightened.
All labels are installed and can be read.
The power cord is with the unit, and it is in good condition: the plug assembly shows no signs of cracking, the plug molding around the blade is not discolored, and the blade is tight in the molding.
Controls and motors
All caregiver controls, along with their indicators, operate correctly; this includes controls on the Allen® Advance Table and on the pendant.
The motors operate correctly: during articulation, there is no friction, abnormal noises, or an audible indication of overload.
Battery
The battery is in good condition: with the unit unplugged, the pendant controls operate all Allen® Advance Table articulations; the battery is charged.
Brakes and casters
The brakes operate correctly.
The casters are in good condition; they do not have cuts, are not worn, and have a good amount of tread.
Electrical checks
The leakage current is less than 100 μ A.
Accessories
The accessories are in good condition.

Storage and transport

Store the Allen® Advance Table in a safe place that is within the recommended environmental conditions on page 55. If you have the optional accessory cart, store the support tops and additional accessories on the cart.

Expected life

The expected life of the Allen® Advance Table 10 years of normal use when correctly maintained. However, certain components have a shorter life cycle and will need to be replaced. They are listed below:

- When correctly maintained, the battery has an expected life of 3 years.

Discard the unit

When the Allen® Advance Table has reached the end of its useful life, discard it in accordance with local and federal standards. Recycle as much as possible.

Parts replacement

WARNING:

Only facility-authorized maintenance persons should do these parts replacement procedures. Otherwise, patient fall, injury, or equipment damage could occur.

WARNING:

Always wear the applicable personal protective equipment when you service the unit. Failure to do so could cause injury.

WARNING:

Before you remove and replace electrical components (power cord, pendant, fuses), make sure to remove power from the Allen® Advance Table. Failure to do so could cause patient fall, injury, or equipment damage.

The parts shown below can be replaced by facility-authorized maintenance persons. If additional parts need to be replaced, contact Allen Medical Technical Support.

Description	Quantity	Part number
Pendant	1	A-71403
Pendant nut	1	C-007239
H-bracket	1	A-71402
Cross pin	1	A-71401
Power cord, EU	1	A-71407
Power cord, UK	1	A-71406
Power cord, US	1	A-71405

Fuse, 5 mm x 20 mm, 6.3 A, 250 V, slow blow (A-71100)	2	C-007318
Fuse, 5 mm x 20 mm, 3.15 A, 250 V, slow blow (A-71101)	2	C-007358

Pendant

Removal

1. Make sure the brakes are locked (the **locked** indicator on the Allen® Advance Table display is **green**).

WARNING:

Failure to remove power from the Allen® Advance Table could cause patient fall, injury, or equipment damage.

2. Turn the Allen® Advance Table **off**, and unplug it from the facility's power source.
3. Unscrew the nut from the Allen® Advance Table connector.

4. Remove the nut from pendant cord, and install the nut on the new pendant cable. Make sure the ribbed end of the nut is toward the connector end of the pendant cable.

5. Disconnect the pendant cable from the Allen® Advance Table connector.

Replacement


1. Make sure the O-ring is installed in the Table connector.

2. Align the tab on the connector end of the new pendant cable with the tab on the inside of the Allen® Advance Table connector.
3. Connect the pendant cable to the Allen® Advance Table connector.
4. Screw the nut on to the Allen® Advance Table connector.
5. Plug the Allen® Advance Table in, and turn it **on**.
6. Make sure all pendant functions operate correctly.

Fuses

Removal

1. Make sure the brakes are locked (the **locked** indicator on the Allen® Advance Table display is **green**).

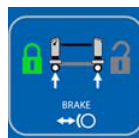⚠ WARNING:

Failure to remove power from the Allen® Advance Table could cause patient fall, injury, or equipment damage.

2. Turn the Allen® Advance Table **off**, and unplug the power cord from the facility's power source and from the Table.

3. Pull the fuse holder from the plug housing.

NOTE:

You may need to use a small screwdriver or a pen to pry the holder from the housing.

4. Remove the fuses from the holder.

Replacement

1. Install the **new** fuses in the holder.
2. Fully insert the holder into the plug housing.
3. Plug the power cord into the Allen® Advance Table and into an applicable power outlet, and turn the Allen® Advance Table **on**.
4. Make sure all Allen® Advance Table functions operate correctly.

Accessories

For installation and use instructions, refer to the manufacturer's documentation included with the accessory.

⚠ WARNING:

Weight ratings may differ for some accessories used with the Allen® Advance Table. Make sure you do not exceed the weight rating that is shown on the accessory. To do so could cause patient fall, injury, or equipment damage.

Part number	Description	Part number	Description
A-40022	Easy Lock® Socket (US)	A-40023	Easy Lock® Socket (EU)
A-40024	Easy Lock® Socket (UK)	A-70310	Comfort Mask Disposables
A-70404	Supine Access Top, 12" (304.8 mm)	A-70434	Supine 12" (304.8 mm) Access Top Section Pad, 2" (50.8 mm) Astrus
A-70436	Advance Operating Top Pad	A-70500-EU	Armboard Rail (EU)
A-70500-UK	Armboard rail (UK)	A-70500-US	Armboard Rail (US)
A-70501	Flex Frame Accessory Rail, 6" (152.4 mm) (US)	A-70501-EU	Flex Frame Accessory Rail, 6" (152.4 mm) (EU)
A-70501-UK	Flex Frame Accessory Rail, 6" (152.4 mm) (UK)	A-70615-EU	Arm Support, (EU)
A-70615-UK	Arm Support, (UK)	A-70615-US	Arm Support, (US)
A-70700	C-Flex® Head Positioner System ^a	A-70701	C-Flex® Base Unit
A-70702	C-Flex® Clean Cape Disposables	A-70710	Starburst Adapter for Standard Mayfield or DORO
A-70715	Prone Mask Head Module	A-70720	Flat Plate Head Module
A-70725	Reverse Starburst Adapter for Standard Mayfield or DORO	A-70800	Bow® Frame ^b
A-70801	Bow® Frame with pads	A-70815	Bow® Frame Cart (US)
A-70815-EU	Bow® Frame cart (EU)	A-70815-UK	Bow® Frame Cart (UK)
A-70821	Crank Arm	A-70862	Foam Head Blocks (case/6)

Part number	Description	Part number	Description
A-70863	Foam Arm Supports (case/6)	O-RSX2	Large Security Strap
Allen® Advance Table Accessories			
A-71261	Allen® Ultra Comfort™ Covers, small (case/6)	A-71262	Allen® Ultra Comfort™ Covers, medium (case/6)
A-71263	Allen® Ultra Comfort™ Covers, large (case/6)	A-71301	Chest Support with pad
A-71302	Chest Support Pad	A-71303	Hip Support, small, with pad
A-71304	Hip Support Pad, small	A-71305	Hip Support, medium, with pad
A-71306	Hip Support Pad, medium	A-71307	Hip Support, large, with pad
A-71308	Hip Support Pad, large	A-71309	Thigh, Support with pad, left
A-71310	Thigh Support Pad, large	A-71311	Thigh, Support with pad, right
A-71312	Sling Support	A-71316	Head Plate
A-71400	Rail Top	A-71404	Cervical Traction Bar
A-71410	Flat Top	A-71411	Safety Straps
A-71500-EU	Adjustable Accessory Rail with clamp (EU)	A-71500-UK	Adjustable Accessory Rail with clamp (UK)
A-71500-US	Crank Arm	A-71501-EU	Adjustable Accessory Rail (EU)
A-71501-UK	Adjustable Accessory Rail (UK)	A-71501-US	Adjustable Accessory Rail (US)
A-71502	Adjustable Accessory Rail Clamp (US)	A-71502-EU	Adjustable Accessory Rail Clamp (EU)
A-71502-UK	Adjustable Accessory Rail Clamp (UK)	A-71600	Accessory Cart (US)
A-71600-EU	Accessory Cart (EU)	A-71600-UK	Accessory Cart (UK)

Specifications

- a. A-70700 includes one each: A-70701, A-70702, A-70310, A-70710, A-70715, and A-70720.
- b. A-70800 includes one each: A-70801, A-70821, and A-70815.

Specifications

Product Identification

Product number	Description
A-71100	Allen® Advance Table, 120 V
A-71101	Allen® Advance Table, 230 V

Dimensions

Feature	Dimension
Width	82.6 cm (32.5")
Column height (measured from the support top to the floor)	67.3 cm to 127 cm (26.5" to 50") (electronically adjustable)
Length (retracted to fully extended)	208.9 to 276.9 cm (82.25" to 109")
Table weight (without a support top)	Approximately 147 kg (325 lbs)
Rail top weight	Approximately 9.5 kg (21 lbs)
Flat top weight	Approximately 16.5 kg (36.5 lbs)
Electronic rotation	25° clockwise and counterclockwise (can be locked at any point within this range)
Manual rotation	360° (can be locked at horizontal and vertical orientations)
Maximum Trendelenburg and Reverse Trendelenburg	10°
Battery life	3 years ^a
Maximum patient weight ^b	272 kg (600 lbs)

- a. If the batteries are maintained correctly, the expected life is 3 years.
- b. Although the Allen® Advance Table and the patient support tops have a weight rating of 272 kg (600 lbs), accessories used with the table may not have that same weight rating. Make sure you do not exceed the weight rating that is shown on the accessory.

Environmental conditions for use

Condition	Range
Ambient temperature	10 °C to 40 °C (50 °F to 104 °F)

Condition	Range
Relative humidity	15% to 70% (non-condensing)
Atmospheric pressure	700 hPa to 1060 hPa

Environmental conditions for transport and storage

Condition	Range
Ambient temperature	-14 °C to 70 °C (7 °F to 158 °F)
Relative humidity	15% to 93% (non-condensing)
Atmospheric pressure	700 hPa to 1060 hPa

Power requirements

Power distribution

Model	Range
120 V	60 Hz
230 V	50 Hz

Maximum current draw

Model	Range
120 V	6 A/m
230 V	3 A

Fuse specification

Model	Type
120 V	2 each, 5 mm x 20 mm, 6.3 A, 250 V, slow blow
230 V	2 each, 5 mm x 20 mm, 3.15 A, 250 V, slow blow

Classification and Standards

Classification	Standard
Technical and Quality Assurance Standards	IEC 60601 Ed. 3 for Class 1 Type B instrument IEC 606012-2-46 Ed. 2.0 relevant sections IEC 60601-1-2 relevant sections IEC 60601-1-6

Classification	Standard
IEC 60601-1 Classification According to Degree Protection against Electrical Shock	Class I equipment
IEC 60601-1 Classification According to Degree Protection against Electrical Shock	Type B applied part
IEC 60601-2-46 Classification According to Degree of Protection against Dust and Fluid Intrusion	IPX 4
Mode of Operation	2 minutes ON/18 minutes OFF

Applied parts (in accordance with IEC 60601-1)

Rail top	Flat top
Positioning supports (such as pads, C-Flex® Head Positioning System, and disposables)	Flat Top pad

Electromagnetic Emissions Guidance

WARNING:

If the Allen® Advance Table is to be used near high-frequency surgical equipment, cardiac defibrillators, and cardiac defibrillator-monitors, refer to the product manufacturer's instructions to make sure the equipment will operate as intended.

This device meets all requirements for electromagnetic compatibility per IEC 60601-1-2. Other products that are used near the Allen® Advance Table should comply with this standard also. Otherwise, electromagnetic interference between the products could cause the products to operate incorrectly. If problems do occur, contact the product manufacturer(s).

WARNING:

The A-71100/A-71101 should not be used adjacent to other electrical equipment. If adjacent use is necessary, observe the A-71100/A-71101 and the other electrical equipment to make sure they operate as intended.

CAUTION:

This device meets all requirements for electromagnetic compatibility per IEC 60601-1-2. It is unlikely that the user will encounter problems with this device because of inadequate electromagnetic immunity.

However, electromagnetic immunity is always relative, and standards are based on anticipated environments of use. If the user observes unusual device behavior, particularly if such behavior is intermittent and associated with nearby use of radio or TV transmitters, cell phones, or electro-surgical equipment, this could be an indication of electromagnetic interference. If such behavior occurs, the user should try to move the interfering equipment further from this device.

Make sure the A-71100/A-71101 operate correctly when they are used near other electronic devices. Portable and mobile radio frequency (RF) communications equipment can affect electrical equipment.

Medical equipment needs special precautions in regard to electromagnetic compatibility (EMC) and needs to be installed and put into service according to the EMC information supplied in the tables that follow.

Guidance and Manufacturer's Declaration—Electromagnetic Emissions		
The A-71100/A-71101 are intended for use in the electromagnetic environment specified below. The customer or the user of the A-71100/A-71101 should make sure it is used in such an environment.		
Emissions Test	Compliance	Electromagnetic Environment—Guidance
RF Emissions CISPR 11	Group 1	The A-71100/A-71101 use RF energy only for its internal functions. Therefore, its RF emissions are low and are not likely to cause any interference in nearby electronic equipment.
RF Emissions CISPR 11	Class A	
Harmonic Emissions IEC 61000-3-2	Class A	
Voltage Fluctuations/ Flicker Emissions IEC 61000-3-3	Complies	

Electromagnetic Immunity Guidance

Guidance and Manufacturer's Declaration - Electromagnetic Immunity			
The A-71100/A-71101 are intended for use in the electromagnetic environment specified below. The customer or the user of the A-71100/A-71101 should make sure it is used in such an environment.			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment—Guidance
Electrostatic Discharge (ESD) IEC 61000-4-2	± 6 kV Contact ± 8 kV Air	± 6 kV Contact ± 8 kV Air	Floors should be wood, concrete, or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical Fast Transient/Burst IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input/output lines	± 2 kV for power supply lines ± 1 kV for input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 1 kV line(s) to line(s) ± 2 kV line(s) to earth	± 1 kV line(s) to line(s) ± 2 kV line(s) to earth	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions, and voltage variations on power supply input lines IEC 61000-4-11	$< 5\% U_T$ ($> 95\%$ dip in U_T) for 0.5 cycles $40\% U_T$ (60% dip in U_T) for 5 cycles $70\% U_T$ (30% dip in U_T) for 25 cycles $< 5\% U_T$ ($> 95\%$ dip in U_T) for 5 seconds	$< 5\% U_T$ ($> 95\%$ dip in U_T) for 0.5 cycles $40\% U_T$ (60% dip in U_T) for 5 cycles $70\% U_T$ (30% dip in U_T) for 25 cycles $< 5\% U_T$ ($> 95\%$ dip in U_T) for 5 seconds	Mains power quality should be of a typical commercial or hospital environment. If the user of the A-71100/A-71101 requires continued operation during power mains interruption, it is recommended to power the equipment from an uninterruptedly power supply or a battery.
Power frequency(50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
Note: U_T is the AC mains voltage prior to application of the test level.			

Guidance and Manufacturer's Declaration - Electromagnetic Immunity			
The A-71100/A-71101 are intended for use in the electromagnetic environment specified below. The customer or the user of the A-71100/A-71101 should make sure it is used in such an environment.			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment—Guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	3 Vrms	Portable and mobile RF communications equipment should be used no closer to any part of the A-71100/A-71101, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance: $d = 1, 17\sqrt{P}$150 kHz to 80 MHz $d = 1, 17\sqrt{P}$80 MHz to 800 MHz $d = 2, 33\sqrt{P}$800 MHz to 2.5 GHz Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and <i>d</i> is the recommended separation in metres (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey^a, should be less than the compliance level in each frequency range^b. Interference may occur in the vicinity of equipment marked with the following symbol. 
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	3 V/m	
Note 1: At 80 MHz and 800 MHz, the higher the frequency range applies. Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

- a. Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the A-71100/A-71101 is used exceeds the applicable RF compliance level above, the equipment should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the equipment.
- b. Over the frequency range of 150 kHz to 80 MHz, field strength should be less than 3 V/m.

Recommended separation distances between portable and mobile RF communications equipment and the A-71100/A-71101			
The A-71100/A-71101 are intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the A-71100/A-71101 can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the A-71100/A-71101 as recommended below, according to the maximum output power of the communications equipment.			
Rated maximum output power of transmitter, W	Separation distance according to frequency of transmitter,m		
	150 kHz to 80 MHz $d = 1,17\sqrt{P}$	80 MHz to 800 MHz $d = 1,17\sqrt{P}$	800 MHz to 2.5 GHz $d = 2,33\sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.37	0.37	0.74
1	1.17	1.17	2.33
10	3.69	3.69	7.38
100	11.67	11.67	23.33
For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer. NOTE 1: At 80 MHz and 800 MHz, the higher the frequency range applies. NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

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部件号	描述	部件号	描述
Allen™ Advance Table 附件			
A-71261	Allen® Ultra Comfort™ Covers, small (case/6)	A-71262	Allen® Ultra Comfort™ Covers, medium (case/6)
A-71263	Allen® Ultra Comfort™ Covers, large (case/6)	A-71301	Chest Support with Pad
A-71302	Chest Support Pad	A-71303	Hip Support, small, with pad
A-71304	Hip Support Pad, small	A-71305	Hip Support, medium, with pad
A-71306	Hip Support Pad, medium	A-71307	Hip Support, large, with pad
A-71308	Hip Support Pad, large	A-71309	Thigh, Support with Pad, left
A-71310	Thigh Support Pad, large	A-71311	Thigh, Support with Pad, right
A-71312	Sling Support	A-71316	Head Plate
A-71400	Rail Top	A-71404	Cervical Traction Bar
A-71410	Flat Top	A-71411	Safety Straps
A-71500-EU	Adjustable Accessory Rail with Clamp (EU)	A-71500-UK	Adjustable Accessory Rail with Clamp (UK)
A-71500-US	Crank Arm	A-71501-EU	Adjustable Accessory Rail (EU)
A-71501-UK	Adjustable Accessory Rail (UK)	A-71501-US	Adjustable Accessory Rail (US)
A-71502	Adjustable Accessory Rail Clamp (US)	A-71502-EU	Adjustable Accessory Rail Clamp (EU)
A-71502-UK	Adjustable Accessory Rail Clamp (UK)	A-71600	Accessory Cart (US)
A-71600-EU	Accessory Cart (EU)	A-71600-UK	Accessory Cart (UK)

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部品番号	説明	部品番号	説明
Allen™ Advance Table Accessories			
A-71261	Allen™ Ultra Comfort™ Covers, Small (case/6)	A-71262	Allen™ Ultra Comfort™ Covers, Medium (case/6)
A-71263	Allen™ Ultra Comfort™ Covers, Large (case/6)	A-71301	Chest Support With Pad
A-71302	Chest Support Pad	A-71303	Hip Support, Small, With Pad
A-71304	Hip Support Pad, Small	A-71305	Hip Support, Medium, With Pad
A-71306	Hip Support Pad, Medium	A-71307	Hip Support, Large, With Pad
A-71308	Hip Support Pad, Large	A-71309	Left Thigh Support With Pad, Left
A-71310	Thigh Support Pad, Large	A-71311	Right Thigh Support With Pad, Right
A-71312	Sling Support	A-71316	Head Plate
A-71400	Rail Top	A-71404	Cervical Traction Bar
A-71410	Flat Top	A-71411	Safety Straps
A-71500-EU	Adjustable Accessory Rail With Clamp (EU)	A-71500-UK	Adjustable Accessory Rail With Clamp (UK)
A-71500-US	Crank Arm	A-71501-EU	Adjustable Accessory Rail (EU)
A-71501-UK	Adjustable Accessory Rail (UK)	A-71501-US	Adjustable Accessory Rail (US)
A-71502	Adjustable Accessory Rail Clamp (US)	A-71502-EU	Adjustable Accessory Rail Clamp (EU)
A-71502-UK	Adjustable Accessory Rail Clamp (UK)	A-71600	Accessory Cart (US)
A-71600-EU	Accessory Cart (EU)	A-71600-UK	Accessory Cart (UK)

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Dalāš numurs	Apraksts	Dalāš numurs	Apraksts
Allen® Advance Table piederumi			
A-71261	Allen® Ultra Comfort™ Covers, small (case/6)	A-71262	Allen® Ultra Comfort™ Covers, medium (case/6)
A-71263	Allen® Ultra Comfort™ Covers, large (case/6)	A-71301	Chest Support with Pad
A-71302	Chest Support Pad	A-71303	Hip Support, small, with pad
A-71304	Hip Support pad, small	A-71305	Hip Support, medium, with pad
A-71306	Hip Support Pad, medium	A-71307	Hip Support, large, with pad
A-71308	Hip Support Pad, large	A-71309	Thigh, Support with Pad, left
A-71310	Thigh Support Pad, large	A-71311	Thigh, Support with Pad, right
A-71312	Sling Support	A-71316	Head Plate
A-71400	Rail Top	A-71404	Cervical Traction Bar
A-71410	Flat Top	A-71411	Safety Straps
A-71500-EU	Adjustable Accessory Rail with Clamp (EU)	A-71500-UK	Adjustable Accessory Rail with Clamp (UK)
A-71500-US	Crank Arm	A-71501-EU	Adjustable Accessory Rail (EU)
A-71501-UK	Adjustable Accessory Rail (UK)	A-71501-US	Adjustable Accessory Rail (US)
A-71502	Adjustable Accessory Rail Clamp (US)	A-71502-EU	Adjustable Accessory Rail Clamp (EU)
A-71502-UK	Adjustable Accessory Rail Clamp (UK)	A-71600	Accessory Cart (US)
A-71600-EU	Accessory Cart (EU)	A-71600-UK	Accessory Cart (UK)

a. A-70700 ietver vienu: A-70701, A-70702, A-70310, A-70710, A-70715, un A-70720.

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Broj dela	Opis	Broj dela	Opis
A-71263	Allen® Ultra Comfort™ Covers, large (case/6)	A-71301	Chest Support with pad
A-71302	Chest Support Pad	A-71303	Hip Support, small, with pad
A-71304	Hip Support Pad, small	A-71305	Hip Support, medium, with pad
A-71306	Hip Support Pad, medium	A-71307	Hip Support, large, with pad
A-71308	Hip Support Pad, large	A-71309	Thigh, Support with pad, left
A-71310	Thigh Support Pad, large	A-71311	Thigh, Support with pad, right
A-71312	Sling Support	A-71316	Head Plate
A-71400	Rail Top	A-71404	Cervical Traction Bar
A-71410	Flat top	A-71411	Safety Straps
A-71500-EU	Adjustable Accessory Rail with clamp (EU)	A-71500-UK	Adjustable Accessory Rail with clamp (UK)
A-71500-US	Crank Arm	A-71501-EU	Adjustable Accessory Rail (EU)
A-71501-UK	Adjustable Accessory Rail (UK)	A-71501-US	Adjustable Accessory Rail (US)
A-71502	Adjustable Accessory Rail Clamp (US)	A-71502-EU	Adjustable Accessory Rail Clamp (EU)
A-71502-UK	Adjustable Accessory Rail Clamp (UK)	A-71600	Accessory Cart (US)
A-71600-EU	Accessory Cart (EU)	A-71600-UK	Accessory Cart (UK)

a. A-70700 obuhvata po jedan: A-70701, A-70702, A-70310, A-70710, A-70715 i A-70720.

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Číslo súčasti	Popis	Číslo súčasti	Popis
A-70501	Flex Frame Accessory Rail, 6" (152.4 mm) (US)	A-70501-EU	Flex Frame Accessory Rail, 6" (152.4 mm) (EU)
A-70501-UK	Flex Frame Accessory Rail, 6" (152.4 mm) (UK)	A-70615-EU	Arm Support, (EU)
A-70615-UK	Arm Support, (UK)	A-70615-US	Arm Support, (US)
A-70700	C-Flex® Head Positioner System ^a	A-70701	C-Flex® base unit
A-70702	C-Flex® Clean Cape Disposables	A-70710	Starburst Adapter for Standard Mayfield or DORO
A-70715	Prone Mask Head Module	A-70720	Flat Plate Head Module
A-70725	Reverse Starburst Adapter for Standard Mayfield or DORO	A-70800	Bow® Frame ^b
A-70801	Bow® Frame with pads	A-70815	Bow® Frame Cart (US)
A-70815-EU	Bow® Frame Cart (EU)	A-70815-UK	Bow® Frame Cart (UK)
A-70821	Crank Arm	A-70862	Foam Head Blocks (case/6)
A-70863	Foam Arm Supports (case/6)	O-RSX2	Large Security Strap
Príslušenstvo stola Allen® Advance Table			
A-71261	Allen® Ultra Comfort™ Covers, small (case/6)	A-71262	Allen® Ultra Comfort™ Covers, medium (case/6)
A-71263	Allen® Ultra Comfort™ Covers, large (case/6)	A-71301	Chest Support with pad
A-71302	Chest Support Pad	A-71303	Hip Support, small, with pad
A-71304	Hip Support Pad, small	A-71305	Hip Support, medium, with pad
A-71306	Hip Support Pad, medium	A-71307	Hip Support, large, with pad

Číslo súčasti	Popis	Číslo súčasti	Popis
A-71308	Hip Support Pad, large	A-71309	Thigh, Support with pad, left
A-71310	Thigh Support Pad, large	A-71311	Thigh, Support with pad, right
A-71312	Sling Support	A-71316	Head Plate
A-71400	Rail Top	A-71404	Cervical Traction Bar
A-71410	Flat Top	A-71411	Safety Straps
A-71500-EU	Adjustable Accessory Rail with clamp (EU)	A-71500-UK	Adjustable Accessory Rail with clamp (UK)
A-71500-US	Crank Arm	A-71501-EU	Adjustable Accessory Rail (EU)
A-71501-UK	Adjustable Accessory Rail (UK)	A-71501-US	Adjustable Accessory Rail (US)
A-71502	Adjustable Accessory Rail Clamp (US)	A-71502-EU	Adjustable Accessory Rail Clamp (EU)
A-71502-UK	Adjustable Accessory Rail Clamp (UK)	A-71600	Accessory Cart (US)
A-71600-EU	Accessory Cart (EU)	A-71600-UK	Accessory Cart (UK)

a. A-70700 obsahuje po jednom kuse z: A-70701, A-70702, A-70310, A-70710, A-70715 a A-70720.

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Številka dela	Opis	Številka dela	Opis
A-70702	C-Flex® Clean Cape Disposables	A-70710	Starburst Adapter for Standard Mayfield or DORO
A-70715	Prone Mask Head Module	A-70720	Flat Plate Head Module
A-70725	Reverse Starburst Adapter for Standard Mayfield or DORO	A-70800	Bow® Frame ^b
A-70801	Bow® Frame with pads	A-70815	Bow® Frame Cart (US)
A-70815-EU	Bow® Frame Cart (EU)	A-70815-UK	Bow® Frame Cart (UK)
A-70821	Crank Arm	A-70862	Foam Head Blocks (case/6)
A-70863	Foam Arm Supports (case/6)	O-RSX2	Large Security Strap
Dodatki Allen® Advance Table			
A-71261	Allen® Ultra Comfort™ Covers, small (case/6)	A-71262	Allen® Ultra Comfort™ Covers, medium (case/6)
A-71263	Allen® Ultra Comfort™ Covers, large (case/6)	A-71301	Chest Support with pad
A-71302	Chest Support Pad	A-71303	Hip Support, small, with pad
A-71304	Hip Support Pad, small	A-71305	Hip Support, medium, with pad
A-71306	Hip Support Pad, medium	A-71307	Hip Support, large, with pad
A-71308	Hip Support Pad, large	A-71309	Thigh, Support with pad, left
A-71310	Thigh Support Pad, large	A-71311	Thigh, Support with pad, right
A-71312	Sling Support	A-71316	Head Plate
A-71400	Rail Top	A-71404	Cervical Traction Bar
A-71410	Flat Top	A-71411	Safety Straps

Številka dela	Opis	Številka dela	Opis
A-71500-EU	Adjustable Accessory Rail with clamp (EU)	A-71500-UK	Adjustable Accessory Rail with clamp (UK)
A-71500-US	Crank Arm	A-71501-EU	Adjustable Accessory Rail (EU)
A-71501-UK	Adjustable Accessory Rail (UK)	A-71501-US	Adjustable Accessory Rail (US)
A-71502	Adjustable Accessory Rail Clamp (US)	A-71502-EU	Adjustable Accessory Rail Clamp (EU)
A-71502-UK	Adjustable Accessory Rail Clamp (UK)	A-71600	Accessory Cart (US)
A-71600-EU	Accessory Cart (EU)	A-71600-UK	Accessory Cart (UK)

- a. A-70700 vključuje po enega: A-70701, A-70702, A-70310, A-70710, A-70715 in A-70720.
- b. A-70800 vključuje po enega: A-70801, A-70821 in A-70815.

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Número de pieza	Descripción	Número de pieza	Descripción
A-70501	Flex Frame Accessory Rail 152,4 mm (6") (US)	A-70501-EU	Flex frame accessory rail, 6" (152.4 mm) (EU)
A-70501-UK	Flex Frame Accessory Rail, 6" (152.4 mm) (UK)	A-70615-EU	Arm Support (EU)
A-70615-UK	Arm Support, (UK)	A-70615-US	Arm Support, (US)
A-70700	C-Flex® Head Positioner System ^a	A-70701	C-Flex® Base Unit
A-70702	C-Flex® Clean Cape Disposables	A-70710	Starburst Adapter for Standard Mayfield o DORO
A-70715	Prone Mask Head Module	A-70720	Flat Plate Head Module
A-70725	Reverse Starburst Adapter for Standard Mayfield o DORO	A-70800	Bow® Frame ^b
A-70801	Bow® Frame with pads	A-70815	Bow® Frame Cart (US)
A-70815-EU	Bow® Frame Cart (EU)	A-70815-UK	Bow® Frame Cart (UK)
A-70821	Crank Arm	A-70862	Foam Head Blocks (caja de 6)
A-70863	Foam Arm Supports (caja de 6)	O-RSX2	Large Security Strap
Accesorios de la Allen® Advance Table			
A-71261	Allen® Ultra Comfort™ Covers, small (caja de 6)	A-71262	Allen® Ultra Comfort™ Covers, medium (caja de 6)
A-71263	Allen® Ultra Comfort™ Covers, large (caja de 6)	A-71301	Chest Support with pad
A-71302	Chest Support Pad	A-71303	Hip Support, small, with pad
A-71304	Hip Support Pad, small	A-71305	Hip Support, medium, with pad
A-71306	Hip Support Pad, medium	A-71307	Hip Support, large, with pad

Número de pieza	Descripción	Número de pieza	Descripción
A-71308	Hip Support Pad, large	A-71309	Thigh, Support with pad, left
A-71310	Thigh Support Pad, large	A-71311	Thigh, Support with pad, right
A-71312	Sling Support	A-71316	Head Plate
A-71400	Rail Top	A-71404	Cervical Traction Bar
A-71410	Flat Top	A-71411	Safety Straps
A-71500-EU	Adjustable Accessory Rail with clamp (EU)	A-71500-UK	Adjustable Accessory Rail with clamp (UK)
A-71500-US	Crank Arm	A-71501-EU	Adjustable Accessory Rail (EU)
A-71501-UK	Adjustable Accessory Rail (UK)	A-71501-US	Adjustable Accessory Rail (US)
A-71502	Adjustable Accessory Rail Clamp (US)	A-71502-EU	Adjustable Accessory Rail Clamp (EU)
A-71502-UK	Adjustable Accessory Rail Clamp (UK)	A-71600	Accessory Cart (US)
A-71600-EU	Accessory Cart (EU)	A-71600-UK	Accessory Cart (UK)

a. El modelo A-70700 incluye uno de cada: A-70701, A-70702, A-70310, A-70710, A-70715 y A-70720.

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© 2014 af Allen Medical Systems, Inc. ALLE RETTIGHEDER