

### Key Features

- Measure copy number and copy-neutral changes, all on one microarray
- No set up fees required
- Wizards provided to make design even easier
- Use probe scores to optimize designs *in silico*

[www.agilent.com/genomics](http://www.agilent.com/genomics)

This item is not approved for use in diagnostic procedures. User is responsible for obtaining regulatory approval or clearance from the appropriate authorities prior to diagnostic use.

Agilent Technologies shall not be liable for errors contained herein or for incidental or consequential damages in connection with the furnishing, performance, or use of this material.

© Agilent Technologies, Inc. 2010  
Published in USA, September 22, 2010  
5990-6423EN

## eArray 7.0

### Create custom SurePrint G3 CGH+SNP arrays


eArray 7.0 has now been extended to include custom design support for Agilent's new SurePrint G3 CGH+SNP microarrays so you can detect LOH, UPD and other copy-neutral changes from a single microarray. This online application offers unprecedented flexibility and functionality of custom array design — at no charge!

Create designs from your own probes, utilize Agilent probe design tools, reference public databases such as NCBI's GenBank, or use Agilent's Human CGH/CNV Probe database, which contains over 28 million prequalified CGH probes and 65,000 SNPs.

- Select content to match experimental needs with 60K, 180K, 400K, or 1M feature array formats
- Supports reference samples NA12878, NA12891, NA18507, NA18517 and NA18579, making it easy to avoid SNPs with zero reference copies

### CGH+SNP: more probes, same high-quality design

The SurePrint CGH+SNP array uses the same Agilent CGH workflow but employs a novel algorithm contained within Agilent's Genomic Workbench software to measure copy-neutral changes based on SNP genotyping data.

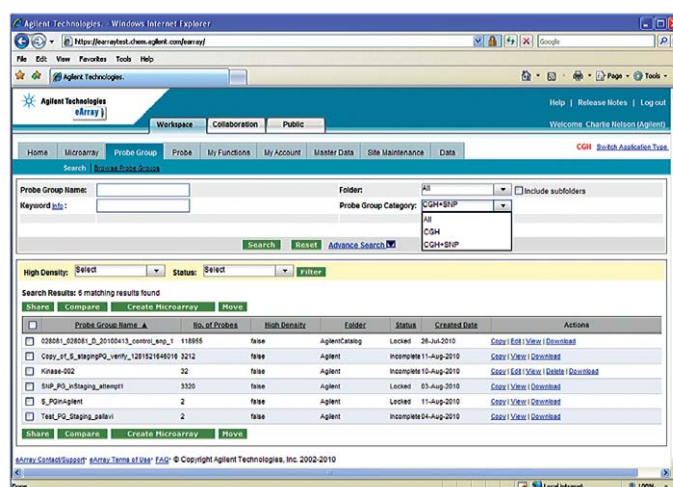

With eArray 7.0, SNP probes are at your fingertips

The Measure of Confidence



Agilent Technologies