EnteraLite® Infinity®
ENTERAL FEEDING PUMP

QUICK REFERENCE INSTRUCTIONS
OPEN SYSTEM USE


PRIMING AND LOADING THE DISPOSABLE SET

THESE INSTRUCTIONS ARE NOT A SUBSTITUTE FOR THE OPERATOR'S MANUAL.

Be sure to read the operator's manual thoroughly prior to operating the pump.
See other side for pump operation instructions.


STEP 1

If you are using the 500 mL or 1200 mL disposable set, hold bag upright and pour in feeding solution. Close cap securely and hang bag.

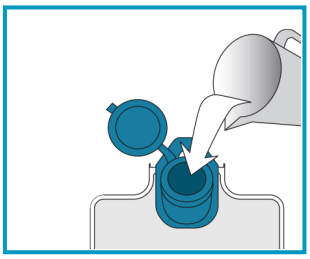


STEP 2


Remove protective cover from the distal end connector.

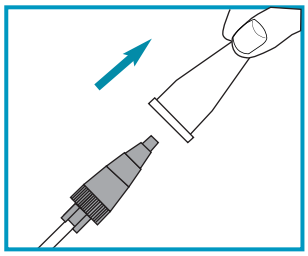

STEP 3 – PUMP PRIMING


Loop silicone tubing around pump wheel stretching lightly.
Seat cassette into pump. Close pump door.
Press and hold **ON/OFF** key for 1.5 seconds to turn pump on.
Press and hold the **PRIME** key until formula reaches the distal end connector.

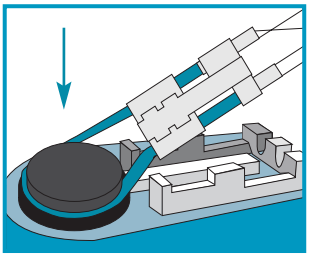


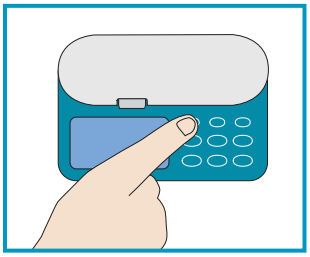

STEP 3 – MANUAL PRIMING

Gently pinch teal colored tubing below the "💧" symbol and squeeze bag at the same time to assist fluid flow.
Hold this position until formula reaches the distal end connector.
Loop silicone tubing around pump wheel stretching lightly.
Seat cassette into pump. Close pump door.

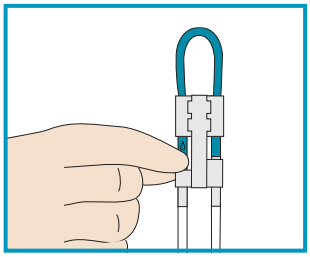

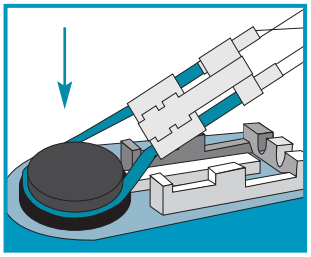

Change set every 24 hours

EnteraLite® and Infinity® are registered trademarks of ZEVEX, Inc.
Except as noted, all trademarks are owned by Société des Produits Nestlé S.A., Vevey, Switzerland.
© 2015 Nestlé HealthCare Nutrition inc. All rights reserved.
201174 Rev A

NestléHealthScience


PROGRAMMING AND RUNNING THE PUMP

THESE INSTRUCTIONS ARE NOT A SUBSTITUTE FOR THE OPERATOR'S MANUAL.

Be sure to read the operator's manual thoroughly prior to operating the pump.
See other side for delivery set priming and loading instructions.

STEP 1

Hang feeding container so that the bottom of the bag is at or above the level of the pump door.

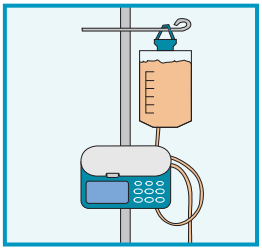

STEP 2

Press and hold **ON/OFF** key for 1.5 seconds to turn pump on.
After running the self test, display will show last programmed rate.




STEP 3 – Program Rate

Program the prescribed rate using the **(+) and (–)** keys.




STEP 4 – Dose Feeding

Press **RATE/DOSE** key to display dose. Program the prescribed dose using the **(+) and (–)** keys. If *FEED INTERVAL* is on display, clear feeding interval by pressing **FEED INT** to display feeding interval, then press **CLEAR** key.

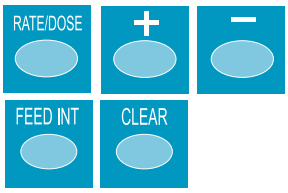


STEP 4 – Continuous Feeding/Infinite Dose

If *FEED INTERVAL* is on display, clear feeding interval by pressing **FEED INT** to display feeding interval, then press **CLEAR** key.To program an infinite dose, press **RATE/DOSE** key twice to display **DOSE** and hold the and hold the **(+)** key until display reads INF.

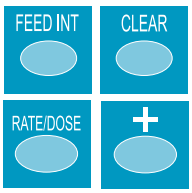

STEP 4 – Intermittent/Interval Feeding

Press **RATE/DOSE** key to display dose. Program the prescribed dose using the **(+) and (–)** keys. Press the **FEED INT** key to display feeding interval. Program the desired number of hours and minutes using the **(+) and (–)** keys.

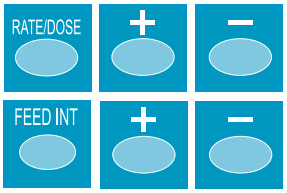

The Infinite Dose or Interval Feeding programming option listed in the operators manual should NOT be used for seriously compromised patients and children, or patients extremely sensitive to receiving air in the stomach with an enteral feeding. It is recomended that a single feeding regimen for specific dose be used for this particular user group.

STEP 5

Press the **RUN/PAUSE** key to start the pump.



STEP 6

When feeding is complete, press the **ON/OFF** key for 1.5 seconds to turn the pump off.



Clearing settings

Press **CLEAR** key to reset displayed function to its minimum value. Press and hold for three seconds to clear rate, dose, feed interval and volume delivered to minimum values at the same time.
Minimum value for all settings is "0" except for rate which is 0.1 mL/hour.



Viewing volume delivered

Press **VOL/TOTAL** key once to display volume delivered for current feeding. Press key again to display total cumulative volume delivered in all feed cycles since total volume was last cleared.



EnteraLite® Infinity® ENTERAL FEEDING PUMP

TROUBLESHOOTING GUIDE

NO FOOD Alarm

If re-priming of the pump delivery set is required, it is advisable to disconnect the enteral adapter from the patient's feeding tube while re-priming.

1. Is bag of pump delivery set empty?

- **YES** → Refill bag, re-prime delivery set and resume feeding.
- **NO**

2. Is air present in the tubing?

- **YES** → Has air been removed from **bag** of pump delivery set?
- **NO**

IF NO: Manually prime tubing by turning bag upside down, gently squeezing lower half of bag while pinching the teal colored tubing just below the drop symbol. Hold this position until all air is removed from bag and tubing.

IF YES: Press and hold the PRIME key to move the air in the tubing past pump delivery set cassette and resume feeding.

3. Is formula blenderized or aggressively mixed?

- **YES** → Let formula sit for 10-15 minutes before preparing pump delivery set for feeding. If foam is present in tubing, flush by priming pump delivery set either manually or with the pump and resume feeding.
- **NO**

4. Are air bubbles trapped inside the cassette?

- **YES** → Remove tubing from pump and manually prime tubing by turning bag upside down, gently squeezing lower half of bag while pinching the teal colored tubing just below the drop symbol. Hold this position until all air has passed through the teal segment fo the tubing.
- **NO** → If alarm continues, contact Nestlé HealthCare Nutrition, (877) 338-4238 for assistance.

ERRA, ERRD or ERRF Alarms

1. Was pump door open while pump was turned on or while running?

- **YES/ UNSURE** → Turn pump off. Unload and reload delivery set cassette. Make sure door is closed securely, then turn pump back on.
- **NO** → If alarm continues, contact Nestlé HealthCare Nutrition, (877) 338-4238 for assistance.

LOAD SET Alarm

1. Is door closed securely?

- **NO/ UNSURE** → Turn pump off. Unload and reload delivery set cassette. Make sure door is closed securely, then turn pump back on.
- **YES**

2. Is door cracked or tab on inside of door broken?

- **YES** → Contact MOOG Customer Advocacy, (800) 970-2337 to order a new door.
- **NO** → Contact Nestlé HealthCare Nutrition, (877) 338-4238 for assistance.

Except as noted, all trademarks are owned by Société des Produits Nestlé S.A., Vevey, Switzerland. EnteraLite® and Infinity® are trademarks of MOOG, Inc. 26084-001 Rev. B
© 2013 Nestlé HealthCare Nutrition inc. All rights reserved.
200765 Rev A

NestléHealthScience
NOURISHING PERSONAL HEALTH

NO FLOW IN or NO FLOW OUT Alarm

- **NO FLOW IN** – Occlusion between delivery set and pump.
- **NO FLOW OUT** – Occlusion between pump and patient.

1. Is delivery set tubing pinched, kinked or clogged?

- **YES** → Check delivery set for obstructions or kinks in tubing. Correct blockage and resume feeding.
- **NO**

2. Is tab on inside of door broken?

- **YES** → **Contact MOOG Customer Advocacy, (800) 970-2337 to order a new door.**
- **NO**

3. Is pressure sensor region of cassette receptacle clean?

- **NO** → Clean pressure sensor area with a cotton swab, soft cloth or dampened sponge, or wash entire pump under running water. Do not use abrasive materials or harsh chemicals.
- **YES**

4. Are their visible signs of damage to pressure sensors area?

- **YES** → **Contact MOOG Customer Advocacy, (800) 970-2337 to return pump for service.**
- **NO** → Disconnect from patient and fill a new delivery set with food. Prime set and start pump. **If alarm continues, contact Nestlé HealthCare Nutrition, (877) 338-4238 for assistance.**

CHARGING PROBLEMS

1. is A/C adapter charger plugged into wall outlet and pump properly?

- **NO/ UNSURE** → Check connections for both wall plug and pump adapter. Green indicator light on charger wall plug should be illuminated and pump display should have plug symbol with scrolling bars between E and F of fuel gauge. When pump is fully charged, bars will stop scrolling.
- **YES**

2. Is wall outlet functioning properly?

- **NO/ UNSURE** → Plug another device into outlet to verify outlet is functioning properly.
- **YES**

3. Does charger port appear to have any damage including bent or missing pins?

- **YES** → **Contact MOOG Customer Advocacy, (800) 970-2337 to order a new charger.**
- **NO**

4. Does charger adapter appear to have any damage?

- **YES** → **Contact MOOG Customer Advocacy, (800) 970-2337 to return pump for service.**
- **NO**

5. Is Battery Symbol and E and F of fuel gauge flashing?

- **YES** → **Contact MOOG Customer Advocacy, (800) 970-2337 to return pump for service.**
- **NO** → **If condition continues, contact MOOG Customer Advocacy, (800) 970-2337 to return pump for service.**