Instruction Manual 1

OMRON

IM1-HEM-7143T2-ESL-01-10/2020
5674648-4A

Automatic Upper Arm Blood Pressure Monitor

X2 Smart (HEM-7143T2-ESL)


English

Read Instruction manual (1) and (2) before use.

For symbols information, refer to "Symbols Description" on the back side of the another "Instruction manual (1)".

All for Healthcare

CE0197

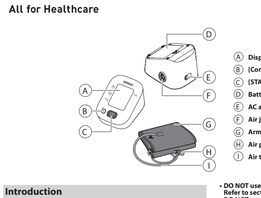A **Display**
B **[Connection] (Memory) button**
C **[START/STOP] button**
D **Battery compartment**
E **AC adapter jack**
F **Air jack**
G **Arm cuff**
H **Air plug**
I **Air tube**

Introduction

Thank you for purchasing the X2 Smart Automatic Upper Arm Blood Pressure Monitor. This blood pressure monitor uses the oscillometric method of blood pressure measurement. It means this monitor detects your blood movement through your brachial artery and converts the movements into a digital reading.

Safety Instructions

This instruction manual provides you with important information about the X2 Smart Automatic Upper Arm Blood Pressure Monitor. To ensure the safe and proper use of this monitor, READ and UNDERSTAND all of these instructions. **If you do not understand these instructions or have any questions, contact your OMRON retail outlet or distributor before attempting to use this monitor. For specific information about your own blood pressure, consult with your physician.**

Intended Use

This device is a digital monitor intended for use in measuring blood pressure and pulse rate in adult patients with upper arm. The monitor detects the appearance of irregular heartbeats during measurement and indicates this via a symbol with readings. It is mainly designed for general household use.

Receiving and Inspection

Remove this monitor from its packaging and inspect for damage. If this monitor is damaged, DO NOT USE and consult with your OMRON retail outlet or distributor.

Important Safety Information

Read the Important Safety Information in this instruction manual before using this monitor. Follow this instruction manual thoroughly for your safety.

Keep this manual for reference. For specific information about your own blood pressure, CONSULT WITH YOUR PHYSICIAN.

⚠Warning **Indicates a potentially hazardous situation which, if not avoided, could result in death or serious injury.**

- **DO NOT use this monitor on infants, toddlers, children or persons who cannot express themselves.**
- **DO NOT adjust medication based on readings from this blood pressure monitor. Take medication as prescribed by your physician. OMRON medication is qualified to diagnose and treat high blood pressure.**
- **DO NOT use this monitor on an injured arm or an arm under medical treatment.**
- **DO NOT apply the arm cuff on your arm while on an intravenous drip or blood transfusion.**
- **DO NOT use this monitor in areas containing high frequency (HF) surgical equipment, magnetic resonance imaging (MRI) equipment, computerized tomography (CT) scanners. This may result in incorrect operation of the monitor and/or cause an inaccurate reading.**
- **DO NOT use this monitor in oxygen rich environments or near flammable gas.**
- **Consult with your physician before using this monitor if you have common arrhythmias such as atrial or ventricular premature beats or atrial fibrillation; arterial sclerosis; poor perfusion; kidney disease; pregnancy; pre-eclampsia; liver disease. NOTE that any of these conditions in addition to patient motion, trembling, or shivering may affect the measurement reading.**
- **Consult your physician if you have diabetes.**
- **DO NOT use this monitor after the durable period has ended. Refer to section 3.**
- **DO NOT crease the arm cuff or the air tube excessively.**
- **DO NOT fold or kink the air tube while taking a measurement. This may cause an injury by interrupting blood flow.**
- **To unplug the air plug, pull on the plastic air plug at the base of the tube, not the tube itself.**
- **DO NOT use this monitor, arm cuff, batteries and accessories specified for this monitor. Use of unsupported AC adapters, arm cuffs and batteries may damage and/or may be hazardous to this monitor.**
- **ONLY use the approved arm cuff for this monitor. Use of other arm cuffs may result in incorrect readings.**
- **Inflating to a higher pressure than necessary may result in bruising of the arm where the cuff is applied. NOTE: refer to "If your systolic pressure is more than 210 mmHg" in section 8 of instruction manual (2) for additional information.**

Data Transmission

- **DO NOT replace batteries or unplug the AC adapter while your readings are being transferred to your smart device. This may result in incorrect operation of this monitor and failure to transfer your blood pressure data.**

AC Adapter (optional accessory) Handling and Usage


- **Fully insert the AC adapter into the outlet.**
- **When unplugging the AC adapter from the outlet, be sure to safely pull from the AC adapter. DO NOT pull from the AC adapter cable.**
- **When handling the AC adapter cable:**
- **DO NOT damage it.**
- **DO NOT break it.**
- **DO NOT twist it.**
- **DO NOT pinch it.**
- **DO NOT forcibly bend or pull it.**
- **DO NOT heat it.**
- **DO NOT use it if it is gathered in a bundle.**
- **DO NOT place heavy objects.**
- **Wipe any dust off of the AC adapter.**
- **Unplug the AC adapter when not in use.**
- **Unplug the AC adapter before cleaning this monitor.**

Battery Handling and Usage

- **DO NOT insert batteries with their polarities incorrectly aligned.**
- **ONLY use 4 "AA" alkaline or manganese batteries with this monitor. Do NOT use other types of batteries. DO NOT use new and used batteries together. DO NOT use different brands of batteries together.**
- **Remove batteries if this monitor will not be used for a long period of time.**
- **If battery fluid should get in your eyes, immediately rinse with plenty of clean water. Consult with your physician immediately.**
- **If battery fluid should get on your skin, wash your skin immediately with plenty of clean, lukewarm water. If irritation, injury or pain persists, consult with your physician.**
- **DO NOT use batteries after their expiration date.**
- **Periodically check batteries to ensure they are in good working condition.**

■ General Precautions

- When you take a measurement on the right arm, the air tube should be at the side of your elbow. Be careful not to rest your arm on the air tube.

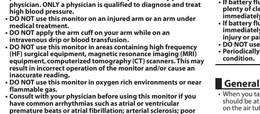- Blood pressure may differ between the right and left arm, and may result in a different measurement value. Always use the same arm for measurement. If the difference between both arms differ substantially, check with your physician on which arm to use for your measurements.
- When using an optional AC adapter, make sure not to place your monitor in a location where it is difficult to plug and unplug the AC adapter.
- To measure in a pressurized room, press the [START/STOP] button while taking a measurement.

Battery Handling and Usage

- Disposal of used batteries must be carried out in accordance with local regulations.
- The supplied batteries may have a shorter life span than new batteries.

1. Error Messages and Troubleshooting

If any of the below problems occur during measurement, check to make sure that the blood pressure device is within 30 cm. If the problem persists, please refer to the following table.

Display/ Problem	Possible Cause	Solution
E1 appears		Press the [START/STOP] button again to turn the monitor off. After inserting the air plug securely and applying the arm cuff correctly, press the [START/STOP] button.
The [START/STOP] button was pressed while the arm cuff is not applied.		
E1 appears or the arm cuff does not inflate.	Air plug is not completely plugged into the monitor.	Insert the air plug securely.
	The arm cuff is not applied correctly.	Apply the arm cuff correctly, then take measurement. Refer to sections 4 and 5 in instruction manual (2).
E2 appears or a measurement cannot be completed after the arm cuff inflates.	You move or talk during a measurement and the arm cuff does not inflate sufficiently.	Remain still and do not talk during a measurement. If measurement appears repeatedly, inflate the arm cuff manually until the systolic pressure is 30 to 40 mmHg above your previous readings. Refer to section 8 of instruction manual (2).
	Due to the systolic pressure is above 210 mmHg.	
E3 appears	The arm cuff is inflated exceeding the maximum allowable pressure.	Do not touch the arm cuff and/or bend the air tube while taking a measurement. If inflating the arm cuff manually, refer to section 8 of instruction manual (2).
E4 appears	You move or talk during a measurement. Vibrations disrupt a measurement.	Remain still and do not talk during a measurement.
The pulse rate is not detected correctly.		Apply the arm cuff correctly, then take another measurement. Refer to section 6 of instruction manual (2).
does not flash during a measurement		Remain still and sit correctly during a measurement. If the "♡" symbol continues to appear, we recommend you to consult with your physician.
Er appears	The monitor has malfunctioned.	Press the [START/STOP] button again. If "Er" still appears, your OMRON retail outlet or distributor.
Er appears	The monitor cannot connect to a smart device or transmit data correctly.	Follow the instructions shown in the "OMRON connect" app. If the "Er" symbol still appears after checking the app, contact your OMRON retail outlet or distributor.
Er appears	The monitor is waiting for pairing with the smart device.	Refer to section 5 of instruction manual (2) for pairing your monitor with your smart device, or press the [START/STOP] button to cancel pairing and turn your monitor off.
P flashes		
flashes	More than 24 readings are not transferred.	
flashes	Your monitor is not paired or not connected with your smart device.	Pair or transfer your readings to the "OMRON connect" app so you can keep your measurements in the app, and this error symbol disappears.
flashes	There are 30 readings in memory to be transferred	
appears		
flashes	Batteries are low.	Replacing all batteries with new ones is recommended. Refer to section 3 of instruction manual (2).
A and appears	Batteries are depleted.	Immediately replace all 4 batteries with new ones. Refer to section 3 of instruction manual (2).
Nothing appears on the display of the monitor.	Batteries are depleted.	Check the battery installation for proper polarities. Refer to section 3 of instruction manual (2).
Readings appear too high or too low.	Batteries are not properly aligned.	
	Blood pressure varies constantly. Many factors including stress, time of day, and how you apply the arm cuff, may affect your blood pressure. Review section 2 of instruction manual (2).	
Any other problem occurs.	Follow the instructions shown in the smart device, or visit the "Help" section in the "OMRON connect" app for further help. If the problem still persists, contact your OMRON retail outlet or distributor.	
Any other problems occur.	Press the [START/STOP] button to turn the monitor off, then press it again to take a measurement. If the problem continues, remove all batteries and wait for 30 seconds. Then re-install batteries. If the problem still persists, contact your OMRON retail outlet or distributor.	

2. Maintenance

2.1 Maintenance

To protect your monitor from damage, follow the directions below:

- Changes or modifications not approved by the manufacturer will void the user warranty.

3. Specifications

Product category	Electronic Sphygmomanometers
Product description	Automatic Upper Arm Blood Pressure Monitor
Model (Code)	X2 Smart (HEM-7143T2-ESL)
Display	LCD digital display
Cuff pressure range	0 to 299 mmHg
Blood pressure measurement range	SYS: 60 to 260 mmHg DIA: 40 to 215 mmHg
Pulse measurement range	40 to 180 beats / min.
Accuracy	Pressure: ±3 mmHg Pulse: ±5% of display reading
Inflation	Automatic by electric pump
Deflation	Automatic pressure release valve
Measurement method	Oscillometric method
Transmission method	**Bluetooth®** Low Energy
Wireless communication	Frequency range: 2.4 GHz (2400 - 2483.5 MHz) Modulation: GFSK Effective radiated power: < 20 dBm
Operation mode	Continuous operation
IP classification	Monitor: IP20 Optional AC adapter: IP21 for HEM-CR401, IP22 for HHP-BFH01
Rating	DC6 V 4 W
Power source	4 "AA" batteries 1.5 V or optional AC adapter (INPUT: AC 100 - 240 V 50 - 60 Hz 0.12 - 0.065 A)
Battery life	Approximately 900 measurements (Using new alkaline batteries and included arm cuff. Depending on the type of battery and arm cuff.)
Durable period (Service life)	Monitor: 5 years Arm cuff: 5 years Optional AC adapter: 5 years
Operating conditions	+10 to +40°C / 15 to 90% RH (non-condensing) / 800 to 1060 hPa
Storage / Transport conditions	-20 to +60°C / 10 to 90% RH (non-condensing)
Weight	Monitor, approximately 250 g (not including batteries) Arm cuff: approximately 110 g
Dimensions (approximately value)	Monitor: 146 mm (W) × 82 mm (H) × 140 mm (L) Arm cuff: 146 mm × 466 mm (air tube: 610 mm)
Cuff circumference applicable to the monitor	17 to 42 cm (included arm cuff: 22 to 32 cm)
Memory	Stores up to 30 readings
Contents	Monitor, arm cuff (HEM-CR24), 4 "AA" batteries, storage case, Instruction Manual (1) and (2)
Protection against electric shock	Internally powered BF equipment (when using only batteries) Class II ME equipment (optional AC adapter)
Applied part	Type BF (arm cuff)

Note

- These specifications are subject to change without notice.
- This monitor is clinically validated according to the requirements of EN ISO 81060-2:2014 and complies with EN ISO 81060-2:2014 and EN 81060-2:2019+A1:2020 (excluding pregnant and pre-eclampsia patients). In the clinical validation study, KS was used on 85 subjects for determination of diastolic blood pressure.
- IP classification is degrees of protection provided by enclosures in accordance with IEC 60529. The monitor and optional AC adapter are protected against solid foreign objects of 12.5 mm diameter and greater such as a finger. The optional AC adapter HHP-CM01 is protected against vertically falling water drops which may cause issues during a normal operation. The optional AC adapter HHP-BFH01 is protected against oblique falling water drops which may cause issues during a normal operation.
- Operation mode classification complies with IEC 60601-1.

Important information regarding Electromagnetic Compatibility (EMC)

HEM-7143T2-ESL manufactured by OMRON HEALTHCARE Co., Ltd. conforms to EN60601-1-2:2015 Electromagnetic Compatibility (EMC) standard.
Further documentation in accordance with this EMC standard is available at: www.omron-healthcare.com
Refer to the EMC information for HEM-7143T2-ESL on the website.

Correct Disposal of This Product (Waste Electrical & Electronic Equipment)

This marking shown on the product or its literature, indicates that it should not be disposed of, with other household waste, at the end of its working life. To prevent possible harm to the environment or human health from uncontrolled waste disposal, please separate this from other types of waste and recycle it responsibly to promote the sustainable reuse of material resources.

Household users should contact either the retailer where they purchased this product, or their local government office, for details of where and how they can return this item for environmentally safe recycling.

Business users should contact their supplier and refer to the terms and conditions of the purchase contract. This product should not be mixed with other commercial waste for disposal.

The **Bluetooth®** word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by OMRON HEALTHCARE Co., Ltd. is under license. Other trademarks and trade names are those of their respective owners. Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries and regions. App Store is a service mark of Apple Inc. Google Play and Google Play logo are trademarks of Google LLC.

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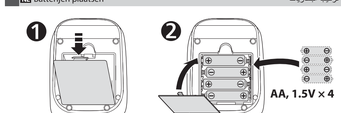
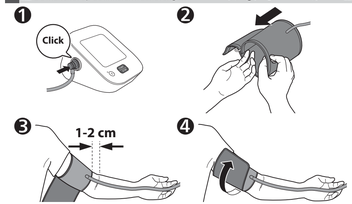
3 Tube side of the cuff should be 1 - 2 cm above the inside elbow.

- **FR** Le côté tuyau du brassard doit être positionné 1 à 2 cm au-dessus de l'intérieur du coude.
- **DE** Das Manschettenstück mit dem Schlauch muss 1 bis 2 cm oberhalb des Ellbogens liegen.
- **IT** Il lato del bracciale con il tubo deve trovarsi al di sopra del gomito interno, a una distanza di circa 1 o 2 cm.
- **ES** El lado del tubo del manguito deberá quedar 1 o 2 cm por encima de la parte interna del codo.
- **NL** De kant met de slang van de manchet moet 1 - 2 cm boven de binnenkant van de elleboog liggen.
- **TR** Kolluğun boru tarafı, dirsek içinin 1 ila 2 cm üstünde olmalıdır.
- **PL** Koniec mankietu z podłączonym przewodem powietrza powinien znajdować się 1–2 cm powyżej zgięcia łokcia.
- **PT** O lado do tubo da braçadeira deve estar 1 - 2 cm acima do interior do cotovelo.
- **SV** Slangsidan av manschetten ska vara 1–2 cm ovanför insidan av armbågen.
- **AR** يجب أن يرتفع جانب الخرطوم من الشريط الضاغط من ١ إلى ٢ سم عن المرفق.

4 Make sure that air tube is on the inside of your arm and wrap the cuff securely so it can no longer slip round.

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If taking measurements on the right arm, refer to:

- **FR** Pour la prise de mesures au bras droit, voir :
- **DE** Bei Messungen am rechten Arm siehe:
- **IT** Se la misurazione viene eseguita al braccio destro, fare riferimento a:
- **ES** Si va a realizar mediciones en el brazo derecho, consulte:
- **NL** Als u metingen aan de rechterarm uitvoert, raadpleeg dan:
- **TR** Sağ koldan ölçüm yapıyorsanız aşağıdakilere bakın:
- **PL** Pomiar ciśnienia na prawym ramieniu, patrz:
- **PT** Se fizer as medições no braço direito, consulte:
- **SV** Om du mäter på höger arm, se:
- **AR** في حال إجراء القياس على الذراع الأيمن، ارجع إلى:

Instruction Manual 1
☞ General Precautions

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1 Sit comfortably with your back and arm supported.

- **FR** S'asseoir confortablement de manière à ce que le dos et le bras soient bien soutenus.
- **DE** Bequem hinsetzen, mit dem Rücken anlehnen und den Arm auf eine Unterlage legen.
- **IT** Sedere comodamente con la schiena e il braccio ben sostenuti.
- **ES** Siéntese cómodamente de modo que la espalda y el brazo estén bien apoyados.
- **NL** Neem een comfortabele zitpositie in, met uw rug en arm ondersteund.
- **TR** Sırtınız ve kolunuz desteklenerek rahat bir şekilde oturun.
- **PL** Usiąść wygodnie, opierając plecy i ramię.
- **PT** Sente-se confortavelmente apoiando as costas e o braço.
- **SV** Sitt bekvämt med stöd för ryggen och armen.
- **AR** اجلس بشكل مريح بحيث يكون الظهر والذراع في وضع استناد.

2 Place the arm cuff at the same level as your heart.

- **FR** Le brassard doit se trouver au même niveau que votre cœur.
- **DE** Die Manschette auf Herzhöhe platzieren.
- **IT** Posizionare il bracciale allo stesso livello del cuore.
- **ES** Coloque el manguito al mismo nivel que el corazón.
- **NL** Plaats de armmanchet op hetzelfde niveau als uw hart.
- **TR** Kolluğu kalbinizle aynı düzeye getirin.
- **PL** Mankiet powinien znajdować się na wysokości serca.
- **PT** Coloque a braçadeira ao mesmo nível do coração.
- **SV** Armmanschetten ska vara på samma nivå som hjärtat.
- **AR** ضع الشريط الضاغط للذراع في مستوى القلب نفسه.

3 Keep feet flat, legs uncrossed, remain still and do not talk.

- **FR** Garder les pieds à plat, les jambes non croisées, ne pas bouger et ne pas parler.
- **DE** Füße flach auf den Boden stellen, die Beine nicht kreuzen, still halten und nicht sprechen.
- **IT** Tenere i piedi ben poggiati, non incrociare le gambe, rimanere fermi e non parlare.
- **ES** Mantenga los pies planos, las piernas sin cruzar, quédese quieto y no hable.
- **NL** Houd de voeten plat, benen niet gekruist, zit stil en praat niet.
- **TR** Ayaklarınızı düz, bacaklarınızı açık tutun, hareketsiz kalın ve konuşmayın.
- **PL** Stopy ułożyć płasko na podłodze, nie krzyżować nóg, siedzieć spokojnie, bez ruchu i nie rozmawiać.
- **PT** Mantenha os pés no chão, não cruze as pernas, mantenha-se imóvel e não fale.
- **SV** Håll fötterna plant mot golvet, benen får inte vara korsade, var stilla och prata inte.
- **AR** اجعل قدميك في وضع مسطح ورساقيك في وضع متباعد وحافظ على ثباتك ولا تتحدث.

8 Taking a Measurement

FR Réalisation d'une mesure
DE Vornehmen einer Messung
IT Misurazione
ES Obtención de una lectura
NL Een meting doen
TR Ölçüm Yapma
PL Wykonywanie pomiaru
PT Realização de uma medição
SV Göra en mätning
AR إجراء القياس

When the [START/STOP] button is pressed, the measurement is taken and saved automatically. Open the app to transfer the reading.


FR Après une pression sur le bouton [START/STOP], la mesure démarre et le résultat est enregistré automatiquement. Ouvrir l'application pour transférer le résultat.
DE Beim Drücken der Taste [START/STOP] erfolgt die Messung. Die Messwerte werden automatisch gespeichert. Öffnen Sie die App, um die Messwerte zu übertragen.
IT Quando si preme il pulsante [START/STOP], la misurazione viene eseguita e salvata automaticamente. Aprire la app per trasferire i risultati.
ES Al pulsar el botón [START/STOP], se obtiene la medición, que se guarda automáticamente. Abra la aplicación para transferir la lectura.
NL Wanneer op de knop [START/STOP] wordt gedrukt, wordt de meting automatisch uitgevoerd en opgeslagen. Open de app om de meting te versturen.
TR [START/STOP] düğmesine basıldığında ölçüm yapılır ve otomatik olarak kaydedilir. Ölçüm değerini aktarmak için uygulamayı açın.
PL Po naciśnięciu przycisku [START/STOP] następuje pomiar i jego automatyczny zapis. Otworzyć aplikację, aby przesłać odczyt.
PT Quando o botão [START/STOP] é premido, a medição é feita e gravada automaticamente. Abra a aplicação para transferir a leitura.
SV När knappen [START/STOP] trycks in tas mätningen och sparas automatiskt. Öppna appen för att överföra mätvärdena.
AR عند الضغط على زر [START/STOP]، يتم إجراء القياس وحفظه تلقائيًا. افتح التطبيق لنقل آخر بيانة.

If your systolic pressure is more than 210 mmHg:
After the arm cuff starts to inflate, press and hold the [START/STOP] button until the monitor inflates 30 to 40 mmHg higher than your expected systolic pressure.

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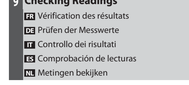❶ Appears if "SYS" is 135 mmHg or above and/or "DIA" is 85 mmHg* or above.

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❷ Appears when your body moves during a measurement. Remove the arm cuff, wait 2 - 3 minutes and try again.

FR S'affiche si vous bougez pendant une mesure. Retirer le brassard, attendre 2 à 3 minutes et essayer à nouveau.
DE Wird angezeigt, wenn Sie sich während der Messung bewegen. Nehmen Sie die Manschette ab, warten Sie 2–3 Minuten und versuchen Sie es erneut.
IT Viene visualizzato se l'utilizzatore si muove durante la misurazione. Rimuovere il bracciale, attendere 2-3 minuti e riprovare.
ES Aparece cuando se mueve el cuerpo durante una medición. Retire el manguito, espere unos 2 o 3 minutos e inténtelo de nuevo.
NL Verschijnt wanneer uw lichaam tijdens een meting beweegt. Verwijder de armmanchet, wacht 2-3 minuten en probeer het opnieuw.
TR Bir ölçüm esnasında vücudunuz hareket ettiğinde görünür. Kolluğu çıkarın, 2 - 3 dakika bekleyip tekrar deneyin.
PL Pojawia się w przypadku wykrycia ruchów ciała podczas pomiaru. Zdjąć mankiet, poczekać 2–3 minuty i spróbować ponownie.
PT Aparece quando o corpo se move durante uma medição. Retire a braçadeira, aguarde 2 - 3 minutos e tente de novo.
SV Visas när du rör dig under en mätning. Ta av manschetten, vänta 2–3 minuter och försök igen.
AR يظهر عند تحركك جسمك في أثناء عملية القياس. قم بإزالة الشريط الخاص بلف الذراع وانتظر من ٢ إلى ٣ دقائق وحاول مرة أخرى.

❸ Appears when an irregular rhythm** is detected during a measurement. If it appears repeatedly, OMRON recommends to consult your physician.

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***The high blood pressure definition is based on the 2018 ESH/ESC Guidelines.**

FR *La définition de l'hypertension est basée sur les recommandations ESH/ESC 2018.
DE *Die Definition für Bluthochdruck basiert auf den 2018 ESH/ESC Richtlinien.
IT *La definizione di alta pressione arteriosa si basa sulle linee guida 2018 di ESH/ESC.
ES *La definición de presión arterial alta se basa en las guías 2018 ESH/ESC.
NL *De definitie van hoge bloeddruk is gebaseerd op de 2018 ESH/ESC-richtlijnen uit 2018.
TR *Yüksek kan basıncı tanımı, 2018 ESH/ESC Kılavuzlarını temel alır.
PL *Definicja wysokiego ciśnienia krwi oparta na wytycznych ESH/ESC z 2018 r.
PT *A definição de tensão arterial alta baseia-se nas Orientações (Guidelines) de 2018 da ESH/ESC.
SV *Definitionen för högt blodtryck baseras på ESH/ESC-riktlinjerna från 2018.
AR *تعريف ضغط الدم المرتفع يستند إلى إرشادات الجمعية الأوروبية لفرط ضغط الدم(ESH)/الجمعية الأوروبية لطب القلب (ESC) لعام 2018.

****An irregular heart beat rhythm is defined as a rhythm that is 25 % less or 25 % more than the average rhythm detected during a measurement.**

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
Error messages or other problems? Refer to:

FR Messages d'erreur ou autres problèmes ? Voir :
DE Weitere Fehlermeldungen oder Probleme siehe:
IT Messaggi di errore o altri problemi? Fare riferimento a:
ES ¿Hay mensajes de error u otros problemas? Consulte:
NL Foutmeldingen of andere problemen? Raadpleeg:
TR Hata mesajları veya başka sorunlar mı var? Bkz:
PL Komunikaty o błędzie lub inny problem? Przejdź do:
PT Mensagens de erro ou outros problemas? Consulte:
SV Felmeddelanden eller andra problem? Se:
AR رسائل الخطأ أو المشكلات الأخرى؟ ارجع إلى:

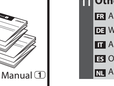Instruction Manual ①

11 Other Settings

FR Autres réglages
DE Weitere Einstellungen
IT Altre impostazioni
ES Otros ajustes
NL Andere instellingen
TR Diğer Ayarlar
PL Inne ustawienia
PT Outras definições
SV Övriga inställningar
AR الإعدادات الأخرى

11.1 Disabling/Enabling Bluetooth®

FR Désactivation/activation de la fonction Bluetooth®
DE Deaktivieren/Aktivieren der Bluetooth®-Funktion
IT Disattivazione/attivazione della funzione Bluetooth®
ES Activar/desactivar el Bluetooth®
NL Bluetooth® uitschakelen/inschakelen
TR Bluetooth® Özelliğini Devre Dışı Bırakma/Etkinleştirme
PL Wyłączanie/włączanie funkcji Bluetooth®
PT Desativação/ativação do Bluetooth®
SV Aktivera/avaktivera Bluetooth®
AR تمكين/تعطيل ميزة Bluetooth®

Bluetooth is enabled by default.

FR Bluetooth est activé par défaut.
DE Die Bluetooth-Funktion ist standardmäßig aktiviert.
IT La funzione Bluetooth è attiva per impostazione predefinita.
ES El Bluetooth está activado por defecto.
NL Bluetooth is standaard ingeschakeld.
TR Bluetooth varsayılan olarak etkindir.
PL Funkcja Bluetooth jest domyślnie włączona.
PT O Bluetooth está ativado por predefinição.
SV Bluetooth är aktiverat som standard.
AR يتم تمكين خاصية Bluetooth بشكل افتراضي.


11.2 Restoring to the Default Settings

FR Réinitialisation aux réglages par défaut
DE Wiederherstellen der Standardeinstellungen
IT Ripristino delle impostazioni predefinite
ES Restablecimiento a los ajustes de fábrica
NL De standaardinstellingen herstellen
TR Varsayılan Ayarları Geri Yükleme
PL Przywracanie ustawień domyślnych
PT Restauração das predefinições
SV Återställa till förvalda inställningar
AR استعادة الإعدادات الافتراضية


12 Optional Medical Accessories

FR Accessoires médicaux optionnels
DE Optionales medizinisches Zubehör
IT Accessori medicali opzionali
ES Accesorios médicos opcionales
NL Optionele medische accessoires
TR Opsiyonel Tıbbi Aksesuarlar
PL Opcjonalne akcesoria medyczne
PT Acessórios médicos opcionais
SV Valfria medicinska tillbehör
AR الملحقات الطبية الاختيارية

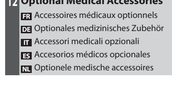Medium Cuff
HEM-CR24
22 - 32 cm

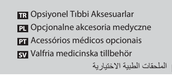Wide Range Soft Cuff
HEM-RML31
22 - 42 cm


Small Cuff
HEM-CS24
17 - 22 cm

Do not throw the air plug away. The air plug can be applicable to the optional cuff.

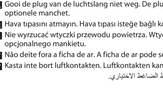FR Ne pas jeter la prise de gonflage. La prise de gonflage peut être utilisée pour le brassard en option.
DE Entsorgen Sie den Luftschlauchstecker nicht. Der Luftschlauchstecker wird für die optionale Manschette verwendet.
IT Non gettare via l'attacco del tubo dell'aria. L'attacco del tubo dell'aria può essere applicato al bracciale opzionale.
ES No tire el conector para tubo de aire. El conector para tubo de aire puede ser utilizado con el manguito opcional.
NL Gooi de plug van de luchtslang niet weg. De plug van de luchtslang kan worden gebruikt op de optionele manchet.
TR Hava tapası atmayın. Hava tapası isteğe bağlı kolluğa uygulanabilir.
PL Nie wyrzucać wtyczki przewodu powietrza. Wtyczkę przewodu powietrza można podłączyć do opcjonalnego mankietu.
PT Não deite fora a ficha de ar. A ficha de ar pode ser aplicável à braçadeira opcional.
SV Kasta inte bort luftkontakten. Luftkontakten kan användas på den valfria manschetten.
AR لا تلتخلص من مدادة الهواء. يمكن استعمال مدادة الهواء مع الشريط المناسط الاختياري.

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OMRON SANTÉ FRANCE SAS www.omron-healthcare.com/distributors