


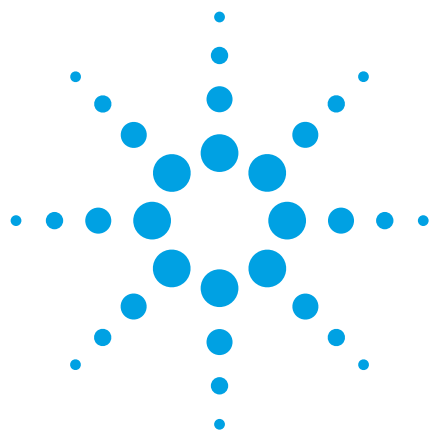

Analysis of Organic Acids using Agilent J&W FactorFour VF-WAXms GC Columns

Application Note

Author

Helena Jacobse
Agilent Technologies, Inc.


Conditions

Technique: GC-FID
Column: VF-WAXms 30 m x 0.53 mm df =
2.0 μ m (part number CP9216)
Temperature: Isotherm, 150 °C
Carrier Gas: Helium, 0.3 bar

Injector: 250 °C, split 75 mL/min
Detector: 275 °C, FID
Sample Size: 0.1 μ L
Sample: 0.1 % (cyclohexane)

Peak Identification:

1. *Acetic acid*
2. *Propionic acid*
3. *Iso-butyric acid*
4. *Butyric acid*
5. *Iso-valeric acid*

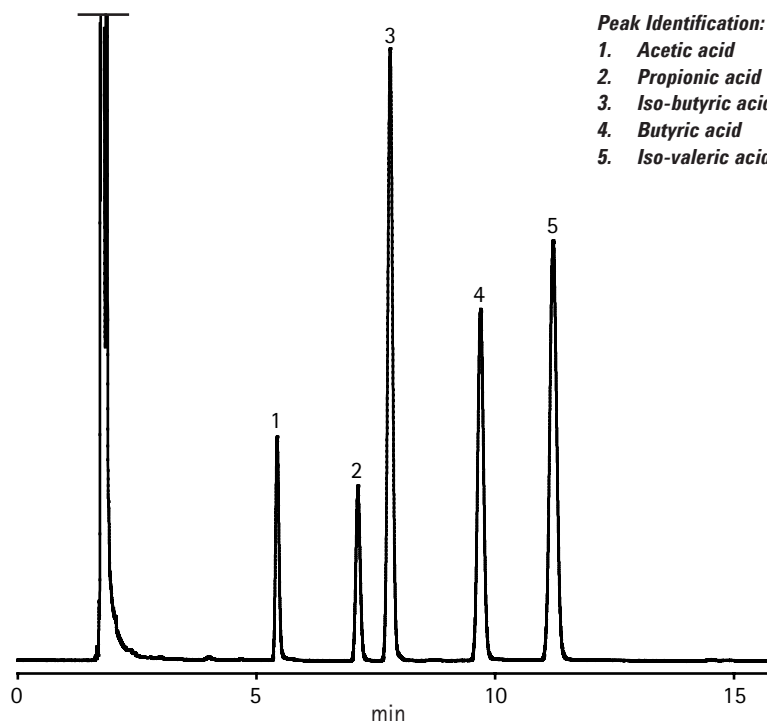

Analysis of organic acids using VF-WAXms GC columns

www.agilent.com/chem

This information is subject to change without notice.

© Agilent Technologies, Inc. 2010

Published in UK, October 08, 2010

SI-00937



Agilent Technologies