




Antihistamines

Excellent Resolution of Drug Mixtures

Application
Pharmaceutical
Robert Ricker

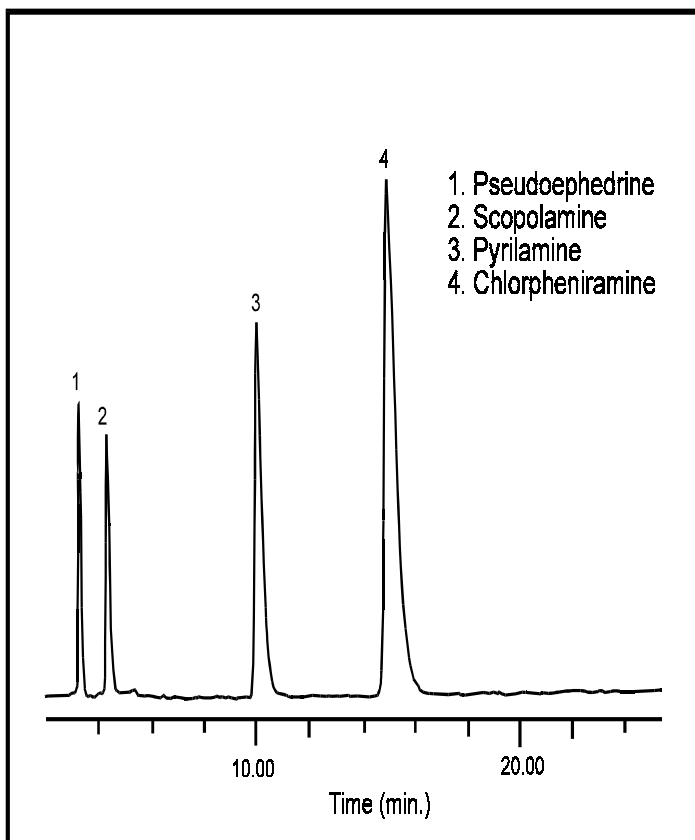

Conditions:
ZORBAX Rx/SB-C8 (4.6 x 150 mm) (Agilent P/N: 883967-901)/(883975-906)
10% ACN, 90% 50mM, KH_2PO_4 , pH 2.5
2.0 mL/min, 25°C, Detect. UV(254 nm)

Highlights

- ZORBAX StableBond C8 has a sterically protected bonded phase that permits stable operation at low mobile-phase pH.
- Basic compounds are separated with good resolution and peak shape on ZORBAX SB-C8 with simple mobile phases.



Agilent Technologies

*Robert Ricker is an application chemist
based at Agilent Technologies, Wilmington,
Delaware.*

For more information on our products and
services, visit our website at:
www.agilent.com/chem

Copyright© 2002 Agilent Technologies, Inc.
All Rights Reserved. Reproduction,
adaptation or translation without prior
written permission is prohibited, except as
allowed under the copyright laws.

Agilent shall not be liable for errors
contained herein or for incidental or
consequential damages in connection with
the furnishing, performance, or use of this
material.

Information, descriptions, and specifications
in this publication are subject to change
without notice.

Printed in the USA
April 25, 2002
5988-6402EN



Agilent Technologies