


Application 265-00

Agilent Monomer Analyzer

Impurities in 1,3-Butadiene

Technical Overview

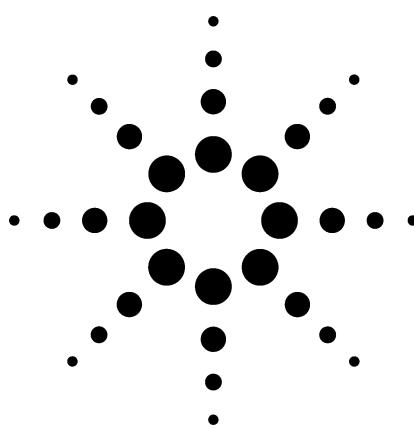


Application Highlights

- Dual Flame Ionization Detectors (FID) to detect the following components to a lower detection limit of 2 ppm, except for those peaks eluting off the tail of the matrix.

Acetylene	Isobutane
Propadiene	N-butane
t-2-Butene	1-butene
Isobutylene	1,2-butadiene
Methyl acetylene	Ethyl acetylene
Acrylonitrile	4-vinyl cyclohexene
Styrene	Ethyl benzene
1,3-pentadiene	Isoprene
1,4-pentadiene	C-1,3-pentadiene

- Analysis time: approximately 25 minutes

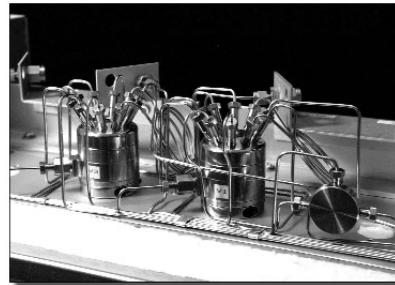


Optional Configurations

- Impurities in isoprene
- Impurities in high purity styrene monomer
- Analysis of vinyl chloride monomer
- Analysis of trace C4 olefins (10 ppb) in polymer grade propylene
- Analysis of 30 different trace oxygenates in polymer grade propylene by MSD
- ASTM D5769/ASTM D5599/ASTM D3606

For More Information


For more information on our products and services, visit our Web site at www.agilent.com/chem.

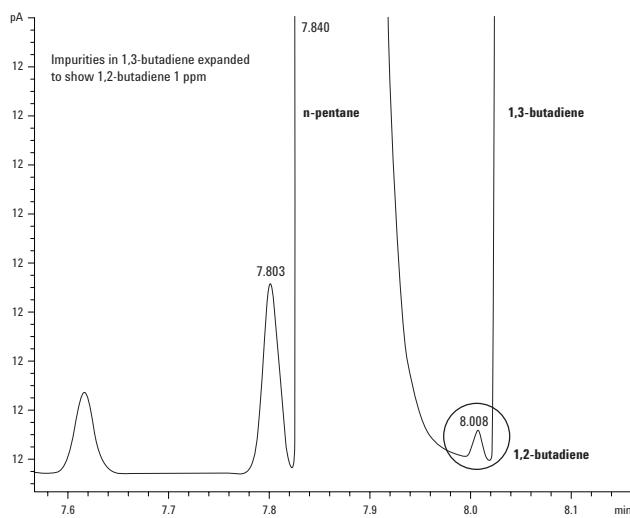

Agilent Technologies



WASSON-ECE
INSTRUMENTATION

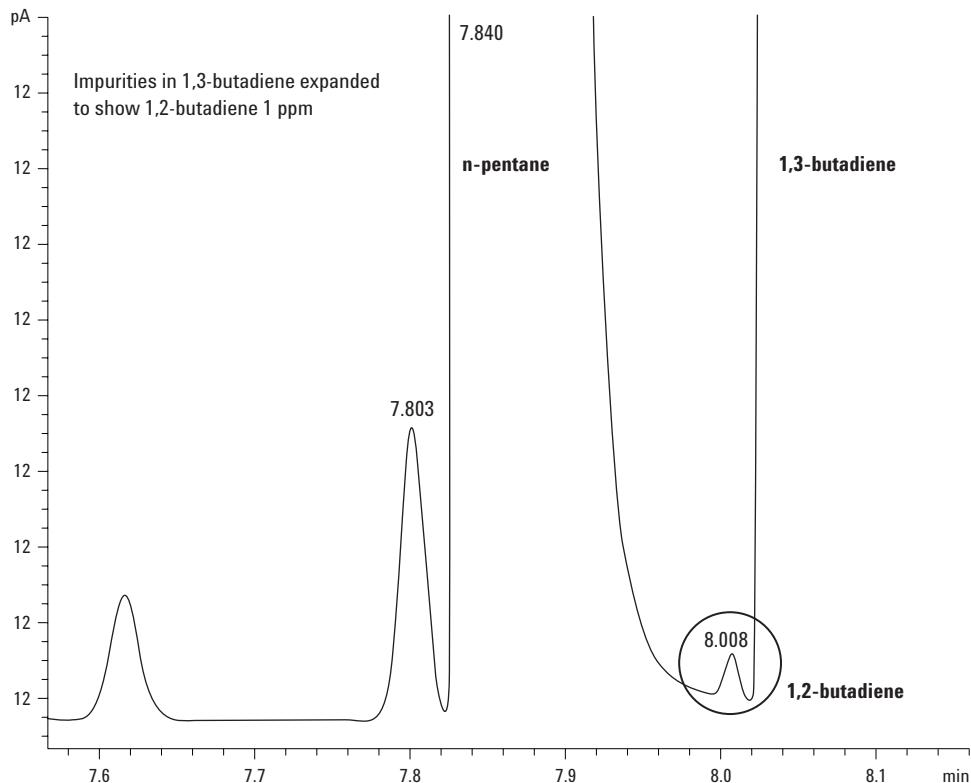

FID output from the Agilent monomer analyzer.

Agilent shall not be liable for errors contained herein or for incidental or consequential damages in connection with the furnishing, performance, or use of this material.

Information, descriptions, and specifications in this publication are subject to change without notice.

© Agilent Technologies, Inc. 2002

Printed in the USA
December 3, 2002
5988-6735EN



Agilent Technologies



WASSON-ECE
INSTRUMENTATION